

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [BEBIRD](#) /

› [BEBIRD Earsight Ultra Visual Ear Wax Removal Tool Camera - Instruction Manual](#)

BEBIRD Bebird Earsight Ultra

BEBIRD Earsight Ultra Visual Ear Wax Removal Tool Camera - Instruction Manual

Model: Bebird Earsight Ultra

INTRODUCTION

The BEBIRD Earsight Ultra is an advanced visual ear cleaning tool designed for safe and effective ear wax removal and ear canal observation. Featuring a 2K HD camera, 5G Wi-Fi connectivity, and a range of specialized tips, it provides a clear view of the ear canal directly on your smartphone, enhancing precision and safety during use. This manual provides essential information for setting up, operating, maintaining, and troubleshooting your device.

IMPORTANT SAFETY INFORMATION

WARNING: Always use this device with extreme caution. Improper use can lead to ear injury. Consult a medical professional if you experience pain, discomfort, or have pre-existing ear conditions before using this product. Keep out of reach of small children without adult supervision.

- Do not insert the device too deeply into the ear canal.
- Always ensure the tip is securely attached before use.
- Clean all accessories thoroughly before and after each use to prevent infection.
- If you experience any discomfort or pain, stop using the device immediately.
- The device is not a substitute for professional medical advice or treatment.

SETUP GUIDE

Follow these steps to prepare your BEBIRD Earsight Ultra for first use:

1. **Charge the Device:** Ensure the main device is fully charged using the provided USB-C charging cable. The charging case also holds a charge for convenience.
2. **Download the App:** Search for "Bebird" in your smartphone's app store (iOS or Android) and download the official application.
3. **Power On:** Press and hold the power button on the main device until the indicator light turns on.
4. **Connect to Wi-Fi:** Open your phone's Wi-Fi settings. Locate and connect to the Wi-Fi network named "bebird-".

XXXXXX" (where XXXXXX is a unique identifier for your device). This connection is local and does not require internet access.

5. **Launch the App:** Open the Bebird app. The live camera feed from the device should appear automatically.

6. **Attach a Tip:** Select the appropriate ear scoop or tweezer tip and securely screw it onto the camera end of the device.

Video: Bebird Earsight Ultra: Unboxing the Ultimate Ear Solution. This video demonstrates the unboxing and initial setup of the Bebird Earsight Ultra, showing the components and how to prepare the device for use.



Image: Step-by-step guide on how to connect the Bebird device to your smartphone's Wi-Fi network using the dedicated app.

OPERATING INSTRUCTIONS

Once set up, the BEBIRD Earsight Ultra offers various functionalities for comprehensive ear care:

1. Visual Ear Canal Examination

- With the device connected and the app open, gently insert the camera tip into your ear canal.
- Observe the live feed on your smartphone to navigate and examine your ear canal.
- The 2K HD camera and bright lighting provide clear, detailed images of your ear.

Video: Bebird Earsight Ultra 5G Ear Wax Removal Tool Camera. This video showcases the device's 5G capabilities and its use for visual ear wax removal, demonstrating the clarity of the camera feed.

2. Ear Wax Removal with Scoop Tips

- Choose from the five different silicone scoop tips based on the type of earwax (wet or dry) and your preference.
- Gently use the scoop tip to clean, scoop, grip, scrape, or pick earwax while monitoring the process on your phone.
- The safety stopper on some tips helps prevent over-insertion, ensuring safe use for all ages.



Image: Illustrations detailing various earwax conditions (dry, oily, flaky) and the recommended specialized cleaning tips for each, ensuring thorough and gentle removal.

3. Using Precision Tweezers

- Attach the patented detachable precision tweezers for removing stubborn or dry wax.
- The tweezers are designed to gently grasp wax without scraping the ear canal, ideal for sensitive ears.
- Maintain hand-eye coordination by watching the live view on your phone.



Image: A detailed view of the Bebird Earsight Ultra with its precision tweezers attachment, highlighting its design for gentle and effective earwax removal.

4. Visual Liquid Drops Delivery

- Utilize the visual dropper accessory for precise application of ear drops.
- The real-time visual guidance ensures accurate dosing and placement of medication.

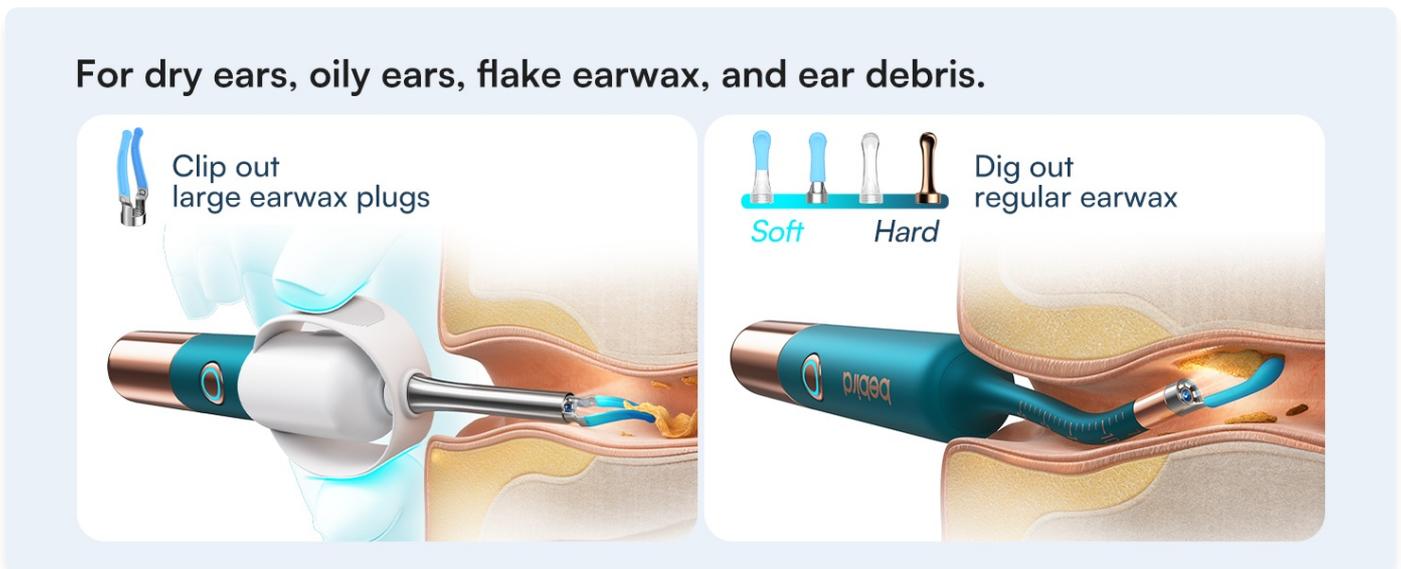


Image: A demonstration of the precise liquid delivery feature using the patented visual ear dropper, allowing for accurate application of ear drops with real-time visual feedback.

MAINTENANCE

Cleaning the Device and Tips

- After each use, detach the ear tip and clean it with a soft cloth or alcohol wipe.
- Wipe the camera lens gently with a lens cleaning cloth to ensure clear images for future use.
- Store all components in the provided storage case to keep them clean and organized.

Charging

- When the battery is low, connect the device to the USB-C charging cable.
- The indicator light will show charging status. A full charge provides approximately 80 minutes of working time.

TROUBLESHOOTING

- **Device not connecting to Wi-Fi:** Ensure the device is powered on and your phone's Wi-Fi is enabled. Try restarting both the device and your phone. Confirm you are connecting to the correct "bebird-XXXXXX" network.
- **Blurry image:** Clean the camera lens carefully with a soft, dry cloth. Ensure there is no debris on the lens or attached tip.
- **Short battery life:** Ensure the device is fully charged before use. If battery life significantly degrades over time, contact customer support.
- **App not displaying live feed:** Close and reopen the Bebird app. Verify the Wi-Fi connection is stable.

SPECIFICATIONS

Feature	Detail
Model Number	Bebird Earsight Ultra
Camera Resolution	2K HD (20MP)
Lens Diameter	3.6 mm
Connectivity	5G Wi-Fi Chip
Endoscope Technology	Snake-shaped flexible endoscope
Working Time	Approximately 80 minutes
Product Dimensions	5.59 x 0.51 x 5.12 inches
Item Weight	1.41 Pounds
Manufacturer	Heifeng Zhizao (Shenzhen) Technology Co., Ltd

WARRANTY AND SUPPORT

Your BEBIRD Earsight Ultra comes with a reliable **24-month warranty** from the date of purchase. For any questions, guidance, or support, our professional customer service team is available via online chat.

How to Contact Online Chat Support:

- **Mobile:**

1. Click "Your Orders"

2. Click "Purchase History" (locate your Bebird order)
3. Select "Chat" from "Get help from BEBIRD"

- **PC:**

1. Go to your Amazon Returns & Orders (locate your Bebird order)
2. Click "Get product support"
3. Click the "Chat" button

Service Hours: Monday-Friday, 11:00-24:00 (Pacific Standard Time, PST).

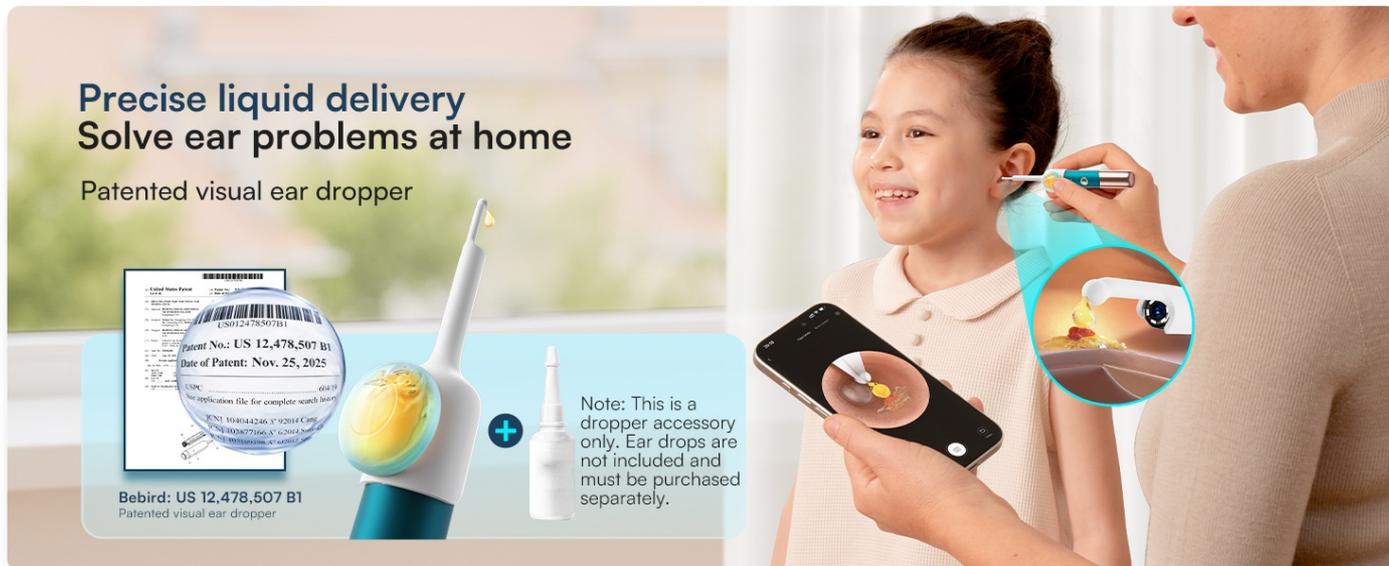


Image: Visual guide detailing the steps to access online chat support for Bebird products on both mobile devices and personal computers.

