

LEIPUT FIND A Pro

LEIPUT Ear Wax Removal Tool with 1080P Camera

Model: FIND A Pro

INTRODUCTION

This instruction manual provides essential information for the safe and effective use of your LEIPUT Ear Wax Removal Tool. This device integrates a 1080P HD camera with 6 LED lights and a 360° wide-angle lens, offering a clear view for inspecting and cleaning the ear canal, nose, mouth, and throat. Please read this manual thoroughly before first use and retain it for future reference.

SAFETY INFORMATION

- Always use the device gently. Avoid inserting the ear scoop too deeply into the ear canal to prevent injury.
- Do not immerse the entire otoscope body in water. The lens is waterproof and can be cleaned, but the main unit is not.
- Keep out of reach of children. This is not a toy.
- Consult a medical professional if you experience pain, discomfort, or have a pre-existing ear condition before using this device.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This device is for personal hygiene and inspection.

PACKAGE CONTENTS

Verify that all items are present in your package:

- Otoscope (Ear Camera) x 1
- Charging Cable (Type-C) x 1
- Silicone Ear Picks x 8 (various sizes/arcs)
- Storage Base x 1
- Ear Wax Removal 8-piece set (black metal tools) x 1

- Pimple Push Tool x 1 (for alternative use)
- User Manual x 1



Image: Assortment of silicone ear picks and metal ear wax removal tools included in the kit.

SETUP

1. Charge the Otoscope

Before first use, fully charge the otoscope. Connect the provided Type-C charging cable to the otoscope's charging port and a suitable USB power source. The built-in 350mAh rechargeable battery provides up to 90 minutes of continuous use on a full charge and up to 30 days of standby time.



Image: Close-up of the otoscope showing the Type-C charging port.

2. Download the App

The earwax removal tool is compatible with both Android and iOS devices. Search for the dedicated application in your device's app store (e.g., Apple App Store for iOS, Google Play Store for Android) and install it. Refer to the packaging or quick start guide for the specific app name.

3. Connect via WiFi

1. Power on the otoscope by pressing and holding the power button.
2. On your smartphone or tablet, go to your WiFi settings.
3. Locate and connect to the otoscope's offline WiFi network. This network will typically have a name related to 'LEIPUT' or 'Otoscope' followed by a series of numbers or letters. No internet connection is required for this

step.

4. Open the downloaded application. The app should automatically detect the connected otoscope and display a real-time image from the camera.



Image: The LEIPUT otoscope with its charging base, demonstrating compatibility with iOS (iPhone, iPad) and Android devices.

4. Attach an Ear Scoop

Select a silicone ear scoop (low, medium, or high arc) that is appropriate for your ear canal size and comfort. Firmly place the chosen silicone ear scoop onto the metal tip of the otoscope. Ensure it is securely tightened to prevent it from falling off during use.

Silicone Ear Scoop With 3 Arcs

Soft and elastic without scratching the ear canal

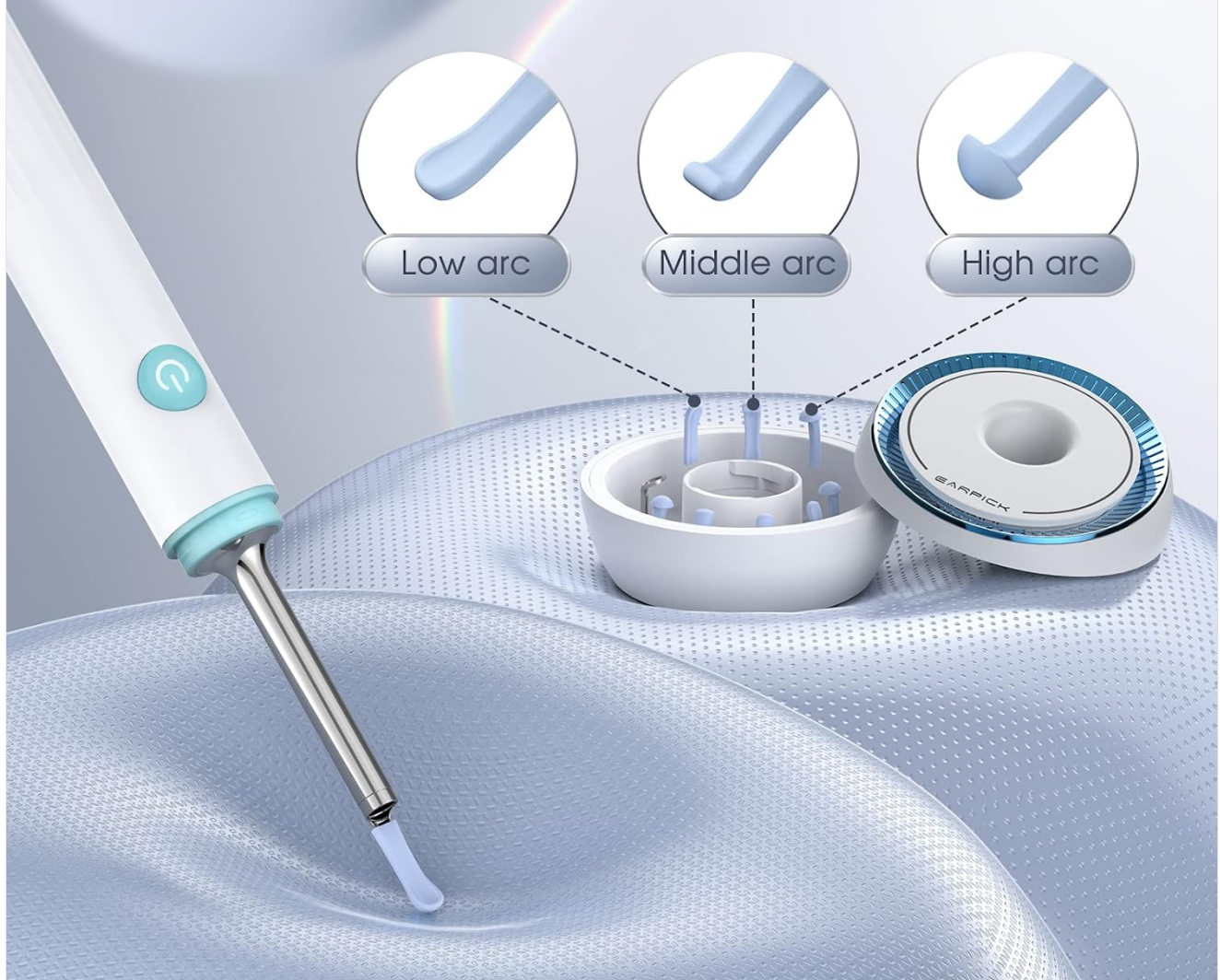


Image: Illustration of the three types of silicone ear scoops (low, middle, high arc) and how they attach to the otoscope.

OPERATING INSTRUCTIONS

1. Power On and View

Press the power button to turn on the otoscope. Open the app on your connected device. The 1080P HD camera with 6 LED lights will provide a clear, real-time view of the area you are examining.

1080P HD Resolution

Provide clear images and videos for your ear care



Ours **Others**



90 Minutes

6 LED Lights

Wifi Connection **360°** Wide Angle

Image: A visual comparison showing the superior clarity of 1080P HD resolution from the LEIPUT otoscope on a smartphone screen, highlighting the ear canal view.

2. Ear Wax Removal

Carefully insert the otoscope with the attached silicone ear scoop into the ear canal while observing the live feed on your device. Gently maneuver the scoop to dislodge and remove earwax. Always proceed with caution and avoid any forceful movements.

3. Multiple Applications

Beyond ear cleaning, the otoscope can be used for daily check-ups of the nose, mouth, and throat. The camera provides a clear view for general inspection of these areas.

Multiple Applications



Image: Demonstrates various applications of the otoscope, including ear inspection, nose inspection, and oral cavity examination.

4. Power Off

After use, press and hold the power button to turn off the otoscope.

MAINTENANCE

Cleaning the Lens and Ear Scoops

The otoscope lens is waterproof. After each use, remove the silicone ear scoop and clean both the scoop and the lens with water or wipe them with alcohol. Ensure the lens is dry before storage. Do not submerge the entire otoscope body.

Lens waterproof

Otoscope body is not waterproof



Image: Demonstrates that the otoscope lens is waterproof and can be cleaned with water, while emphasizing that the main body should not be submerged.

Storage

Store the otoscope and its accessories in the provided storage base or a clean, dry place when not in use. This helps protect the device from dust and damage.

TROUBLESHOOTING

- **No Image on App:** Ensure the otoscope is powered on and your device is connected to its specific WiFi network (not your home internet). Restart both the otoscope and the app if necessary.
- **Poor Image Quality:** Check if the lens is clean. Ensure there is sufficient light (the otoscope has built-in LED lights).
- **Device Not Charging:** Verify the Type-C cable is securely connected to both the otoscope and the power source. Try a different charging cable or power adapter.
- **Ear Scoop Falls Off:** Ensure the silicone ear scoop is firmly pressed onto the metal tip and tightened before use.

- **App Crashing/Freezing:** Close and reopen the app. Ensure your device's operating system is up to date.

SPECIFICATIONS

Model Number	FIND A Pro
Camera Resolution	1080P HD
LED Lights	6
Field of View	360° Wide-Angle
Connectivity	WiFi (Offline Network)
Battery Capacity	350mAh Lithium Polymer
Battery Life	Up to 90 minutes (continuous use)
Standby Time	Up to 30 days
Lens Waterproof Rating	Yes (main body not waterproof)
Compatibility	iOS, Android
Package Dimensions	6.85 x 4.25 x 0.87 inches
Item Weight	4.66 ounces

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided on the product packaging or the official LEIPUT website. Keep your purchase receipt as proof of purchase.

© 2024 LEIPUT. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question