

[manuals.plus](#) /

› [HOMCOM](#) /

› [HOMCOM EVA Padded Shower Chair with Arms and Back \(Model 713-116V00GY\) - Instruction Manual](#)

**HOMCOM 713-116V00GY**

# HOMCOM EVA Padded Shower Chair with Arms and Back (Model 713-116V00GY) - Instruction Manual



## 1. INTRODUCTION

This manual provides essential information for the safe and effective use of your HOMCOM EVA Padded Shower Chair. Please read these instructions thoroughly before assembly and operation. Retain this manual for future reference.



*Image 1.1: The HOMCOM EVA Padded Shower Chair, featuring a blue padded seat and backrest, with adjustable aluminum legs and armrests.*

## 2. SAFETY INFORMATION

- **Weight Capacity:** Do not exceed the maximum weight capacity of 299 lbs (135 kg).

- **Stability:** Ensure all four legs are evenly adjusted and the suction foot pads are firmly secured to a clean, dry, non-porous surface before use.
- **Inspection:** Regularly inspect the chair for loose parts, damage, or wear. Do not use if any components are compromised.
- **Assembly:** Follow all assembly instructions carefully. Ensure all connections are secure before use.
- **Wet Surfaces:** Exercise caution on wet surfaces. While the chair is designed with anti-slip features, wet environments can still pose a risk.
- **Children:** Keep children away from the shower chair during assembly and unsupervised use.

### 3. PACKAGE CONTENTS

---

Verify that all components are present before beginning assembly:

- Padded Seat Panel with Backrest
- Aluminum Leg Frames (4 pieces)
- Armrests (2 pieces)
- Suction Foot Pads (4 pieces, pre-attached to legs)
- Assembly Hardware (e.g., screws, knobs for tool-free assembly)

### 4. ASSEMBLY INSTRUCTIONS

---

The HOMCOM Shower Chair is designed for tool-free assembly. Please follow these steps:

1. **Unpack Components:** Carefully remove all parts from the packaging and lay them out on a clean, flat surface.
2. **Attach Legs to Seat:** Insert the top ends of the aluminum leg frames into the designated slots on the underside of the seat panel. Ensure the spring pins snap securely into place.
3. **Secure Legs:** Use the provided hand-tightened knobs to further secure the leg frames to the seat. Turn clockwise until snug.
4. **Attach Armrests:** Slide the armrests into the side openings of the seat frame until they click into position.
5. **Adjust Height:** Refer to Section 5 for detailed instructions on adjusting the chair's height.
6. **Final Check:** Before first use, gently shake the chair to ensure all connections are firm and stable.

# TOOL FREE ASSEMBLY

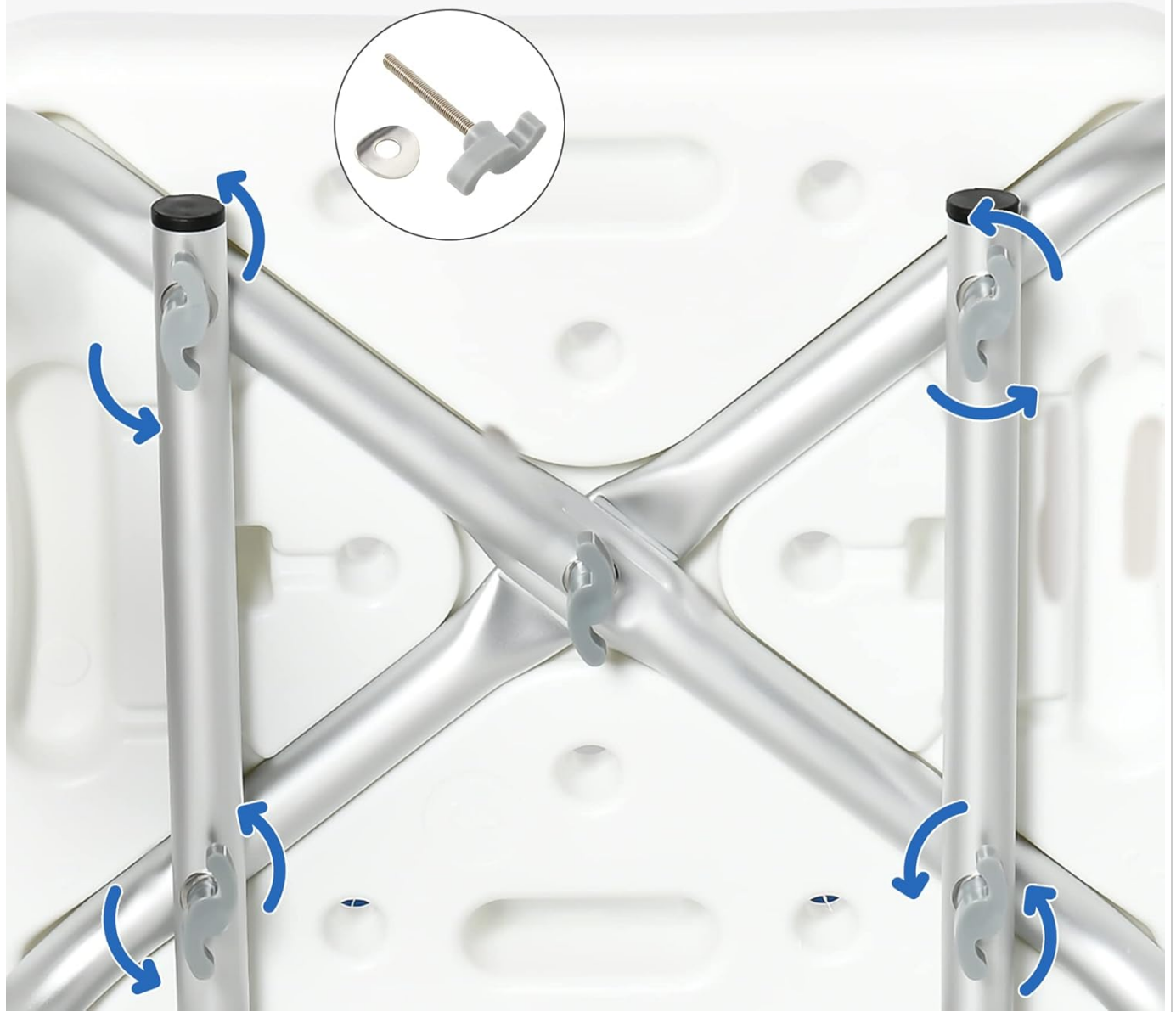


Image 4.1: Diagram illustrating the tool-free assembly process, showing how the legs attach to the seat with secure clips and knobs.

## 5. OPERATING INSTRUCTIONS

### 5.1. Height Adjustment

The shower chair features 6 levels of height adjustment to suit individual user needs. To adjust the height:

1. Ensure the chair is empty and on a stable surface.
2. Locate the spring pins on each leg.
3. Press the spring pin inward and slide the leg extension up or down to the desired height setting.
4. Ensure the spring pin fully protrudes through the corresponding hole on the leg extension for all four legs.
5. Verify that all four legs are set to the same height to maintain stability.



Image 5.1: Close-up view of the adjustable leg with numbered height settings and the spring pin mechanism.

## 5.2. Using the Chair

- Position the chair securely in the shower or bath area. Ensure the suction foot pads are firmly attached to the floor.
- Use the armrests for support when sitting down or standing up.
- The EVA padded seat and backrest provide comfort and help prevent slipping.
- Drainage holes in the seat prevent water accumulation.



Image 5.2: The shower chair positioned in a bathroom next to a bathtub, demonstrating its intended use.

Video 5.3: An official product video demonstrating the features and use of the HOMCOM EVA Padded Shower Chair, including its adjustable height and padded surfaces.

## 6. MAINTENANCE AND CARE

- **Cleaning:** Clean the chair regularly with mild soap and water. Rinse thoroughly and wipe dry to prevent water spots and mineral buildup.
- **Do Not Use Abrasives:** Avoid using abrasive cleaners, scouring pads, or harsh chemicals, as these can damage the surfaces.
- **Storage:** Store the chair in a dry place when not in use for extended periods.
- **Inspection:** Periodically check all screws, pins, and connections to ensure they remain tight and secure.

## 7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Chair feels unstable or wobbly.	Legs not adjusted to the same height; suction foot pads not properly secured; loose connections.	Ensure all legs are at the same height. Press down firmly on the seat to engage suction cups. Check and tighten all assembly connections.
Difficulty adjusting leg height.	Spring pin stuck or misaligned.	Ensure the spring pin is fully depressed. Wiggle the leg gently while adjusting. Clean any debris from the pin mechanism.
Water accumulates on the seat.	Drainage holes are blocked.	Clear any debris from the drainage holes on the seat.

## 8. SPECIFICATIONS

Feature	Detail
Model Number	713-116V00GY
Brand	HOMCOM
Material	Aluminum frame, HDPE panel, EVA padding
Color	Blue
Overall Dimensions	21"D x 22.25"W x 29.25"H (53.34 cm D x 56.52 cm W x 74.3 cm H)
Adjustable Height Range	15.75" - 20.75" (40 cm - 52.7 cm) (Seat Height)
Seat Dimensions	19.5"W x 13"D (49.5 cm W x 33 cm D)
Weight Capacity	299 lbs (135 kg)
Item Weight	9.02 lbs (4.09 kg)





Image 8.1: Detailed diagram showing the key dimensions of the shower chair, including overall height, width, and seat measurements.

## 9. WARRANTY AND SUPPORT

For warranty information or product support, please contact the manufacturer, HOMCOM (Aosom LLC), directly through their official website or the retailer where the product was purchased. Please have your model number (713-116V00GY) and proof of purchase available.

**Manufacturer:** Aosom LLC

**UPC:** 196393349199