

SKINACT 164242

SKINACT Bellage Electric Treatment Table (Chair) User Manual

Model: 164242

INTRODUCTION

This manual provides essential information for the safe and effective operation and maintenance of your SKINACT Bellage Electric Treatment Table. Please read this manual thoroughly before initial use and keep it for future reference.

The SKINACT Bellage Electric Treatment Table is a versatile, electrically adjustable chair designed for professional use in salons, spas, and medical offices. It features four independent motors for precise positioning of the backrest, leg rest, seat tilt, and overall height, ensuring optimal comfort for clients and ergonomic convenience for practitioners.



Figure 1: SKINACT Bellage Electric Treatment Table (Beige)

SAFETY INFORMATION

Adherence to these safety guidelines is crucial for preventing injury and damage to the equipment.

- Do not use near water or in excessively damp environments.
- Keep out of reach of children.
- Avoid using harsh chemicals or abrasive cleaners on any part of the table.

- Always unplug the device from the power outlet before performing any adjustments, cleaning, or maintenance.
- Do not exceed the maximum weight limit of 300 lbs (136 kg).
- Ensure all electrical connections are secure and undamaged before operation.
- Do not operate the chair if any part is damaged or malfunctioning. Contact customer support immediately.

SETUP

The SKINACT Bellage Electric Treatment Table typically arrives pre-assembled or requires minimal assembly. Follow these general steps for initial setup:

1. Carefully remove all packaging materials. Inspect the table for any signs of shipping damage.
2. Place the table on a stable, level surface. Ensure adequate space around the table for full range of motion and safe operation.
3. Connect the power cord to a grounded electrical outlet (100-120 volts, 50-60 hertz).
4. Connect the wired remote control and any optional foot controls to their respective ports on the table's base.
5. Adjust the leveling feet on the base to ensure the table is perfectly stable and does not wobble.
6. Attach the desired headrest (with or without face hole) and armrests as needed.



Figure 2: Key Features and Electric Adjustments

OPERATING INSTRUCTIONS

The SKINACT Bellage Electric Treatment Table offers multiple points of adjustment for optimal client comfort and practitioner convenience. All adjustments are controlled via the wired remote control or optional foot controls.

Using the Remote Control

The wired remote control features intuitive buttons for each adjustable section of the chair. Familiarize yourself with the button layout:

- **Height Adjustment:** Use the up/down arrows to raise or lower the entire table.
- **Backrest Angle:** Adjust the backrest from 0° to 80° incline.
- **Leg Rest Angle:** Adjust the leg rest from 0° to 85° decline.
- **Seat Tilt Angle:** Adjust the bed tilt angle from 0° to 10°.



Figure 3: Wired Remote Control

Adjusting Table Sections

- **Backrest Adjustment:** Use the corresponding buttons on the remote to raise or lower the backrest.



Figure 4: Headrest without Face Hole

- **Leg Rest Adjustment:** Use the corresponding buttons on the remote to adjust the leg rest angle. The leg rest also features an extension for added length.



Figure 5: Leg Rest Adjustment

- **Height Adjustment:** Use the remote to raise or lower the entire table to your desired working height.



Figure 6: Height Adjustment

- **Armrest and Headrest Removal:** Armrests and headrests are removable for versatility.



Figure 7: Armrest Removal



Figure 8: Headrest Removal

MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your treatment table.

- **Cleaning:** Wipe the surface with a dry cloth. For deeper cleaning, use a mild, non-abrasive cleaner specifically designed for PU upholstery. Avoid harsh chemicals.
- **Inspection:** Regularly inspect all moving parts, electrical cords, and connections for wear or damage.
- **Lubrication:** Moving parts may require periodic lubrication. Consult a qualified technician if you notice any stiffness or unusual noises during operation.
- **Storage:** When not in use for extended periods, store the table in a clean, dry environment, away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

If you encounter issues with your SKINACT Bellage Electric Treatment Table, try these basic troubleshooting steps before contacting support:

- **Table Not Responding:**
 - Check if the power cord is securely plugged into both the table and the wall outlet.
 - Ensure the wall outlet is functional by testing it with another device.
 - Verify that the remote control and foot controls are properly connected.
- **Unusual Noises During Operation:**
 - Ensure no objects are obstructing the moving parts of the table.
 - If noises persist, contact customer support; do not attempt to disassemble the motors.
- **Intermittent Functionality:**
 - Check all cable connections for looseness or damage.
 - Avoid overloading the table beyond its specified weight capacity.

For issues not resolved by these steps, please refer to the Warranty and Support section.

SPECIFICATIONS

Feature	Detail
Brand	SKINACT
Model Number	164242
Item Dimensions (L x W x H)	68 x 28 x 32 inches
Length (with headrest)	76.8 inches
Height Range	25.5 - 31.5 inches
Width (with armrests)	35.4 inches
Backrest Angle	0° to 80° incline
Bed Tilt Angle	0° to 10°
Leg Rest Angle	0° to 85° decline
Base Dimensions	36" L x 23" W

Power Requirements	100-120 volts, 50-60 hertz, 200 watts, 2 amp
Material	Metal frame, plastic cover, foam, PU upholstery
Item Weight	200 pounds
Maximum Capacity	350 lbs
Color	Beige
Country of Origin	China

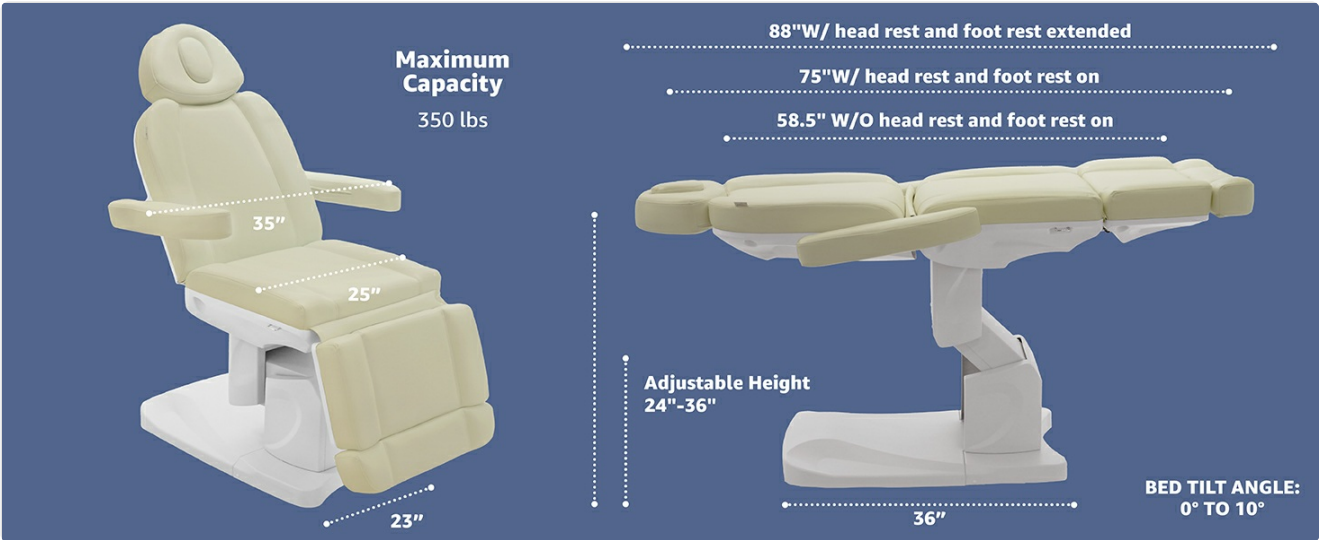


Figure 9: Detailed Dimensions and Capacity

WARRANTY AND SUPPORT

Your SKINACT Bellage Electric Treatment Table comes with a**ONE YEAR WARRANTY** covering manufacturing defects and malfunctions under normal use.

For warranty claims, technical assistance, or any questions regarding your product, please contact SKINACT Customer Support. You can find contact information on the official SKINACT website or through your purchase platform.

Additional protection plans may be available for purchase to extend coverage beyond the standard warranty period. Please refer to your purchase details for more information on available protection plans.



Toll-free: 1.800.418.0630	www.gilbert.com www.gasandequipment.com	International: 818.824.6640
---------------------------	--	-----------------------------

Detailed overview and specifications for the SkinAct Slate Metal Cart, a durable all-metal utility trolley designed for medical, spa, and salon environments. Features include ample storage, smooth-rolling casters, and a sleek design.