

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [HORIGEN](#) /

› [HORIGEN USB Rechargeable Electric Nasal Aspirator User Manual](#)

HORIGEN Electric Nasal Aspirator

HORIGEN USB Rechargeable Electric Nasal Aspirator User Manual

Model: Electric Nasal Aspirator

INTRODUCTION

Thank you for choosing the HORIGEN Electric Nasal Aspirator. This device is designed to safely and effectively clear nasal passages for newborns, infants, and toddlers. Please read this manual thoroughly before first use to ensure proper operation, maintenance, and safety.

SAFETY INFORMATION

- This product is intended for nasal aspiration only. Do not use for other purposes.
- Keep out of reach of children when not in use. Adult supervision is required during use.
- Ensure all parts are correctly assembled before use.
- Do not immerse the main unit in water. The main unit is not waterproof, only the collection cup and silicone tips are washable.
- Do not use if any part of the device is damaged.
- Use only the provided USB charging cable.
- Avoid using harsh chemicals or abrasive cleaners on any part of the device.
- If your child experiences discomfort or bleeding, discontinue use immediately and consult a medical professional.
- The device is made from BPA-free and food-grade materials.

PACKAGE CONTENTS

The package includes:

- 1 x Electric Nasal Aspirator Main Unit
- 3 x Silicone Suction Tips (various sizes/shapes)
- 1 x USB Charging Cable

- User Manual (this document)

PRODUCT OVERVIEW

Familiarize yourself with the components of your HORIZEN Electric Nasal Aspirator.



Figure 1: Main unit with labeled controls and display. This image shows the LED screen, power display, light button, music button, and power button, along with indicators for the three suction levels.

Protection Systems for Prevent Backflow

Easier to Suck the Mucus, Booger and Snot



Figure 2: Disassembled view of the suction and collection components, including the silicone tip, suction bin, storage bin, and anti-reflux cover. These parts are designed to prevent backflow and are easy to clean.

Soft and Food Grade Silicone



Figure 3: Selection of soft, food-grade silicone tips. Different tip designs are provided for various needs, such as sucking deepest, deep, or shallow nasal mucus, and for picking out hard boogers.

SETUP

1. Charging the Device

1. Before first use, fully charge the nasal aspirator.
2. Connect the USB charging cable to the charging port on the device and to a standard USB power adapter (DC 5V 1A, not included).
3. The power display on the LED screen will indicate charging status. A full charge typically takes approximately 2-3 hours.
4. Disconnect the charger once fully charged.

2. Assembly

1. Ensure all parts are clean before assembly (refer to "Maintenance and Cleaning" section).
2. Attach the anti-reflux cover to the main unit.
3. Place the storage bin onto the anti-reflux cover.
4. Place the suction bin onto the storage bin.

5. Select the appropriate silicone tip for your child's needs and securely attach it to the suction bin.

OPERATING INSTRUCTIONS

1. Power On/Off

- To power on, press and hold the **Power button** (marked with a circle and vertical line) for 2 seconds.
- To power off, press and hold the **Power button** for 2 seconds.

2. Adjusting Suction Levels

The aspirator features 3 suction levels to accommodate different needs.

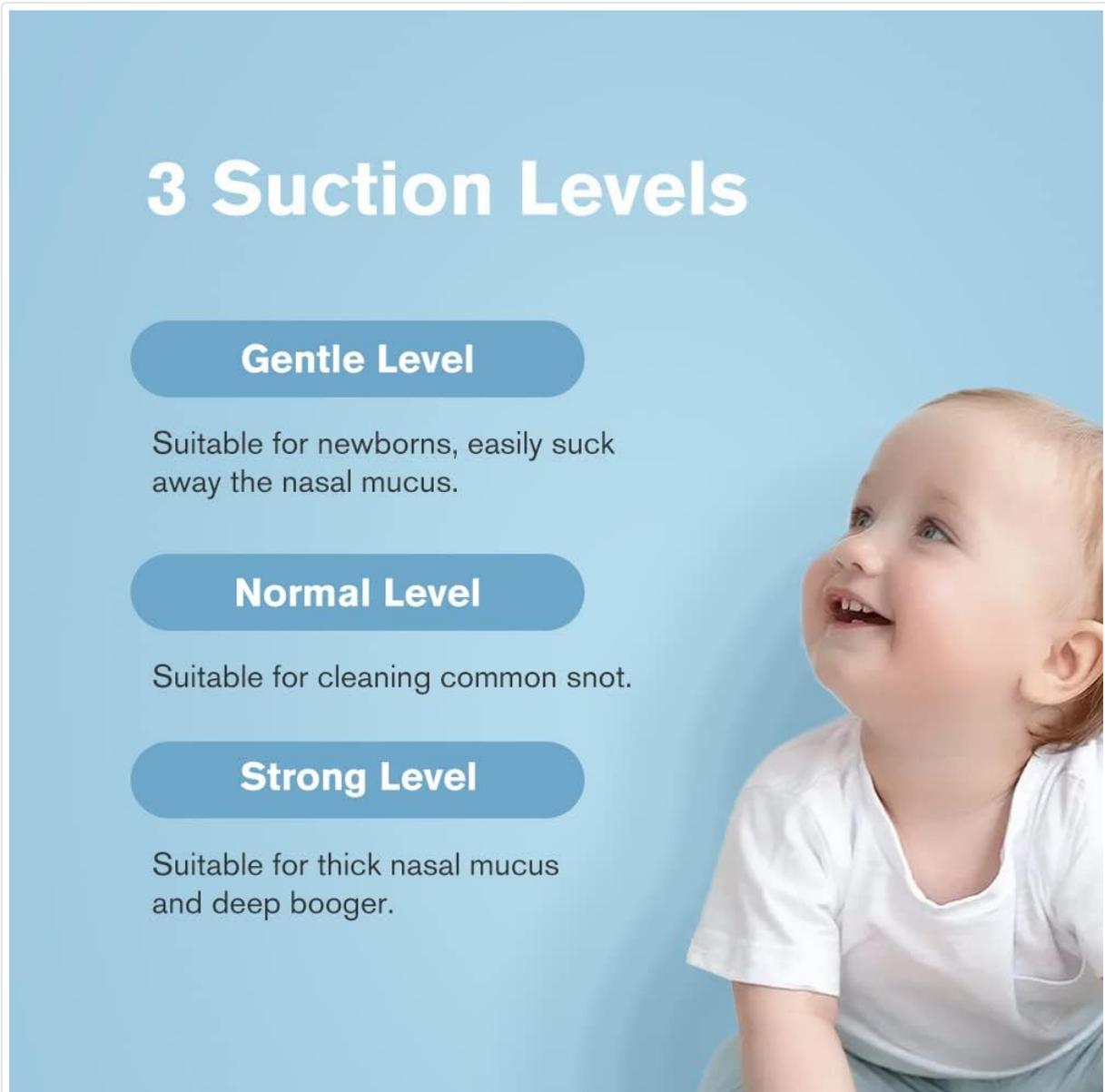


Figure 4: Explanation of the three suction levels. SP-1 (Gentle) is suitable for newborns, SP-2 (Normal) for common snot, and SP-3 (Strong) for thick nasal mucus.

- Press the **Mode button** (marked 'M') to cycle through the suction levels: SP-1 (Gentle), SP-2 (Normal), SP-3 (Strong).
- The current suction level will be displayed on the LED screen.
- Choose the appropriate level based on the consistency of the mucus and your child's comfort.

3. Music and Light Function

The aspirator includes a music and light function to help relax your baby during use.



Figure 5: The device features a light and music function to help soothe the baby. It offers colorful light and 10 gentle music tracks.

- Press the **Music button** (marked with a music note) to play gentle melodies. Press again to cycle through or turn off.
- Press the **Light button** (marked with a light bulb) to activate the colorful light. Press again to turn off.

4. Usage Steps

1. Ensure your baby is in a comfortable position, preferably lying down or slightly reclined.
2. Turn on the aspirator and select the desired suction level.
3. Gently insert the silicone tip into one of your baby's nostrils, ensuring a snug but not tight fit.
4. Allow the aspirator to clear the mucus. Move the tip slightly if needed to reach different areas.
5. Repeat for the other nostril.
6. Turn off the device after use.
7. Clean the removable parts immediately after each use.

MAINTENANCE AND CLEANING

Proper cleaning ensures hygiene and extends the life of your aspirator. All removable parts are waterproof and easy to clean.



Figure 6: The image demonstrates the ease of removing and cleaning the suction nozzle and collector parts under running water.

1. Turn off the device and disconnect it from any power source.
2. Carefully detach the silicone tip, suction bin, storage bin, and anti-reflux cover from the main unit.
3. Rinse these removable parts thoroughly under warm running water. You may use mild dish soap if necessary.
4. Use a small brush (not included) to clean any stubborn residue inside the bins or tips.
5. Ensure all parts are completely dry before reassembling or storing. Air dry or use a clean cloth.
6. Wipe the main unit with a soft, damp cloth. **Do not immerse the main unit in water.**
7. Store the aspirator in a clean, dry place.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully.
Device does not turn on.	Power button not pressed correctly.	Press and hold the Power button for 2 seconds.
Weak suction.	Silicone tip not properly attached or clogged.	Ensure the silicone tip is securely attached. Clean the tip and collection bins.
Weak suction.	Low battery.	Charge the device.
Weak suction.	Incorrect suction level selected.	Increase the suction level using the Mode button.
Device makes unusual noise.	Parts not assembled correctly.	Disassemble and reassemble all parts, ensuring they are securely fitted.

SPECIFICATIONS

- **Model:** Electric Nasal Aspirator
- **Brand:** HORIGEN
- **Dimensions:** 188 x 40 x 75 mm (Product), 18.49 x 11.81 x 4.7 cm (Package)
- **Weight:** 240 grams (Package)
- **Material:** Plastic and Silicone (BPA-free, food-grade)
- **Operating Voltage:** 3.7 V
- **Charging Voltage:** DC 5 V 1 A
- **Maximum Power:** 1.5 W
- **Battery Capacity:** 700 mAh (Lithium-ion)
- **Suction Levels:** 3 (Gentle, Normal, Strong)
- **Water Resistance:** Removable parts are waterproof for cleaning. Main unit is not waterproof.
- **Features:** Music function (10 melodies), Colorful light function

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the retailer where you purchased the product or visit the official HORIGEN website. Keep your purchase receipt as proof of purchase.