

Omron 7143T1

Omron HEM 7143T1 Digital Bluetooth Blood Pressure Monitor User Manual

Model: HEM 7143T1

INTRODUCTION

This user manual provides essential information for the safe and effective operation of your Omron HEM 7143T1 Digital Bluetooth Blood Pressure Monitor. Please read this manual thoroughly before using the device and keep it for future reference. This monitor is designed for home use to measure blood pressure and pulse rate.



Figure 1: Omron HEM 7143T1 Digital Bluetooth Blood Pressure Monitor and arm cuff.

IMPORTANT SAFETY INFORMATION

- Consult your physician before using this device if you have any medical conditions.
- Do not self-diagnose or self-treat based on measurement results. Always consult a healthcare professional.
- Keep the device out of reach of infants, toddlers, and children.
- Use only the authorized cuff for this monitor.
- Do not use the device in a moving vehicle or airplane.
- Avoid using the device in areas with strong electromagnetic fields.

PRODUCT COMPONENTS

The Omron HEM 7143T1 Blood Pressure Monitor package typically includes the following items:

- Main Unit (Blood Pressure Monitor)
- Arm Cuff (Standard size, suitable for adult arm circumference)

- Instruction Manual (This document)
- Storage Case (May vary by region)
- Batteries (Typically 4 x AA batteries)

ADVANTAGES & FEATURES



Bluetooth
Connectivity



Hypertension
Indicator



Irregular
Heart Beat
Detection



Body
Movement
Detection



Cuff
Wrapping
Guide



IntelliSense
Technology



30 Memory



Pulse



TECHNOLOGY
BY OMRON HEALTHCARE
IN JAPAN



Figure 2: Key features of the Omron HEM 7143T1 monitor, including Bluetooth Connectivity, Hypertension Indicator, Irregular Heartbeat Detection, Body Movement Detection, Cuff Wrapping Guide, IntelliSense Technology, 30 Memory, and Pulse measurement.

SETUP GUIDE

1. Installing Batteries

1. Open the battery compartment cover on the back of the monitor.
2. Insert four "AA" alkaline batteries according to the polarity indications (+/-) inside the compartment.
3. Close the battery compartment cover securely.

Note: If the battery indicator appears on the display, replace all four batteries with new ones. Do not mix old and new batteries or different types of batteries.

2. Connecting the Air Plug

1. Locate the air jack on the side of the main unit.

2. Firmly insert the air plug from the arm cuff into the air jack until it clicks into place.

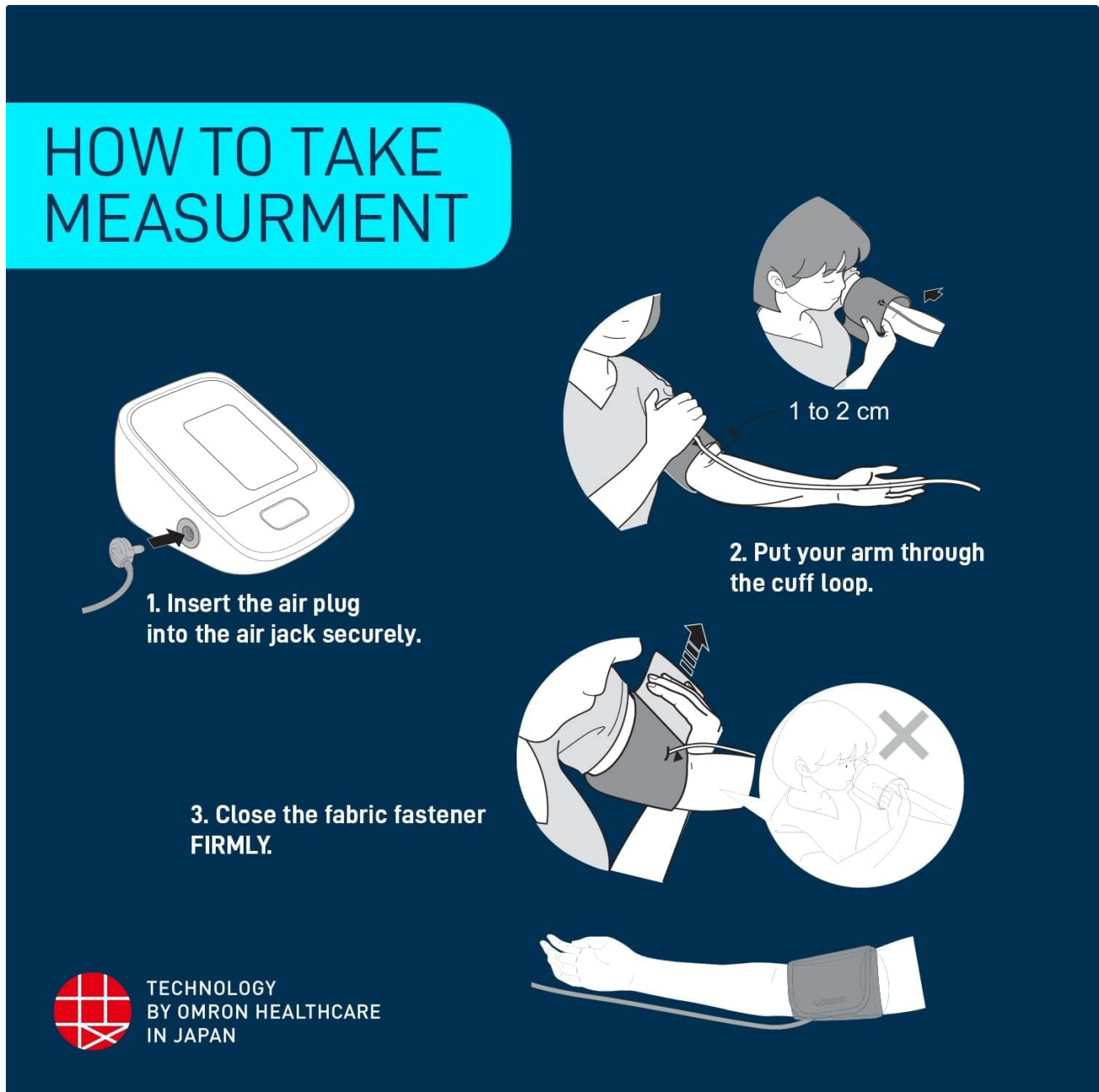


Figure 3: Step-by-step guide for connecting the air plug and correctly positioning the arm cuff for measurement.

OPERATING INSTRUCTIONS

1. Applying the Arm Cuff

1. Sit in a chair with your feet flat on the floor. Rest your arm on a table so the cuff is at heart level.
2. Slide your arm through the cuff loop. The air tube should run down the inside of your arm and be aligned with your middle finger.
3. The bottom edge of the cuff should be approximately 1 to 2 cm above your elbow.
4. Wrap the cuff firmly around your upper arm, ensuring it is snug but not too tight. The "OK" symbol for cuff wrapping will appear on the display if correctly applied.

2. Taking a Measurement

1. Ensure you are relaxed and have not exercised, eaten, or smoked for at least 30 minutes prior to measurement.
2. Press the **START/STOP** button. The cuff will begin to inflate automatically.
3. Remain still and silent during the measurement.

4. Once the measurement is complete, the cuff will deflate, and your systolic pressure, diastolic pressure, and pulse rate will be displayed.
5. Press the **START/STOP** button again to turn off the monitor. The measurement will be automatically stored in the device's memory.

3. Using Bluetooth Connectivity and Omron Connect App

The Omron HEM 7143T1 features Bluetooth connectivity, allowing you to sync your blood pressure readings with the Omron Connect app on your smartphone or tablet. This enables easy tracking, analysis, and sharing of your health data.


1. Download the "OMRON connect" app from your device's app store (Google Play Store or Apple App Store).
2. Follow the in-app instructions to pair your Omron HEM 7143T1 monitor with your smartphone via Bluetooth.
3. After each measurement, the data will automatically transfer to the app, providing a comprehensive history of your readings.
4. The app allows you to view trends, set goals, and share reports with your healthcare provider.



Omron Connect App

Your Smart Health Partner for sharing Accurate & Reliable results from OMRON *Devices to your Doctor


The Smart App to manage your blood pressure effectively

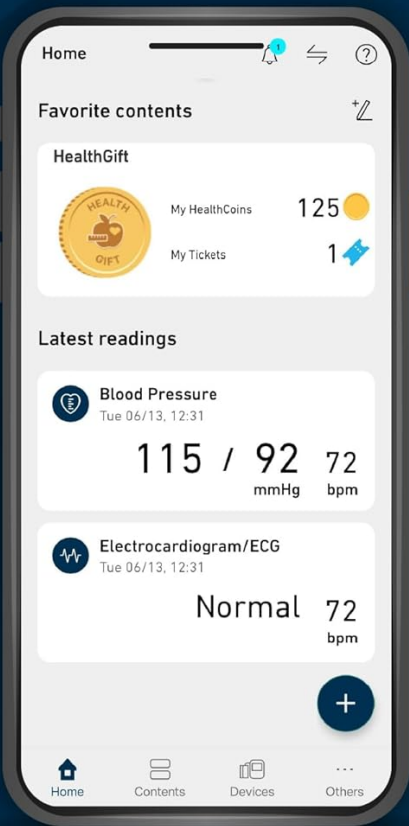



Instant Health Overview to keep track of your health data

Discover patterns, track your progress, and set goals to improve your health with intuitive dashboard

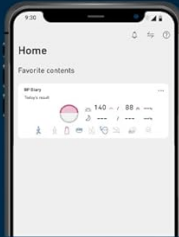

BP diary to set goals & manage your lifestyle

Track and optimize your blood pressure with ease using Omron Connect app's BP diary



Connect Your Fitness Apps

Seamlessly integrate fitness & wellness apps for unified health data in one place



HealthGift to win amazing prizes

Earn HealthCoins & Tickets for each measurement using the regular tracking features & redeem them for exciting gifts & offers

*Applicable on Blood Pressure monitors & selected Body Composition monitors



Figure 4: The Omron Connect App interface, showing features like health data tracking, BP diary, fitness app integration, and rewards.

MAINTENANCE AND STORAGE

Cleaning the Monitor

- Wipe the main unit and cuff with a soft, dry cloth.
- If necessary, use a cloth dampened with water and a mild soap, then wipe dry.
- Do not use abrasive cleaners or immerse the device in water.

Storage

- Store the monitor and cuff in the provided storage case or a clean, dry place.
- Avoid extreme temperatures, humidity, direct sunlight, and dust.
- Remove batteries if the unit will not be used for an extended period to prevent leakage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No display when pressing START/STOP.	Batteries are dead or incorrectly inserted.	Replace all batteries with new ones, ensuring correct polarity.
Cuff does not inflate or inflates too little.	Air plug is not securely connected or cuff is not properly applied.	Ensure air plug is fully inserted. Reapply cuff snugly.
Error message (e.g., "E1", "E2") displayed.	Movement during measurement, irregular heartbeat, or cuff too loose/tight.	Remain still and silent. Reapply cuff. If error persists, consult manual for specific error codes or contact customer support.
Cannot connect to Omron Connect App.	Bluetooth is off on phone, app not updated, or pairing issue.	Ensure Bluetooth is enabled on your phone. Update the app. Try re-pairing the device through the app settings.

PRODUCT SPECIFICATIONS

Feature	Detail
Model	HEM 7143T1
Display Type	Digital LCD
Measurement Method	Oscillometric method
Measurement Range	Pressure: 0 to 299 mmHg, Pulse: 40 to 180 beats/min
Accuracy	Pressure: ± 3 mmHg, Pulse: $\pm 5\%$ of display reading
Memory	30 measurements with date and time
Power Source	4 "AA" batteries 1.5V or optional AC adapter (not included)
Battery Life	Approximately 1000 measurements with new alkaline batteries
Product Dimensions	14L x 10.3W x 8.2Th cm (approx.)
Product Weight	Approx. 250 g (excluding batteries)
Cuff Circumference	22 to 32 cm (standard cuff)
Connectivity	Bluetooth® Low Energy
Operating Temperature/Humidity	10°C to 40°C (50°F to 104°F) / 15% to 90% RH (non-condensing)

WARRANTY AND SUPPORT



Omron products are manufactured to the highest quality standards. Specific warranty terms and conditions for your Omron HEM 7143T1 Blood Pressure Monitor are typically provided with the product packaging or can be found on the official Omron website for your region.


For technical support, troubleshooting assistance beyond this manual, or warranty claims, please contact Omron customer service in your country. Contact information can usually be found on the Omron website or on the product packaging.

You can also visit the official Omron website for FAQs and additional resources: www.omron-healthcare.com

© 2024 Omron Healthcare. All rights reserved.

Related Documents - 7143T1

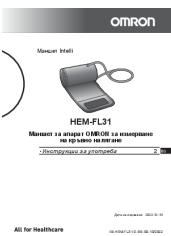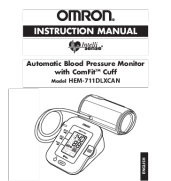	Omron RS7 Intelli IT and RS2 Blood Pressure Monitors Information on Omron's RS7 Intelli IT and RS2 blood pressure monitors, including features, specifications, and usage for accurate blood pressure monitoring.
	OMRON HEM-FL31 Intelli Cuff User Manual and Instructions Comprehensive user manual for the OMRON HEM-FL31 Intelli blood pressure monitor cuff, providing instructions for use, maintenance, and technical specifications.
	OMRON HEM-711DLXCAN IntelliSense Automatic Blood Pressure Monitor Instruction Manual User guide for the OMRON HEM-711DLXCAN IntelliSense Automatic Blood Pressure Monitor with ComFit Cuff. Learn how to operate, maintain, and troubleshoot your device for accurate home blood pressure readings.
	OMRON HEM-7143T1 Automatic Blood Pressure Monitor Instruction Manual Comprehensive instruction manual for the OMRON HEM-7143T1 Automatic Blood Pressure Monitor. This guide provides detailed information on setup, proper usage, measurement procedures, understanding readings, troubleshooting common errors, and device specifications to ensure accurate and reliable blood pressure monitoring.
	OMRON Intelli Wrap Cuff HEM-FL31 User Manual and Specifications Detailed user manual and specifications for the OMRON Intelli Wrap Cuff (HEM-FL31), including instructions for use, maintenance, and product features. Compatible with OMRON blood pressure monitors.

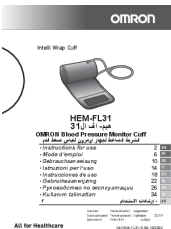

OMRON

Intelli Wrap Cuff

HEM-FL31

312

OMRON Blood Pressure Monitor Cuff and Intelli Wrap Cuff User Manual

Model's language

- Deutschsprachige Anleitung
- Instructions de l'utilisateur
- Gebruiksaanwijzing
- Instrucciones de uso
- Gebruiksaanwijzing
- Paradokshet për përdorim
- Käyttöohje

All for Healthcare

[OMRON Intelli Wrap Cuff HEM-FL31 User Manual & Specifications](#)

User manual and technical specifications for the OMRON Intelli Wrap Cuff HEM-FL31 blood pressure monitor cuff. Includes instructions for use, maintenance, safety precautions, and device specifications.

Documents - Omron – 7143T1




[\[pdf\]](#) User Manual Instructions

IM1 HEM 7143T1 EBK E1 02 04 2022 Manual de instruções Instruction im1 2 web im hem 7143t1 ebk e1 omron healthcare storage files manuals |||

Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK IM1-HEM-**7143T1**-EBK-E1-02-04/2022 English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . D A Display...


lang:en **score:52** filesize: 1.47 M page_count: 6 document date: 2022-09-05

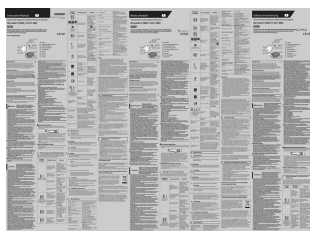

[\[pdf\]](#) User Manual Instructions Catalog

IM1 2 web HEM 7143T1 EBK 01 10 2020 M2 Intelli IT IM TR omron healthcare on demandware static Sites master catalog default dwa2742035 s ||| ||| Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-7143T1-EBK English Read Instruction manual and before use. For symbols information, refer to Symbols Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . IM1-HEM-**7143T1**-EBK-01-10/2020 5673055-3A D A ...

lang:en **score:49** filesize: 5.17 M page_count: 6 document date: 2021-05-20

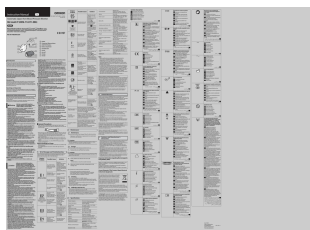

[\[pdf\]](#) User Manual Instructions

IM1 2 web HEM 7143T1 EBK 01 10 2020 Instruction Manual M49903EN gimaitaly DocumentiGIMA Manuali EN ||| Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . IM1-HEM-**7143T1**-EBK-01-10/2020 5673055-3A D A ...

lang:en **score:47** filesize: 3.69 M page_count: 3 document date: 2023-09-01




[\[pdf\]](#) User Manual

210312 HEM 7143T1 AP IM2 OH IM Multi fm DS 01 041 Instruction Manual 2 e4791dfc6f2e4afe9f081f2d4687dba1 omronhealthcare ap Content products |||

Instruction Manual 2 Automatic Blood Pressure Monitor Model HEM**7143T1** Read Instruction manual and before use. f . Bahasa indonesia Baca buku petunjuk dan sebelum penggunaan. Ting Vit c hng dn s dng v trc khi s dng. Bahasa Melayu Baca Manual arahan dan sebe...


lang:id **score:41** filesize: 18.97 M page_count: 2 document date: 2023-04-26

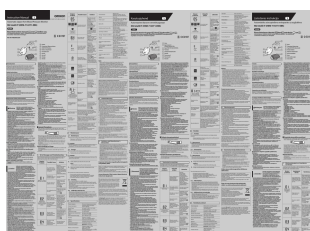


[\[pdf\]](#) User Manual Instructions

IM1 HEM 7143T1 EBK E3 02 04 2022 im1 2 web im hem 7143t1 ebk e3 omron healthcare storage files manuals ||| Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK IM1-HEM-**7143T1**-EBK-E3-02-04/2022 English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . D A Display...

lang:en **score:39** filesize: 1.95 M page_count: 6 document date: 2022-09-05

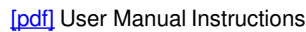


lang:en score:38 filesize: 1.79 M page_count: 6 document date: 2022-09-05

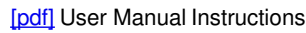

lang:fr score:38 filesize: 4.02 M page_count: 3 document date: 2023-09-01




lang:id score:35 filesize: 19.02 M page_count: 2 document date: 2023-03-02

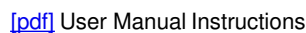


lang:en score:34 filesize: 5.07 M page_count: 6 document date: 2022-06-08



lang:i-klingon score:29 filesize: 3.67 M page count: 3 document date: 2023-09-01

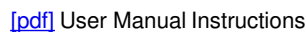

lang:en score:29 filesize: 5.1 M page count: 6 document date: 2022-06-08

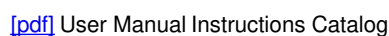

lang:en score:28 filesize: 1.45 M page count: 6 document date: 2021-07-29



[\[pdf\]](#) User Manual Instructions Catalog

IM1 HEM 7143T1 EBK E3 01 10 2020 M2 Intelli IT web IM RO omron healthcare on demandware static Sites master catalog default dw71a52c4b s |||

Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK IM1-HEM-**7143T1**-EBK-E3-01-10/2020 English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . D A Display...

lang:en **score:28** filesize: 1.92 M page_count: 6 document date: 2021-07-29




[\[pdf\]](#) User Manual

230302 HEM 7143T1 AP IM2 OH IM Multi fm DS 01 041 Omron HEM7143T1 Standard Blood Pressure Monitor Instruction Manual 2 5653123 2E 230308 smartwellness au 2023 10 |||

...

lang:en **score:27** filesize: 22.49 M page_count: 2 document date: 2023-09-28

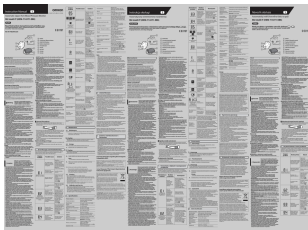


[\[pdf\]](#) User Manual Instructions Catalog

IM1 HEM 7143T1 EBK E2 01 10 2020 M2 Intelli IT web IM HU omron healthcare on demandware static Sites master catalog default dwfbab1a62 s |||

Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK IM1-HEM-**7143T1**-EBK-E2-01-10/2020 English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . D A Display...

lang:en **score:27** filesize: 1.76 M page_count: 6 document date: 2021-07-29

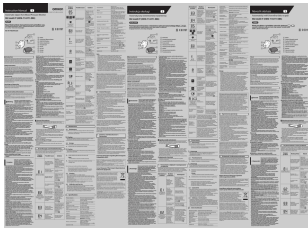

[\[pdf\]](#) User Manual Instructions

IM1 HEM 7143T1 EBK E2 02 04 2022 im1 2 web im hem 7143t1 ebk e2sklep stanmed24 pl 2023 11 e2 |||

Instruction Manual Automatic Upper Arm Blood Pressure Monitor M2 Intelli IT HEM-**7143T1**-EBK IM1-HEM-**7143T1**-EBK-E2-02-04/2022 English Read Instruction manual and before use. For symbols information, refer to Symbols Description on the back side of the another Instruction manual . D A Display...

lang:en **score:27** filesize: 1.79 M page_count: 6 document date: 2022-09-05



[\[pdf\]](#) User Manual Specifications Guide

demhel OMRON Europe B V Model code M300 Intelli IT HEM 7143T1 D PM 2070 01 102021 NPIS EN mindtecstore mediafiles Sonstiges Shop |||

New Product Information Sheet Specifications subject to change OMRON HEALTHCARE EUROPE B.V. Model code : M300 Intelli IT HEM-**7143T1**-D Versions 1 Remark 1st version, without key market visual Date 13/10/2021 PM-2070-01-10/2021 M300 Intelli IT Accurate measurements in one touch Product Name ...

lang:en **score:26** filesize: 1.13 M page_count: 19 document date: 2021-10-13




[pdf] User Manual Declaration of Conformity

Prateek 01 10 cdr Lenovo 1mile healthcare 2 days ago — Bluetooth BP HEM 7142 7143T1 PAN INDIA Freight Paid
Krups Manual Scales Complete Range Available At 1Mile Healthcare dr morepen glucose bg03 25 test strips 5
iming data5 SELLER Doc 2023 4 298684072 GC OG PI 49602298 |||

We Are Authorised Distributors Of: R 1MILE HEALTHCARE R I N V E S T I N Y O U R H E A L
T H POINT ... t 1Mile Healthcare BP HEM-7121J BP HEM-7124 BP Bluetooth BP Bluetooth BP
HEM-8712 HEM-7142 HEM-**7143T1** Bluetooth BP BP BP Bluetooth BP HEM-**7143T1**-AP
HEM-7156AP HEM-7156AAP HEM-6161AP HEM-...

lang:en score:26 filesize: 18.7 M page_count: 18 document date: 2023-03-29

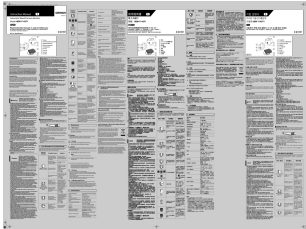


[pdf] User Manual

230207 HEM 7143T1 AP IM1 OH IM en zhtw ko my hi fm DS 01 041 Omron HEM7143T1 Standard Blood Pressure
Monitor Instruction Manual 1 5653098 8F 230308 smartwellness au 2023 10 |||

...

lang:en score:26 filesize: 32.51 M page_count: 2 document date: 2023-09-28

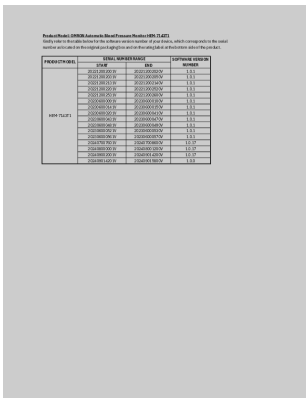


[pdf] Label

Jia Xin Koh OHS OMRON Automatic Blood Pressure Monitor HEM 7143T1 SERIAL NUMBER RANGE SOFTWARE
VERSION Product Model PRODUCT MODEL Kindly refer to the AP Software Version AU omronhealthcare ap
Content software version |||

Product Model: OMRON Automatic Blood Pressure Monitor HEM**7143T1** Kindly refer to the
table below for the software version number of your device, which corresponds to the serial
number as located on the original packaging box and on the rating label at the bottom side of
the product. PRODUCT MODEL H...

lang:i-klngon score:25 filesize: 36.47 K page_count: 1 document date: 2024-11-12

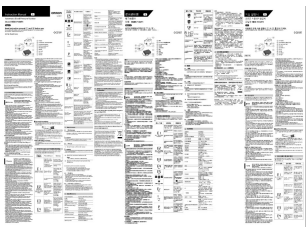

OMRON HEM-7143T1 Automatic Blood Pressure Monitor Instruction Manual

Comprehensive instruction manual for the OMRON HEM-7143T1 Automatic Blood Pressure
Monitor. This guide provides detailed information on setup, proper usage, measurement
procedures, understanding readings, troubleshooting common errors, and device
specifications to ensure accurate and reliable blood pressure monitoring.

lang:i-klngon score:25 filesize: 3.45 M page_count: 4 document date: 2023-01-17




[pdf] User Manual

220125 HEM 7143T1 AP IM2 OH IM Multi fm DS 01 041 Instruction Manual 2 NEW STAR COMPANY LTD OMRON
Upper Arm Blood Pressure Monitor drive google file d 1b2fCSmT yqt4ZL4sl40Qyg3b3QSos TA view usp sharing |||

...

lang:i-klngon score:25 filesize: 22.49 M page_count: 2 document date: 2023-06-12

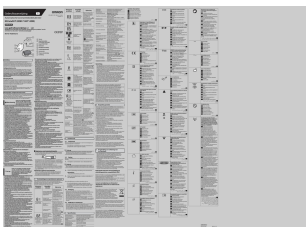


[pdf]

IM1 2 web HEM 7143T1 EBK 01 10 2020 M49903NL gimaitaly DocumentiGIMA Manuali NL |||

Gebruiksaanwijzing Automatische bovenarmbloeddrukmeter M2 Intelli IT HEM-**7143T1**-EBK
Nederlands Lees gebruiksaanwijzing Raadpleeg voor informatie over symbolen Beschrijving
van symbolen aan de achterkant van de andere Gebruiksaanwijzing . IM1-HEM-**7143T1**-EBK-
01-10/2020 5673056-1A D A Scherm ...

lang:nl score:25 filesize: 3.73 M page_count: 3 document date: 2023-09-01

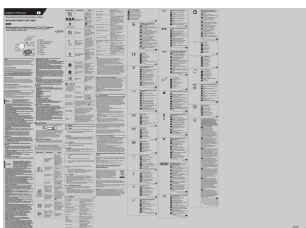


[pdf]

IM1 2 web HEM 7143T1 EBK 01 10 2020 M49903TR gimaitaly DocumentiGIMA Manuali TR |||

Kullanım Kilavuzu Otomatik st Kol Kan Basıncı Im Cihazı M2 Intelli IT HEM-**7143T1**-EBK Trke
Kullanmadan nce kullanim kilavuzu ve yi okuyun. Simgelerle ilgili bilgi iin, bir dier Kullanım
kilavuzu belgesinin arka yzndeki Simgelerin Aiklamasi blmne bakın. D A Grnt B Balantı Hafıza
dmesi ...

lang:id score:24 filesize: 3.99 M page_count: 3 document date: 2023-09-01

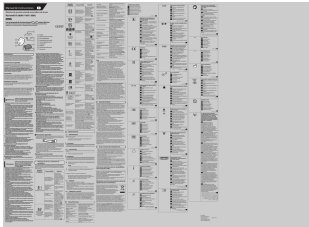


[pdf] User Manual

IM1 2 web HEM 7143T1 EBK 01 10 2020 M49903ES gimaitaly DocumentiGIMA Manuali ES |||

Manual de instrucciones Monitor de presión arterial automático de brazo M2 Intelli IT

HEM-7143T1-EBK Español Lea el manual de instrucciones y antes del uso. Para obtener información acerca de los símbolos, consulte la Descripción de los símbolos en el anverso del otro Manual de instrucciones . D A ...

lang:es score:24 filesize: 3.72 M page_count: 3 document date: 2023-09-01

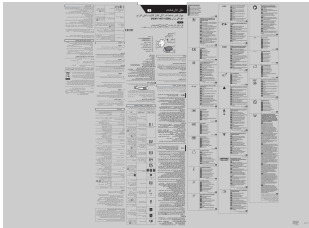

[pdf]


IM1 2 web HEM 7143T1 EBK 01 10 2020 May 17 2021 · Pièce appliquée Type BF Degré de protection contre les

chocs électriques courant fuite Anwendungsteil Typ M49903SA gimaitaly DocumentiGIMA Manuali SA |||

. .OMRON . . OMRON . OMRON . - ... EN1060 : : . .OMRON HEALTHCARE Co., Ltd

HEM-7143T1-EBK .2014/53/1EU www.omron-healthcare.com : OMRON . OMR...

lang:i-klnglon score:24 filesize: 3.69 M page_count: 3 document date: 2023-09-01

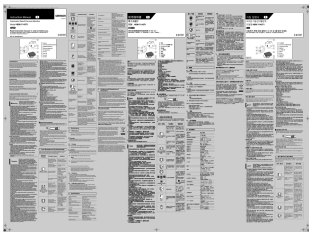


[pdf] User Manual

220125 HEM 7143T1 AP IM1 OH IM en zhtw ko my hi fm DS 01 041 Instruction Manual 1 NEW STAR COMPANY

LTD OMRON Upper Arm Blood Pressure Monitor drive google file d 1jX d2rQFmwzGy YmmNb2DSCaYBnlAcy view
usp sharing |||

...

lang:i-klnglon score:23 filesize: 32.34 M page_count: 2 document date: 2023-06-01

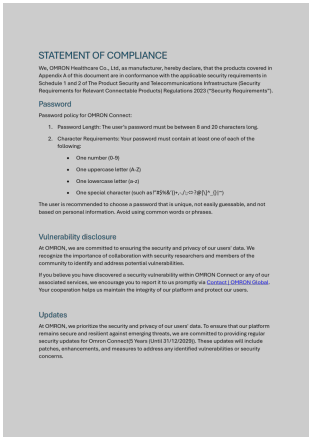


[pdf]

Jeremy Costa Pires HCB PTN STATEMENT OF COMPLIANCE Omron healthcare co uk compliance statement 2024

12 20 — We OMRON Healthcare Co Ltd as manufacturer hereby declare that the products covered in Appendix A of
this document are conformance with PSTI Act EMEA psti omron storage files policies |||

STATEMENT OF COMPLIANCE We, OMRON Healthcare Co., Ltd, as manufacturer, hereby
declare, that the pro ... d HEM-7381T1-EBK/EOSL OMRON Healthcare Co., Ltd HEM-7155T-
EBK/ESL OMRON Healthcare Co., Ltd HEM-7143T1-EBK /HEM 7143T2-ESL OMRON
Healthcare Co., Ltd HEM-7146T2-EBK/ESL OMRON Healthcare Co., Ltd...

lang:en score:22 filesize: 212.8 K page_count: 3 document date: 2024-12-11

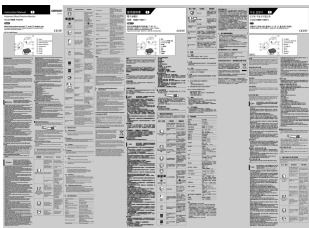

[pdf] User Manual

231020 HEM 7143T1 AP IM1 OH IM en zhtw ko my hi fm DS 01 041 ManualInstruction Manual Omron Healthcare

Malaysia89036650d9da4dc29fb9f34d8d72ce5domronhealthcare ap Content products

89036650d9da4dc29fb9f34d8d72ce5d ManualHEM MalaysiaInstruction Blood Pressure Monitor

Malaysiaomronhealthcare 89036650d9da4dc29fb9f34d8d72ce5domronhealthcare 443 Instruction manualInstruction
Monitors Upper Arm Asia PacificProduct Support PacificInstruction

Malaysia4448966057f446a68c5ae18388d687234448966057f446a68c5ae18388d68723omronhealthcare

4448966057f446a68c5ae18388d68723omronhealthcare ||| |||

...

lang:en score:22 filesize: 16.54 M page_count: 2 document date: 2024-10-31




[pdf] Declaration of Conformity
SKM C550i24011514530 903PP6215 EU Declaration of Conformity DoCProduct Support Omron Healthcare Asia PacificHEM 7143T1 AP OHQ CS DoC RED 5636391Bomronhealthcare ap Content document HEM 5636391B |||



[pdf]
View MoreUpper Arm Automatic Blood Pressure Monitor HEM 7143T1 OMRON Healthcare Brand Shop11zon JPEG to 13cdn shopify s files 1 0262 4130 2577 11zon 13 v 1706692466 |||



[pdf]
View MoreUpper Arm Automatic Blood Pressure Monitor HEM 7143T1 OMRON Healthcare Brand Shop11zon JPEG to 17cdn shopify s files 1 0262 4130 2577 11zon 17 v 1706692655 |||

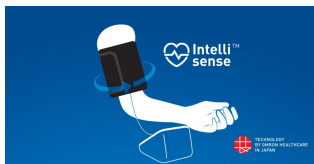

[pdf]
View MoreUpper Arm Automatic Blood Pressure Monitor HEM 7143T1 OMRON Healthcare Brand Shop11zon JPEG to 19cdn shopify s files 1 0262 4130 2577 11zon 19 v 1706692706 |||



[pdf]
View MoreUpper Arm Automatic Blood Pressure Monitor HEM 7143T1 OMRON Healthcare Brand Shop11zon JPEG to 15cdn shopify s files 1 0262 4130 2577 11zon 15 v 1706692552 |||



[pdf]
View MoreUpper Arm Automatic Blood Pressure Monitor HEM 7143T1 OMRON Healthcare Brand Shop11zon JPEG to 16cdn shopify s files 1 0262 4130 2577 11zon 16 v 1706692607 |||

