

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [MCombo](#) /

› [MCombo Dual Motor Power Lift Recliner Chair Model 6160-7815PU Instruction Manual](#)

MCombo 6160-7815PU

MCombo Dual Motor Power Lift Recliner Chair

MODEL: 6160-7815PU



Introduction

This manual provides detailed instructions for the assembly, operation, and maintenance of your MCombo Dual Motor Power Lift Recliner Chair, Model 6160-7815PU. Please read this manual thoroughly before use to ensure proper function and safety. Keep this manual for future reference.



Image: The MCombo Dual Motor Power Lift Recliner Chair in Cream White, showcasing its design and features.

Important Safety Information

Always use the HOME button on the remote control to lift from the prone position. The backrest will move first to ensure your safety. Allow up to 16 seconds for the chair to fully lift.

- Ensure all connections are secure before operating the chair.
- Keep children and pets away from the chair's moving parts.
- Do not allow anyone to sit on the footrest when it is extended.
- Avoid placing objects under the chair that could obstruct its movement.
- Do not attempt to disassemble or repair the chair yourself. Contact customer support for assistance.
- Maximum weight recommendation: 350 pounds.

Package Contents

Before assembly, verify that all components are present:

- Backrest
- Lumbar Pillow
- Accessory Box (contains remote, power adapter, cables)
- Instruction Manual (this document)
- Left Armrest
- Right Armrest
- Seat Base

The chair is typically shipped in 3 boxes. Ensure all boxes are received before beginning assembly.

Setup and Assembly

Assembly of the MCombo Power Lift Recliner Chair is designed to be straightforward and requires no tools. It can typically be completed by a novice in approximately 20 minutes.

1. **Unpack Components:** Carefully remove all parts from their packaging. Retain packaging materials until assembly is complete and the chair is functioning correctly.
2. **Attach Backrest:** Slide the backrest onto the metal posts located on the seat base until it clicks securely into place.
3. **Attach Armrests:** Slide the left and right armrests onto their respective slots on the seat base. Ensure they are firmly seated.
4. **Connect Wiring:** Locate the wiring harnesses extending from the backrest, armrests, and seat base. Connect the corresponding wires. Ensure all connections are tight and secure. This step may require getting down to floor level.
5. **Connect Power:** Plug the power adapter into the chair's power input, then plug the adapter into a wall outlet.
6. **Initial Test:** Use the remote control to test all functions (lift, recline, massage, heat) to ensure proper operation.

Operating Instructions

Remote Control Functions




Lay-flat



Reclined



TV

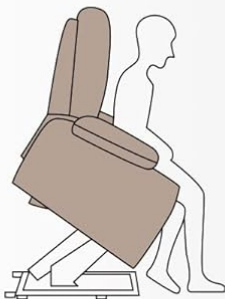

Lift

Home Button

(Back to upright position)

Infinite Position



Image: The primary remote control for adjusting chair positions, including lift, recline, TV, lay-flat, and a central HOME button.

- **Dual Motor System:** The chair features a dual motor system, allowing independent adjustment of the backrest and footrest. This enables a wide range of positions, including lift, TV viewing, and lay-flat. The maximum reclining angle is 165 degrees.
- **HOME Button:** Pressing the HOME button on the remote control will return the chair to its upright position from any reclined or lifted state. For safety, when lifting from a prone position, the backrest will move first. This process takes approximately 16 seconds.
- **Position Adjustment:** Use the dedicated buttons on the remote to adjust the backrest and footrest to your desired position.

Massage and Heating Functions

Massage and Heating

Heating function works separately from vibration.



8 Vibration Nodes

3 Massage Modes

3 Intensity Levels



• Lumbar Heat



Image: A separate remote control for the massage and heating features, detailing controls for 8 vibration nodes, 3 massage modes, 3 intensity levels, and lumbar heat.

- The heating function operates independently from the vibrating function.
- **Vibration:** Select from 8 vibration nodes, 3 massage modes, and 3 intensity levels to target specific areas and customize your massage experience. A timer function is also available.
- **Lumbar Heat:** Activate the lumbar heating function for soothing warmth in the lower back area.

Convenience Features

Everything is right at hand.

- 2 Anti-falling Ropes
 - 2 USB & Type C Ports
 - 2 Cup Holders
 - 2 Side & Front Pockets
- ## Overstuffed Armrests



Image: Detail of the chair's armrest, highlighting the integrated dual USB and Type-C charging ports, cup holders, and convenient side pockets.

- **USB Ports & Cup Holders:** The chair includes dual USB ports for charging mobile devices and cup holders for beverages.
- **Extended Footrest:** The 4.7-inch extended footrest provides additional support for taller individuals (recommended height: 5'7"-6'2"), promoting better circulation.
- **Side & Front Pockets:** Convenient pockets are integrated into the design for storing remotes, magazines, or other small items.

Maintenance

- **Cleaning:** Wipe clean with a damp cloth. Avoid harsh chemicals or abrasive cleaners that may damage the faux leather upholstery.
- **Regular Inspection:** Periodically check all connections and moving parts for wear or looseness. Tighten any loose fasteners.

- **Power Cord Care:** Ensure the power cord is not pinched or placed under heavy objects. Avoid tripping hazards.

Troubleshooting

If you encounter any issues with your MCombo Power Lift Recliner Chair, please refer to the following common solutions:

- **Chair Not Moving:**
 - Check if the power cord is securely plugged into both the chair and the wall outlet.
 - Ensure all internal wiring connections are firm.
 - Verify that the wall outlet is functional by plugging in another device.
- **Noise During Operation:**
 - Ensure all assembly bolts and connections are tightened.
 - Check for any obstructions under or around the chair's mechanism.
- **Massage/Heat Not Working:**
 - Confirm the massage/heat remote is properly connected and has power.
 - Ensure the functions are activated on the remote.

If the problem persists after attempting these steps, please contact MCombo customer support.

Specifications



Image: Detailed dimensional diagram of the recliner chair, illustrating measurements for height, width, depth, seat height, and recline angles.

Feature	Detail
Brand	MCombo
Model Name	6160-7815PU
Color	Cream White
Product Dimensions (Upright)	37.4"D x 36.6"W x 46"H
Seat Height	19.7 Inches
Seat Depth	22 Inches
Maximum Weight Recommendation	350 Pounds

Special Features	Dual Motor, Power Lift, Heat, Massage, Extended Footrest, USB Ports
Upholstery Material	Faux Leather
Fill Material	Foam
Frame Material	Wood and Metal (FSC certified wooden frame)
Item Weight	142.2 pounds



Image: Certifications indicating the use of CertiPUR-US certified foam, OEKO-TEX Standard 100 certified fabric, and FSC certified wood in the chair's construction.

Warranty

The MCombo Dual Motor Power Lift Recliner Chair comes with a**2-Year Limited Warranty**. This warranty covers:

- Reclining and Lift Mechanism
- Electrical Parts
- Upholstery Material
- Cushions
- Wood and Metal Frame







Please retain your proof of purchase for warranty claims.

Customer Support

For any questions, concerns, or assistance with your MCombo Power Lift Recliner Chair, please contact MCombo customer support. Refer to the contact information provided at the point of purchase or on the official MCombo website.

Online Resources: For additional information and support, visit the official MCombo store:[MCombo Store on Amazon](#)

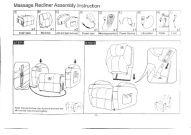

	<p>Mcombo Lift Recliner Owner's Manual: Assembly, Operation, and Care</p> <p>Comprehensive owner's manual for the Mcombo Lift Recliner, covering assembly instructions, operating procedures, care and maintenance, troubleshooting, warranty, and return policy. Includes model numbers 6160-7893, 6160-7679, 6160-7680.</p>
	<p>MCombo Lift Recliner Owner's Manual: Assembly, Operation, and Care Guide</p> <p>Comprehensive owner's manual for the MCombo Lift Recliner, detailing safety guidelines, assembly instructions, operating procedures for lift and massage functions, care and maintenance tips, troubleshooting common issues, lift motor replacement, disassembly, warranty information, and return policy. Includes model numbers 6160-7893, 6160-7890, and 6160-7896.</p>
	<p>Mcombo Lift Chair 6160-7040DB Assembly and Operation Guide</p> <p>Comprehensive assembly and operation guide for the Mcombo Lift Chair 6160-7040DB, detailing parts, assembly steps, remote control functions, and troubleshooting tips.</p>
	<p>Massage Recliner Assembly Instructions and User Guide</p> <p>Comprehensive assembly instructions and FAQ for the Mcombo Massage Recliner, detailing setup, operation of the remote control, reclining functions, and troubleshooting common issues.</p>
	<p>MCombo Lift Recliner Owner's Manual: Assembly, Operation, and Care</p> <p>Comprehensive owner's manual for the MCombo Lift Recliner (Large-Regular), model 6160-7539/7549. Includes assembly instructions, operating procedures, safety guidelines, maintenance tips, troubleshooting, warranty, and return policy.</p>
	<p>MCombo Lift Recliner Owner's Manual and Assembly Guide</p> <p>Comprehensive owner's manual for the MCombo Lift Recliner, covering safety guidelines, assembly instructions, operating procedures, care and maintenance, troubleshooting, warranty, and return policy. Includes detailed diagrams and instructions for model numbers 6160-7631, 6160-7634, and 6160-7888.</p>