

## Lionelo LOC-FIDI FREE WHITE

# Fidi Free Electric Breast Pump User Manual

**Brand:** Lionelo | **Model:** LOC-FIDI FREE WHITE

## 1. INTRODUCTION

The Lionelo Fidi Free Electric Breast Pump is designed for convenience and efficiency, offering a hands-free and discreet pumping experience. Its quiet operation and portable design allow for use anytime, anywhere, without disturbing others. This manual provides essential information for the safe and effective use of your breast pump.



Image: The Fidi Free Electric Breast Pump, showcasing its compact design and transparent milk collection cup.

# Fidi FREE



## VERY QUIET OPERATION



Image: A chart illustrating the quiet operation of the Fidi Free breast pump, with a noise level of 50dB, comparable to a whisper.

## 2. SAFETY INFORMATION

Please read all safety instructions carefully before using the device. Failure to follow these instructions may result in injury or damage to the product.

- Ensure all components are clean and sterilized before each use.
- The funnels are made of skin-friendly silicone, intended for contact with food.
- Do not immerse the main pump unit in water or other liquids.
- Keep out of reach of children.
- Use only the provided USB cable for charging.
- Do not use the device if any part is damaged.
- This product meets the required European standards for safety and hygiene.

### 3. PRODUCT COMPONENTS

The Fidi Free Electric Breast Pump package includes the following items:

- Breast pump main unit
- Milk bottle (1 pc, 150 ml) with matching teat
- Milk storage bags (2x)
- Cotton bag for storing the breast pump (20x30 cm)
- USB charging cable
- Funnels (1x)
- Funnel lid (1x)



Image: The three included funnel sizes (17mm, 21mm, 25mm) for optimal comfort and fit.

### 4. SETUP

#### 4.1 Assembly

Follow these steps to assemble your breast pump before first use:

1. Ensure all parts are clean and sterilized according to the Maintenance section.
2. Attach the appropriate size funnel to the main pump unit. Ensure it is securely fitted.
3. Connect the milk collection bottle to the bottom of the pump unit.
4. Place the funnel lid on top of the funnel if not in use or for storage.

## 4.2 Charging the Device

Before first use, fully charge the breast pump. The integrated LEDs will indicate the charging status.

- Connect the USB cable to the charging port on the pump unit.
- Plug the other end of the USB cable into a compatible USB power adapter (not included).
- The LEDs will illuminate to show charging progress.
- A full charge provides a full day's use.

### ● OPERATION INDICATOR IN THE FORM OF LED DIODES



### ● BATTERY CHARGE INDICATOR

### ● SUCTION POWER: UP TO 0.031 KPA (5 WORK INTENSITY LEVELS)

Image: The control panel of the breast pump, highlighting the LED operation and battery charge indicators.

## 5. OPERATING INSTRUCTIONS

The Fidi Free Electric Breast Pump features two working modes and multiple intensity levels for comfortable and efficient milk expression.

## 5.1 Powering On/Off

- To turn on, press and hold the power button. The LEDs will light up.
- To turn off, press and hold the power button again.

## 5.2 Working Modes and Intensity Levels

The device offers two distinct working modes, each with 5 intensity levels:

- **Massage Mode:** This mode mimics a baby's initial rapid suckling to stimulate milk let-down. Use this mode to initiate pumping.
- **Milk Expression Mode:** Once milk flow begins, switch to this mode for steady and efficient milk extraction. This mode simulates a baby's deeper, slower suckling.

Use the '+' and '-' buttons to adjust the intensity level within each mode. Start with a lower intensity and gradually increase to a comfortable level.

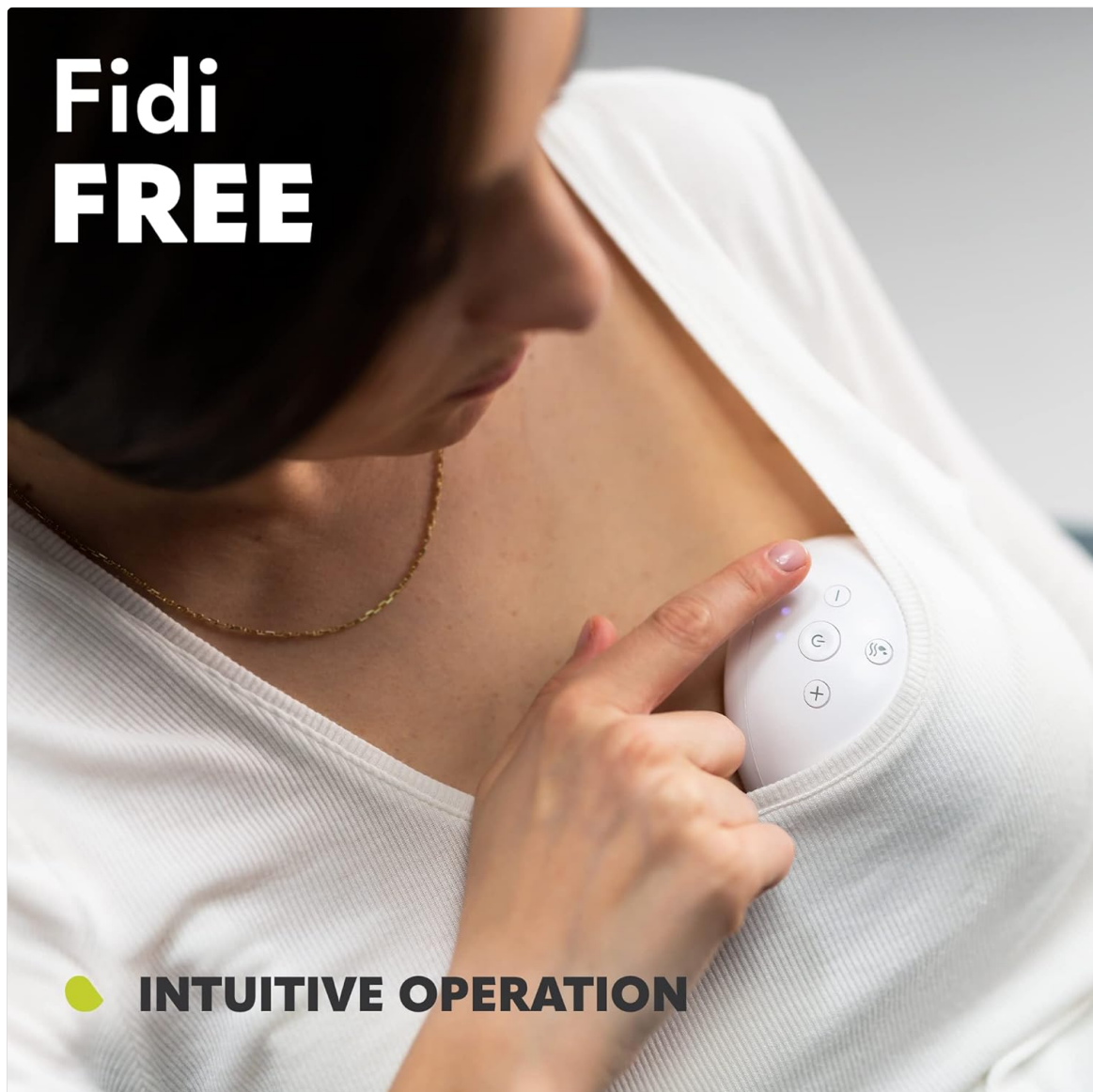


Image: The breast pump being used hands-free and discreetly, showing the intuitive button controls.

## 5.3 Timer Function

The breast pump is equipped with a timer and will automatically turn off after 30 minutes of continuous



operation. This feature helps prevent over-pumping and conserves battery life.

## **6. MAINTENANCE AND CLEANING**

Proper cleaning and sterilization are crucial for hygiene and the longevity of your breast pump.

### **6.1 Disassembly for Cleaning**

Before cleaning, ensure the device is turned off and unplugged. Carefully disassemble all parts that come into contact with breast milk: the funnel, milk collection bottle, and any other removable silicone components.

### **6.2 Cleaning Instructions**

- Wash all disassembled parts (except the main pump unit) with warm soapy water immediately after each use.
- Rinse thoroughly with clean water.
- Do not use abrasive cleaners or harsh chemicals.
- The main pump unit should be wiped clean with a damp cloth only. Do not submerge it in water.

### **6.3 Sterilization**

Silicone parts and the milk container can be sterilized at a temperature of 100°C (212°F). This can be done using a steam sterilizer, boiling water, or a microwave sterilizer according to the sterilizer's instructions.



Image: An exploded view demonstrating the sterilizable components of the breast pump, including the funnel and milk collection parts.

## 7. TROUBLESHOOTING

If you encounter any issues with your breast pump, refer to the table below for common problems and solutions.

Problem	Possible Cause	Solution
No suction or weak suction	Improper assembly; Damaged silicone parts; Low battery	Check all connections are secure; Inspect silicone parts for tears and replace if necessary; Charge the battery fully.
Pump not turning on	Battery is depleted; Device malfunction	Charge the breast pump; If still not working, contact customer support.
Uncomfortable pumping experience	Incorrect funnel size; Intensity level too high	Ensure you are using the correct funnel size (17mm, 21mm, or 25mm); Reduce the intensity level.



Problem	Possible Cause	Solution
Milk leaking	Improper assembly; Damaged bottle or seal	Ensure bottle is securely attached and all seals are in place; Inspect bottle and seals for damage.

## 8. SPECIFICATIONS

Feature	Detail
Model Number	LOC-FIDI FREE WHITE
Dimensions	4.33 x 2.76 x 5.12 inches
Item Weight	12.3 ounces
Material	Silicone, PP
Power Source	Battery operated (1 Lithium Ion battery required)
Working Modes	2 (Massage, Milk Expression)
Intensity Levels	5 per mode
Noise Level	Less than 50dB
Funnel Sizes Included	17 mm, 21 mm, 25 mm
Milk Container Capacity	150 ML
Automatic Shut-off	After 30 minutes

## 9. WARRANTY AND SUPPORT


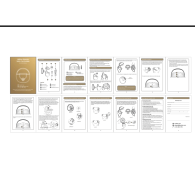
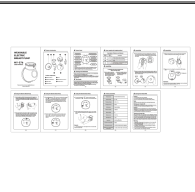
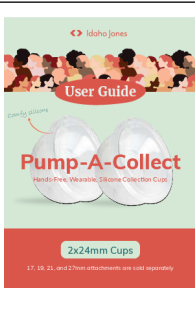


The Lionelo Fidi Free Electric Breast Pump is designed and manufactured to high-quality standards. This product conforms to the standards required in the EU.

For any questions, technical assistance, or warranty claims, please contact the seller or manufacturer directly. Keep your purchase receipt as proof of purchase.

**Manufacturer:** Lionelo

**UPC:** 5903771703659

## Related Documents - LOC-FIDI FREE WHITE

	<p><a href="#">Lionelo Fidi Go Breast Pump User Manual</a></p> <p>This document serves as the official user manual for the Lionelo Fidi Go breast pump. It provides comprehensive guidance on the setup, operation, maintenance, and safety of the Fidi Go breast pump. Users will find detailed instructions on assembly, usage modes, milk storage, troubleshooting, and cleaning procedures.</p>
	<p><a href="#">Mykitor Wearable Breast Pump User Manual</a></p> <p>User manual for the Mykitor Wearable Breast Pump, detailing its overview, operation, assembly, cleaning, and storage instructions.</p>
	<p><a href="#">MY-376 Wearable Electric Breast Pump User Manual</a></p> <p>User manual for the MY-376 Wearable Electric Breast Pump, detailing product introduction, components, assembly, control panel, usage instructions, cleaning, disinfection, product parameters, and breast milk storage.</p>
	<p><a href="#">Idaho Jones Pump-A-Collect User Guide: Hands-Free Wearable Silicone Collection Cups</a></p> <p>Comprehensive user guide for Idaho Jones Pump-A-Collect hands-free, wearable silicone collection cups. Learn about assembly, cleaning, usage, troubleshooting, and compatible pumps.</p>
	<p><a href="#">MomMed S21 Wearable Breast Pump User Manual</a></p> <p>Comprehensive user manual for the MomMed S21 Wearable Breast Pump. Provides detailed instructions on assembly, operation modes, cleaning, maintenance, and breast milk storage. Essential guide for users of MomMed's hands-free electric breast pump.</p>
	<p><a href="#">MomMed S21 Wearable Breast Pump: Installation and User Guide</a></p> <p>Comprehensive guide for the MomMed S21 wearable breast pump, covering installation, assembly, operation, cleaning, and storage of breast milk. Learn how to use your hands-free breast pump effectively.</p>