

Momcozy M6

Momcozy M6 Hands-Free Wearable Breast Pump Instruction Manual

Model: M6 (BP311)

1. IMPORTANT SAFETY INFORMATION

Warning: Do NOT use this device if you have a cardiac pacemaker, implanted defibrillator, or any other metallic or electronic implant, as it may cause electromagnetic interference or serious injury.

Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

2. PRODUCT OVERVIEW

The Momcozy M6 Hands-Free Wearable Breast Pump is designed for discreet and efficient milk expression. It features a 2nd Gen motor for strong suction, 3 modes, and 9 adjustable levels to optimize comfort and milk production. The DoubleFit Flange ensures a secure and comfortable fit, while the stand-up milk collector simplifies pouring. The pump is designed for easy cleaning with only 4 parts.

Key Features:

- **Ideal Rhythm:** Optimized pumping rhythm for enhanced comfort and increased milk production.
- **Powerful Suction:** Advanced 2nd Gen motor with 295 mmHg suction.
- **Multiple Modes & Levels:** 3 modes (Stimulation, Expression, Mixed) and 9 adjustable levels.
- **Discreet Design:** Contours to breast shape for seamless fit within a bra.
- **DoubleFit Flange:** Soft and hard silicone dual-seal for comfort and safety.
- **Stand-Up Milk Collector:** Simplifies milk transfer without spills.
- **Custom Fit:** Available in 4 sizes (17/19/21/24mm) for personalized comfort.
- **Long-Lasting Battery:** Up to 180 minutes (6 sessions) on one charge.
- **User-Friendly Display:** Large LED display shows time, battery, and intensity.
- **Easy Maintenance:** Only 4 easy-to-clean parts.

3. PRODUCT COMPONENTS

The Momcozy M6 Breast Pump consists of the following main components:

- **Pump Motor:** The main unit that generates suction and houses controls.
- **Flange (24mm):** The part that comes into contact with the breast.
- **Flange Inserts (17/19/21mm):** For personalized nipple sizing.
- **Milk Collector:** Gathers expressed milk.
- **Flange Cover:** Protects the flange.
- **Diaphragm:** Prevents backflow of milk.
- **Silicone Valve:** Directs milk flow into the collector.
- **USB Type-C Cable:** For charging the pump.
- **Storage Case:** For convenient transport and storage.



The infographic features a photograph of two women sitting at a table outdoors, smiling. In front of them are various pump components. To the right, a dark red background contains white text and icons. Below the photo, three white boxes with icons and text are displayed: a calendar icon for '365 Days', a speech bubble icon for '15 Minutes Fast Reply', and a truck icon for 'Fast Delivery'. On the right side, the text reads: '60+ Countries and regions', '4,000,000+ Users (Since brand's inception)', and '#1 The First Pump Inspired by Bra Design'.

Image: Detailed view of all components included in the Momcozy M6 two-pack box, showing pump motors, flanges, milk collectors, covers, diaphragms, valves, flange inserts, USB cables, user manual, flange size ruler, operating instructions, and storage case.

4. SETUP AND ASSEMBLY

Before first use, sterilize all pump parts that come into contact with breast milk (flange, milk collector, flange cover, diaphragm, silicone valve) according to the cleaning instructions.

Assembly Steps:

1. Insert the silicone valve firmly into the bottom of the milk collector. Ensure it is securely seated to prevent leaks.
2. Place the diaphragm onto the top of the milk collector, ensuring a tight seal.
3. Attach the flange (with any necessary insert) to the milk collector, aligning the pouring spout with the groove. Ensure the lip of the flange is securely over the edge of the collector.
4. Connect the pump motor to the assembled flange and collector unit. Ensure the electrical contacts align and the motor clicks securely into place.

Your browser does not support the video tag.

Video: Visual guide on assembling the Momcozy M6 Wearable Pump components.

Flange Sizing:

Correct flange sizing is crucial for comfort and effective pumping. Use the provided nipple measurement card to determine your size. If your nipple diameter requires a size smaller than 24mm, use the appropriate flange insert (17/19/21mm) with the 24mm flange.



Image: Nipple Measurement Card to help determine the correct flange or insert size for optimal pumping.

5. OPERATING INSTRUCTIONS

The Momcozy M6 pump features an LED display and intuitive controls for easy operation.

Control Panel:

210+
Patents Authorized

150+
Trademarks Registered

60+
Copyrights Registered

20+
Countries & Regions of Trademarks Deloyment.

USD1009251S1

Image: Diagram of the pump motor's control panel, indicating running time/suction level display, battery display, On/Off/Pause button, Mode Selection button, Increase button, Decrease button, and Milk Outlet.

Modes and Levels:

The pump offers 3 modes and 9 adjustable suction levels:

- **Stimulation Mode:** Quick and short sucking pattern to initiate milk flow.
- **Expression Mode:** Slow sucking pattern for efficient milk expression.
- **Mixed Mode:** Combines quick and short sucking with a slow and long sucking pattern.

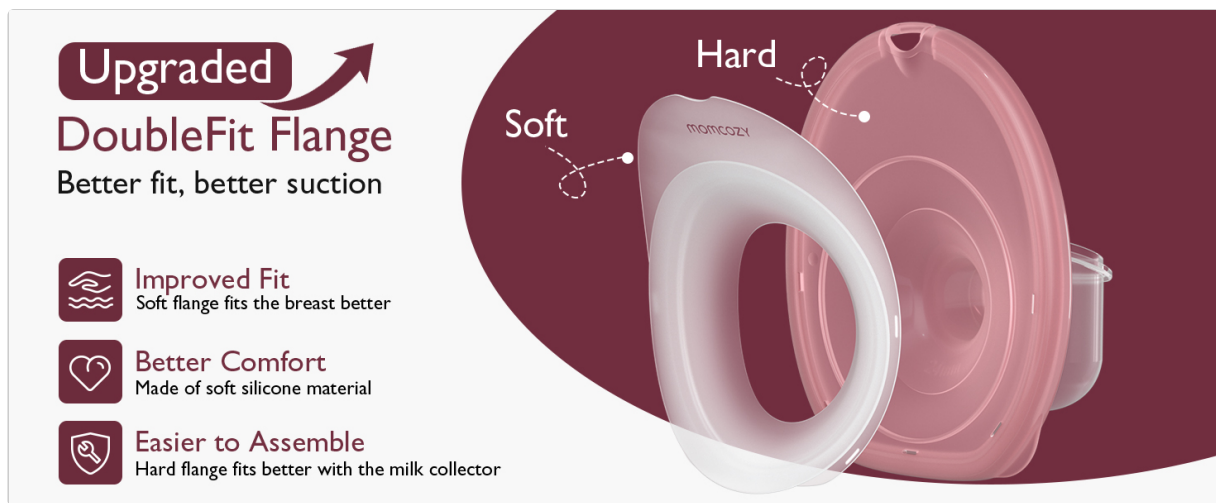


Image: Visual representation of the 3 modes (Stimulation, Expression, Mixed) and 9 suction levels of the Momcozy M6 pump.

Powering On/Off:

- Press and hold the Power button to turn the pump On or Off.
- A short press of the Power button will pause or resume pumping.

Adjusting Settings:

- Use the Mode Selection button to switch between Stimulation, Expression, and Mixed modes.
- Use the Increase (+) and Decrease (-) buttons to adjust the suction level (1-9) within each mode.
- It is recommended to start with a lower suction level and gradually increase to a comfortable setting that provides effective milk expression.

Your browser does not support the video tag.

Video: Demonstration on how to use the Momcozy M6 Mobile Style pump, including assembly and operation.

6. MAINTENANCE AND CLEANING

Proper cleaning and maintenance are essential for hygiene and optimal pump performance. The Momcozy M6 is designed with only 4 easy-to-clean parts.

Cleaning Instructions:

1. Disassemble all pump parts that come into contact with breast milk (flange, milk collector, flange cover, diaphragm, silicone valve).
2. Wash these parts with warm soapy water.
3. Rinse thoroughly with clean water.
4. Air dry completely on a clean surface or use a dedicated drying rack.
5. The pump motor should NOT be washed or submerged in water. Wipe it clean with a damp cloth if necessary.

Sterilization:

Sterilize all milk-contacting parts at least once a day. These parts are dishwasher safe and can be sterilized in a cooking pot (boiling water) or a steam sterilizer. Do NOT use a microwave for sterilization.



Image: The four main parts of the Momcozy M6 pump shown disassembled and submerged in water, highlighting their easy-to-clean design.

7. TROUBLESHOOTING

- **Low Suction:** Ensure all pump parts are correctly assembled and securely sealed. Check that the silicone valve and diaphragm are not damaged or worn.
- **No Power:** Check the battery indicator on the LED display. Recharge the pump using the provided USB Type-C cable if the battery is low.
- **Discomfort/Pain:** Verify that you are using the correct flange size. Refer to the nipple measurement card. Adjust suction levels to a comfortable setting. If discomfort persists, consult a lactation professional.
- **Milk Leakage:** Ensure all parts are correctly assembled and tightly sealed. Check for any cracks or damage to the milk collector or flange.
- **Pump Not Turning On:** Ensure the pump is fully charged. Press and hold the power button for a few seconds.

8. SPECIFICATIONS

- **Product Dimensions:** 2 x 3 x 3 inches
- **Item Weight:** 3.15 Pounds
- **Item Model Number:** BP311
- **Manufacturer:** Momcozy
- **Special Feature:** Portable
- **Number of Items:** 1 (per pump unit, product is a 2-pack)
- **Color:** Cozy Red
- **Size:** 2 Count (referring to the pack quantity)

9. WARRANTY AND SUPPORT

For warranty information, please refer to the warranty card included in your product packaging or visit the official Momcozy website. For customer support, troubleshooting assistance, or to order replacement parts, please contact Momcozy customer service through their official channels.