

Fruhdi MLSK13-B

Fruhdi Lift Chair Remote Control (MLSK13-B) Instruction Manual

Model: MLSK13-B

INTRODUCTION

This manual provides instructions for the Fruhdi MLSK13-B 4-button, 5-pin prong remote control, designed as a replacement for lift chairs and power recliners. This remote features a backlight for improved visibility during use.





Image: Fruhdi MLSK13-B remote control with illuminated buttons.

SETUP

Before installation, ensure your lift chair or power recliner is compatible with a 4-button, 5-pin round plug remote control. Disconnect power to the chair before proceeding.

1. **Identify the Connection:** Locate the existing remote control cable connection on your lift chair. It should be a 5-pin round male plug.
2. **Connect the New Remote:** Align the 5-pin round female connector of the Fruhdi MLSK13-B remote with the male plug on your chair. Gently push until it clicks securely into place. Ensure the connection is firm.



Image: Example of a 90-degree 5-pin round connector. Verify your chair's connection type for proper fit.

Note: Some chairs may have a straight plug connection, while this remote may feature a 90-degree plug. Confirm compatibility with your chair's design to ensure the plug can be securely latched.

- **Restore Power:** Once the remote is securely connected, restore power to your lift chair.

OPERATING INSTRUCTIONS

The Fruhdi MLSK13-B remote control features four buttons to operate two independent motors, typically controlling the headrest and power recline/footrest functions of your chair. The buttons are backlit for ease of use in low light conditions.

- **Headrest Adjustment:**

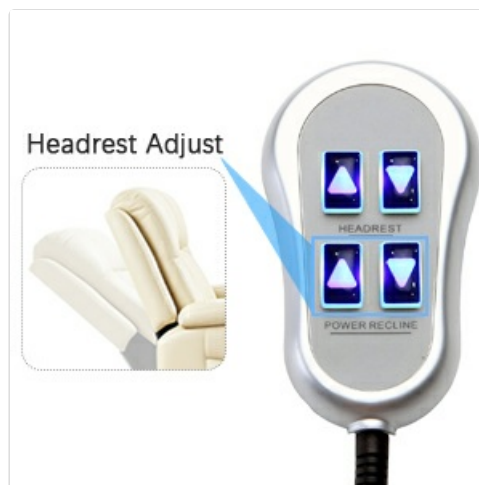


Image: Remote control indicating headrest adjustment buttons.

- Press the **UP arrow button (top-left)** under 'HEADREST' to raise the headrest.
- Press the **DOWN arrow button (top-right)** under 'HEADREST' to lower the headrest.
- **Power Recline / Footrest Adjustment:**



Image: Remote control indicating power recline/footrest adjustment buttons.

- Press the **UP arrow button (bottom-left)** under 'POWER RECLINE' to extend the footrest or recline the chair.
- Press the **DOWN arrow button (bottom-right)** under 'POWER RECLINE' to retract the footrest or return the chair to an upright position.

Hold the desired button until the chair reaches the desired position. Release the button to stop movement.

MAINTENANCE

- **Cleaning:** Wipe the remote control with a soft, dry cloth. Avoid using abrasive cleaners or solvents.
- **Liquid Exposure:** Keep the remote control away from liquids to prevent damage to internal components.
- **Cable Care:** Avoid bending or twisting the cable excessively. Ensure the cable is not pinched or caught under the chair, which could lead to damage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Remote is not responding.	Loose or disconnected cable. No power to the chair.	Ensure the 5-pin plug is securely connected to the chair's assembly. Verify the chair is plugged into a working power outlet.
Buttons are reversed or not performing expected function.	Compatibility issue with chair's wiring.	This remote is designed for specific wiring configurations. If buttons are reversed, contact customer support for assistance.
Remote plug does not stay securely connected.	Incompatible plug angle (e.g., 90-degree plug on a chair designed for a straight plug).	Ensure the plug type (straight vs. 90-degree) is compatible with your chair's connection point and allows for secure latching. If not, contact customer support.

SPECIFICATIONS

- **Brand:** Fruhdi
- **Model:** MLSK13-B
- **Number of Buttons:** 4
- **Pin Configuration:** 5-pin round plug
- **Special Feature:** Backlight
- **Compatible Devices:** Lift Chairs, Power Recliners
- **Manufacturer:** Mulin

SUPPORT

For further assistance, troubleshooting not covered in this manual, or warranty inquiries, please contact Fruhdi customer support through the retailer where the product was purchased. Provide your product model number (MLSK13-B) and a detailed description of the issue for efficient service.