

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [WEUANY](#) /

› [WEUANY Ear Wax Removal Kit Instruction Manual](#)

WEUANY W10

WEUANY Ear Wax Removal Kit Instruction Manual

MODEL: W10

Introduction

This manual provides comprehensive instructions for the safe and effective use of your WEUANY Ear Wax Removal Kit. Designed for gentle and efficient ear cleaning at home, this manual irrigation system helps remove earwax buildup, dirt, and debris without causing discomfort or damage to the sensitive skin of your ears. It is FSA/HSA eligible.

Safety Information

WARNING:

- Use this ear cleaner with **warm water only** (98.6-104°F / 37-40°C), as room temperature and cold water can cause dizziness.
- Consult your physician before use if you have a pre-existing ear condition.
- DO NOT USE if you have ear tubes (tympanostomy tubes), have had ear surgery, if you suspect/have a perforated eardrum, or have an ear infection. Please consult your ENT doctor.
- Children need to be accompanied by their parents to use this device.

Package Contents

Your WEUANY Ear Wax Removal Kit includes the following components:

- 1x Ergonomic Water Bottle (250ml capacity)
- 1x Water Outlet Button/Pump Head
- 1x High-Quality Transparent Hose with Ear Tip Connector
- 10x Single-Hole Disposable Ear Tips
- 1x Five-Hole Disposable Ear Tip
- 1x Ear Catch Basin (with measurements)
- 1x Microfiber Towel (for drying)
- Instruction Manual



Figure 1: All components included in the WEUANY Ear Wax Removal Kit.

Setup Instructions

1. **Prepare the Washing Solution:** Fill the ergonomic water bottle with lukewarm water (98.6-104°F / 37-40°C). For enhanced earwax softening, you may use a solution of 4 parts warm water to 1 part hydrogen peroxide.
2. **Assemble the Device:** Screw the pump head onto the filled water bottle. Then, attach the transparent hose to the pump head. Select either a single-hole or five-hole disposable ear tip and screw it onto the end of the hose.

Easy to use, no need to charge.



Figure 2: Assembly steps for the ear cleaning system.

Video: Product Assembly and Usage Overview

Your browser does not support the video tag.

Video 1: A comprehensive guide to assembling and using the WEUANY Ear Wax Removal Kit, demonstrating each step from filling the bottle to the ear irrigation process.

Operating Instructions

- Positioning:** Place the ear catch basin under your ear to collect draining water and earwax. You may also place the microfiber towel over your shoulder to prevent splashes.
- Irrigation:** Gently insert the ear tip into your ear canal. Press the water outlet button on the pump head to begin irrigation. Apply gentle, continuous pressure. For stubborn earwax, you may need to irrigate multiple times.
- Drainage:** After irrigation, tilt your head towards the affected ear to allow any remaining water and loosened earwax to drain out.
- Drying:** Use the microfiber towel to gently dry your ear.



No need to pay more expensive doctor bills.

This Ear Irrigation Kit removes
earwax quickly, effectively
and safely at home.



Figure 3: Proper positioning for ear irrigation using the catch basin.

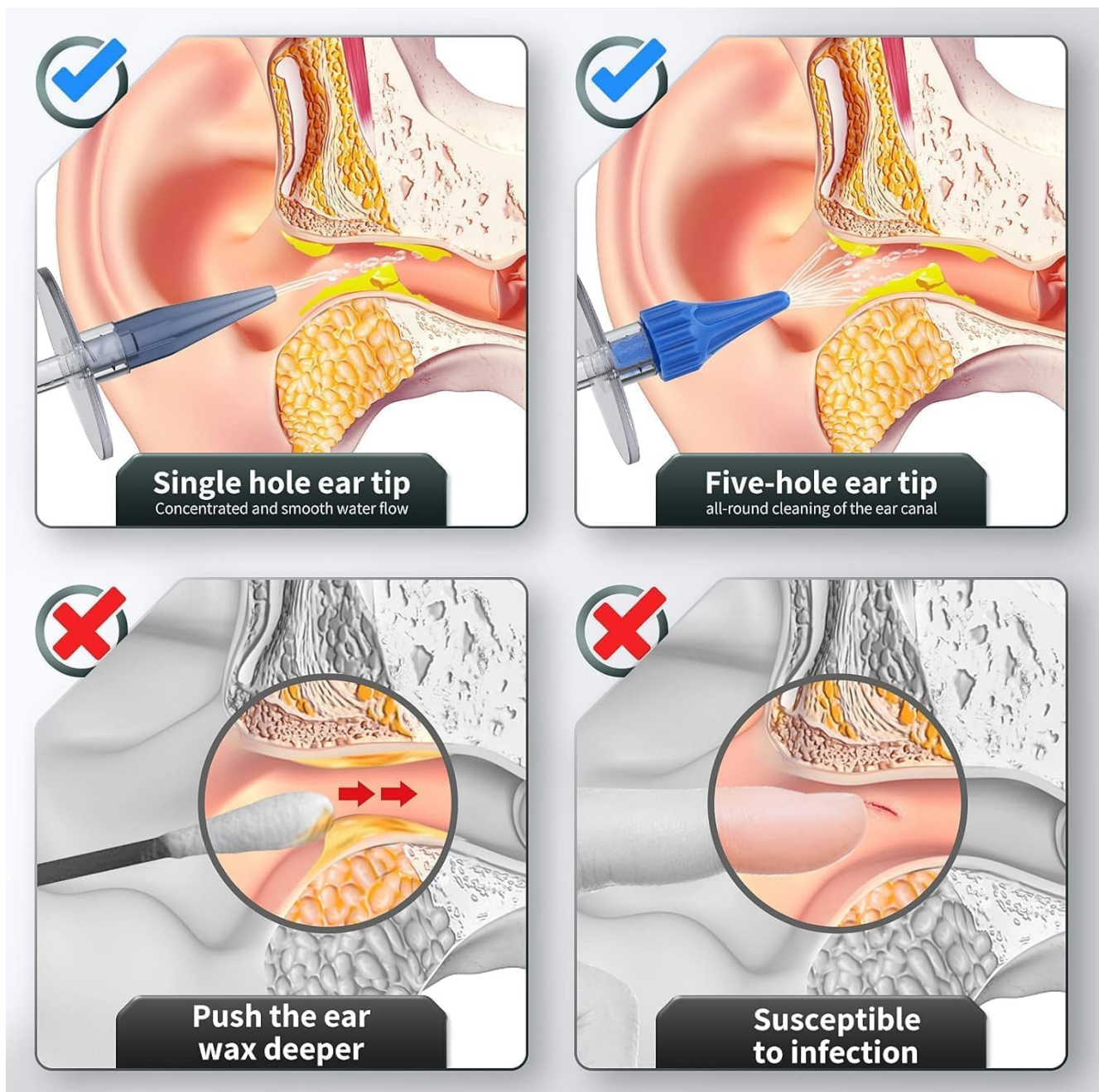


Figure 4: Different ear tip types for varied water flow and a warning against using cotton swabs.

Video: User-Friendly & Gentle Ear Wax Removal Flush Kit in 3 Steps

Your browser does not support the video tag.

Video 2: A quick demonstration of the ear wax removal process, highlighting its user-friendly and gentle nature.

Maintenance and Storage

- After each use, unscrew the pump head from the bottle and discard the used ear tip.
- Rinse all reusable components (bottle, pump head, hose, catch basin) with clean water and allow them to air dry completely before storing.
- Store the kit in a cool, dry place away from direct sunlight.

Full set ear wax removal kit.



Figure 5: The complete ear wax removal kit, ready for use or storage.

Troubleshooting

- **No water flow:** Ensure the hose is securely attached to the pump head and the ear tip is firmly screwed on. Check if the water bottle contains enough liquid.
- **Dizziness after use:** This can occur if cold water is used. Always use lukewarm water (98.6-104°F / 37-40°C) to prevent dizziness.
- **Earwax not loosening:** Pre-soften stubborn earwax with earwax softening drops or a hydrogen peroxide solution before irrigation for better results.

Product Specifications

Brand	WEUANY
Model Number	W10
Package Dimensions	6.65 x 3.19 x 2.83 inches
Item Weight	5.28 ounces

Water Bottle Capacity	250ml
------------------------------	-------

Warranty and Support

WEUANY is committed to providing high-quality products and excellent customer service. If you encounter any issues with your ear cleaning tool wax removal system, please reach out to us. We offer 100% quality product and service to our customers. For support, please visit the official WEUANY Store on Amazon or contact customer service through your purchase platform.