

[Manuals.plus](#) /

› [AFDEAL](#) /

› AFDEAL Heated Shoulder Brace with Vibration User Manual

AFDEAL Heated Shoulder Brace B0BLMYC37G

AFDEAL Heated Shoulder Brace with Vibration User Manual

Model: B0BLMYC37G

INTRODUCTION

The AFDEAL Heated Shoulder Brace with Vibration is designed to provide targeted heat therapy and soothing vibration massage to the shoulder area. This device is suitable for both left and right shoulders and features adjustable settings for personalized comfort and relief from muscle pain, joint stiffness, and for general relaxation. It is designed for wireless use, allowing for mobility during operation.

SAFETY INFORMATION

- Do not use the device while sleeping or leave it unattended.
- Avoid using on irritated, broken, or injured skin.
- Keep the device out of reach of children.
- Store in a dry and cool place when not in use.
- Ensure the device is fully charged before first use to optimize performance.

INDICATIONS FOR USE

This heated shoulder brace is indicated for temporary relief of muscle pain, joint stiffness, and for general relaxation of the shoulder area.

PACKAGE CONTENTS

- 1 x Wireless Heated Shoulder Brace

SETUP AND ADJUSTMENT

1. **Connect Strap:** Connect the main strap to the buckle on the brace.

2. **Position Brace:** Place the brace on your shoulder.
3. **Secure Strap:** Pull the long strap from under your other arm, thread it through the buckle, and pull it from the other side of the buckle to secure.
4. **Adjust Fit:** Adjust the straps to ensure a firm yet comfortable fit. The brace should not immobilize your shoulder, allowing for a specific range of motion. An extension strap is included if additional length is needed.

Heated Shoulder Brace

Easy to Use



Step 1:

Connect strap to the buckle that is on the brace. Once strap is in place, put the brace on your shoulder.



Step 2:

Stick it on your arm and make sure that it is in a comfortable and desirable position.



Step 3:

Pull the long strap from under the other arm, put it through the buckle, and now pull it from the other side of the buckle. Use the extend strap included to prolong the whole strap if it is too short for your need.



Step 4:

Secure the strap and make sure to have a firm, then press the heating and the massage button to the proper level to enjoy vibration combined with heating massage.

Figure 1: Step-by-step guide for wearing the shoulder brace.



Figure 2: Key features of the AFDEAL Heated Shoulder Brace.

OPERATING INSTRUCTIONS

Power On/Off and Mode Selection

- **Power:** Press and hold the power button on the control panel to turn the device on or off.
- **Heat Modes:** Short press the heat button to cycle through 3 heating modes. The LED indicator will show: **Red Light** for High heat, **Blue Light** for Medium heat, and **Green Light** for Low heat.
- **Vibration Modes:** Short press the vibration button to cycle through 3 vibration modes. The LED indicator will show: **Red Light** for Strong vibration, **Blue Light** for Medium vibration, and **Green Light** for Low vibration.
- You can use heat and vibration modes independently or combined.

2 EN 1 Chaleur et Massage

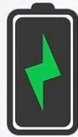


Figure 3: Heat and Vibration Mode Controls.

Battery Information

- The device features a rechargeable battery. Charge fully before first use.
- Charging time is approximately 3 hours.
- Expected battery life: up to 14 hours for massage only, 2 hours for heating only, and 2 hours for combined heat and massage.

AUTONOMIE DE BATTERIE INCROYABLE



3 HEURES
Temps de charge



14 HEURES
Massage



2 HEURES
Chauffage



2 HEURES
Les deux



Figure 4: Battery autonomy details.

MAINTENANCE

- **Cleaning:** Clean the product after each use. Use a soft, damp cloth to wipe the surface. Do not immerse the device in water.
- **Storage:** Store the device in a dry and cool place, away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No heat or vibration felt	Incorrect strap adjustment; Low battery; Incorrect mode selected.	Adjust the straps to ensure proper contact with the skin. Charge the device fully. Check the LED indicators and cycle through heat/vibration modes.
Battery drains quickly	Using high heat/vibration settings; Incomplete charging.	Reduce heat/vibration intensity. Ensure the device is fully charged before each use.

Problem	Possible Cause	Solution
Device does not turn on	Low or no battery charge.	Connect the device to a power source and charge it fully.

SPECIFICATIONS

- **Brand:** AFDEAL
- **Model:** B0BLMYC37G
- **Color:** Dark Grey
- **Material:** Neoprene
- **Special Feature:** Adjustable
- **Number of Settings:** 6 (3 Heat, 3 Vibration)
- **Heating Elements:** 2
- **Adjustable Chest Circumference:** 27-55 inches (68.6-139.7 cm)
- **Adjustable Arm Circumference:** 9-16 inches (22.9-40.6 cm)
- **Components Included:** Wireless Heated Shoulder Brace

WARRANTY INFORMATION

This product comes with a standard manufacturer's warranty. Please refer to the product packaging or contact customer support for specific warranty terms and conditions.

CUSTOMER SUPPORT

For any questions, concerns, or technical assistance, please contact AFDEAL customer support through the retailer's platform or visit the official AFDEAL brand store.



Figure 5: AFDEAL Heated Shoulder Brace in use, targeting various pain points.

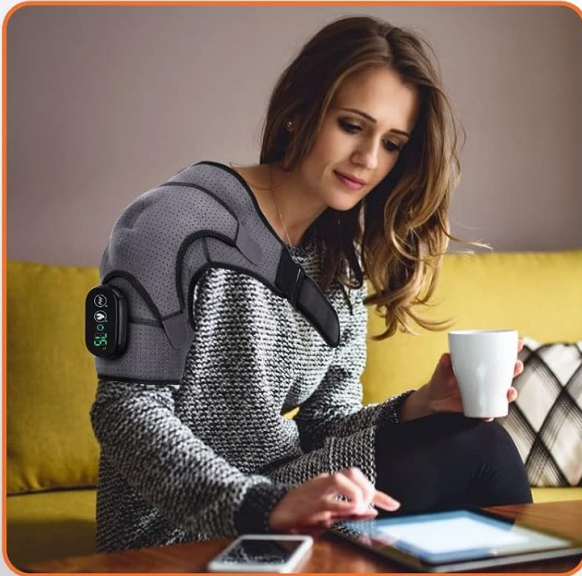


Figure 6: Versatile use of the AFDEAL Heated Shoulder Brace in daily activities.

UNISEXE ET DEUX ÉPAULES



Figure 7: Unisex design and adaptability for both shoulders.