

Metro Mobility Max Sport

Metro Mobility Max Sport Lightweight Folding Mobility Scooter User Manual

Your comprehensive guide to setting up, operating, and maintaining your Metro Mobility Max Sport Scooter.

1. INTRODUCTION

The Metro Mobility Max Sport Scooter is designed to provide convenient and reliable personal mobility. Its lightweight and folding design make it ideal for travel and everyday use. This manual provides essential information to ensure safe and efficient operation of your scooter.

Product Overview

- **2 in 1 Cup & Phone Holder:** Integrated accessory for convenience.
- **Front & Rear LED Light:** Enhances visibility and safety.
- **Rear Suspension:** Provides a smoother ride.
- **Intelligent Electromagnetic brake:** Ensures reliable stopping power.
- **Adjustable Seat:** For personalized comfort.
- **Max Drive Range:** Up to 10 miles on a full charge.
- **Max Speed:** Up to 4.97 mph.
- **Weight Capacity:** Supports up to 300 lbs.
- **Product Weight:** Approximately 101 lbs.
- **Tires:** 9"x3" solid front and rear tires.



Figure 1: Metro Mobility Max Sport Red Lightweight Folding Mobility Scooter. This image shows the scooter from a front-right angle, highlighting its compact design, front basket, and comfortable seating.

2. SAFETY INFORMATION

Always read and understand all instructions and warnings before operating your scooter. Failure to do so may result in injury or damage to the scooter.

- Do not operate the scooter under the influence of alcohol or medication.
- Always wear appropriate safety gear, if recommended by local regulations.
- Ensure the battery is fully charged before long trips.
- Avoid sudden stops or turns, especially on uneven terrain.
- Do not exceed the maximum weight capacity of 300 lbs.
- Keep hands and feet clear of moving parts.
- Consult a healthcare professional if you have any medical conditions that may affect your ability to operate the scooter safely.

3. SETUP

Your Metro Mobility Max Sport Scooter comes largely pre-assembled. Follow these steps for initial setup:

1. **Unpacking:** Carefully remove the scooter and all components from the packaging. Retain packaging for future transport or storage.
2. **Unfolding:** Place the scooter on a flat, stable surface. Release any locking mechanisms and gently unfold the scooter until it clicks into its operational position. Ensure all parts are securely locked.
3. **Seat Installation:** Insert the seat post into the receptacle on the scooter base. Adjust the seat height to a comfortable position and secure it.
4. **Handlebar Adjustment:** Adjust the tiller (handlebar column) to a comfortable height and angle. Secure the adjustment lever.

5. **Battery Charging:** Locate the charging port. Connect the charger to the scooter and then to a standard wall outlet. Allow the battery to fully charge before first use (typically 8-12 hours). The indicator light on the charger will change when charging is complete.
6. **Basket and Holder Attachment:** Attach the front basket and the 2-in-1 cup and phone holder to the tiller as desired.



Figure 2: The Metro Mobility Max Sport Scooter in a compact, folded configuration. This view demonstrates how the scooter can be collapsed for easier transport and storage, highlighting its portability.

4. OPERATING INSTRUCTIONS

Familiarize yourself with the controls before operating the scooter.

Controls Overview



Figure 3: Close-up view of the Metro Mobility Max Sport Scooter's control panel. This image displays the key ignition, battery indicator, speed adjustment dial, and horn button, which are essential for operation.

- **Key Ignition:** Insert the key and turn to the 'ON' position to power on the scooter.
- **Battery Indicator:** Displays the current battery charge level.
- **Speed Adjustment Dial:** Controls the maximum speed of the scooter. Turn clockwise for faster, counter-clockwise for slower.
- **Forward/Reverse Lever:** Push one side to move forward, the other side to move in reverse. Release to

engage the electromagnetic brake.

- **Horn Button:** Press to sound the horn.
- **Headlight Switch:** Activates the front LED light.

Driving the Scooter

1. Ensure the scooter is on a level surface and you are seated comfortably.
2. Turn the key to the 'ON' position.
3. Set the speed adjustment dial to a low setting for initial operation.
4. Gently press the forward lever to begin moving. Release the lever to stop.
5. To turn, gently steer the tiller in the desired direction.
6. When finished, turn the key to the 'OFF' position and remove it.

5. MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your scooter.

- **Battery Care:** Always charge the battery fully after each use. Do not allow the battery to completely discharge. If storing for an extended period, charge the battery at least once a month.
- **Cleaning:** Wipe down the scooter with a damp cloth. Do not use abrasive cleaners or spray water directly onto electrical components.
- **Tire Inspection:** Regularly check the solid tires for wear or damage. While solid tires do not require air pressure, ensure they are free from embedded debris.
- **Fasteners:** Periodically check all nuts, bolts, and fasteners to ensure they are tight. Do not overtighten.
- **General Inspection:** Before each use, quickly inspect the scooter for any loose parts, unusual noises, or signs of damage.

6. TROUBLESHOOTING

This section addresses common issues you might encounter.

Problem	Possible Cause	Solution
Scooter does not power on.	Battery is low or dead; Key not fully inserted or turned; Circuit breaker tripped.	Charge the battery; Ensure key is fully inserted and turned to 'ON'; Check and reset the circuit breaker (if applicable).
Scooter moves slowly or loses power.	Low battery charge; Speed dial set too low; Overload (exceeding weight capacity).	Recharge battery; Adjust speed dial to a higher setting; Reduce load on the scooter.
Brake not engaging properly.	Loose connection; Brake mechanism issue.	Ensure all connections are secure. If issue persists, contact customer support.
Unusual noises during operation.	Loose components; Debris in wheels; Motor issue.	Inspect for loose parts and tighten; Clear any debris from wheels; If noise persists, seek professional assistance.

7. SPECIFICATIONS

Detailed technical specifications for the Metro Mobility Max Sport Scooter.

Feature	Detail
Product Dimensions	44 x 20 x 37 inches
Item Model Number	Max Sport
Item Weight	101 Pounds
Manufacturer	Metro Mobility
Brand	Metro Mobility
Color	Red
Wheel Size	9 Inches
Max Drive Range	10 miles
Max Speed	4.97 mph
Weight Capacity	300 lbs

8. WARRANTY AND SUPPORT


For warranty information, please refer to the specific warranty card included with your product or visit the official Metro Mobility website. Keep your purchase receipt as proof of purchase.





Customer Support


If you require further assistance, have questions about operation, maintenance, or need to report an issue, please contact Metro Mobility customer support. Contact details can typically be found on the manufacturer's website or on your product packaging.

For the most up-to-date support information, please visit:[www.metromobility.com/support](#) (Placeholder link, refer to actual manufacturer website)

Related Documents - Max Sport

 Owner's Manual (M1)	<p>Metro Mobility M1 Owner's Manual: Safe Operation and Maintenance Guide</p> <p>Comprehensive owner's manual for the Metro Mobility M1 mobility scooter, detailing safe operation, assembly, battery charging, controls, troubleshooting, and maintenance for optimal performance and longevity.</p>
---	---

 <p>METRO MOBILITY</p> <p>USER MANUAL Model: M1 Lite</p> <p>See model at Metro mobility.com or call 1-800-828-8888 for more information. ©2018 Metro Mobility, Inc.</p>	<p>Metro Mobility M1 Lite User Manual: Operation, Safety, and Maintenance Guide</p> <p>Comprehensive user manual for the Metro Mobility M1 Lite mobility scooter, covering operation, safety guidelines, maintenance, troubleshooting, and specifications. Learn how to safely use and care for your M1 Lite.</p>
	<p>Metro Mobility M1 Lite User Manual: Operation, Safety, and Maintenance Guide</p> <p>Comprehensive user manual for the Metro Mobility M1 Lite mobility scooter, covering operation, safety guidelines, battery charging, maintenance, troubleshooting, and warranty information.</p>
 <p>METRO MOBILITY</p> <p>Owner's Manual (S800)</p>	<p>Metro Mobility S800 Owner's Manual: Operation, Maintenance, and Safety</p> <p>Comprehensive owner's manual for the Metro Mobility S800 mobility scooter, covering safety instructions, technical specifications, operation, charging, maintenance, troubleshooting, and warranty information.</p>
 <p>METRO MOBILITY</p> <p>Owner's Manual (Patriot)</p>	<p>Metro Mobility Patriot Owner's Manual: Safe Operation and Maintenance</p> <p>Comprehensive guide for the Metro Mobility Patriot mobility scooter, covering safe operation, features, battery charging, maintenance, and troubleshooting for optimal user experience and device longevity.</p>
 <p>METRO MOBILITY</p> <p>USER MANUAL Model: S500</p> <p>See model at Metro mobility.com or call 1-800-828-8888 for more information. ©2018 Metro Mobility, Inc.</p>	<p>Metro Mobility S500 User Manual: Operation, Safety, and Maintenance Guide</p> <p>Comprehensive user manual for the Metro Mobility S500 mobility scooter, covering operation, safety guidelines, battery charging, maintenance, troubleshooting, and specifications. Learn how to safely use and care for your S500 scooter.</p>

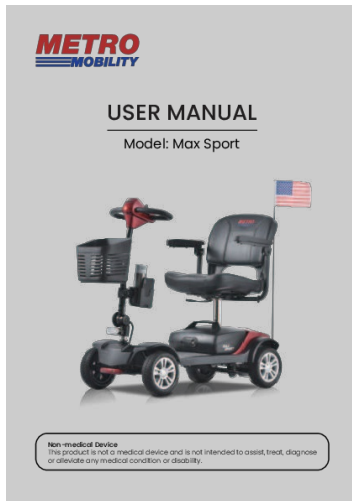

[\[pdf\]](#) User Manual

MAX SPORT Download User Manual Metro Mobility Max Sport Wheelchair Scooter USA USER MANUAL v 1722147297 cdn.shopify.com/files/1/0810/4889/5795/

USER MANUAL Model: **Max Sport** Non-medical Device This product is not a medical device and is not intended to assist, treat, diagnose or alleviate any medical condition or disability. How to Use This Manual Please read the user manual carefully before taking the product into use. This manual applies...

lang:en score:33 filesize: 4.99 M page_count: 29 document date: 2024-07-28

[\[pdf\]](#) User Manual

MAX SPORT USER MANUAL Model Max Sport Rehabmart To ensure that the scooter is OFF remove key Charging Socket The socket to connect charger located on battery pack see figure 5 and 6 max sport user manual 1 rehabmart s max 1 srsltid AfmBOoqcQpUnAw8zWEsrqucD4tPHEv9pMv76A bMyFS97mo mKWykuH Rehabmart USER Rehabmart Make sure switch in position before getting on or off To The To 1max AfmBOorXBIBL120NrYLCv 2wnW0jAUtzg5S0G5WR6 qxVyoZAc8Z9m2lrehabmart 6 This covers operations assembly maintenance fault solutions for It includes safety operation seat adjustment and max mKWykuHrehabmart AfmBOorg7BcTuyWDPa1oHMPSIm2R KNr4SqZnZIHGudBEsSDr2liLTrh and This KNr4SqZnZIHGudBEsSDr2liLTrhrehabmart AfmBOooH76HpoEiLPQlxMlnR5ProZk BSCebWDtE L0gWaKSkHsjzyld and The mobility intended comfortably transport those with walking difficulties · This designed person only How This AfmBOoq0z jo3EkfV4cm0UYkNPYsfqKxKVWd2bEOq8R

RhX9shxeQzKrehabmart How The

AfmBOoqQQLtHMg8OgU0qNU3Tt3nAkK9y5Tbr8xSvoSZrurnINPnuweXrehabmart

AfmBOooaPeiNfDrBks YXiKhITBeZQlebsL6cOvmqz x7pfRirrKogrzrehabmart

L0gWaKSkHsjzyldrehabmart AfmBOoo3m0Y2E8EZJeK6nJrYxk MMfribuFimgRe4oWGluzCQ2F009wx

MMfribuFimgRe4oWGluzCQ2F009wxrehabmart AfmBOopkmSqeF2ujmQysZ1HRhv22b

JOojJGxYv97w4pSSh8zPjdu4tG AfmBOopR EDYn8oC2WFAOFdu k8NsMliQnhMNfkcYt2QqhYPLuChA

eSrehabmart AfmBOorhk0wgsFmA31nIX3ONqEukgNSTRmwJnZO0rGyGBiYvMM1Qn VDrehabmart

AfmBOoo9TWDgullL6RmnqXq7 aA4kPZEyusJdEgPkaCsaQgUgPd8JPe0rehabmart

JOojJGxYv97w4pSSh8zPjdu4tGrehabmart

AfmBOord10XyD9EX682ez56qrltomouNjKOOAUex4aaecMNpzW ha5QE ha5QEREhabmart

AfmBOooEHaQ qZO4lx7cbq07PNCkm53tPLBk5qecqGdFnNKDZSQBhO1o

AfmBOorlcN7FsE9NZrIFKyBLP5aaFDFdkR1SJa3bNyUbmVuT4dLSOtDWrehabmart

qZO4lx7cbq07PNCkm53tPLBk5qecqGdFnNKDZSQBhO1orehabmart AfmBOooZjCVwclwa

KA3MyeXy3mWqgGHYicBW2PokqE5mwvwhcBLcn8N ||| ||| USER MANUAL Model: Max Sport Non-

medical Device This product is not a medical device and is not intended to assist, treat, diagnose or

alleviate any medical condition or disability. How to Use This Manual Please read the user manual careful

||| USER MANUAL Model: Max Sport Non-medical Device This product is not a medical device and is not

intended to assist, treat, diagnose or alleviate any medical condition or disability. How to Use This Manual

Please read the user manual careful ||| USER MANUAL Model: Max Sport Non-medical Device This

product is not a medical device and is not intended to assist, treat, diagnose or alleviate any medical

condition or disability. How to Use This Manual Please read the user manual careful ||| USER MANUAL

Model: Max Sport Non-medical Device This product is not a medical device and is not intended to assist,

treat, diagnose or alleviate any medical condition or disability. How to Use This Manual Please read the

user manual careful ||| USER MANUAL Model: Max Sport Non-medical Device This product is not a



Owner's Manual

(Max Sport)



[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide Specifications

Metro Mobility Max Sport Manual Scooter Free Shippingcdn shopify s files 1 0764 6098 2607 v

1715708252 ShippingMax 1715708252cdn ShippingLumex Liberty Half No Gap Bed Rails GF6590BH

American Quality Health Productscdn 0656 5081 9239 1757098022 ||| ||| Owner s Manual Max Sport How

to use this manual How to use this manual Your scooter will reach you in excellent condition having been personally inspected before leaving the factory. By following the guidelines for maintenance your scooter

Owner s Manual **Max Sport** How to use this manual How to use this manual Your scooter will reach you in excellent condition having been personally inspected before leaving the factory. By following the guidelines for maintenance your scooter will maintain its first class condition and give you ...

lang:en score:27 filesize: 5.53 M page_count: 58 document date: 2023-06-26

USER MANUAL

Model: Max Sport



Non-medical Device
This product is not a medical device and is not intended to assist, treat, diagnose or alleviate any medical condition or disability.

[\[pdf\]](#) User Manual

MAX SPORT USER MANUAL Model Max Sport 15 mag 2024 — Please read the user manual

carefully before taking product into use This applies to our model contains operations 29 pagine max sport 1 rehabmart s |||

USER MANUAL Model: **Max Sport** Non-medical Device This product is not a medical device and is not intended to assist, treat, diagnose or alleviate any medical condition or disability. How to Use This Manual Please read the user manual carefully before taking the product into use. cThis manual applie...

lang:en score:27 filesize: 29.8 M page_count: 29 document date: 2023-08-07



Owner's Manual

(Max Sport)



[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide Specifications

Metro Mobility Max Sport Manual Owner s American Express If you need more information on the

specifications installation etc Please check out v 1715197520 cdn shopify files 1 0565 3115 7057 |||

Owner s Manual **Max Sport** How to use this manual How to use this manual Your scooter will reach you in excellent condition having been personally inspected before leaving the factory. By following the guidelines for maintenance your scooter will maintain its first class condition and give you ...

lang:en score:25 filesize: 5.53 M page_count: 58 document date: 2023-06-26



[\[pdf\]](#) User Manual

MAX SPORT MetroMobility Max Sport Mobility 4 Wheel Scooter — MobilityActive Manual v 1712775507 cdn shopify s files 1 0685 3627 1068 |||

USER MANUAL Model: **Max Sport** Non-medical Device This product is not a medical device and is not intended to assist, treat, diagnose or alleviate any medical condition or disability. How to Use This Manual Please read the user manual carefully before taking the product into use. cThis manual applie...

lang:en score:24 filesize: 5.41 M page_count: 29 document date: 2023-09-28

METRO
MOBILITY
Owner's Manual
(Max Plus)



[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide Specifications

Metro Mobility Max Plus Manual Lumex Liberty Half No Gap Bed Rails GF6590BH 1 American Quality Health Productscdn shopify s files 0656 5081 9239 v 1757088903 ProductsMetro 4 Wheel Scooter Free Shippingcdn 1757088903cdn 0764 6098 2607 1715719607 ||| ||| Owner s Manual Max Plus How to use this manual How to use this manual Your scooter will reach you in excellent condition having been personally inspected before leaving the factory. By following the guidelines for maintenance your scooter w

Owner s Manual Max Plus How to use this manual How to use this manual Your scooter will reach ... o malfunction electrically and mechanically and may cause rusting. Maximum User Weight Limitations. **Max Sport/** Max Plus = 120Kg 264lbs . Running noisy70db. DISPOSAL The symbol above means that in acc...

lang:en score:21 filesize: 5.89 M page_count: 58 document date: 2023-06-26

METRO
MOBILITY
Owner's Manual
(Max Plus)



[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide Specifications

Metro Mobility Max Plus Manual Owner s American Express If you need more information on the specifications installation etc Please check out v 1715202580 cdn shopify files 1 0565 3115 7057 |||

Owner s Manual Max Plus How to use this manual How to use this manual Your scooter will reach ... o malfunction electrically and mechanically and may cause rusting. Maximum User Weight Limitations. **Max Sport/** Max Plus = 120Kg 264lbs . Running noisy70db. DISPOSAL The symbol above means that in acc...

lang:en score:20 filesize: 5.89 M page_count: 58 document date: 2023-06-26

METRO
MOBILITY



Owner's Manual
(Max Sport)

[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide

leoz MAX SPORT

USER S MANUAL content toolots s3 amazonaws attachment file M A |||

Owner s Manual **Max Sport** How to use this manual How to use this manual Your scooter will reach you in excellent condition having been personally inspected before leaving the factory. By following the guidelines for maintenance your scooter will maintain its first class condition and give you years...

lang:en score:12 filesize: 1.04 M page_count: 52 document date: 2023-07-17

METRO
MOBILITY



Owner's Manual
(Max Plus)

[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide

leoz Metro

Mobility 4 Wheel Foldable Scooter for Adults 300 lbs Capacity Powered Scooters Seniors Travel USB Charging Port 2 Cup Holders Red User Manual USA LLC A1hQvLudlXL m media amazon images I |||

Owner s Manual Max Plus How to use this manual How to use this manual Your scooter will reach you ... o malfunction electrically and mechanically and may cause rusting. Maximum User Weight Limitations. **Max Sport/ Max Plus** = 120Kg 264lbs . Running noisy70db. DISPOSAL The symbol above means that in acc...

lang:en score:11 filesize: 1.07 M page_count: 52 document date: 2021-11-24

METRO
MOBILITY



Owner's Manual
(Max Sport)

[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide

leoz User Manual

Metro Mobility 4 Wheel Foldable Scooter for Adults 300 lbs Capacity Powered Scooters Senior Travel Long Range Power Extended Battery with Comfortable Larger Legroom Blue Health Household A1q87DRJW L m media amazon images I |||

Owner s Manual **Max Sport** How to use this manual How to use this manual Your scooter will reach you in excellent condition having been personally inspected before leaving the factory. By following the guidelines for maintenance your scooter will maintain its first class condition and give you years...

lang:en score:11 filesize: 1.04 M page_count: 52 document date: 2022-10-26

METRO
MOBILITY



Owner's Manual
(Max Plus)

[\[pdf\]](#) User Manual Owner's Manual Quick Start Guide

leoz MAX PLUS

USER S MANUAL content toolots s3 amazonaws attachment file M A |||

Owner s Manual Max Plus How to use this manual How to use this manual Your scooter will reach you ... o malfunction electrically and mechanically and may cause rusting. Maximum User Weight Limitations. **Max Sport/** Max Plus = 136kg 300lbs . Running noisy70db. DISPOSAL The symbol above means that in acc...

lang:en score:10 filesize: 1.07 M page_count: 52 document date: 2023-07-17