



Manuals.plus /

- › Pokjet /
- › Pokjet Oticon ProWax MiniFit Hearing Aid Wax Guards Instruction Manual

Pokjet ProWax MiniFit

Pokjet Oticon ProWax MiniFit Hearing Aid Wax Guards Instruction Manual

Model: ProWax MiniFit


1. INTRODUCTION AND PRODUCT OVERVIEW

The Pokjet Oticon ProWax MiniFit Hearing Aid Wax Guards are designed to protect your hearing aid's receiver from earwax accumulation and debris. Regular replacement of these wax guards helps maintain optimal sound quality and extends the lifespan of your hearing aid by preventing internal damage. Each pack contains 6 filters, and this product includes 5 packs, totaling 30 filters.

2. COMPATIBILITY

These wax guards are compatible with most Oticon hearing aids that utilize ProWax MiniFit filters. This includes Oticon standard speakers 60, 85, and 100, as well as Oticon ProWax MiniFit receiver hearing aids such as Oticon Nera, Oticon Alta, and Alta Pro models with new MiniFit speakers.

**Fits most Oticon
Hearing Aids that
require Wax
Guards. Including
Oticon standard
speaker 60 85 100
Oticon ProWax
minifit receiver
hearing aid.**

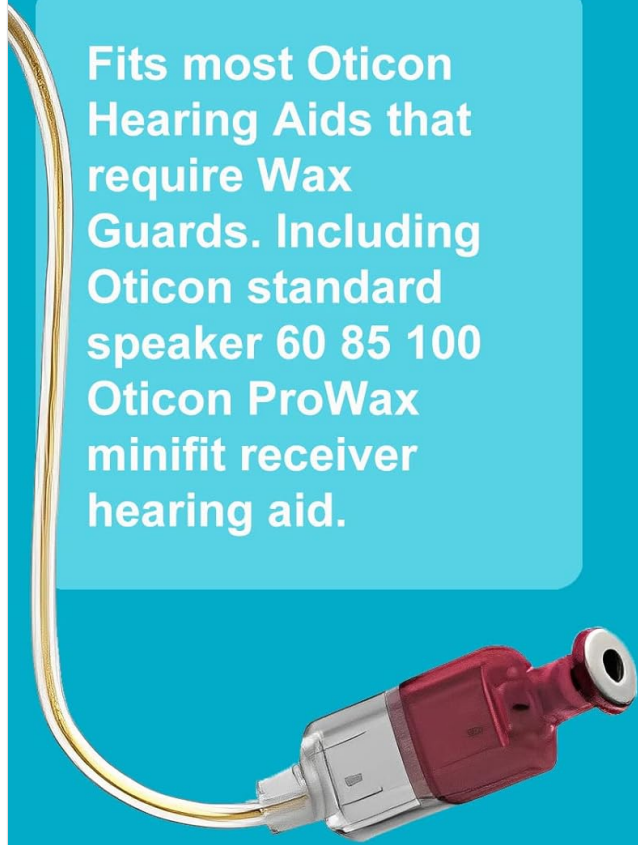

This image illustrates the compatibility of the wax guards with Oticon standard speakers 60, 85, and 100, showing the different speaker types.



This image indicates that the wax guards fit Oticon MiniFit hearing aids, including Oticon Nera, Oticon Alta, and Alta Pro models.

3. INSTALLATION AND REPLACEMENT INSTRUCTIONS

Follow these steps to replace the wax guards in your Oticon hearing aid:

1. **Prepare the Tool:** Locate the wax guard replacement tool. It is typically integrated into the wax guard dispenser.
2. **Remove the Old Wax Guard:** Insert the empty end of the replacement tool into the center of the old wax guard on your hearing aid's receiver. Gently pull to remove the used wax guard.
3. **Insert the New Wax Guard:** Rotate the tool to align the new wax guard with the receiver opening. Push the new wax guard firmly into place until it is securely seated.
4. **Discard the Used Tool:** Once the new wax guard is installed, remove the tool and dispose of it properly.



Replace ProWax miniFit filter


1. Tool



2. Remove



3. Insert

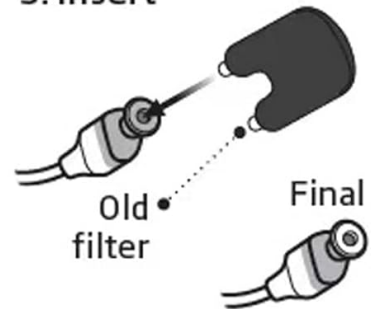

This diagram illustrates the four-step process for replacing a ProWax MiniFit filter: using the tool, removing the old filter, inserting the new filter, and discarding the tool.



Step 1: Illustration showing how to remove the old wax guard from the hearing aid receiver using the tool.



Step 2: Illustration showing the process of replacing the old wax guard with a new one from the dispenser.



Step 3: Illustration showing how to install the new wax guard securely into the hearing aid receiver.



Step 4: Illustration showing the disposal of the used wax guard and tool into a trash bin.

4. MAINTENANCE

Regularly inspect your hearing aid's wax guard for signs of blockage or discoloration. Replace the wax guard as soon as you notice any accumulation of earwax or if sound quality diminishes. This proactive maintenance ensures consistent performance and hygiene.

5. IMPORTANT SAFETY INFORMATION

- **Perforated Eardrum:** Do not use these wax guards if you have a perforated eardrum. Consult with a medical professional for appropriate hearing aid maintenance solutions in such cases.
- **Choking Hazard:** Keep wax guards and replacement tools out of reach of children and pets due to their small size.

***General Disclaimer:** Statements regarding hearing aid accessories have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. Always consult with a qualified audiologist or medical professional for advice regarding your hearing health and hearing aid maintenance.*

6. PRODUCT SPECIFICATIONS

- **Product Dimensions:** 2 x 1.4 x 2.4 inches
- **Item Weight:** 1.76 ounces
- **Manufacturer:** LIKODA TEK
- **Quantity:** 5 packs (30 pieces total)
- **Material:** High-quality rubber

7. WARRANTY AND SUPPORT

Specific warranty information for Pokjet products is not provided in this manual. For any product-related inquiries, technical support, or warranty claims, please refer to the contact information provided by your retailer or the manufacturer, Pokjet.