

Lansinoh 53351

Lansinoh Smartpump 3.0 Double Electric Breast Pump User Manual

Model: 53351

1. PRODUCT OVERVIEW

The Lansinoh Smartpump 3.0 Double Electric Breast Pump is designed to support breastfeeding mothers with efficient and comfortable milk expression. This portable and rechargeable pump connects to the Lansinoh Baby App for tracking sessions and offers hospital-strength suction with customizable settings.



Image: Lansinoh Smartpump 3.0 with bottles and app on phone, highlighting its comfortable and cordless features.

Key Features:

- **Comfortable:** Includes 3 sizes of ComfortFit flanges with soft, flexible rims for ideal fit and suction.

- **Cordless:** Rechargeable battery allows pumping anywhere without needing to be plugged into a wall.
- **Connected:** Automatically tracks pumping sessions and allows setting adjustments via the Lansinoh Baby app.
- **Customizable Settings:** Features 3 pumping styles to mimic baby's natural feeding patterns and 2-phase technology (let-down and expression).
- **Maximize Milk Flow:** Offers 8 suction levels for adjustable hospital-strength suction.

Your browser does not support the video tag.

Video: An overview of the Lansinoh Double Electric Breast Pump, demonstrating its features and ease of use.

2. WHAT'S IN THE BOX

Your Lansinoh Smartpump 3.0 package includes the following components:

- (1) Pump base unit
- (4) 5 oz Bottles
- (2) Breast flange body
- (4) Diaphragms
- (2) Diaphragm caps
- (4) Bottle storage lids
- (2) ComfortFit™ Flanges (25 mm)
- (2) ComfortFit™ Flanges (30.5 mm)
- (2) ComfortFit™ Flanges (28.5 mm)
- (2) 300 mm Connection tubes
- (1) 600 mm Connection tube
- (1) Tubing strap
- (1) Y tubing connector
- (1) Pump connector
- (4) White valve
- (1) AC Adapter/power cord
- (1) Bottle collar
- (1) Natural Wave Nipple Slow Flow Size 2
- (1) Nipple cover
- (1) Rechargeable Lithium-ion Battery
- (1) Cooler bag

US

THEM



VS



- ✓ 3 breast flange sizes included ✗
- ✓ Pump unit only weighs 1.4 lbs ✗
- ✓ 4 bottles and cooler included ✗

Image: All components of the Lansinoh Smartpump 3.0 laid out for easy identification.


Your browser does not support the video tag.

Video: A detailed look at the contents of the Lansinoh Smartpump 3.0 box, including the pump, bottles, flanges, and accessories.

3. SETUP

1. **Clean Components:** Before first use, wash all parts that come into contact with breast milk (flanges, bottles, valves, diaphragms, bottle storage lids, bottle collar, nipple) in warm, soapy water. Rinse thoroughly and air dry.
2. **Assemble Flanges and Bottles:** Attach the appropriate size ComfortFit flange to the breast flange body. Insert the white valve into the bottom of the breast flange body. Screw a 5 oz bottle onto the breast flange body.
3. **Connect Tubing:** Connect the Y-tubing connector to the pump connector. Attach the 300mm connection tubes to the breast flange bodies and the Y-tubing connector. Connect the 600mm connection tube to the Y-tubing connector and the pump base unit.
4. **Charge Pump Unit:** Connect the AC adapter/power cord to the pump base unit and plug it into a power outlet. Ensure the rechargeable Lithium-ion battery is fully charged before cordless use.

ComfortFit[®] flanges *hug* your curves

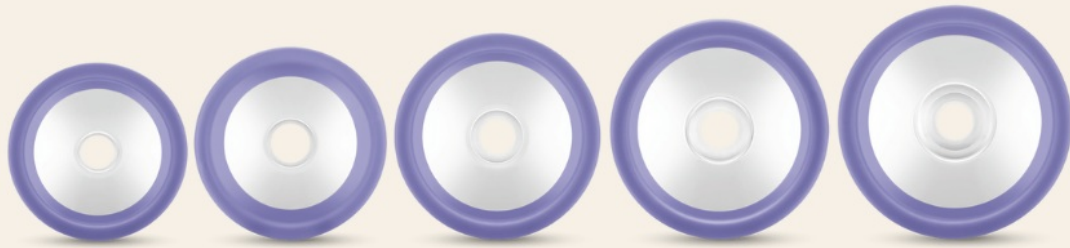

5 sizes available:
21, 25 (included) 28 (included), 30.5 (included),
and 36mm

Image: Different sizes of ComfortFit flanges, designed to hug curves for comfort and effective milk expression.

4. OPERATING INSTRUCTIONS

1. **Power On:** Press the power button on the pump base unit to turn it on. The display will illuminate.
2. **Select Pumping Style:** The pump offers 3 distinct pumping styles. Use the mode button to cycle through these styles to find the one that best mimics your baby's natural feeding pattern.
3. **Adjust Suction and Speed:** The pump features 2-phase technology (let-down and expression) and 8 adjustable suction levels. Use the '+' and '-' buttons to increase or decrease suction strength and speed to your comfort and milk flow needs.
4. **Start Pumping:** Place the assembled flanges onto your breasts, ensuring a secure seal. The pump will begin its cycle.
5. **Monitor and Track:** Connect the pump to the Lansinoh Baby app via Bluetooth to automatically track your pumping sessions, milk volume, and other feeding data.
6. **Cordless Pumping:** Utilize the rechargeable battery for pumping on the go, without the need for a power outlet.



Cordless pumping

Pump anywhere without plugging into a wall using rechargeable battery

Image: A woman using the Lansinoh Smartpump 3.0, demonstrating its cordless and portable design.

ComfortFit[®] Flanges in 3 Sizes Included



For a comfortable, secure fit and
maximum milk expression
(25mm, 28mm, & 30.5mm size flanges included)

Image: The Lansinoh Smartpump 3.0 connected to the Lansinoh Baby app, illustrating its smart tracking capabilities.

Your browser does not support the video tag.

Video: A user demonstrating the operation of the Lansinoh Smartpump 3.0, highlighting its quietness and ease of use.

5. MAINTENANCE AND CLEANING

Proper cleaning and maintenance are essential for the hygiene and longevity of your breast pump. All parts that come into contact with breast milk are BPA and BPS free and are designed for easy cleaning.

Cleaning Instructions:

- **After Each Use:** Disassemble all parts that come into contact with breast milk. Wash immediately in warm, soapy water. Rinse thoroughly with clean water.
- **Sterilization:** Sterilize components regularly by boiling for 5 minutes or using a steam sterilizer. Refer to the product manual for specific sterilization guidelines for each part.
- **Dishwasher Safe:** Many components are dishwasher safe (top rack only). Check individual part instructions for confirmation.
- **Tubing:** The connection tubes should not be washed unless milk enters them. If milk enters the tubing,

wash and air dry thoroughly. Ensure tubing is completely dry before reattaching to prevent mold growth.

- **Pump Unit:** Wipe the pump base unit with a clean, damp cloth. Do not immerse the pump unit in water.

6. TROUBLESHOOTING

- **No Suction or Weak Suction:**

- Check that all pump parts are correctly assembled and securely connected.
- Ensure the white valve is clean, undamaged, and properly seated. Replace if necessary.
- Verify that the tubing is not kinked, pinched, or damaged.
- Confirm the diaphragm is correctly placed and not torn.

- **Pump Not Turning On:**

- Ensure the pump unit is fully charged or properly connected to the AC adapter.
- Check the power outlet for functionality.

- **Milk Backflow into Tubing:**

- This indicates a potential issue with the diaphragm or valve. Inspect these parts for damage or incorrect assembly.
- Clean and dry tubing thoroughly if milk has entered it.

- **Pain or Discomfort During Pumping:**

- Ensure you are using the correct flange size. The ComfortFit flanges come in multiple sizes (25mm, 28.5mm, 30.5mm) to help achieve a comfortable fit.
- Adjust suction levels to a comfortable setting. Pumping should not be painful.
- If discomfort persists, consult a lactation consultant.

7. SPECIFICATIONS

Model Number	53351
Brand Name	Lansinoh
Operation Mode	Electric
Power Source	Battery Powered (Rechargeable Lithium-ion Battery included)
Material Type	Plastic
Material Type Free	BPA Free, Bisphenol-S (BPS) Free
Color	White
Package Dimensions	11.52 x 11.5 x 8.37 inches
Item Weight	4.98 Pounds
Number of Items	39
Manufacturer Recommended Age	1 month - 3 years
Date First Available	May 16, 2023

8. WARRANTY AND SUPPORT

Lansinoh is committed to supporting mothers throughout their breast pump journey. For any questions or warranty assistance, please reach out to Lansinoh's dedicated support team. Refer to the product packaging or the official Lansinoh website for contact information and detailed warranty terms.

Your browser does not support the video tag.

Video: A user review discussing the long-term reliability and performance of the Lansinoh Smartpump 3.0 over two years of use.