



[Manuals.plus](#) /

> [CANMOV](#) /

> CANMOV Large Power Lift Recliner Chair with Massage and Heat (Model CD011) - Instruction Manual

CANMOV CD011

CANMOV Large Power Lift Recliner Chair with Massage and Heat

MODEL: CD011

Instruction Manual

1. Introduction

Thank you for choosing the CANMOV Large Power Lift Recliner Chair. This manual provides essential information for the safe and efficient operation, assembly, and maintenance of your new recliner. Designed for comfort and ease of use, this chair features power lift and recline functions, massage and heat therapy, integrated USB charging, and convenient storage options. Please read this manual thoroughly before use and retain it for future reference.

2. Safety Instructions

- Ensure the chair is used by individuals weighing no more than 450 pounds (204.12 kg) for safe and optimal use.
- Keep children and pets away from the chair's moving parts.
- Do not allow anyone to sit on the footrest when it is extended.
- Avoid placing hands or feet near the lifting mechanism to prevent injury.
- Always ensure the power cord is not pinched or damaged.
- Unplug the chair from the power outlet before cleaning or maintenance.
- Do not operate the chair if any parts are damaged or missing. Contact customer support immediately.

3. Package Contents

Your CANMOV Power Lift Recliner Chair is packaged in several parts for safe transport. Please ensure all components are present before beginning assembly.

- Chair Backrest

- Chair Seat (with integrated armrests, side pockets, cup holders, and control panel)
- Power Adapter and Cord
- Massage/Heat Remote Control
- Lift/Recline Remote Control
- Instruction Manual

4. Setup and Assembly

This power lift chair requires minimal assembly with no tools needed. Follow these steps to set up your recliner:

1. Carefully remove all components from the packaging.
2. Align the metal grooves on the chair backrest with the corresponding brackets on the chair seat. Slide the backrest down firmly until it locks into place.
3. Connect the power cords from the backrest and seat to the main power adapter. Ensure all connections are secure.
4. Plug the power adapter into a wall outlet.

For a visual guide on assembly and features, please refer to the following video:

Video: Large Power Lift Recliner Chair with Massage and Heat - Assembly and Feature Overview. This video demonstrates the assembly process and highlights key features of the recliner.

5. Operating Instructions

5.1 Power Lift and Recline Functions

The chair features a power lift mechanism and a reclining function, controlled by a dedicated remote. This allows for smooth transitions between sitting, reclining, and standing positions.

- **To Recline:** Use the recline remote to adjust the backrest and footrest to your desired position. The chair offers multiple reclining positions for reading, watching TV, or napping.
- **To Lift:** Press and hold the lift button on the remote. The chair will gently tilt forward, assisting you to a standing position without straining your back or knees.


Product Function



Reading



Watching TV

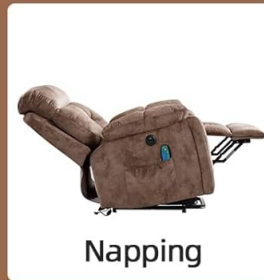

Napping



Stand Up

Image: Product Function - Illustrates the various positions of the recliner, from upright to fully reclined and the lift-assist function.

Video: CANMOV Large Power Lift Chair with Massage for Elderly - Demonstrates the power lift, recline, massage, heat, and USB charging features.

5.2 Massage and Heating Functions

The recliner is equipped with 5 massage modes and 8 massage points, along with a lumbar heating function, all controlled by a separate remote.

- **To Activate Massage/Heat:** Use the massage remote to select your desired massage mode, intensity, and specific massage points (back, lumbar, thigh, leg). Press the 'HEAT' button to activate lumbar heating.
- **Timer:** Set a timer for your massage and heat sessions using the remote.



Image: Massage & Heating - Highlights the 8 massage points and lumbar heating area on the chair, along with the massage remote.

5.3 USB Port and Cup Holders

For added convenience, the chair includes a USB charging port and hidden cup holders.

- **USB Charging:** Locate the USB port on the side of the chair. Plug in your device's charging cable to power up your electronics while you relax.
- **Cup Holders:** The armrests feature hidden cup holders. Gently push the armrest cover to reveal the cup holder.
- **Side Pocket:** A side pocket is available for storing remote controls, magazines, or other small items.



Pullable Cup Holder

Cup holders bring you an immersive cinema experience. And when you don't need it, you can hide it.

Image: Pullable Cup Holder - Shows the concealed cup holder in the armrest, designed for convenience.

6. Maintenance

To maintain the longevity and appearance of your recliner, follow these simple care instructions:

- **Cleaning:** Wipe clean with a damp cloth. Avoid harsh chemicals or abrasive cleaners.
- **Spills:** Blot spills immediately with a clean, dry cloth.
- **Regular Inspection:** Periodically check all connections and moving parts to ensure they are secure and functioning correctly.

7. Troubleshooting

If you encounter any issues with your recliner, please refer to the following common troubleshooting tips:

- **Chair Not Responding:** Ensure the power cord is securely plugged into both the chair and the wall outlet. Check if the wall outlet is functioning by plugging in another device.
- **Motor Noise:** A slight motor hum during operation is normal. If the noise is excessive or unusual, discontinue use and contact customer support.

- **Massage/Heat Not Working:** Verify the massage/heat remote is properly connected and powered on. Check the intensity settings.
- **USB Port Not Charging:** Ensure the chair is powered on and the USB cable is securely connected to both the port and your device.

If the problem persists, please contact CANMOV customer support.

8. Specifications

Feature	Detail
Brand	CANMOV
Model Name	Large Power Lift (CD011)
Color	Camel
Product Dimensions	38"D x 39"W x 43"H
Seat Dimensions	21"W x 22"D x 20.8"H
Backrest Height	26.8 Inches
Maximum Weight Recommendation	450 pounds (204.12 Kilograms)
Material	Textile (Spandex upholstery, Foam fill)
Special Features	Adjustable Backrest, Cup Holders, Ergonomic, Heated, Side Pocket, Massage
Reclining Position Count	3
Item Weight	100 pounds

9. Warranty and Support

CANMOV provides a 1-year warranty guarantee for this product. This warranty covers installation problems, damaged, and missing parts. For any purchase-related issues or support, please feel free to contact CANMOV customer service. We are committed to providing an unmatched customer experience and respond to every message within 24 hours.

Additional protection plans may be available through your retailer.