

PEPE Mobility P10020

PEPE Mobility P10020 Folding Transport Wheelchair User Manual

Model: P10020

INTRODUCTION

This manual provides comprehensive instructions for the safe and effective use of your PEPE Mobility P10020 Folding Transport Wheelchair. Please read all instructions carefully before assembly and operation to ensure user safety and product longevity. This transport wheelchair is designed to assist individuals with mobility challenges, allowing them to be transported by a capable attendant.

SAFETY INFORMATION

Your safety and the safety of the user are paramount. Always adhere to the following safety guidelines:

- Ensure the wheelchair is fully unfolded and locked into position before use.
- Always engage the parking brakes before assisting a user into or out of the wheelchair.
- The safety belt must be fastened securely around the user at all times during transport.
- Do not exceed the maximum user weight of 250 lbs (113 kg).
- Only a third party with sufficient strength should push the wheelchair.
- Regularly check all nuts, bolts, and moving parts for tightness and wear.
- Be cautious when navigating uneven terrain, inclines, or obstacles. Use the step lift lever to safely overcome small obstacles.
- Avoid leaving the parking brake engaged for extended periods when not in use, as this may affect wheel integrity.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a medical device and should be used as directed.

PRODUCT FEATURES

The PEPE Mobility P10020 Transport Wheelchair is designed for comfort, safety, and ease of use.

- **Large Wheels:** Equipped with solid polyurethane wheels featuring a wrinkled design for improved grip and reduced slipping. Front wheels can be tightened with the integrated spanner if they become loose over time.
- **Foldable Design:** Features folding armrests for easier user transfer and a quick folding system for compact storage and transport.

- **Comfort:** Padded armrests enhance user comfort. A small rear pocket provides convenient storage for personal belongings.
- **Safety Features:** Includes a locking brake system for complete immobilization and a safety belt for user support. The steel structure ensures durability and resistance.
- **Adjustable Footrests:** Removable and height-adjustable footrests accommodate various user needs.



Image: The wheelchair features flip-back armrests, allowing for easier side transfers for the user.



Image: Close-up of the padded armrests, designed for enhanced user comfort during transport.



Image: Detail of the large polyurethane wheels, engineered for improved grip and stability on various surfaces.

SETUP

Follow these steps to set up your PEPE Mobility Transport Wheelchair:

1. **Unfolding the Wheelchair:** Gently pull the armrests outwards to unfold the wheelchair until it clicks into a fully open and locked position. Ensure the seat is flat and secure.
2. **Attaching Footrests:** Align the footrest posts with the corresponding receivers on the wheelchair frame. Push down firmly until they lock into place.
3. **Adjusting Footrests:** Press the button on the footrest post and slide the footrest up or down to the desired height for the user's comfort. Release the button to lock it in place.
4. **Checking Brakes:** Verify that the parking brakes are functional by pushing the brake levers down to engage and pulling them up to release.

Your browser does not support the video tag.

Video: This video demonstrates the assembly process of the PEPE Mobility Transport Wheelchair, including unfolding and attaching components.



Image: Illustration of the swing-away and height-adjustable footrests, showing how they can be moved for user convenience.

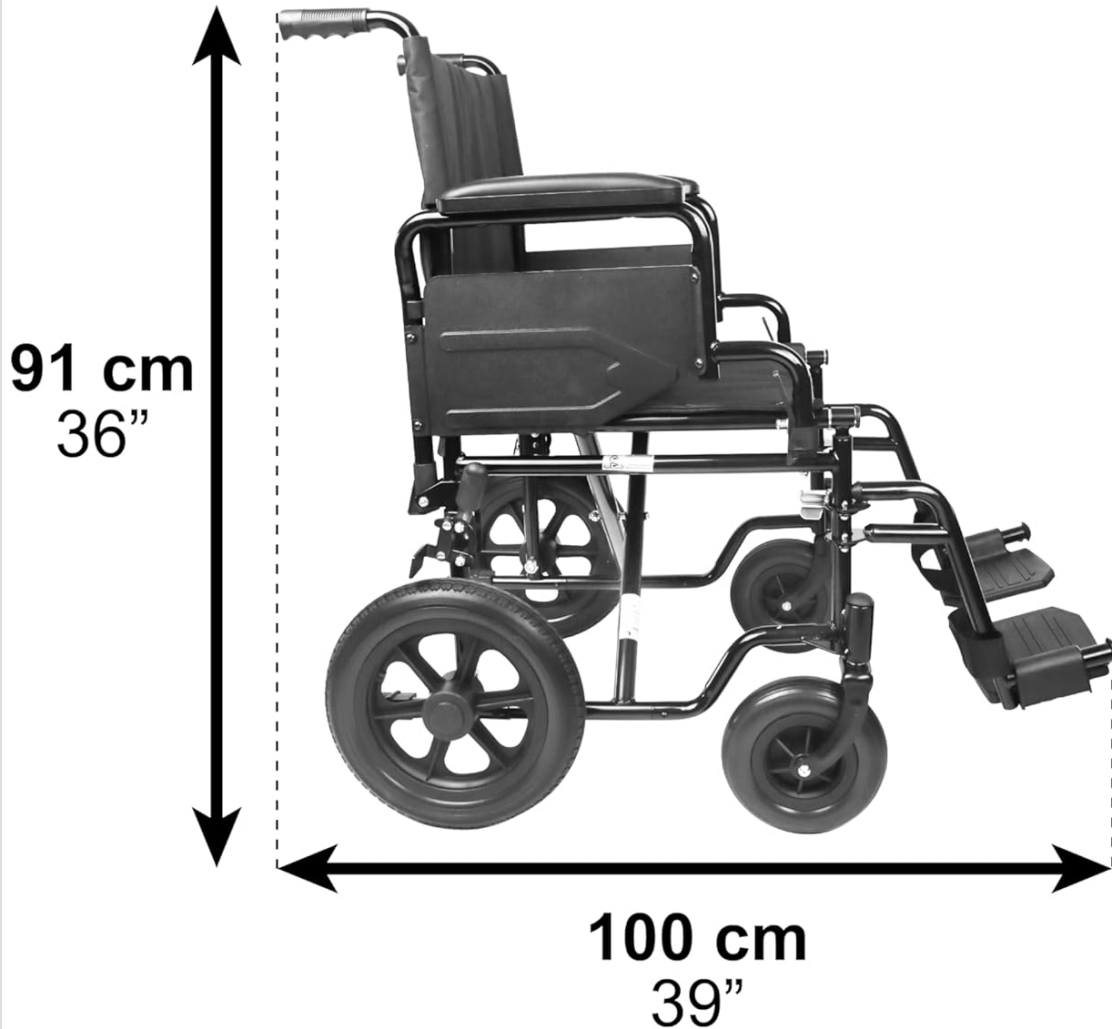


Image: The wheelchair in its folded state, demonstrating its compact design for storage and transport.

OPERATING INSTRUCTIONS

To operate the PEPE Mobility Transport Wheelchair safely:

1. **Seating the User:** Ensure the parking brakes are engaged. Assist the user to sit down, making sure they are centered and comfortable.
2. **Fasten Safety Belt:** Securely fasten the safety belt around the user's waist.
3. **Adjust Footrests:** Position the user's feet comfortably on the footrests.
4. **Releasing Brakes:** Disengage the parking brakes by pulling the levers upwards.
5. **Maneuvering:** Push the wheelchair using the handles. Use the step lift lever on the rear wheels to help overcome small obstacles like curbs or door thresholds.
6. **Stopping:** To stop, push the brake levers downwards to engage the parking brakes.
7. **Folding for Transport:** To fold, engage the brakes, lift the seat upwards, and push the sides of the frame together until it collapses. Remove footrests if desired for a more compact fold.

Your browser does not support the video tag.

Video: This video demonstrates the proper use of the PEPE Mobility Transport Wheelchair, including seating, securing, and maneuvering.



Image: An attendant pushing a user in the PEPE Mobility Transport Wheelchair, demonstrating typical use.



Image: Close-up of hands fastening the safety belt, an essential step for user security.



Image: Detail of a hand engaging the parking brake lever on the wheelchair handle.

MAINTENANCE

Proper maintenance ensures the longevity and safe operation of your wheelchair:

- **Cleaning:** Wipe down the frame and seat fabric with a damp cloth and mild soap. Avoid abrasive cleaners.

- **Wheel Inspection:** Regularly check wheels for wear and tear. Ensure the front wheels are securely tightened using the provided spanner if needed.
- **Brake Check:** Test the brakes regularly to ensure they engage and disengage properly.
- **Fastener Check:** Periodically inspect all nuts, bolts, and screws for tightness. Tighten any loose fasteners.
- **Storage:** Store the wheelchair in a dry place, away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

Here are solutions to common issues you might encounter:

- **Wheelchair does not unfold completely:** Ensure no fabric or components are snagged. Apply gentle pressure to the armrests until the frame locks into place.
- **Footrests are loose:** Re-attach them firmly, ensuring they click into place. If they wobble, check for any obstructions or damage to the locking mechanism.
- **Brakes are not holding:** Inspect the brake pads and cables. If worn or loose, they may need adjustment or replacement by a qualified technician.
- **Front wheels feel loose:** Use the integrated spanner to tighten the nuts on the front wheels as demonstrated in the assembly video.

SPECIFICATIONS

Feature	Detail
Model Number	P10020
Product Dimensions (Unfolded)	35.5"D x 24"W x 39"H (90.17 cm D x 60.96 cm W x 99.06 cm H)
Product Dimensions (Folded)	10.5" x 27.5" x 36" (approx. 26.67 cm x 69.85 cm x 91.44 cm)
Seat Width	18 inches (46 cm)
Seat Height	19 inches (49 cm)
Item Weight	33.07 Pounds (15 Kilograms)
Maximum User Weight	250 lbs (113 kg)
Frame Material	Steel
Wheel Material	Polyurethane
Color	Black

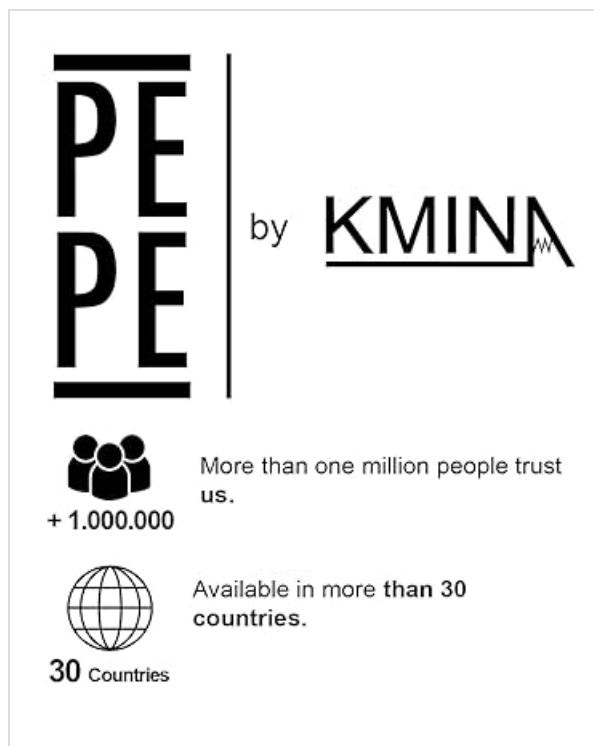


Image: Front view of the wheelchair with key dimensions including seat width (18 inches / 46 cm) and folded width (10.5 inches / 27 cm).

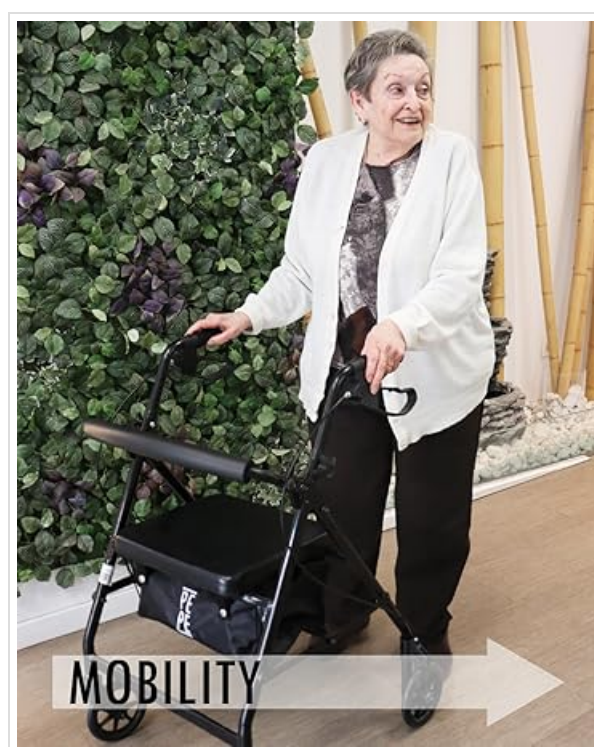


Image: Side view of the wheelchair with overall dimensions (39 inches / 100 cm length, 36 inches / 91 cm height) and maximum weight capacity (250 lbs / 113 kg).

WARRANTY AND SUPPORT

For warranty information or technical support regarding your PEPE Mobility P10020 Transport Wheelchair, please refer to the documentation included with your purchase or contact PEPE Mobility customer service. The manufacturer is IMD SL.