



[Manuals.plus](#) /

> [Vive](#) /

> [Vive Bariatric Heavy-Duty Wheelchair \(Model B0C76S5MKD\) - User Manual](#)

Vive B0C76S5MKD

Vive Bariatric Heavy-Duty Wheelchair User Manual

Model: B0C76S5MKD | Brand: Vive

[Information](#) [Setup](#) [Operation](#) [Introduction](#) [Safety](#) [Maintenance](#) [Troubleshooting](#) [Specifications](#) [Warranty & Support](#)

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your Vive Bariatric Heavy-Duty Wheelchair. Designed for strength, comfort, and convenience, this wheelchair supports individuals up to 400 pounds. It features a durable steel frame, an extra-wide 20-inch seat, padded armrests, and adjustable footrests. Please read this manual thoroughly before initial use and retain it for future reference.



Image 1: The Vive Bariatric Heavy-Duty Wheelchair, showcasing its robust design.

2. SAFETY INFORMATION

Always prioritize safety when using the wheelchair. Failure to follow these guidelines may result in injury or damage to the product.

- **Weight Capacity:** Do not exceed the maximum weight capacity of 400 pounds.
- **Brakes:** Always engage the wheel locks (brakes) when transferring in or out of the wheelchair, or when the wheelchair is stationary.
- **Stability:** Avoid leaning excessively to one side or reaching too far, as this may cause the wheelchair to tip.
- **Footrests:** Ensure footrests are properly adjusted and secured before use. Do not stand on the footrests.
- **Maneuvering:** Be aware of your surroundings, especially when navigating through doorways or tight spaces. The total width of the wheelchair is 27 inches.

- **Assistance:** If you require assistance, ensure the attendant is trained in safe wheelchair operation and transfers.
- **Regular Inspection:** Periodically check all components for wear, damage, or loose parts. Do not use the wheelchair if any part is damaged.

3. SETUP

The Vive Bariatric Heavy-Duty Wheelchair comes largely pre-assembled, requiring minimal setup before first use.

3.1 Unfolding the Wheelchair

1. Place the folded wheelchair on a flat, stable surface.
2. Grasp the top of the backrest and the front of the seat.
3. Pull the backrest and seat apart until the frame fully expands and locks into place.
4. Press down firmly on the seat rails to ensure the seat is fully open and flat.

Folds Quickly & Easily
For transports or compact storage



3.2 Attaching and Adjusting Footrests

1. Align the footrest attachment pins with the corresponding holes on the wheelchair frame.
2. Push the footrest into place until it securely locks.
3. To adjust the length, locate the adjustment lever on the footrest assembly.
4. Pull the lever, adjust the footrest to the desired height, and release the lever to lock it. Ensure both footrests are at an equal height.
5. To swing the footrests away for easier transfers, pull the release lever and pivot the footrest outwards.

Swing-Away Footrests

Adjust with the pull of a lever



Image 3: Illustration of the swing-away and adjustable footrests.

3.3 Overview Video

Your browser does not support the video tag.

Video 1: An overview of the Vive Heavy-Duty Wheelchair, demonstrating its features and ease of use.

4. OPERATION

Operating your Vive Bariatric Heavy-Duty Wheelchair is straightforward, whether you are self-propelling or being assisted.

4.1 Self-Propelling

To move the wheelchair independently, grasp the handrims on the large rear wheels and push forward or backward. Use both hands to steer by applying more force to one handrim than the other.

4.2 Attendant Assistance

For assisted movement, an attendant can push the wheelchair using the handles located at the top of the backrest. Ensure the user's feet are securely on the footrests before moving.

4.3 Engaging and Disengaging Brakes

The wheelchair is equipped with push-to-lock brakes on the rear wheels. To engage, push the brake levers forward until they lock the wheels. To disengage, pull the levers backward.

4.4 Transfers

For easier transfers, the armrests are designed to swing back. Pull the release mechanism on the armrest to swing it out of the way. This allows for closer access to beds, chairs, or other surfaces.

Swingback Armrests

Let you get into & exit your wheelchair more easily



Image 4: Swingback armrests facilitate easier transfers into and out of the wheelchair.

4.5 Portability and Storage

The wheelchair folds quickly and easily for transport or compact storage. Remove the footrests, then pull up on the center of the seat fabric to collapse the frame. The folded width is approximately 18 inches, making it suitable for car trunks.

Your browser does not support the video tag.

Video 2: Demonstration of the Vive Bariatric Wheelchair's features, including folding for transport.

5. MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your wheelchair.

- **Cleaning:** Wipe down the frame and seat with a damp cloth and mild soap. Avoid abrasive cleaners. The material is easy to clean for a hygienic experience.

- **Wheel Inspection:** Regularly check the large rear wheels and front casters for debris, wear, or damage. Ensure the tires are in good condition.
- **Brake Check:** Test the brakes regularly to ensure they engage and disengage properly and hold the wheels securely.
- **Fasteners:** Periodically check all nuts, bolts, and screws to ensure they are tight. Do not overtighten.
- **Storage:** Store the wheelchair in a dry, clean environment when not in use.

6. TROUBLESHOOTING

If you encounter any issues with your wheelchair, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Wheelchair does not roll smoothly	Debris in wheels, loose axle, worn bearings	Clean wheels, check axle for tightness, contact support if bearings are worn.
Brakes not holding	Brake pads worn, brake mechanism loose	Inspect brake pads for wear, adjust brake mechanism if loose. If worn, contact support for replacement.
Difficulty folding/unfolding	Obstruction, frame not aligned	Check for any obstructions. Ensure the seat is fully pulled up when folding, or pushed down when unfolding.
Footrests are loose or difficult to attach	Pins not aligned, release lever stuck	Ensure pins are correctly aligned with the holes. Check if the release lever is fully engaged or disengaged.

If you cannot resolve an issue using the table above, please contact Vive customer support.

7. SPECIFICATIONS

Feature	Detail
Product Dimensions (D x W x H)	42" x 27" x 36"
Item Weight	45 Pounds
Maximum Weight Recommendation	400 Pounds
Seat Width	20 Inches
Seat Height	19.5 Inches
Seat Length	16 Inches
Wheel Size	24 Inches (rear)
Turn Radius	27 Inches

Feature	Detail
Frame Material	Alloy Steel
Material Type	Heavy-Duty Steel
Brake Style	Push-to-Lock
Included Components	1 Standard Wheelchair

Heavy Duty Wheelchair for Adults

Maximum Stability

Supports up to 400 pounds

Item Weight:
45 Pounds

Seat Width: **20"**

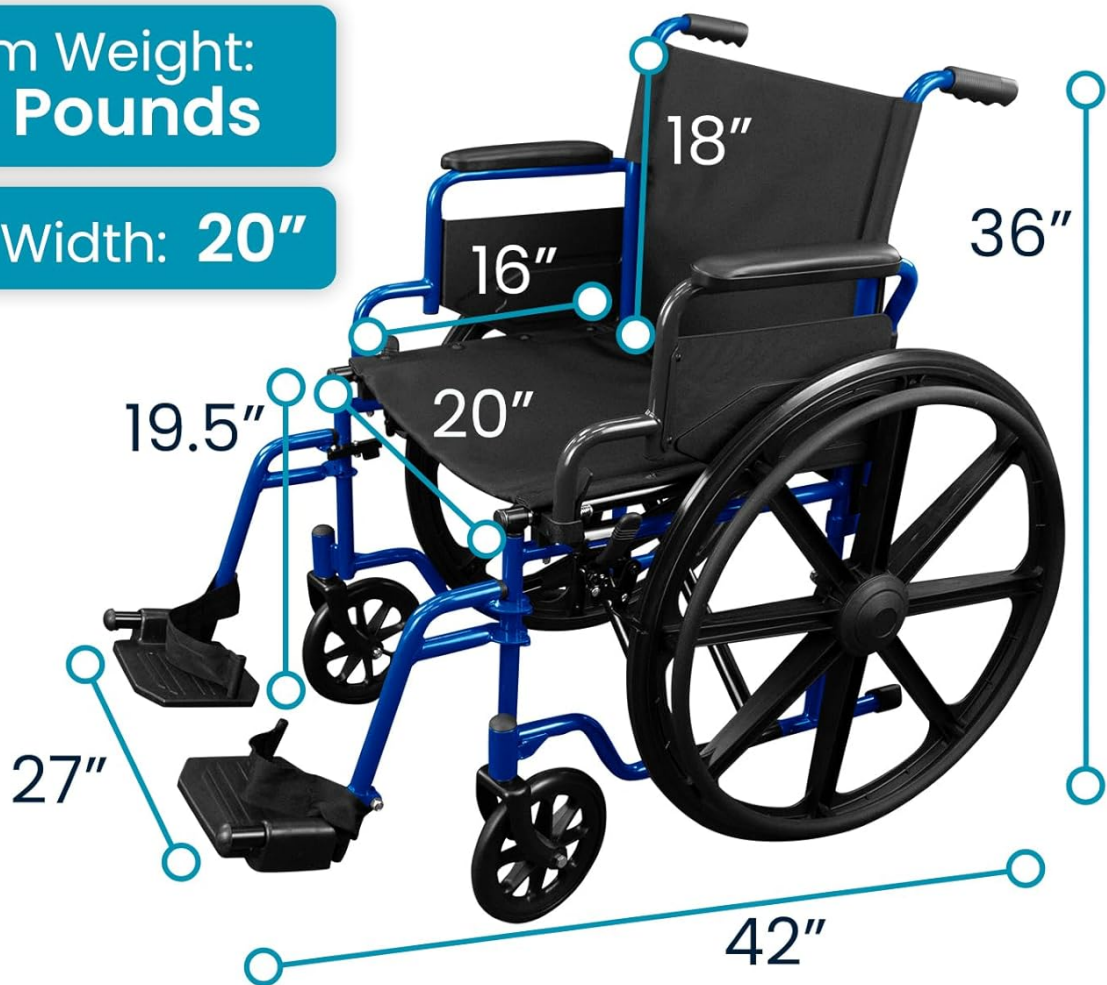


Image 5: Key dimensions of the wheelchair for reference.

8. WARRANTY AND SUPPORT

Your Vive Bariatric Heavy-Duty Wheelchair is backed by a **1-Year Guarantee**. This warranty covers manufacturing defects and ensures your confidence in the product's quality.

Additionally, the product is **FSA/HSA Approved**, making it eligible for purchase with health savings and flexible

savings accounts.

For customer support, returns, or warranty claims, please contact Vive Health Direct. They offer a 30-day free return policy for refunds or replacements.

Vive One Year Guarantee

Purchase now with confidence



Image 6: Vive Health 1-Year Guarantee for peace of mind.