



Manuals.plus /

› vasdaren /

› vasdaren Ear Wax Removal Kit with 1296P Camera - Model S8PRO User Manual

vasdaren S8PRO

vasdaren Ear Wax Removal Kit with 1296P Camera - Model S8PRO User Manual

Comprehensive instructions for safe and effective ear cleaning.

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your vasdaren Ear Wax Removal Kit, Model S8PRO. This device features a 1296P HD camera, Wi-Fi connectivity, and a range of tools designed to assist in ear hygiene. Please read this manual thoroughly before initial use.

PACKAGE CONTENTS

Verify that all items listed below are included in your package:

- Earwax Removal Kit with 1296P Camera and 6 LED lights: 1 unit
- Silicone Ear Scoops: 9 pieces
- Traditional Ear Care Tools: 8 pieces
- Earbud Storage Bottle: 1 unit
- USB Type-C Charging Cable: 1 unit
- Instruction Manual: 1 unit



Image: All components of the vasdaren Ear Wax Removal Kit, including the camera device, silicone ear scoops, traditional tools, charging cable, and storage bottle.

PRODUCT OVERVIEW

The vasdaren S8PRO Ear Wax Removal Kit is designed for detailed visual inspection and gentle removal of earwax. Key features include:

- **1296P HD Camera:** Ultra-fine lens with 360° wide-angle view and 6 cold light LED lamps for clear visualization.
- **Built-in Wi-Fi:** Connects directly to your smartphone or tablet for real-time viewing.
- **6-Axis Gyroscope:** Provides a stable and clear image during use.
- **Rechargeable Battery:** 130mAh battery, USB Type-C fast charging, up to 70 minutes working time, 30 days standby.
- **IP67 Waterproof Lens:** Allows for easy cleaning of the camera lens with water.
- **Multi-purpose Use:** Suitable for earwax removal, oral examinations, nasal and auditory canal examinations, and

preliminary ear canal checks.

- **Compatibility:** Works with both Android and iOS devices.

Ultra HD Camera 1296P

With 6 LED lights and 360 wide angle



OURS 1296P

OTHERS 720P

- 1296P
- WiFi Connection
- IP67 Waterproof
- multi-axis gyroscope
- System Compatible
- Photos and Videos

Image: Detailed view of the 1296P HD camera lens and integrated LED lights on the device tip.

SETUP INSTRUCTIONS

1. **Charge the Device:** Before first use, fully charge the ear cleaner using the provided USB Type-C cable. The charging port is located at the bottom of the device. A blue indicator light will show charging status.
2. **Download the App:** Scan the QR code provided in the instruction manual or search for "EndscopeTool" in your device's app store (compatible with Android and iOS).
3. **Power On:** Press and hold the power button on the ear cleaner for a few seconds to turn it on. A blue indicator light will illuminate.
4. **Connect via Wi-Fi:**
 - Go to your smartphone or tablet's Wi-Fi settings.
 - Locate and select the Wi-Fi network named "softish-xxx" (where xxx is a unique identifier).

- Connect to this network. No password is required.
- The built-in Wi-Fi does not support internet access or data transmission, ensuring privacy.

5. **Launch the App:** Open the "EndscopeTool" app. The device should automatically connect, and you will see a live view from the camera.



Image: A smartphone screen showing the live feed from the ear camera, demonstrating the Wi-Fi connection and app interface.

OPERATING INSTRUCTIONS

Always use the device with caution and ensure a clear view on your connected device.

1. **Attach a Silicone Scoop:** Gently attach one of the soft silicone ear scoops to the tip of the ear cleaner. Ensure it is securely fastened.
2. **Insert into Ear Canal:** Carefully insert the ear cleaner into your ear canal while observing the live feed on your smartphone. The 6 cold light LEDs will illuminate the area, and the 1296P HD camera will provide a detailed view.
3. **Navigate and Clean:** Use the live feed to guide the ear scoop. Gently scrape or scoop out earwax. Avoid pushing

wax deeper into the ear. The 6-axis gyroscope helps maintain a stable image.

4. **Capture Media:** The app allows you to take photos and record videos of your ear canal for monitoring or sharing with a healthcare professional.
5. **Other Applications:** The device can also be used for oral examinations, nasal and auditory canal examinations, and preliminary checks before a doctor's visit. Always use appropriate caution and ensure the tip is clean for each application.

Smart WIFI Connection

Connect effortlessly to your device for real-time viewing

The image shows a woman using the ear cleaner on a child's ear. The device is black with a blue tip. A smartphone is held in front of the ear, displaying a live video feed of the ear canal. The device is connected to the smartphone via WiFi. The bottom banner features three icons: a WiFi symbol labeled 'Wifi Connection', a smartphone icon labeled 'APP Smart Control', and a '1296P HD Not Blurry' label. The device itself is shown on the right side of the banner.

Image: A user demonstrating the ear cleaner in action, with the live camera feed visible on a connected smartphone.

Your browser does not support the video tag.

Video: An official demonstration of the vasdaren S8PRO ear cleaner's operation guide, showing how to use the device with the app.

MAINTENANCE

- **Cleaning the Lens:** The camera lens is IP67 waterproof. After each use, remove the silicone scoop and gently rinse the lens with water. Use the small cleaning brush provided to remove any debris. Do not immerse the entire device in water.

- **Replacing Scoops:** The silicone ear scoops are disposable. Replace them after each use or if they show signs of wear.
- **Charging:** Recharge the device when the battery indicator light is low. A full charge takes approximately 30 minutes and provides up to 70 minutes of continuous use or 30 days of standby time.
- **Storage:** Store the ear cleaner and its accessories in the provided storage bottle or a clean, dry place away from direct sunlight and extreme temperatures.

6-Axis Upgrader Gyroscope

Provide a clear and stable image



Support Pictures

Record the key parts

Switch between left and right earphone modes

Video recording function

Monitor ear health regularly

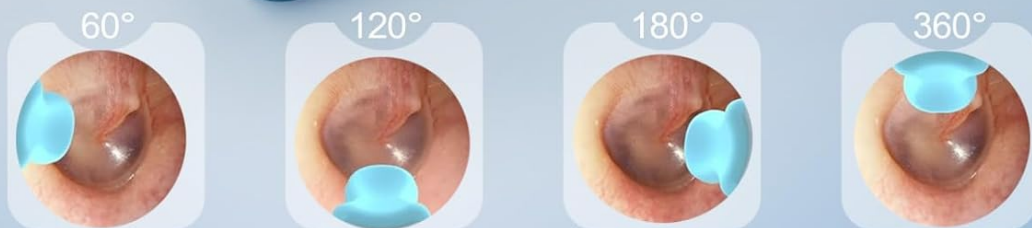


Image: The ear cleaner's waterproof lens being cleaned under running water.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully using the provided USB Type-C cable.
Cannot connect to Wi-Fi.	Incorrect Wi-Fi network selected or device not powered on.	Ensure the device is powered on and select the "softish-xxx" network in your phone's Wi-Fi settings. Restart both the device and your phone if necessary.
Blurred image.	Dirty camera lens.	Clean the camera lens with water and the provided cleaning brush.
App not displaying live feed.	App permissions not granted or Wi-Fi not connected.	Check app permissions (especially local network access) in your phone settings. Confirm Wi-Fi connection to "softish-xxx". Restart the app.

SPECIFICATIONS

- **Model Number:** S8PRO
- **Brand:** vasdaren
- **Camera Resolution:** 1296P HD
- **Field of View:** 360° wide-angle
- **Lighting:** 6 cold light LED lamps
- **Connectivity:** Built-in Wi-Fi (direct connection to device)
- **Battery Capacity:** 130mAh
- **Charging Time:** Approximately 30 minutes
- **Working Time:** Up to 70 minutes
- **Standby Time:** Up to 30 days
- **Charging Port:** USB Type-C
- **Waterproof Rating:** IP67 (lens only)
- **Compatibility:** Android and iOS devices
- **Product Dimensions:** 6.46 x 2.68 x 0.71 inches
- **Item Weight:** 2.08 ounces

WARRANTY AND SUPPORT

The vasdaren Ear Wax Removal Kit comes with a **2-year warranty**. If you encounter any issues or require assistance, please contact our customer support team. Refer to the contact information provided in the original product packaging or on the official vasdaren website.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.