

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Osim](#) /
- › [OSIM uVision 3 Eye Massager Instruction Manual](#)

Osim OS-180

OSIM uVision 3 Eye Massager Instruction Manual

Brand: Osim | Model: OS-180

1. INTRODUCTION

The OSIM uVision 3 Eye Massager is designed to provide relief from digital eye strain through a combination of vibration massage and soothing warmth. This portable device offers multiple massage programs to address various eye discomforts, promoting relaxation and revitalization.



Image 1.1: A person wearing the OSIM uVision 3 Eye Massager. This image illustrates the product in use, highlighting its design and how it covers the eye area.

2. SAFETY INFORMATION

Please read all instructions carefully before using the OSIM uVision 3 Eye Massager. Retain this manual for future reference.

- Do not use this device if you have any eye conditions, injuries, or have recently undergone eye surgery. Consult a medical professional if you have concerns.
- Keep the device away from water and moisture. Do not use it in wet environments.
- Do not disassemble or attempt to repair the device yourself. Contact authorized service personnel for assistance.
- Ensure the device is charged using the provided Micro USB cable and a compatible power source.
- Discontinue use immediately if you experience any discomfort, pain, or adverse reactions.
- This device is not intended for medical use.

3. PRODUCT OVERVIEW

The OSIM uVision 3 Eye Massager features a compact and portable design, incorporating vibration

massage and soothing warmth to alleviate eye strain. It is equipped with a one-touch button for easy operation and a see-through lens for convenience.

Key Features:

- **Vibration Massage:** Stimulates circulation around the eye area to relieve tension.
- **Soothing Warmth:** Promotes blood circulation and enhances relaxation.
- **One-Touch Button:** Simplifies power control and program selection.
- **See-Through Lens:** Allows for daily activities while using the massager.
- **Micro USB Powered:** Convenient charging and operation.
- **Light and Portable:** Weighs approximately 200g for easy transport.

See-through lens

Equipped with see-through lens, you can go on with your daily activities while relaxing your eyes.



Image 3.1: Close-up view of the see-through lens on the OSIM uVision 3 Eye Massager. This feature allows users to maintain some visibility during use.

One-touch button

You can simply switch on and toggle between the different massage programs with just one button.



Image 3.2: Detail of the one-touch button on the side of the massager, indicating its simple control interface.

4. SETUP

4.1 Charging the Device

The OSIM uVision 3 Eye Massager is powered via Micro USB. Connect the provided Micro USB cable to the device and a suitable power source (e.g., power bank, USB wall adapter, computer USB port). The device does not contain a built-in battery, ensuring continuous operation when connected to power.



Image 4.1: The Micro USB port on the OSIM uVision 3 Eye Massager, indicating its power input method.

4.2 Wearing the Eye Massager

Place the massager over your eyes, ensuring a comfortable fit. Adjust the strap at the back of your head to secure the device. The adjustable strap accommodates different head sizes.

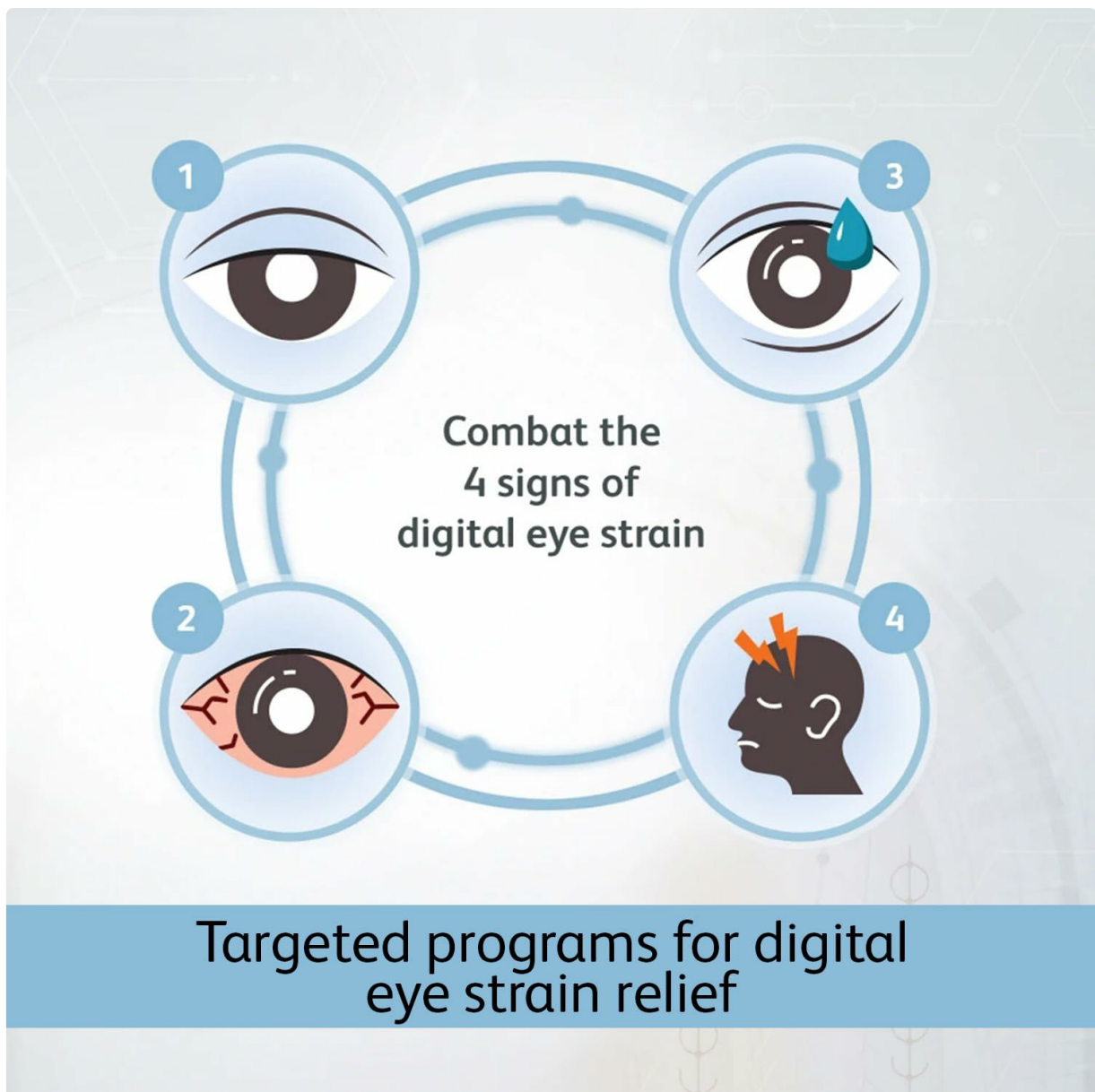


Image 4.2: A person holding the OSIM uVision 3 Eye Massager, demonstrating its compact size and the adjustable strap.

5. OPERATING INSTRUCTIONS

The OSIM uVision 3 Eye Massager features a single button for intuitive control of its massage programs.

5.1 Power On/Off and Program Selection

1. To power on the device, press and hold the one-touch button for a few seconds.
2. Once powered on, short press the button to cycle through the available massage programs.
3. To power off the device, press and hold the one-touch button again until it switches off.

5.2 Massage Programs

The massager offers four distinct programs, each combining vibration and warmth to target specific eye discomforts:

- **For Tired Eyes (Energize):** This program uses energizing vibration to stimulate circulation around the eye area, effectively helping to relieve tension and revitalize tired eyes.
- **For Sore Irritated Eyes:** The effects of vibration massage are amplified with soothing warmth to relieve discomfort and rejuvenate your eyes.
- **For Dry Eyes:** This program provides a relaxing massage with warmth. The warmth softens the


meibum (eyes' natural oil barrier to prevent tear evaporation), which can improve eye moisture levels.

- **For Headaches:** This program focuses on soothing warmth that envelops the eye area, helping to alleviate tension and pain associated with headaches.

PRODUCT SPECIFICATIONS

Product Dimension

Measurements in CM (*Note: all measurements are subjected to ±5% tolerance)

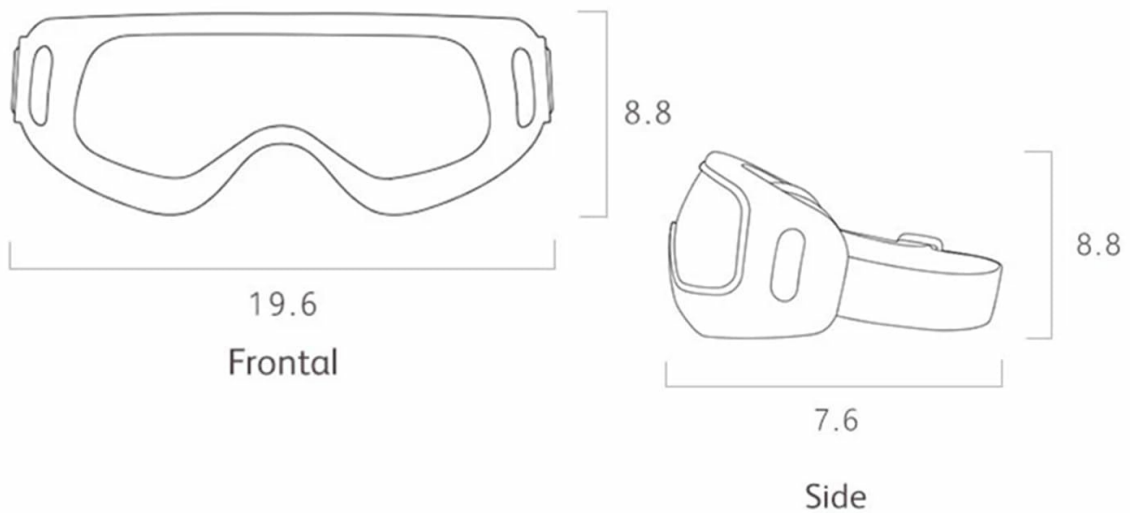

Model Specification

Model	: OSIM Eye Massager	USB Power Cord Length	: Approx. 1M
Model Number	: OS-180	Voltage	: 5V

Weight

Gross Weight	: Approx. 370g	Net Weight	: Approx. 200g
---------------------	----------------	-------------------	----------------

Image 5.1: A visual representation of the vibration massage and soothing warmth functions of the eye massager.

5.3 Product Demonstration

Video 5.1: A short preview video demonstrating the OSIM uVision 3 Eye Massager's features, including its one-touch operation, vibration, warmth, and see-through lens. This video provides a quick visual guide to the product's functionality.

Video 5.2: An informational video detailing the specific benefits of the OSIM uVision 3 Eye Massager's programs for tired eyes, sore eyes, dry eyes, and headaches. This video helps users understand which program to select for their specific needs.

6. MAINTENANCE

6.1 Cleaning

- Ensure the device is powered off and disconnected from any power source before cleaning.
- Wipe the surface of the massager with a soft, slightly damp cloth.

- Do not use abrasive cleaners, solvents, or harsh chemicals, as these may damage the device.
- Allow the device to air dry completely before storing or next use.

6.2 Storage

- Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep it out of reach of children and pets.

7. TROUBLESHOOTING

If you encounter any issues with your OSIM uVision 3 Eye Massager, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not power on.	Not connected to power, or power source issue.	Ensure the Micro USB cable is securely connected to both the device and a working power source. Try a different USB port or power adapter.
No vibration or warmth.	Incorrect program selected, or device malfunction.	Short press the one-touch button to cycle through programs. If the issue persists, power off and restart the device.
Discomfort during use.	Improper fit, or sensitivity.	Adjust the strap for a more comfortable fit. If discomfort continues, discontinue use.

If the problem persists after attempting these solutions, please contact OSIM customer support for further assistance.

8. SPECIFICATIONS

The following are the technical specifications for the OSIM uVision 3 Eye Massager:

Feature	Detail
Model Number	OS-180
Product Dimensions (LxWxH)	22 x 10 x 8 cm
Item Weight	390 g
Power Source	Corded Electric (Micro USB)
Material	Fiber
Use For	Eyes
Special Feature	Relieve Digital Eye Strain
Product Benefits	Stress Relief



Image 8.1: A diagram illustrating the frontal and side dimensions of the OSIM uVision 3 Eye Massager.

9. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the product packaging or contact OSIM directly through their official website or authorized service centers. Keep your proof of purchase for any warranty claims.