

PUMA PU0362O 001 54

PUMA PU0362O 001 54 Eyeglasses User Manual

Model: PU0362O 001 54

1. INTRODUCTION

This manual provides essential information for the proper use, care, and maintenance of your PUMA PU0362O 001 54 eyeglasses. Please read these instructions carefully before using your new eyewear to ensure optimal performance and longevity. Retain this manual for future reference.

2. PRODUCT OVERVIEW

The PUMA PU0362O 001 54 eyeglasses are designed for men, featuring a square frame shape crafted from durable acetate material. These frames are intended for prescription lenses, which must be fitted by a qualified optician.

2.1 Key Components and Dimensions

Familiarize yourself with the key components and measurements of your eyeglasses:

- **Frame Front:** The main part of the eyeglasses that holds the lenses.
- **Lenses:** (To be fitted by an optician)
- **Bridge:** The part connecting the two lens rims, resting on your nose.
- **Temples (Arms):** The side pieces that extend from the frame front over your ears.
- **Hinges:** Connect the temples to the frame front, allowing them to fold.

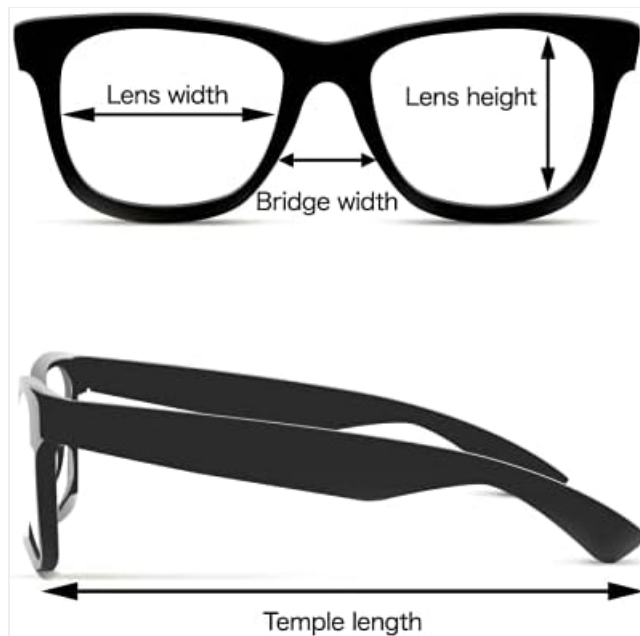


Figure 1: Eyeglasses Dimension Diagram. This image shows the key measurements of the eyeglasses: lens width, bridge width, lens height, and temple length, which are crucial for proper fit and lens ordering.

3. SETUP AND INITIAL ADJUSTMENT

Your PUMA eyeglasses frames require professional lens fitting and initial adjustment for optimal comfort and vision.

1. **Lens Fitting:** Take your new frames to a qualified optician. They will fit your prescription lenses into the frames.
2. **Initial Adjustment:** The optician will also adjust the frames to fit your face comfortably. This includes adjusting the temples for a secure fit behind your ears and ensuring the bridge rests properly on your nose.
3. **Personal Fit:** If the eyeglasses feel too tight, too loose, or cause discomfort, return to your optician for further adjustments. Do not attempt to bend or adjust the frames yourself, as this may cause damage.

4. OPERATING (WEARING INSTRUCTIONS)

Proper handling and wearing will extend the life of your eyeglasses.

- **Putting On/Taking Off:** Always use both hands to put on or take off your eyeglasses. Grasp the temples firmly near the hinges. This prevents bending the frames and loosening the screws.
- **Wearing Position:** Ensure the eyeglasses sit comfortably on your nose bridge and the temples rest gently behind your ears. The lenses should be centered in front of your eyes.
- **Avoid Headwear:** Do not push your eyeglasses onto the top of your head, as this can stretch the frames and cause them to lose their shape.

5. MAINTENANCE AND CARE

Regular cleaning and proper storage are crucial for maintaining the clarity of your lenses and the condition of your frames.

5.1 Cleaning Your Eyeglasses

1. **Rinse:** Hold your eyeglasses under lukewarm running water to rinse off any dust or debris. Avoid hot water, as it can damage lens coatings and frame materials.
2. **Clean:** Apply a small drop of mild dish soap (lotion-free) to each lens. Gently rub both sides of the lenses and all parts of the frame for a few seconds.

3. **Rinse Thoroughly:** Rinse the eyeglasses thoroughly under lukewarm running water to remove all soap residue.
4. **Dry:** Gently shake off excess water. Dry the lenses and frames with a clean, soft microfiber cloth specifically designed for eyeglasses. Avoid using paper towels, tissues, or clothing, as they can scratch the lenses.

5.2 Storage

- Always store your eyeglasses in a protective hard case when not in use.
- Avoid placing eyeglasses face down on any surface to prevent scratching the lenses.
- Keep eyeglasses away from extreme heat (e.g., car dashboards, direct sunlight) as this can warp the frame material.

6. TROUBLESHOOTING

Here are solutions to common issues you might encounter:

| Issue | Solution |
|--|---|
| Eyeglasses feel loose or slide down. | Visit your optician for a professional adjustment of the temples and nose pads. |
| Eyeglasses feel too tight or cause pressure. | Visit your optician for a professional adjustment to relieve pressure points. |
| Lenses are constantly smudged or dirty. | Follow the cleaning instructions in Section 5.1. Ensure you are using a clean microfiber cloth. Avoid touching lenses with fingers. |
| Screws on hinges are loose. | Gently tighten with a small optical screwdriver. Do not overtighten. If unsure, visit your optician. |

7. SPECIFICATIONS

Detailed technical specifications for the PUMA PU0362O 001 54 eyeglasses:

Brand: PUMA

Model: PU0362O 001 54

Department: Men

Frame Material: Acetate

Item Shape: Square

Color: Black

Bridge Width: 17 Millimeters

Arm Length (Temple Length): 140 Millimeters

Product Dimensions: 6 x 2 x 2 inches (Frame only)

Item Weight: 6 Ounces (Frame only)

UPC: 889652366234

Manufacturer: Puma

Date First Available: February 9, 2023

8. WARRANTY AND SUPPORT

PUMA is committed to quality and craftsmanship. While specific warranty details are typically provided with your purchase, PUMA products generally come with a manufacturer's warranty covering defects in materials and workmanship.

For detailed warranty information, claims, or product support, please refer to the documentation included with your eyeglasses or visit the official PUMA eyewear website. Your optician may also be able to assist with warranty-related inquiries.

Important: Unauthorized repairs or modifications may void your warranty.