

[Manuals.plus](#) /

> [Lictin](#) /

> Lictin Reusable Incontinence Pads for Men - Instruction Manual (Model ABLIUI0101)

## Lictin ABLIUI0101

# Lictin Reusable Incontinence Pads for Men

Model: ABLIUI0101

## 1. INTRODUCTION

---

The Lictin Reusable Incontinence Pads for Men are designed to provide effective management for urinary leakage. These washable guards offer a sustainable and cost-effective alternative to disposable products, reducing waste while ensuring comfort and confidence. This manual provides essential information on the product's features, proper usage, and care to ensure optimal performance.



Image 1.1: Two Lictin Reusable Incontinence Pads for Men, showcasing their ergonomic design.

## 2. PRODUCT FEATURES

- **Reusable Design:** These pads are washable and designed for multiple uses, offering an eco-conscious and economical choice.
- **3-Layer Protection:** Constructed with a 100% polyester outer layer, a soft fiber inner layer, and an absorbent core of terry cloth (100% polyester) plus TPU for enhanced reliability.
- **High Absorbency & Odor Control:** Each pad has a maximum absorption capacity of up to 140ml, providing reliable leakage protection and effective odor control for all-day freshness.
- **Comfortable & Discreet Fit:** The pads feature a cup-like, anatomical shape tailored for men's bodies, with special protective sides to ensure both comfort and discretion.



Image 2.1: Illustration detailing the three-layer protection system of the Lictin pad.

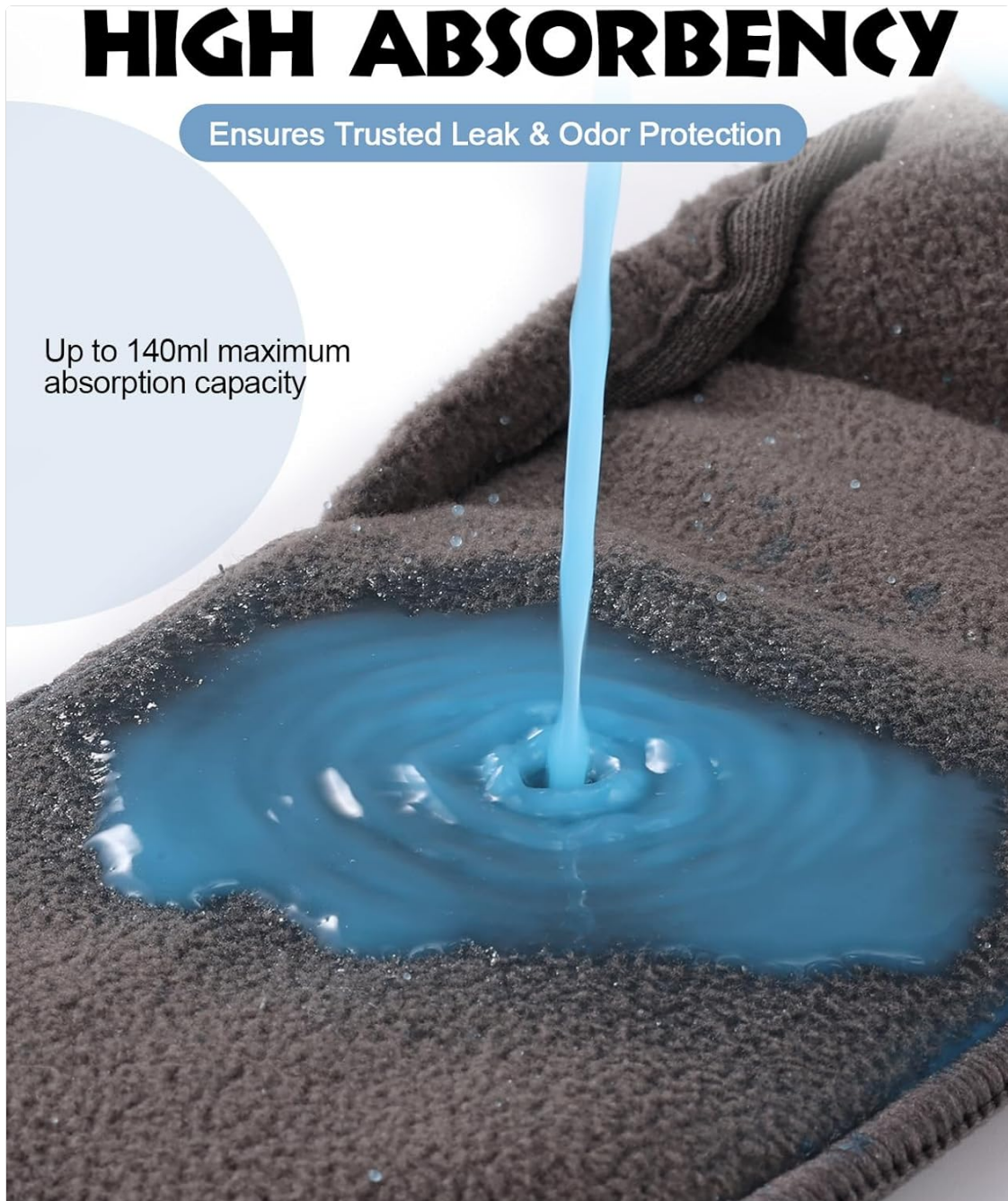


Image 2.2: A demonstration of the pad's high absorbency, showing liquid quickly being absorbed into the material.

### 3. USAGE INSTRUCTIONS

For optimal performance and comfort, follow these steps for using your Lictin Reusable

Incontinence Pads:

1. **Prepare the Pad:** Ensure the pad is clean and dry before use.
2. **Positioning:** Place the pad inside your underwear with the wider section positioned towards the front, covering the area requiring protection. The anatomical shape is designed to conform to the male body.
3. **Secure Fit:** For increased stability and to prevent shifting, it is recommended to wear men's briefs. The snug fit of briefs helps keep the pad securely in place.
4. **Check Comfort:** Adjust the pad as needed to ensure a comfortable and discreet fit.

## COMFORTABLE AND DISCREET FIT



Image 3.1: A man wearing briefs, illustrating how the Lictin pad provides a comfortable and discreet fit under clothing.

Your browser does not support the video tag.

Video 3.1: A demonstration of how to use washable protective incontinence pads for men, showing placement and fit.

## 4. MAINTENANCE

Proper care of your Lictin Reusable Incontinence Pads will extend their lifespan and maintain their effectiveness.

- **Washing:** The pads are easy to maintain and can be washed like ordinary clothes. They support both hand washing and machine washing.
- **Machine Washing:** Use a gentle cycle with cold or warm water. Avoid using bleach or fabric softeners, as these can reduce absorbency.
- **Hand Washing:** Rinse thoroughly with cold or warm water and a mild detergent.
- **Drying:** Air dry the pads or tumble dry on a low heat setting. High heat may damage the materials.

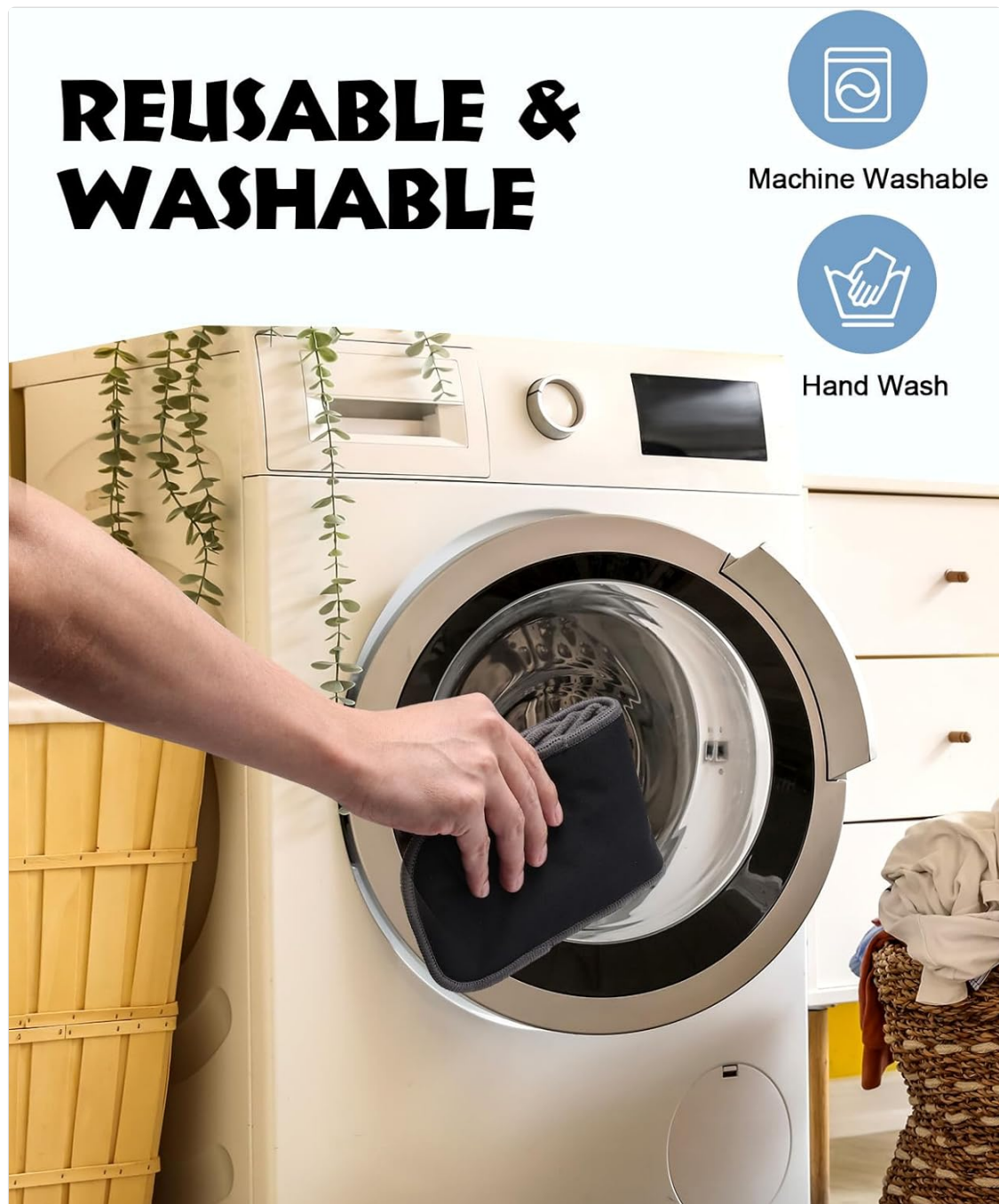


Image 4.1: A hand placing a Lictin incontinence pad into a washing machine, illustrating the ease of cleaning.

## 5. TROUBLESHOOTING

If you encounter any issues with your Lictin Reusable Incontinence Pads, consider the following:

- **Pad Shifting:** Ensure the pad is correctly positioned within your underwear. Wearing snug-fitting briefs is crucial for maintaining stability.
- **Reduced Absorbency:** Check if fabric softeners or harsh detergents were used during washing, as these can affect absorbency. Ensure pads are thoroughly dry before use.
- **Discomfort:** Adjust the pad's placement. If discomfort persists, ensure your underwear is not too tight or too loose, affecting the pad's natural fit.
- **Odor Persistence:** Ensure pads are washed promptly after use and thoroughly dried. Consider an additional rinse cycle if odors persist.

## 6. SPECIFICATIONS

Feature	Detail
Brand	Lictin
Model Number	ABLIUI0101
Material	Polyester, TPU, Terry Cloth
Dimensions (approx.)	32.5 cm x 11 cm (12.8 in x 4.3 in)
Absorbency Capacity	Up to 140ml
Reusability	Reusable
Care Instructions	Hand or Machine Washable

## 7. WARRANTY AND SUPPORT

For specific warranty information or customer support regarding your Lictin Reusable Incontinence Pads, please refer to the product packaging or contact Lictin customer service directly. Contact details are typically provided with your purchase or on the manufacturer's official website.