

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [MCombo](#) /

› [MCombo Power Lift Recliner Chair R7096 User Manual](#)

## MCombo 6160-7095Fabric

# MCombo Power Lift Recliner Chair R7096 User Manual

Model: 6160-7095Fabric

## 1. INTRODUCTION

---

This manual provides essential information for the safe and effective use of your MCombo Power Lift Recliner Chair, model R7096. Please read these instructions thoroughly before assembly and operation to ensure proper function and to prevent injury.

The MCombo Power Lift Recliner Chair is designed to provide comfort and assistance, featuring a power lift mechanism, heat, and massage functions. It is constructed with OEKO-TEX (Standard 100) certified fabric and CertiPUR-US certified sponge for safety and environmental responsibility.



Image: An elderly woman seated in the MCombo Power Lift Recliner Chair, highlighting its comfortable design and suitability for users.

## 2. SAFETY INFORMATION

---

### Important Safety Instructions:

- Do not exceed the maximum weight capacity of 350 lbs (158 kg).
- Use the lift mechanism with caution to avoid injury.
- Keep fingers, clothing, and other body parts away from moving parts of the chair, especially during recline and lift operations.
- Always follow the assembly and operation instructions provided in this manual.
- Ensure the chair is placed on a level surface.
- Keep children and pets away from the chair's moving parts.
- Unplug the chair from the power outlet before cleaning or maintenance.
- Do not operate the chair if the power cord or plug is damaged. Contact customer support for assistance.

- Avoid placing liquids directly on the chair's controls or electrical components.

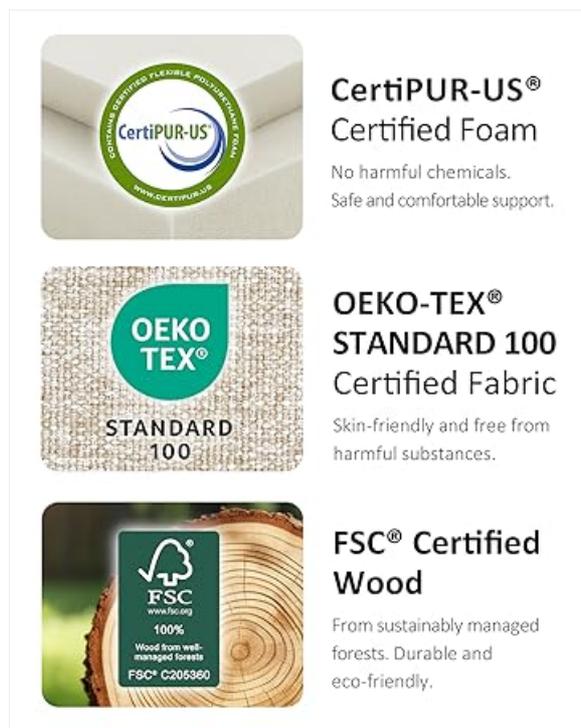


Image: Certifications for the chair's materials, including CertiPUR-US for foam, OEKO-TEX Standard 100 for fabric, and FSC for wood, ensuring safety and sustainability.

### 3. PACKAGE CONTENTS

---

Before beginning assembly, ensure all components are present and undamaged. If any parts are missing or damaged, do not proceed with assembly and contact customer support.

- Armrests (Left and Right)
- Backrest
- Base (Seat Frame with Mechanism)
- Accessories Box (contains power adapter, remote controls, anti-falling ropes)
- User Manual

### 4. ASSEMBLY INSTRUCTIONS

---

Assembly is designed to be straightforward and typically takes about 20 minutes for one person, with no additional tools required. It is recommended to have two people for moving the boxes due to weight.

1. **Unpack Components:** Carefully remove all parts from the packaging. Keep packaging materials until assembly is complete in case of returns.
2. **Attach Armrests:** Slide the left and right armrests onto the corresponding metal brackets on the chair base until they click securely into place.
3. **Connect Backrest:** Align the backrest with the receiving brackets on the chair base and push down firmly until it locks into position.
4. **Connect Wiring:** Locate the electrical cables extending from the backrest and base. Connect them according to the color-coded or labeled connectors. Ensure all connections are secure.
5. **Connect Power:** Plug the power adapter into the chair's power input, then plug the adapter into a wall outlet.
6. **Secure Anti-Falling Ropes (if applicable):** Attach the anti-falling ropes to the designated points on the chair and to a stable wall or furniture piece to prevent tipping during full recline or lift.

For visual guidance, please refer to the assembly video below:

Your browser does not support the video tag.

Video: Detailed assembly instructions for the MCombo Power Lift Recliner Chair, demonstrating how to connect the backrest, armrests, and wiring.

## 5. OPERATING INSTRUCTIONS

---

Your MCombo Power Lift Recliner Chair comes with two remote controls: one for the power lift and recline functions, and another for the heat and massage features.

### 5.1 Power Lift and Recline Functions

Use the 2-button remote control for adjusting the chair's position.


- **To Recline:** Press and hold the **down arrow** button. The chair will smoothly recline backward, and the footrest will extend. Release the button when you reach your desired position. The chair can recline up to 150 degrees.
- **To Lift/Stand:** Press and hold the **up arrow** button. The chair will gently lift and tilt forward, assisting you to a standing position. Release the button when you are ready to stand.
- **To Return to Upright Position:** From a reclined position, press and hold the **up arrow** button. The footrest will retract, and the backrest will return to an upright sitting position.



**Reclined**



**TV**

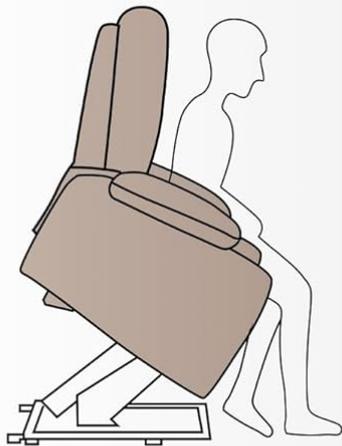

**UP**

**2 Buttons**  
Easy Touch, Easy Life



Image: The simple two-button remote control for operating the chair's lift and recline functions.



# Recline for True Relaxation

10° deeper than usual lift chairs, it can elevate basic comfort into deep, immersive relaxation you can feel.



Image: The chair fully reclined, illustrating the extended footrest and comfortable angle for relaxation.

## 5.2 Heat and Massage Functions

Use the multi-button remote control for activating and adjusting the heat and massage features.

- **Power On/Off:** Press the **POWER** button to turn the massage and heat functions on or off.
- **Heat Function:** Press the **HEAT** button to activate lumbar heating. The heating function operates independently of the vibration massage.
- **Massage Modes:** The chair features 8 vibration nodes and 5 massage modes. Press the **MODE** button to cycle through the different massage patterns.
- **Intensity Levels:** Adjust the massage intensity using the **INTENSITY +/-** buttons. There are 2 intensity levels.
- **Target Zones:** Select specific massage zones (**BACK**, **LUMBAR**, **THIGH**, **LEG**) by pressing the corresponding buttons.
- **Timer:** Set an auto-off timer for 10, 20, or 30 minutes using the **TIMER** button.

# Lumbar Heating Massage

Heating function works separately from vibration.

- 8 Vibration Nodes
- 3 Massage Modes
- 3 Intensity Levels
- Lumbar Heat



Image: The remote control for heat and massage, detailing the various buttons and an overlay showing the 8 vibration nodes and lumbar heating area on the chair.

Your browser does not support the video tag.

Video: Demonstration of the heat and massage functions of the MCombo Power Lift Recliner Chair, showing remote control usage and effects.

## 5.3 Convenient Features

- USB Charging Ports:** The chair includes a USB and a Type-C charging port, allowing you to charge your devices while seated.
- Cup Holders:** Two built-in cup holders are located on the armrests for easy access to beverages.
- Side Pockets:** Side pockets provide storage for remote controls, magazines, or other small items.

## Expect more Get more.

- 2 Anti-falling Ropes
- 2 USB & Type-C Ports
- 2 Cup Holders
- 2 Side & Front Pockets



Image: Detailed view of the chair's convenient features, including USB/Type-C ports, cup holders, and side pockets.

## 6. CARE AND MAINTENANCE

Proper care will extend the life and appearance of your recliner chair.

- **Cleaning:** Spot clean the fabric with a mild detergent and a damp cloth. Do not use harsh chemicals or abrasive cleaners. Test any cleaning solution on an inconspicuous area first.
- **Spills:** Blot spills immediately with a clean, dry cloth.
- **Dusting:** Regularly dust the chair with a soft, dry cloth.
- **Inspection:** Periodically check all connections and fasteners to ensure they are tight.
- **Power Cord:** Ensure the power cord is not pinched or placed under heavy objects.

## 7. TROUBLESHOOTING

If you encounter issues with your MCombo Power Lift Recliner Chair, refer to the following common solutions:

Problem	Possible Cause	Solution
Chair does not operate (no movement, heat, or massage)	<ul style="list-style-type: none"><li>• Power cord unplugged</li><li>• Loose electrical connections</li><li>• Power outlet not functioning</li></ul>	<ul style="list-style-type: none"><li>• Ensure power cord is securely plugged into the chair and a working wall outlet.</li><li>• Check all internal wiring connections for tightness.</li><li>• Test the outlet with another appliance.</li></ul>
Lift or recline mechanism is stuck or moves slowly	<ul style="list-style-type: none"><li>• Obstruction in mechanism</li><li>• Weight exceeding capacity</li><li>• Loose connections</li></ul>	<ul style="list-style-type: none"><li>• Inspect for any objects obstructing the moving parts.</li><li>• Ensure user weight does not exceed 350 lbs.</li><li>• Verify all electrical connections are secure.</li></ul>
Heat or massage not working	<ul style="list-style-type: none"><li>• Remote control battery (if applicable)</li><li>• Loose connections</li><li>• Function not activated</li></ul>	<ul style="list-style-type: none"><li>• Ensure the massage/heat remote is powered on.</li><li>• Check all wiring connections.</li><li>• Confirm the specific function (heat or massage mode) is selected on the remote.</li></ul>

If the problem persists after attempting these solutions, please contact MCombo customer support for further assistance.

## 8. SPECIFICATIONS

Feature	Detail
Model Name	6160-7095Fabric
Product Dimensions (D x W x H)	38.2"D x 41.7"W x 43.7"H
Seat Depth	21 inches
Seat Height	19.1 Inches
Seat Length	27.2 Inches
Chair Backrest Width	38.6 Inches
Arm Height	26 Inches
Maximum Weight Recommendation	350 Pounds (158 kg)
Material	Alloy Steel, Wood frame; Foam fill; Fabric upholstery
Special Features	Adjustable Backrest, Heated, Power Lift, Power Recliner, Wide Seat
Reclining Position Count	3 (Upright, TV, Reclined)
Recommended Height	5'1"-5'9" (for Medium-Wide size)

Ideal Height : 5'1"-5'9" Max Capacity : 350 lbs



Image: Dimensional diagram of the chair, illustrating measurements for height, width, depth, seat dimensions, and recline angles, along with recommended user height and weight capacity.

## 9. WARRANTY AND SUPPORT

Your MCombo Power Lift Recliner Chair is backed by a **2-year limited warranty**. This warranty covers free replacements for:

- Reclining/Lifting Mechanisms
- Electrical Systems
- Upholstery & Cushions
- Wood/Metal Frame Structures

For warranty claims, technical support, or any questions regarding your product, please contact MCombo customer service. Refer to the contact information provided with your purchase documentation or visit the official MCombo website.



