

ERYITRDK B0CKW6TC9L

ERYITRDK Ear Wax Removal Tool with Camera - Model B0CKW6TC9L - User Manual

Comprehensive instructions for safe and effective use.

1. INTRODUCTION

This manual provides detailed instructions for the ERYITRDK Ear Wax Removal Tool with Camera. This device is designed to assist in visual ear cleaning by providing a clear view of the ear canal directly on your smartphone or tablet. It features a high-definition camera, integrated LED lights, and a comfortable design for effective ear wax removal. Please read this manual thoroughly before first use to ensure proper operation and safety.



Figure 1.1: ERYITRDK Ear Wax Removal Tool with Camera and included accessories. This image displays the main device, various silicone covers, and charging cable, alongside a smartphone screen showing the ear canal view.

2. PRODUCT FEATURES

- **High-Definition Camera:** Provides a clear, detailed view of the ear canal for precise cleaning.
- **360° Visual Ear Cleaning:** Built-in gyroscope supports stable, omnidirectional viewing.
- **Wireless Connectivity:** Connects to smartphones or tablets via Wi-Fi hotspot for real-time viewing.
- **Integrated LED Lights:** Six complementary LED lights illuminate the ear canal for enhanced visibility.
- **Soft Silicone Covers:** Includes six soft silicone covers for comfortable and safe use, suitable for various users.
- **Intelligent Temperature Control:** Maintains an optimal device temperature during use to prevent discomfort.
- **Long Battery Life:** One-hour charge provides up to four hours of continuous operation.

1080P HD EAR CAMERA

Full HD resolution can provide superior clear image.



Figure 2.1: Illustration of the 1080P HD Ear Camera's capabilities, highlighting 360° TakeWide view, WiFi wireless connection, and photo capture functions. It also compares the clarity of 1080P resolution to 720P.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- ERYITRDK Ear Wax Removal Tool with Camera
- 6 x Soft Silicone Ear Spoon Covers
- USB Charging Cable (Type-C)
- User Manual (this document)

WIRELESS EAR CLEANER

Help you see clearly inside your ear



Figure 3.1: Detailed view of the wireless ear cleaner components, including the 6 LED chip light, 200W miniature camera, Type-C Charging Port, Power Button, and Silicone Ear Spoons.

4. SETUP INSTRUCTIONS

Follow these steps to set up your ERYITRDK Ear Wax Removal Tool:

1. **Charge the Device:** Before first use, fully charge the device using the provided USB Type-C cable. The charging port is located on the device body. A full charge takes approximately one hour.
2. **Download the Application:** Scan the QR code provided in the product packaging or search for the recommended application (e.g., "360TOOL Smart Ear Collector" or similar as indicated in the product instructions) on your smartphone's app store (iOS or Android). Install the application and grant necessary permissions.
3. **Power On the Device:** Press and hold the power button on the ear wax removal tool until the indicator light illuminates.
4. **Connect to Wi-Fi:**
 - Open the downloaded application on your smartphone.
 - In the app, locate the option to connect a new device (e.g., "Device AP Direct Connection").
 - Navigate to your phone's Wi-Fi settings.

- Find and select the Wi-Fi hotspot signal broadcast by the device. This signal typically begins with "Cooler-XXXX" or a similar identifier.
 - Connect to this Wi-Fi network. No password is usually required for direct device hotspots.
5. **Start Viewing:** Return to the application. Select the connected device (e.g., "360TOOL Smart Ear Collector") to access the live camera feed from the ear wax removal tool.
 6. **Attach Silicone Cover:** Gently attach one of the soft silicone ear spoon covers to the tip of the device. Ensure it is securely fitted.

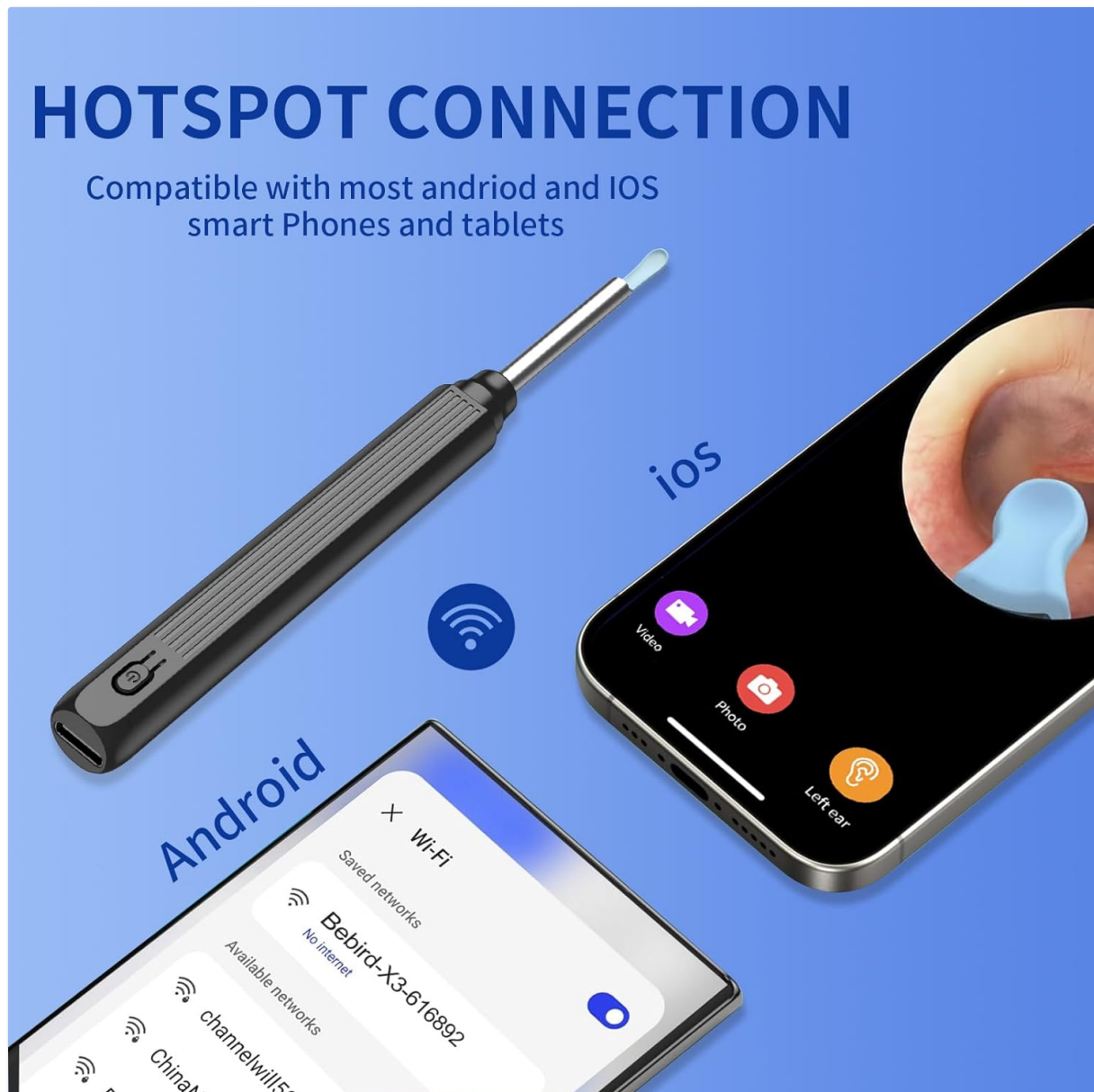


Figure 4.1: Visual guide demonstrating the hotspot connection process for the ear camera with both iOS and Android devices. It shows the device broadcasting a Wi-Fi signal and a smartphone connecting to it.

5. OPERATING INSTRUCTIONS

Once the device is set up and connected to your smartphone, you can begin the ear cleaning process.

1. **Positioning:** Gently insert the tip of the device into your ear canal while observing the live feed on your smartphone. Do not force the device into the ear.
2. **Visual Guidance:** Use the clear, illuminated view from the camera to locate earwax. The 360° gyroscope ensures a stable image regardless of the device's orientation.

- 3. **Gentle Removal:** Carefully use the silicone ear spoon to scoop or dislodge earwax. Always move slowly and deliberately. Avoid pushing wax deeper into the ear canal.
- 4. **Monitor Temperature:** The device features intelligent temperature control. If you feel any discomfort or excessive heat, immediately remove the device from your ear.
- 5. **Photography/Video (Optional):** The application may allow you to capture photos or record videos of your ear canal for personal reference or to share with a healthcare professional.
- 6. **Completion:** After cleaning, gently remove the device from your ear.

Important: Always exercise extreme caution when inserting any object into the ear canal. If you experience pain, discomfort, or hearing changes, discontinue use immediately and consult a medical professional.

6. MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your device.

- **Cleaning the Tip:** After each use, remove the silicone ear spoon cover. Clean the cover and the device tip with a soft cloth dampened with rubbing alcohol or an antiseptic wipe. Ensure no debris remains on the camera lens or LED lights.
- **Replacing Covers:** Replace silicone covers as needed, especially if they show signs of wear or damage.
- **Charging:** Recharge the device when the battery is low. Avoid overcharging for extended periods.
- **Storage:** Store the device in a clean, dry place, away from direct sunlight and extreme temperatures. Keep it out of reach of children.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not power on.	Low battery.	Charge the device for at least one hour using the provided USB cable.
Cannot connect to Wi-Fi hotspot.	Device not powered on; incorrect Wi-Fi selection; app permissions not granted.	Ensure the device is powered on. Verify you are connecting to the correct "Cooler-XXXX" Wi-Fi network. Check app permissions in your phone settings. Restart the device and the app.
No image or blurry image in the app.	Poor Wi-Fi connection; camera lens dirty; app not fully launched.	Ensure a stable Wi-Fi connection. Clean the camera lens gently. Close and reopen the application.
Device becomes uncomfortably warm.	Extended use; environmental factors.	Discontinue use immediately and allow the device to cool down. The device has intelligent temperature control, but prolonged use can still generate heat.
Silicone ear spoon broke.	Wear and tear; improper handling.	Replace with a new silicone cover. Handle covers gently during attachment and removal.

8. SPECIFICATIONS

- **Model:** B0CKW6TC9L
- **Brand:** ERYITRDK
- **Power Source:** 1 Lithium Polymer battery (included)
- **Charging Time:** Approximately 1 hour
- **Operating Time:** Up to 4 hours (continuous use)
- **Connectivity:** Wi-Fi Hotspot (compatible with iOS and Android)
- **Camera Resolution:** High Definition (e.g., 1080P as per images)
- **Dimensions:** 6.57 x 2.01 x 0.67 inches
- **Weight:** 1.13 ounces
- **Manufacturer:** ERYITRDK



Figure 8.1: Diagram illustrating the internal "Powerful Smart Chip" with features such as smart temperature control, 20 FPS transmission rate, and a 3-axis gyroscope for stable operation.

9. SAFETY INFORMATION

Please adhere to the following safety guidelines to prevent injury or damage to the device:

- Do not insert the device too deeply into the ear canal.
- Always use the soft silicone ear spoon covers. Do not use the device without a cover.
- If you experience any pain, discomfort, dizziness, or bleeding, stop using the device immediately and consult a medical professional.
- This device is not a substitute for professional medical advice or treatment. Consult a doctor for persistent ear problems or excessive earwax buildup.
- Keep the device out of reach of children. Adult supervision is required if used on children.
- Do not use the device if it is damaged or malfunctioning.
- Avoid exposing the device to water or excessive moisture.
- Do not attempt to disassemble or repair the device yourself.



Figure 9.1: Examples of multiple applications for the ear camera, including cleaning ears, nose, mouth, and use with pets. While versatile, always prioritize safety and hygiene for each application.

10. LEGAL DISCLAIMER

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended for personal hygiene and visual assistance only.

Documents - ERYITRDK – B0CKW6TC9L

no relevant documents