



[Manuals.plus](#) /

> [Inspire](#) /

> Inspire BoostMAX Incontinence Booster Pads User Manual

Inspire INSBML72

Inspire BoostMAX Incontinence Booster Pads User Manual

Model: INSBML72

1. INTRODUCTION

The Inspire BoostMAX Incontinence Booster Pads are designed to significantly increase the absorbency of adult diapers and pull-ups. These extra-long pads (18.5" x 3.9") provide superior coverage and are engineered to prevent leaks, ensuring dryness and comfort for extended periods, including overnight use. The pads feature a super absorbent polymer that locks in fluid and reduces odor, along with a soft top sheet for enhanced skin comfort.

This manual provides essential information on the safe and effective use of your Inspire BoostMAX Booster Pads.

2. SAFETY INFORMATION

- Always use these booster pads in conjunction with adult diapers or pull-ups. Do not use them independently as a primary incontinence product.
- Dispose of used pads properly and immediately in a waste receptacle.
- Keep out of reach of children and pets to prevent accidental ingestion or misuse.
- Do not flush booster pads down the toilet, as they are not designed to be flushable and can cause plumbing issues.
- Always wash hands thoroughly after handling used incontinence products.

3. SETUP AND INSTALLATION

To properly install the Inspire BoostMAX Booster Pad, follow these steps:

1. Ensure the adult diaper or pull-up is clean and ready for use.
2. Locate the adhesive strip on the back of the booster pad.
3. Carefully peel off the protective backing from the adhesive strip.
4. Position the booster pad inside the adult diaper or pull-up, ensuring it is centered and covers the area requiring additional absorbency. Press firmly to secure the adhesive strip to the inner surface of the diaper.
5. Adjust the booster pad as needed to ensure a comfortable and secure fit within the primary incontinence garment.

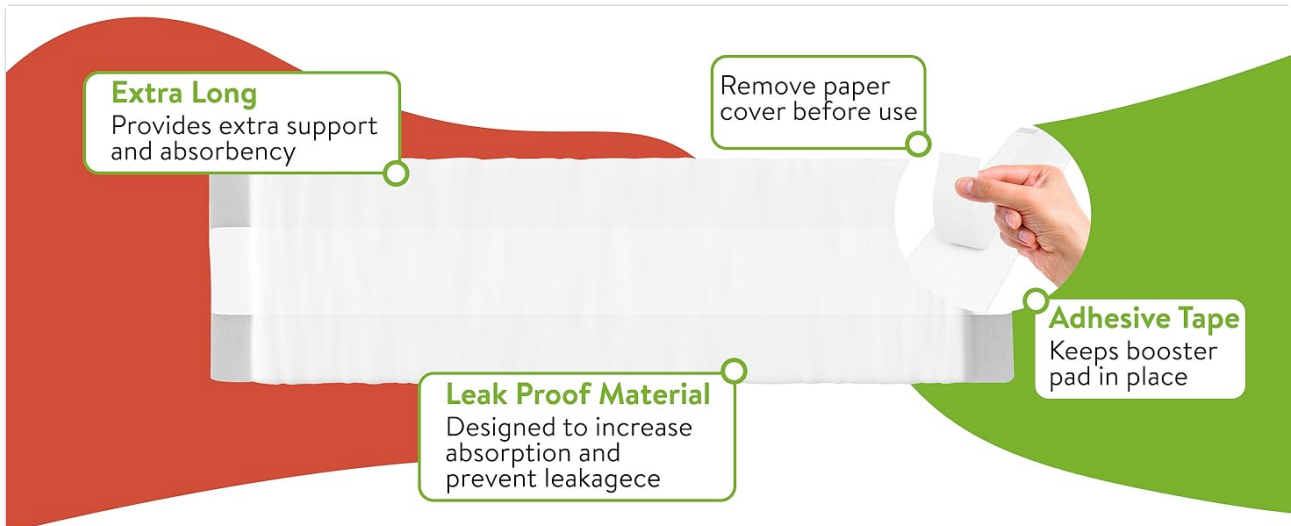


Figure 1: Peeling the adhesive strip for secure placement.

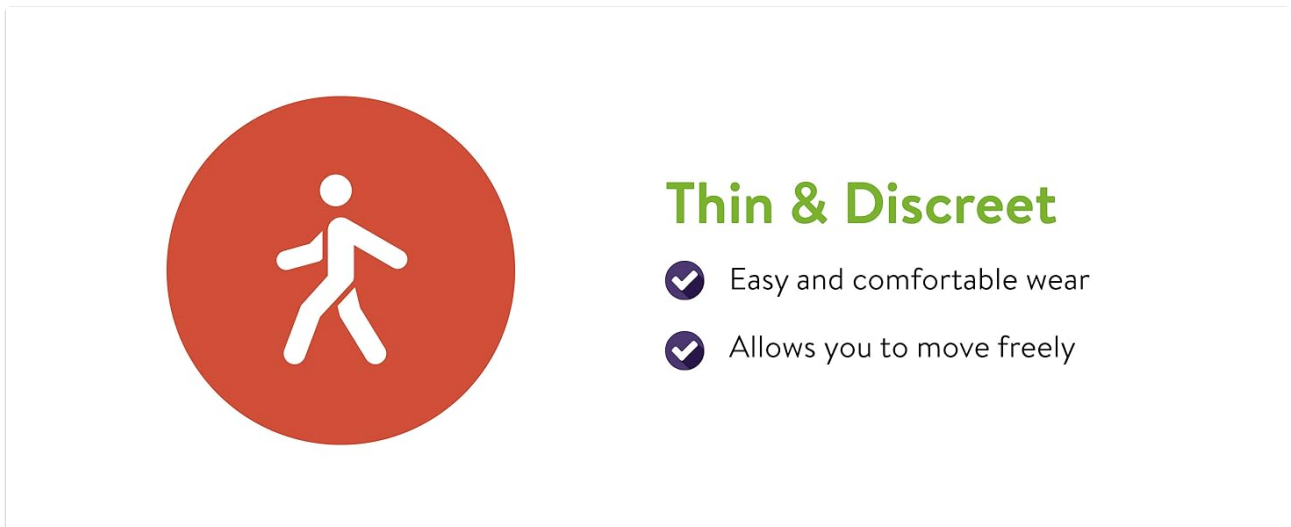


Figure 2: Booster pad placed inside an adult diaper for discreet and thin design.

4. OPERATION

The Inspire BoostMAX Booster Pad works by providing an additional layer of absorbency to your existing adult diaper or pull-up. Its advanced design ensures efficient fluid management:

- **Flow-Through Design:** The pad is engineered to rapidly absorb fluid and then allow any additional flow to pass through to the primary adult diaper, maximizing the overall absorbency capacity of your incontinence system.
- **Super Absorbent Polymer:** The core of the booster pad contains super absorbent polymer (SAP) which quickly gels upon contact with fluid, effectively locking in moisture and reducing odors for extended wear.

- **Soft Top Sheet:** A soft, comfortable top sheet wicks away moisture from the skin, helping to maintain skin dryness and reduce the risk of irritation or breakdown.

*Inspire*TM
BoostMAX
Booster Pads

Figure 3: Illustration of the flow-through design for enhanced absorbency.



Ultra Absorbent & Leakproof

- ✓ Offers an extra layer of protection
- ✓ Provides extra support when you need it most

Figure 4: Demonstration of the pad's maximum absorption capabilities.

5. MAINTENANCE AND DISPOSAL

Inspire BoostMAX Booster Pads are designed for single use and are disposable. Proper disposal is essential for hygiene and environmental responsibility.

- After use, remove the booster pad from the adult diaper or pull-up.
- Roll or fold the used pad and secure it with the adhesive strip or tape from the primary garment, if available, to contain waste.
- Dispose of the used pad in a designated waste bin. Do not flush down the toilet.
- Always wash your hands thoroughly with soap and water after handling used incontinence products.

6. TROUBLESHOOTING

If you experience issues with your Inspire BoostMAX Booster Pads, consider the following:

Problem	Possible Cause	Solution
Leaks still occur	Booster pad not positioned correctly; primary diaper is saturated; absorbency needs exceed product capacity.	Ensure the booster pad is centered and securely attached. Check if the primary diaper is adequately absorbent for your needs. Consider changing the primary diaper and booster pad more frequently.
Pad shifts or bunches	Adhesive strip not fully secured; excessive movement.	Press the adhesive strip firmly onto the primary diaper during installation. Ensure the pad is smooth and flat before wearing.
Skin irritation	Infrequent changes; allergic reaction to material.	Change pads and primary diapers more frequently to maintain skin dryness. If irritation persists, discontinue use and consult a healthcare professional.

7. PRODUCT SPECIFICATIONS

Feature	Detail
Product Name	Inspire BoostMAX Incontinence Booster Pads
Model Number	INSBML72
Brand	Inspire
Protector Type	Booster Pad
Age Range	Adult
Target Audience	Unisex-Adults
Item Length	18.5 inches (Long)
Item Width	3.9 inches
Number of Items	72 Count
Material Type	Fluff, Super Absorbent Polymer
Reusability	Disposable
Water Resistance Level	Waterproof (outer layer)
Color	White
Absorbency	5+ Cups

8. WARRANTY AND SUPPORT

Inspire offers a **Leak Free Guarantee** on its BoostMAX Incontinence Booster Pads. If you experience issues with product performance, please contact Inspire customer support for assistance.

For further inquiries or support, please visit the official Inspire store or contact K2 Health Products, the manufacturer.

Manufacturer: K2 Health Products

Brand Store: [Inspire Store on Amazon](#)

© 2026 Inspire. All rights reserved.

This manual is for informational purposes only. Product specifications are subject to change without notice.