



Manuals.plus /

› AVENT /

› Philips Avent Manual Breast Pump, Easy Pumping with Natural Motion Technology, Silicone Cushion, 4oz Natural Bottle, Sealing Disc, Screw Ring, Transparent, Model SCF430/03 Manual Breast Pump Manual

AVENT SCF430/03

Philips Avent Manual Breast Pump User Manual

Model: SCF430/03

INTRODUCTION

The Philips Avent Manual Breast Pump is designed to provide gentle comfort and efficient milk expression, supporting mothers throughout their breastfeeding journey. It incorporates Natural Motion Technology to mimic a baby's natural suckling, ensuring quick milk flow and comfort. This manual pump is ideal for on-the-go use due to its compact and lightweight design.

WHAT'S INCLUDED

Your Philips Avent Manual Breast Pump set includes the following components:

- 1x Philips Avent Manual Breast Pump
- 1x Silicone Cushion
- 1x 4oz Natural Baby Bottle
- 1x Sealing Disc
- 1x Screw Ring

KEY FEATURES

- **Quick Milk Flow:** Features Natural Motion technology with a soft silicone shield that mimics a baby's unique suction motion, stimulating the breast for fast milk flow.
- **Mom-Friendly Pumping:** The innovative soft silicone cushion stimulates the breast to initiate milk flow faster, ensuring a more comfortable and efficient pumping experience.
- **One Size Fits All:** The silicone shield gently and flexibly adapts to various breast shapes and sizes, fitting 99.98% of all nipple sizes (up to 1.18 inches).
- **Easy to Use and Clean:** Designed with a simple construction consisting of only a few parts, making it easy to set up, use, and clean.

- **Compact and Portable:** Small and lightweight, this breast pump is easy to store and transport, enabling convenient and discreet pumping on the go.
- **Upright Pumping Position:** The angle of the pump allows for a comfortable upright pumping position, eliminating the need to lean forward.

SETUP

1. **Clean and Sterilize:** Before first use, disassemble all parts of the breast pump and wash them thoroughly with warm soapy water. Rinse well. Sterilize all parts by boiling them for 5 minutes or using a steam sterilizer. Ensure all parts are completely dry before assembly.
2. **Assemble the Pump:**
 - Insert the white valve into the pump body from below. Ensure it is securely in place.
 - Attach the silicone diaphragm to the pump body, pressing down firmly around the edges to ensure a perfect seal.
 - Place the pump handle onto the diaphragm stem, pushing down until it clicks into place.
 - Screw the 4oz Natural baby bottle onto the bottom of the pump body.
 - Attach the soft silicone cushion to the funnel of the pump body, ensuring it fits snugly around the rim.
3. **Check Assembly:** Verify that all parts are correctly assembled and securely fitted to prevent leaks and ensure proper suction.



Image: The Philips Avent Manual Breast Pump fully assembled, showing the pump body, handle, silicone cushion, and the attached 4oz Natural bottle with a separate sealing lid.

OPERATING INSTRUCTIONS


1. **Positioning:** Sit comfortably in an upright position. Place the assembled breast pump's silicone cushion over your breast, ensuring your nipple is centered in the funnel. Create a good seal between the cushion and your skin.
2. **Initiating Suction:** Gently press the handle until you feel suction. Release the handle to allow milk to flow. Repeat this motion rhythmically, mimicking your baby's natural suckling pattern.
3. **Pumping Technique:** Begin with a faster, lighter rhythm to stimulate milk let-down. Once milk starts flowing, switch to a slower, deeper rhythm to express milk efficiently. Adjust the rhythm and pressure to what feels most comfortable and effective for you.
4. **Upright Pumping:** The pump's design allows you to maintain an upright posture during pumping, which can enhance comfort and milk flow.
5. **Milk Storage:** Once you have finished pumping, carefully remove the pump from your breast. Unscrew the bottle from the pump body. Use the included sealing disc and screw ring to seal the bottle for storage in the refrigerator or

freezer.

Natural Motion technology: suction and nipple stimulation for quick milk flow



Image: A woman demonstrating the use of the Philips AVENT Manual Breast Pump, illustrating the Natural Motion technology for efficient milk expression.

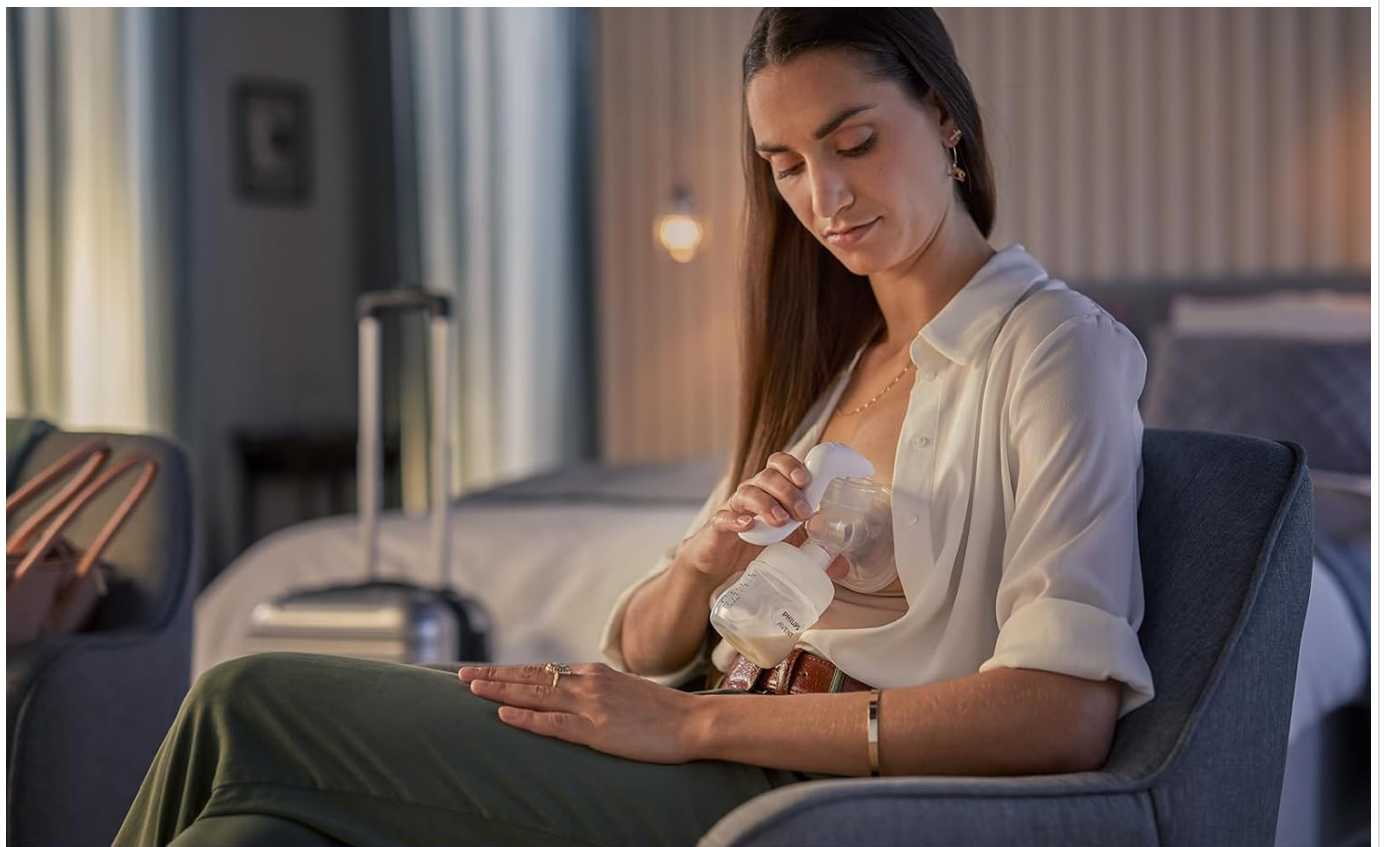

**Upright pumping
position** thanks
to a no lean
forward design

Image: A woman comfortably using the Philips Avent Manual Breast Pump while sitting upright, showcasing the no-lean-forward design.

CLEANING & MAINTENANCE

Proper cleaning and maintenance are crucial for the hygiene and longevity of your breast pump. The Philips Avent Manual Breast Pump is designed for easy cleaning with its minimal parts.

1. **Disassembly:** After each use, fully disassemble all parts of the breast pump.
2. **Washing:** Wash all disassembled parts immediately in warm water with a mild dish soap. Use a soft brush to clean hard-to-reach areas.
3. **Rinsing:** Rinse all parts thoroughly under cold running water to remove any soap residue.
4. **Sterilization:** Sterilize all parts daily. You can do this by:
 - **Boiling:** Place all parts in a large pot of boiling water for 5 minutes. Ensure parts do not touch the bottom or sides of the pot to prevent damage.
 - **Steam Sterilizer:** Use a Philips Avent steam sterilizer or another compatible steam sterilizer according to its instructions.

- Dishwasher: All parts are dishwasher safe (top rack only).

5. **Drying and Storage:** After sterilization, allow parts to air dry completely on a clean drying rack. Store the dry parts in a clean, covered container until next use.



Few parts for
easy set up,
use and cleaning

Image: Hands cleaning the few, easily detachable parts of the Philips Avent Manual Breast Pump, demonstrating its simple cleaning process.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No or low suction	<ul style="list-style-type: none"> • Improper assembly of parts. • Silicone diaphragm or valve not sealed correctly. • Silicone cushion not creating a tight seal on the breast. 	<ul style="list-style-type: none"> • Disassemble and reassemble all parts, ensuring they are securely fitted. • Check the white valve for tears or damage; replace if necessary. • Ensure the silicone cushion is firmly pressed against the breast, with the nipple centered.
Milk leaking from the pump	<ul style="list-style-type: none"> • Bottle not screwed on tightly. • Parts not assembled correctly. • Silicone cushion not properly attached or damaged. 	<ul style="list-style-type: none"> • Ensure the bottle is securely screwed onto the pump body. • Check all connections for proper fit and seal. • Inspect the silicone cushion for any tears or deformities; replace if damaged.
Discomfort during pumping	<ul style="list-style-type: none"> • Incorrect breast positioning. • Excessive pressure or too fast a rhythm. • Silicone cushion not suitable for nipple size. 	<ul style="list-style-type: none"> • Reposition the pump to ensure the nipple is centered and comfortable. • Adjust pumping rhythm and pressure to a comfortable level. • While the cushion is designed to fit most, ensure it feels comfortable. If discomfort persists, consult a lactation specialist.

SPECIFICATIONS

Attribute	Detail
Model Number	SCF430/03
Product Dimensions	3.54 x 5.91 x 6.18 inches
Item Weight	11 ounces
Material Type	Plastic
Material Free	BPA Free
Dishwasher Safe	Yes
Style	Manual Breast Pump

WARRANTY & SUPPORT

Information regarding specific warranty details and customer support for the Philips AVENT Manual Breast Pump (Model SCF430/03) is typically provided with the product packaging or can be found on the official Philips AVENT website. Please refer to these resources for the most accurate and up-to-date information.

For further assistance, you may visit the AVENT Store on Amazon or contact Philips AVENT customer service directly.

