

HEYVALUE Nasensauger Baby Elektrisch

HEYVALUE Electric Nasal Aspirator User Manual

Model: Nasensauger Baby Elektrisch

1. PRODUCT OVERVIEW

The HEYVALUE Electric Nasal Aspirator is designed to provide gentle and effective relief from nasal congestion for newborns, infants, and toddlers. It features adjustable suction levels, a soothing night light, and integrated nursery rhymes to enhance comfort during use. This device is equipped with food-grade silicone tips for safety and is designed for easy cleaning.



Image 1.1: The HEYVALUE Electric Nasal Aspirator, including the main unit, suction chamber, tubing, and two silicone tips, alongside its product packaging.

Key Features:

- **9 Adjustable Suction Levels:** Provides customizable suction power for varying levels of congestion.
- **Integrated Night Light:** Features a 3-level adjustable ambient night light for use in low-light conditions.
- **Soothing Nursery Rhymes:** Includes calming tunes to help distract and soothe the baby during aspiration.
- **Food-Grade Silicone Tips:** Two sizes of soft, hypoallergenic tips ensure comfort and safety for delicate nasal passages.
- **Portable and Rechargeable:** Compact design with USB-C charging for convenient on-the-go use.
- **Easy to Clean:** Detachable and dishwasher-safe components for hygienic maintenance.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the HEYVALUE Electric Nasal Aspirator. Failure to follow these guidelines may result in injury or damage to the device.

- Limit suctioning to no more than four times a day to prevent irritation of the nasal passages.
- Always clean the aspirator thoroughly with warm, soapy water or a sterilizing solution after each use.
- Never insert the aspirator tip too far into the baby's nostril.
- Avoid using the aspirator if the baby is experiencing ear pain.
- Keep the device out of reach of children when not in use.
- Ensure all components are properly assembled before operation.



Image 2.1: The aspirator features ultra-soft, food-grade silicone tips designed to be gentle on a baby's delicate skin.

3. PRODUCT COMPONENTS

The HEYVALUE Electric Nasal Aspirator package includes the following components:

- Main Aspirator Unit
- Transparent Suction Chamber
- Silicone Tubing
- Small Silicone Tip (for newborns/infants)
- Large Silicone Tip (for toddlers)
- USB-C Charging Cable
- User Manual

4. SETUP

Before first use, and after each cleaning, ensure the device is properly set up.

1. **Charge the Device:** Connect the provided USB-C cable to the charging port on the main unit and to a suitable USB power adapter (not included). The battery indicator on the display will show charging status. Ensure the device is fully charged for optimal performance.
2. **Assemble Components:**
 - Attach the transparent suction chamber securely to the main aspirator unit.
 - Connect one end of the silicone tubing to the suction chamber.
 - Select the appropriate food-grade silicone tip size (small for newborns/infants, large for toddlers) and attach it firmly to the other end of the silicone tubing.
3. **Initial Cleaning:** Before the first use, clean all detachable parts (suction chamber, tubing, silicone tips) as described in the 'Maintenance and Cleaning' section. Ensure all parts are completely dry before assembly.

USB-C Rechargeable & Portable

More convenient for on-the-go use



Image 4.1: The compact and USB-C rechargeable design makes the aspirator portable for on-the-go use.

5. OPERATING INSTRUCTIONS

Follow these steps to safely and effectively use the nasal aspirator:

1. **Prepare the Baby:** Position your baby comfortably, ideally in a slightly reclined position. If mucus is dry or thick, consider using a saline spray a few minutes before aspiration to loosen it.
2. **Power On:** Press and hold the power button (U) on the main unit to turn on the device. The display will illuminate.
3. **Adjust Suction Level:** Use the '+' and '-' buttons to select one of the 9 available suction levels. Start with a lower level and gradually increase if needed, observing your baby's comfort.
4. **Activate Night Light (Optional):** Press the light bulb button () to cycle through the 3 adjustable night light

levels. This can provide gentle illumination during nighttime use.

5. **Play Nursery Rhymes (Optional):** Press the music note button () to play soothing nursery rhymes, which can help distract and calm your baby.
6. **Perform Aspiration:** Gently insert the silicone tip into one of your baby's nostrils, ensuring a snug but not tight fit. Avoid pushing the tip too deep. Allow the device to suction for a few seconds until the mucus is cleared. Repeat for the other nostril as necessary.
7. **Power Off:** After use, press and hold the power button (U) to turn off the device.



Image 5.1: The aspirator offers 9 levels of powerful suction to effectively clear nasal passages.



Image 5.2: The device includes soothing nursery rhymes to help calm and distract the baby during use.

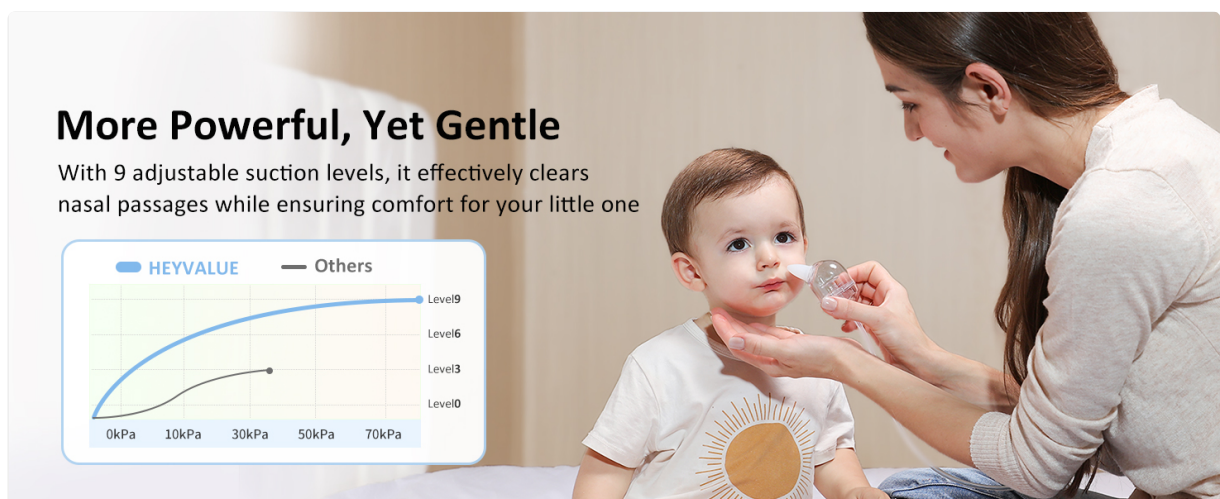


Image 5.3: The cozy night light provides perfect illumination for nighttime nose-cleaning and feeding.

6. MAINTENANCE AND CLEANING

Proper cleaning after each use is essential for hygiene and to maintain the device's performance.

1. **Disassemble:** Turn off the aspirator. Detach the silicone tip, silicone tubing, and transparent suction chamber from the main unit.
2. **Rinse:** Rinse the detached components (silicone tip, tubing, and suction chamber) thoroughly under clean running water.
3. **Clean:** For a more thorough cleaning, wash these components with warm, soapy water or a sterilizing solution. The suction chamber and tubing are also dishwasher safe for convenience.
4. **Dry:** Ensure all cleaned parts are completely dry before reassembling or storing the device. Air drying is recommended. Do not immerse the main aspirator unit in water.
5. **Store:** Store the clean and dry aspirator in a cool, dry place, away from direct sunlight.



Image 6.1: The detachable components are easy to hand wash and are also dishwasher safe for convenient cleaning.

7. TROUBLESHOOTING

If you encounter any issues with your HEYVALUE Electric Nasal Aspirator, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on	Low battery; Device not charged	Charge the device using the provided USB-C cable. Ensure the charging cable is securely connected.
Weak or no suction	Improper assembly; Blockage in tubing or tip; Suction chamber not sealed; Low battery	Check that the suction chamber, tubing, and silicone tip are securely attached and properly sealed. Inspect for any visible blockages and clean if necessary. Ensure the device is fully charged.
Unusual noise during operation	Normal operating sound; Loose component	A slight motor sound is normal. If the noise is excessive or unusual, check that all detachable components are securely fastened.

Problem	Possible Cause	Solution
Night light or music not working	Feature not activated; Low battery	Press the respective button (light bulb for night light, music note for rhymes) to activate. Ensure the device has sufficient battery charge.

8. SPECIFICATIONS

Model Name	NASA603
Item Model Number	Nasensauger Baby Elektrisch
Brand	HEYVALUE
Color	White
Age Range (Description)	Newborns, Infants and Toddler
Material Type	Plastic
Material Type Free	Latex Free, BPA Free
Power Source	Battery Powered
Water Resistance Level	Waterproof (for detachable components)
Product Care Instructions	Dishwasher Safe (for detachable components)
Item Weight	14.9 ounces
Package Dimensions	6.5 x 6.42 x 2.28 inches
UPC	717225339975
Manufacturer	HEYVALUE

9. WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the contact details provided with your product packaging or visit the official HEYVALUE website. Our customer service team is available to assist with any questions or concerns regarding your nasal aspirator.