

Redux Home Dryer

Redux Home Hearing Aid Dryer and Dehumidifier User Manual

Brand: Redux | **Model:** Home Dryer (ASIN: B0CLS9LLLT)

1. INTRODUCTION

The Redux Home Hearing Aid Dryer and Dehumidifier is designed to effectively remove moisture from hearing aids, AirPods, and earbuds. Utilizing active drying technology with a self-recharging desiccant system, it ensures optimal performance and extends the lifespan of your valuable devices. This manual provides essential information for the safe and efficient operation of your Redux Home Dryer.



Redux Home Dryer, a device designed for moisture removal from hearing aids and other small electronics.

2. SAFETY INFORMATION

Please read and adhere to the following safety guidelines to prevent damage to the device or injury:

- Do not immerse the Redux Home Dryer in water or any other liquids.
- This device is intended solely for drying hearing aids and other specified electronic devices.
- Use only the included power supply.
- Avoid bringing excessive liquids into contact with the device.
- The Redux Home Dryer is a tool, not a toy. Keep out of reach of children.
- Damage resulting from misuse is not covered by warranty.

3. PRODUCT OVERVIEW AND FEATURES

The Redux Home Dryer incorporates advanced technology to provide comprehensive moisture removal for your electronic devices. Key features include:

- **100% Moisture Removal & No Refills:** Utilizes active drying technology with smart sensors to eliminate all traces of harmful moisture. The self-recharging desiccant system requires no replacements.
- **Fast 45-Minute Cycles:** Delivers a complete dry in approximately 45 minutes, ideal for daily use.
- **Safe for Rechargeable Devices:** The gentle drying process is suitable for all rechargeable hearing aids, AirPods, and delicate earbuds, unlike harsh UV dryers.
- **Smart App Connectivity:** Connects via Bluetooth to the Redux app, providing drying completion confirmation and detailed moisture removal data.
- **Maintenance-Free:** The internal desiccant automatically recharges after each cycle, ensuring continuous effective operation without ongoing costs.



**Reliable,
Efficient,
and Safe.**



Hearing Aids Implants Hearables

Professional-grade drying for your devices.



**IMPROVES
HEARING
AIDS IN
ANY
CONDITION**

DEVICES THOUGHT TO BE DEAD:	DEVICES THAT SOUND WEAK:	DEVICES THAT SOUND GOOD:
50% REVIVED	82% IMPROVED	80% IMPROVED

Revolutionary desiccant technology.

4. SETUP

Follow these steps for initial setup:


1. **Unpack:** Carefully remove the Redux Home Dryer and all accessories from its packaging.
2. **Power Connection:** Connect the provided power supply to the device and plug it into a suitable electrical outlet.
3. **App Download (Optional):** For enhanced features and monitoring, download the Redux app from your device's app store and connect via Bluetooth.

5. OPERATING INSTRUCTIONS

Using your Redux Home Dryer is simple:

1. **Place Devices:** Open the lid of the Redux Home Dryer. Carefully place your hearing aids, AirPods, or earbuds into the drying chamber.
2. **Close Lid:** Securely close the lid.
3. **Start Cycle:** Press the 'Start' button. The device will display steady blue lights, indicating that the drying cycle is active.
4. **Drying Completion:** The blue lights will turn off when the drying process is complete (typically around 45 minutes). You can then unlock the lid and retrieve your dried devices.
5. **App Monitoring:** If using the Redux app, you will receive confirmation when drying is complete and can view detailed moisture removal data.


REVOLUTIONARY DESICCANT TECHNOLOGY.

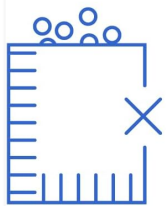

Redux Home Dryer regenerates its desiccant, meaning no refills are required.



Hearing aids in the drying chamber.

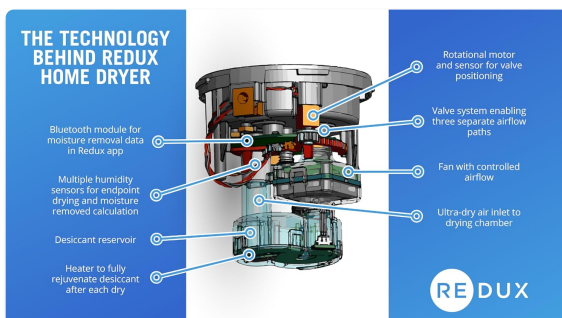

AirPods in the drying chamber.

Official Product Video

Your browser does not support the video tag.

Video: Redux Home: Professional-grade hearing aid drying at home! This video demonstrates the Redux Home Dryer's operation and benefits for maintaining hearing devices.

6. MAINTENANCE

The Redux Home Dryer is designed for minimal maintenance:

- **Self-Recharging Desiccant:** The internal desiccant automatically regenerates after each drying cycle. No manual refills or replacement desiccant bricks are required.
- **Cleaning:** Periodically wipe the exterior of the device with a soft, dry cloth. Ensure the drying chamber is free of debris.

7. TROUBLESHOOTING

Refer to the table below for common issues and their solutions:

Issue	Possible Cause	Solution
Orange light flashing	Temporary error or desiccant regeneration cycle.	Unplug the device, wait 10 seconds, then plug it back in. If the issue persists, contact support.
Device not drying effectively	Lid not properly closed, device not powered on, or cycle not initiated.	Ensure the lid is securely closed. Verify the device is powered on and the drying cycle has started (steady blue light).
App connectivity problems	Bluetooth disabled, app not running, or temporary connection issue.	Ensure Bluetooth is enabled on your mobile device and the Redux app is open. Restart the app or the dryer if necessary.

8. SPECIFICATIONS

- **Product Dimensions:** 4.4 x 4.4 x 5.3 inches
- **Weight:** 1.5 Pounds
- **Manufacturer:** Redux
- **ASIN:** B0CLS9LLLT

9. WARRANTY AND SUPPORT

For product support, warranty information, or to address any concerns not covered in this manual, please refer to the manufacturer's official website or contact their customer service directly. It is recommended to retain your proof of purchase for any warranty claims.

