

Behaesty P3

Behaesty Ear Wax Removal Tool Camera User Manual

Model: P3

INTRODUCTION

Welcome to the user manual for your Behaesty Ear Wax Removal Tool Camera. This device combines an 800W Pixels FHD camera with a soft silicone ear scoop, providing a clear visual for safe and effective ear wax removal. It connects wirelessly to your smartphone or tablet, allowing for real-time viewing, photo capture, and video recording of your ear canal.

WHAT'S INCLUDED

- 1 x Visual Earwax Remover with Lid
- 10 x Silicone Ear Scoops
- 1 x Traditional Earwax Removal 6-piece Set
- 1 x USB Type-C Charging Cable
- 1 x Instruction Manual



The complete package includes the visual earwax remover, multiple silicone ear scoops for hygiene, a set of traditional metal ear picks for various needs, a USB-C charging cable, and this user manual.

KEY FEATURES

High-Resolution Camera

- 800W Pixels FHD Camera: Provides clear and detailed images of the ear canal.
- 360° Wide Angle Lens: Offers a comprehensive view for easier navigation.
- 6 LED Super Lights: Ensures bright and thermostatic illumination inside the ear.

Upgraded 800W Pixel HD Camera

Clearer And Steadier
Transmission



HD Camera



65 minutes



360°Wide Angle



Vido Recording



Take Pictures



WiFi Wireless

The device features an 800W pixel Full HD camera for superior image clarity and six integrated LED lights to illuminate the ear canal, ensuring a bright and stable view.

Smart Connectivity

- WiFi Connection: Seamlessly connects to your smartphone or tablet via WiFi for real-time viewing.
- App Compatibility: Works with most Android (5.0 and above) and iOS (10.0 and above) devices via the "Soulear" or "JesHome" APP.
- Image & Video Capture: Allows users to save images and videos directly to their device.



Wide Compatibility

Work with ios system and android system



ios



android

Search JesHome APP in Apple Store / Google Play / App Store or scan QR code to download and install directly



App store



Google play

Mobilephone system support instructions

- 1.For Android phones: support 5.0 and above
- 2.For iPhone: supports 10.0 and above.

Note: The APP is safe, if any questions, please check the instructions.

The ear camera is compatible with both iOS and Android smartphones and tablets. Simply download the dedicated 'JesHome' or 'Soulear' app to connect and view the live feed. [Download App Here](#)

User-Friendly Design

- Soft Silicone Ear Tip: Ensures comfortable and safe use, protecting the ear canal.
- IP67 Waterproof Lens: Allows for easy cleaning and maintenance of the lens.
- Diamond-shaped Body Design: Provides a secure and comfortable grip.

All- In- One Design Ear Cleaner

High-speed WIFI chip makes your earwax
nowhere to hide




Video Recording



Take Pictures



WiFi Signal

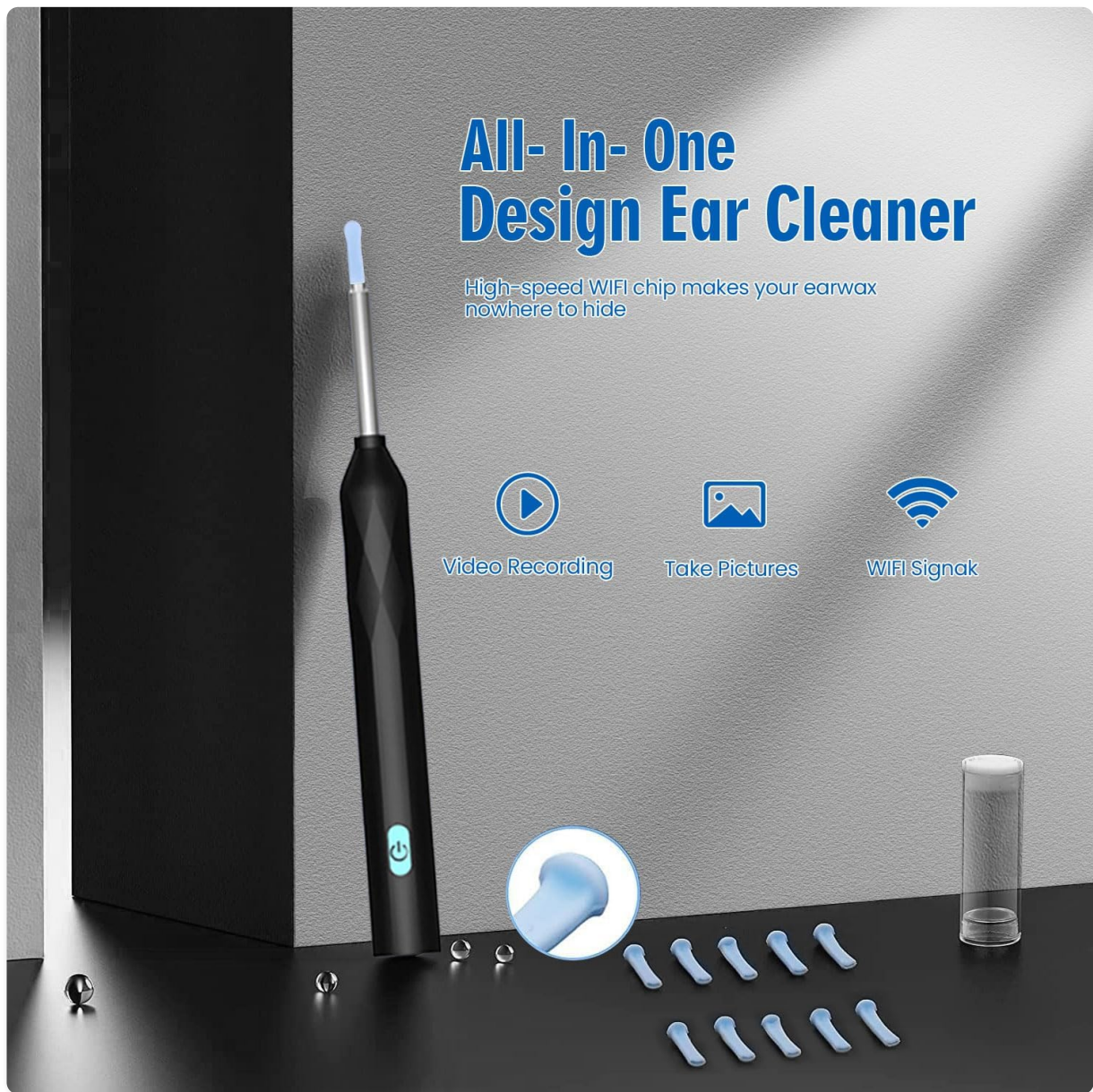

Designed for ease of use, the device features a soft silicone ear tip for gentle insertion and an IP67 waterproof lens for hygienic cleaning. Its ergonomic design ensures a comfortable and secure grip during use.

SETUP GUIDE

Step 1: Download the App

Scan the QR code provided in the manual or search for "Soulear" or "JesHome" APP in your device's app store (Apple App Store or Google Play Store).

[Download App Here](#)

Step 2: Power On the Device

Press and hold the power button on the ear wax removal tool for 1 second until the blue indicator light illuminates.

Step 3: Connect via WiFi

On your smartphone or tablet, go to your WiFi settings. Select the network named "JesHome-PX-XXXX" (where XXXX is a unique identifier for your device) and connect to it. This is a direct device-to-phone connection, not requiring internet access.

Step 4: Start Using the App

Open the "Soulear" or "JesHome" APP. The app should automatically detect the connected device. Click "Start" to begin real-time viewing.



Follow these steps to connect your device: 1. Download the JesHome APP. 2. Hold the power button to boot. 3. Open the APP and click Add. 4. Select 'JesHome-PX-XXXX' WiFi. 5. Return to the APP and click Start.

Watch this video for a visual guide on setting up and operating your Behaesty Ear Wax Removal Tool Camera, including connecting to the app and using the camera for ear inspection.

OPERATING INSTRUCTIONS

General Use

Once connected, the live view of your ear canal will appear on your device's screen. Use the silicone ear scoop attached to the camera to gently remove ear wax while observing the process on your screen. The 360° wide-angle lens ensures the image stays upright even as the camera rotates, providing stable orientation.

Multi-Purpose Application

Beyond ear wax removal, this versatile tool can be used for various inspections, including:

- Kids Oral Health Examination
- Nasal Examination
- Dermatological Observation (e.g., skin pores)
- Scalp Health Check
- Pet Ear Check



The Behaesty visual ear tool is not limited to ear cleaning; it can also be used for examining the mouth, nose, skin, scalp, and even pet ears, making it a versatile health tool for the whole family.

MAINTENANCE AND CARE

Cleaning the Ear Scoops

Always clean the silicone ear scoops and traditional ear picks with alcohol before and after each use to ensure hygiene and prevent contamination.

Cleaning the Lens

The IP67 waterproof lens can be safely rinsed under water. Ensure it is completely dry before storage.

Charging

Recharge the device using the provided USB Type-C cable. A full charge takes approximately 40 minutes. The device can last up to 180 days on standby and up to 90 minutes of continuous use.

Storage

Store the device and accessories in a cool, dry place, away from direct sunlight and high humidity. Do not place heavy objects on top of the device.


TROUBLESHOOTING

S/N	Fault	Solution
1	No response after startup	Please charge with the provided charging cable. The charging duration is 40 minutes.
2	Camera indicator does not light up	Shut down and restart or reconnect to WiFi.
3	APP crashes	Uninstall the existing APP and rescan the QR code to install the latest version of APP.
4	Silicone sleeve is loose	Remove the silicone sleeve and slowly re-insert the earpick, or replace the silicone sleeve.
5	Blue indicator flashes quickly	Low battery, please charge the product with the provided charging cable.
6	Unstable connection, it is easy to disconnect or freeze when using	Restart the earpick and reconnect to the APP.
7	After the product is connected to APP, image cannot be displayed normally	Disconnect the phone from the WiFi of the earpick, restart the earpick, and reconnect.

If the above solutions fail to resolve the problem, please contact our customer service staff. We will further provide you with practical solutions on APP-Personal Center-FAQ.

Troubleshooting

S/N	Fault	Solution	S/N	Fault	Solution
1	No response after startup	Please charge with the provided charging cable. The charging duration is 40 minutes	5	Blue indicator flashes quickly	Low battery, please charge the product with the provided charging cable
2	Camera indicator does not light up	Shut down and restart or reconnect to WIFI	6	Unstable connection, it is easy to disconnect or freeze when using	Restart the earpick and reconnect to the APP
3	APP crashes	Uninstall the existing APP and rescan the QR code to install the latest version of APP	7	After the product is connected to APP, image cannot be displayed normally	Disconnect the phone from the WiFi of the earpick, restart the earpick, and reconnect to WiFi
4	Silicone sleeve is loose	Remove the silicone sleeve and slowly re-insert the earpick, or replace the silicone sleeve	If the above operation still fails to solve the problem, please contact our customer service staff, we will further provide you with practical solutions on APP-Personal Center-FAQ.		

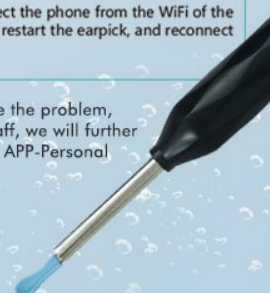

Refer to this table for quick solutions to common issues you might encounter with your device, such as connection problems or indicator lights.

PRODUCT SPECIFICATIONS

- **Model Number:** P3
- **Camera Pixels:** 800W FHD
- **Lens:** 360° Wide Angle, 6 LED Lights
- **Waterproof Rating:** IP67 (Lens)
- **Connectivity:** WiFi (2.4GHZ Operating Frequency)
- **Battery Life:** Up to 90 minutes continuous use, 180 days standby
- **Charging Time:** Approximately 40 minutes (USB Type-C)
- **Compatibility:** Android 5.0 and above, iOS 10.0 and above
- **Dimensions:** 7.36 x 2.95 x 0.75 inches
- **Weight:** 2.4 ounces
- **Manufacturer:** Behaesty

IMPORTANT SAFETY INFORMATION

- It is recommended to carefully clean the lens with cotton swabs.
- Please observe your surroundings before use. Do not use it in an environment where it may come into contact with moving people or animals. Do not use it in a moving vehicle, or in other unstable environments.
- Please keep your body stable when using it.
- This product is not suitable for children under 3.
- Avoid accidental injury; children are prohibited from using this product alone.
- Do not place the machine in a humid environment for a long time, or immerse the whole machine in liquid.
- The product has a built-in polymer battery, which must be charged once a month if not used frequently.
- Please charge the device with a 5V charger. It is forbidden to use this device when it is under charging status.
- Avoid exposing the product to direct sunlight.
- Do not place heavy objects on top of the machine.
- During use, the temperature of the product will rise slightly (up to about 3°C), which is normal.
- Do not put the product in water for clearing to avoid damage.

Legal Disclaimer

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided in the "Personal Center" or "FAQ" section of the "Soulear" or "JesHome" APP, or visit the official Behaesty store on Amazon.

[Visit Behaesty Store on Amazon](#)

© 2025 Behaesty. All rights reserved.