

[Manuals.plus](#) /

› [Behaesty](#) /

› Behaesty Smart Visual Ear Cleaner User Manual

## Behaesty P3

# Behaesty Smart Visual Ear Cleaner

MODEL: P3

User Manual

## INTRODUCTION

---

Thank you for choosing the Behaesty Smart Visual Ear Cleaner. This innovative device combines an 800W Pixels FHD camera with a comprehensive ear wax removal kit, providing a safe and effective solution for personal ear hygiene. This manual will guide you through the setup, operation, maintenance, and troubleshooting of your new device to ensure optimal performance and a clear view of your ear canal.

## PACKAGE CONTENTS

---

Upon opening your Behaesty Smart Visual Ear Cleaner package, please verify that all components listed below are present and in good condition.



Image: The Behaesty Smart Visual Ear Cleaner along with its full set of accessories, including silicone ear scoops, traditional ear picks, and a USB Type-C charging cable.

- 1 x Behaesty Visual Earwax Remover with Lid
- 10 x Silicone Ear Scoops
- 1 x Traditional Earwax Removal 6-piece Set
- 1 x USB Type-C Charging Cable
- 1 x Instruction Manual

Video: An overview of the Behaesty Smart Visual Ear Cleaner and its included accessories, demonstrating the contents of the package.

## KEY FEATURES

### 800W Pixels FHD Camera & 360° Wide Angle Lens

Experience crystal-clear visuals with the 800W Pixels FHD camera and a 360° wide-angle professional-grade lens. This provides a full, clear view of your ear canal, making earwax removal easier and safer. Equipped with 6 LED IP67 waterproof lights, the device ensures bright illumination and thermostatic

conditions for real-time observation.



Image: Detailed view of the 800W Pixel HD camera lens.



Image: Illustration of the 360° wide-angle view.



Image: View of the 6 super lights on the lens.

## Wide Compatibility & WiFi Connectivity

The ear cleaning tool is compatible with most Android and iOS smartphones and tablets that have WiFi functionality. It connects to the ear pick and the accompanying 'Soulear' app via WiFi, offering a video transmission rate of up to 20 FPS. You can easily save images and videos for later review. Operating frequency: 2.4GHZ.



Image: Compatibility with iOS and Android devices.

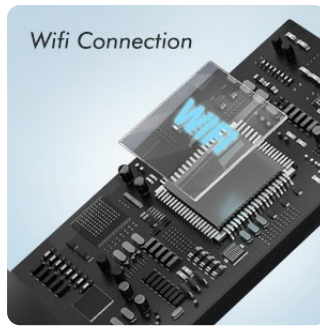


Image: Visual representation of WiFi connection.

## Long-Lasting Battery & Fast Charging

The built-in battery ensures extended use, lasting up to 180 days on standby and up to 90 minutes of continuous use. It can be conveniently recharged anytime using the included USB Type-C cable, making it eco-friendly and practical for daily use.



Image: Built-in battery and Type-C charging port.

## Multi-scene, Multi-purpose Application

Beyond ear wax removal, this versatile tool can be used for various applications, including kids' oral health examination, nasal examination, dermatological observation, scalp health checks, and even pet ear checks, making it a valuable health tool for the entire family.

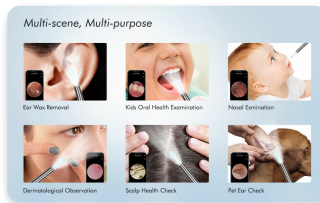


Image: Multi-scene, multi-purpose uses of the device.

## Temperature Control System & Ergonomic Design

The device features an intelligent temperature control system, maintaining the handle at approximately 32°C and the lens at 25°C for comfortable and safe use. Its diamond-shaped body design ensures a secure and comfortable grip during operation.

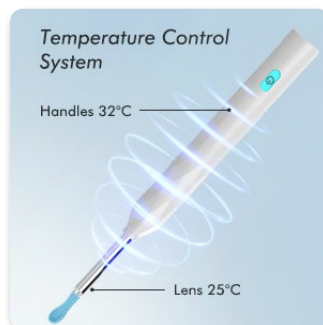


Image: Temperature control system details.



Image: Diamond-shaped body design for ergonomic grip.

## SETUP GUIDE

Follow these steps to set up your Behaesty Smart Visual Ear Cleaner and begin using it.

1. **Download the App:** Scan the QR code below or search for "Soulear" (or "JesHome") in your device's app store (Apple App Store / Google Play Store) to download and install the application. [Download App Here](#)
2. **Power On:** Press and hold the power button on the ear cleaner for 1 second to turn it on. The blue indicator light will illuminate.
3. **Connect to WiFi:** On your smartphone or tablet, go to your WiFi settings. Select the network named "JesHome-PX-XXXX" (where XXXX is a unique identifier) and connect to it. No password is required.
4. **Start Using:** Open the "Soulear" (or "JesHome") app. Click on the "Start" button to begin real-time observation of your ear canal.

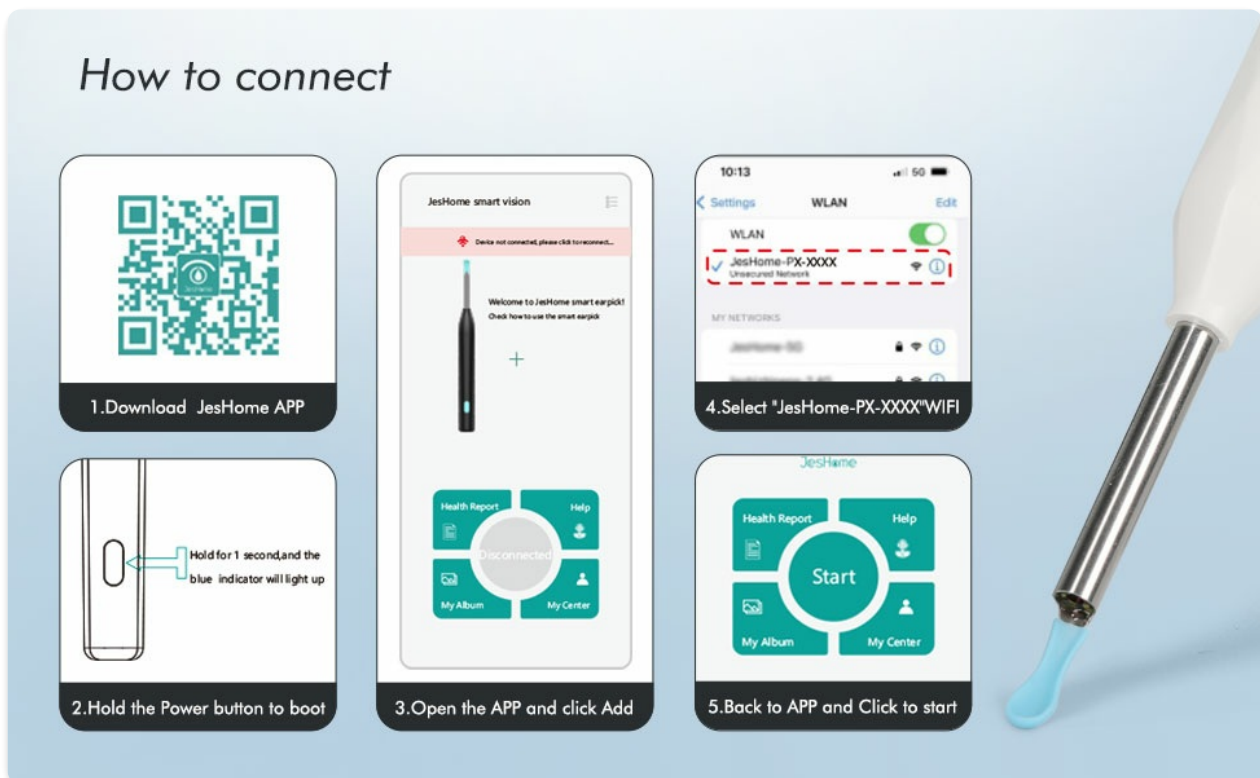


Image: Visual guide for connecting the device to your smartphone via WiFi.

Video: A demonstration of the power-on sequence and WiFi connection process for the Behaesty Smart Visual Ear Cleaner.

## OPERATING INSTRUCTIONS

Once connected, you can safely and effectively use your visual ear cleaner.

1. **Attach Ear Scoop:** Gently attach one of the silicone ear scoops to the tip of the visual ear cleaner. Ensure it is securely fitted.
2. **Position the Device:** Carefully insert the ear cleaner into your ear canal while observing the live feed on your smartphone screen. The 360° wide-angle lens provides a comprehensive view.
3. **Gentle Removal:** Using the silicone ear scoop, gently and slowly maneuver to remove visible earwax. Always move with caution and avoid pushing wax deeper into the ear canal.
4. **Monitor Temperature:** The device maintains a safe temperature. If you feel any discomfort or excessive warmth, discontinue use immediately.
5. **Capture Images/Videos:** Use the app's functions to take photos or record videos of your ear canal for personal record or consultation with a professional.
6. **Proper Technique:** Avoid using sharp objects like pliers or cotton swabs, which can push wax further in or cause injury. The visual ear cleaner with soft silicone tips is designed for safer removal.

**pliers** ❌

**cotton swab** ❌

**finger** ❌

## 6-Axis Gyroscope

Easy to use and provide smooth image HD transmission

- 45 Minutes Charging time
- 6 Axis Gyroscope
- 3.0 Aperture
- 65 Minutes Service time

Photo Video

Right ear Photo album

Image: Illustrates safe and unsafe methods of ear cleaning.

# All- In- One Design Ear Cleaner

High-speed WIFI chip makes your earwax  
nowhere to hide



Video Recording



Take Pictures



WIFI Signak



Image: Demonstrates the soft silicone ear tip and the IP67 waterproof lens.

## MAINTENANCE AND CARE

Proper maintenance ensures the longevity and hygiene of your Behaesty Smart Visual Ear Cleaner.

- **Cleaning Ear Scoops:** Always clean the silicone ear scoops and traditional ear picks with alcohol before and after each use. This prevents the spread of bacteria and ensures hygiene.
- **Cleaning the Lens:** Gently wipe the camera lens with a soft, dry cloth. Avoid using abrasive materials or harsh chemicals that could damage the lens.
- **Storage:** Store the device and its accessories in a clean, dry place, away from direct sunlight and extreme temperatures.
- **Charging:** Only use the provided USB Type-C cable for charging. Do not use damaged cables or chargers.

## TROUBLESHOOTING

If you encounter any issues with your device, please refer to the table below for common problems and their

solutions.

S/N	Fault	Solution
1	No response after startup	Please charge with the provided charging cable. The charging duration is 40 minutes.
2	Camera indicator does not light up	Shut down and restart or reconnect to WiFi.
3	APP crashes	Uninstall the existing APP and rescan the QR code to install the latest version of APP.
4	Silicone sleeve is loose	Remove the silicone sleeve and slowly re-insert the earpick, or replace the silicone sleeve.
5	Blue indicator flashes quickly	Low battery, please charge the product with the provided charging cable.
6	Unstable connection, it is easy to disconnect or freeze when using	Restart the earpick and reconnect to the APP.
7	After the product is connected to APP, image cannot be displayed normally	Disconnect the phone from the WiFi of the earpick, restart the earpick, and reconnect to WiFi.

If the above operation still fails to solve the problem, please contact our customer service staff. We will further provide you with practical solutions on APP-Personal Center-FAQ.

## Troubleshooting

S/N	Fault	Solution	S/N	Fault	Solution
1	No response after startup	Please charge with the provided charging cable. The charging duration is 40 minutes	5	Blue indicator flashes quickly	Low battery, please charge the product with the provided charging cable
2	Camera indicator does not light up	Shut down and restart or reconnect to WIFI	6	Unstable connection, it is easy to disconnect or freeze when using	Restart the earpick and reconnect to the APP
3	APP crashes	Uninstall the existing APP and rescan the QR code to install the latest version of APP	7	After the product is connected to APP, image cannot be displayed normally	Disconnect the phone from the WiFi of the earpick, restart the earpick, and reconnect to WiFi
4	Silicone sleeve is loose	Remove the silicone sleeve and slowly re-insert the earpick, or replace the silicone sleeve	If the above operation still fails to solve the problem, please contact our customer service staff, we will further provide you with practical solutions on APP-Personal Center-FAQ.		

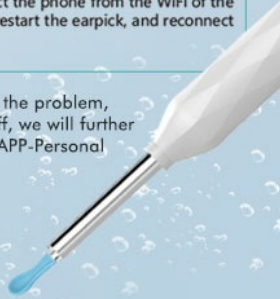


Image: Troubleshooting guide for common device issues.

## SPECIFICATIONS

- **Model Number:** P3

- **Manufacturer:** Behaesty
- **Package Dimensions:** 7.28 x 2.95 x 0.75 inches
- **Item Weight:** 2.4 ounces
- **Camera Resolution:** 800W Pixels FHD
- **Lens:** 360° Wide Angle, 6 LED Lights
- **Waterproof Rating:** IP67 (Lens)
- **Connectivity:** WiFi (2.4GHZ)
- **Battery Life:** Up to 90 minutes continuous use, 180 days standby
- **Charging Port:** USB Type-C

## WARRANTY AND SUPPORT

---

Behaesty products are designed for reliability and performance. This product comes with a standard manufacturer's warranty covering defects in materials and workmanship under normal use. For specific warranty details, please refer to the warranty card included in your package or contact Behaesty customer support.

For technical assistance, troubleshooting beyond this manual, or warranty claims, please contact Behaesty Customer Service through the 'Personal Center-FAQ' section within the 'Soulear' (or 'JesHome') app, or visit the official Behaesty store on Amazon.