



Manuals.plus /

› SKG /

› SKG EYE1 Bluetooth Heated Eye Massager User Manual

SKG EYE1

SKG EYE1 Bluetooth Heated Eye Massager User Manual

Model: EYE1

INTRODUCTION

This manual provides comprehensive instructions for the safe and effective use of your SKG EYE1 Bluetooth Heated Eye Massager. Please read this manual thoroughly before operating the device and retain it for future reference. The SKG EYE1 is designed to provide a soothing eye massage experience with heat therapy, air compression, and optional music playback.

IMPORTANT SAFETY INFORMATION

WARNING: To reduce the risk of electric shock, burns, fire, or injury:

- Do not use this device if you have any eye conditions such as glaucoma, cataracts, retinal detachment, or have undergone eye surgery recently. Consult a physician before use if you have any medical conditions.
- Keep out of reach of children. This device is not a toy.
- Do not immerse the device in water or other liquids.
- Do not use while driving or operating heavy machinery.
- Do not disassemble or attempt to repair the device yourself. Contact authorized service personnel.
- Use only the provided charging cable and adapter.
- Discontinue use immediately if you experience any discomfort or pain.
- Statements regarding this product have not been evaluated by the FDA and are not intended to

diagnose, treat, cure, or prevent any disease or health condition.

PRODUCT OVERVIEW

The SKG EYE1 eye massager features a sleek design with a comfortable fit, incorporating heat, air compression, and vibration massage. It includes a unique vision window for situational awareness.



Figure 1: Front view of the SKG EYE1 Bluetooth Heated Eye Massager. This image displays the device in its blue color, highlighting its ergonomic shape and the SKG logo.

Key Features:

- **Vision Window:** Allows for partial visibility during use.
- **Heat Therapy:** Gentle warmth (104-107°F / 40-42°C) for relaxation.
- **Air Compression Massage:** Targets acupoints around the eyes.
- **Vibration Massage:** 9-zone AirBreath vibration for temple massage.
- **Bluetooth Connectivity:** For personalized music playback.
- **Built-in Music:** Pre-loaded soothing tracks.
- **Foldable Design:** For easy portability and storage.
- **One-Button Operation:** Simplified control.

SETUP

1. Charging the Device

Before first use, fully charge the eye massager. Connect the provided charging cable to the device's charging port and a suitable USB power adapter. The charging indicator light will show the charging status. A full charge typically takes a few hours and provides approximately four to five 15-minute sessions.

Note: Do not leave the device plugged in after it is fully charged to prolong battery life.

2. Wearing the Massager

Ensure the device is powered off. Place the massager over your eyes, aligning the vision window. Adjust the elastic strap for a comfortable and secure fit. The device should feel snug but not overly tight.



Figure 2: Illustration of how to wear the SKG EYE1 eye massager, showing the adjustable headband for a custom fit.

3. Bluetooth Pairing (Optional)

To play your own music, activate Bluetooth on your smartphone or other audio device. Power on the SKG EYE1 massager. Search for "SKG EYE1" (or similar name) in your device's Bluetooth settings and pair. Once connected, you can stream audio to the massager.

OPERATING INSTRUCTIONS

1. Power On/Off

Press and hold the power button (usually located on the side of the device) for 2-3 seconds to turn the massager on or off.

2. Mode Selection

Once powered on, short press the power button to cycle through the available massage modes. The SKG EYE1 typically offers 4 built-in modes:

- **Vitality Mode:** Combines heat, air pressure, and vibration for an invigorating massage.
- **Soothing Mode:** Focuses on gentle air pressure and heat for relaxation.
- **Sleep Mode:** Primarily uses soft air pressure and minimal heat to aid in relaxation before sleep.
- **Eye Exercise Mode:** A specific sequence of massage and heat designed for eye relief.

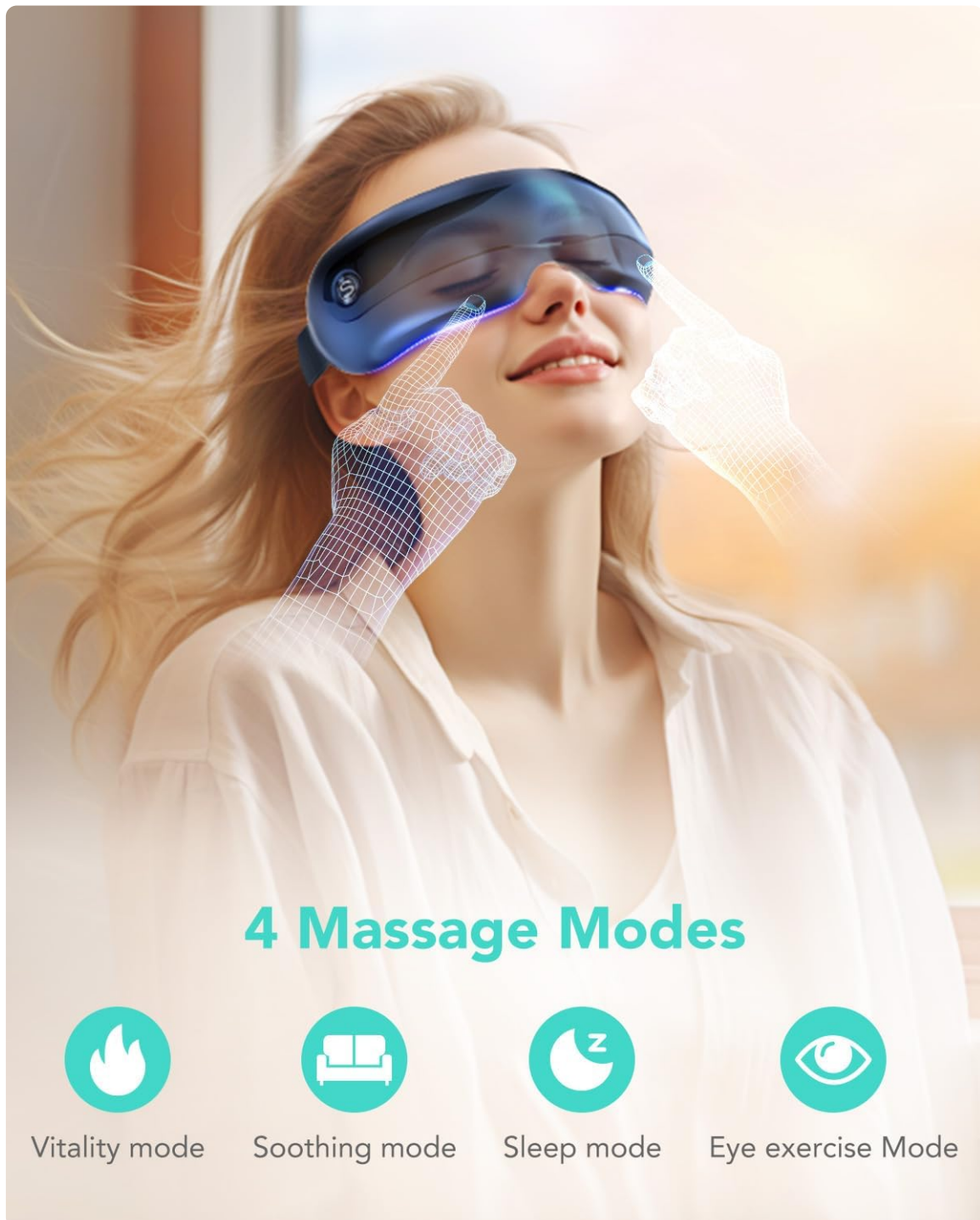


Figure 3: Visual representation of the four distinct massage modes available on the SKG EYE1 eye massager.

3. Heat Function

The heat function is integrated into most massage modes, providing a gentle warmth between 104-107°F (40-42°C). This helps to soothe tired eyes and promote blood circulation.

Constant Temperature



Figure 4: Diagram illustrating the constant temperature heat therapy provided by the SKG EYE1 eye massager, ranging from 104°F to 107°F.

4. Music Playback

The device comes with 3 built-in soothing music tracks. If you have paired your device via Bluetooth, your personal playlist will override the built-in music. To turn off the music, double-click the power button.



Figure 5: A user enjoying music while wearing the SKG EYE1 eye massager, demonstrating the Bluetooth music feature.

5. Vision Window

The integrated vision window allows for limited visibility, enabling you to perform simple tasks or be aware of your surroundings during the massage session.

Figure 6: A user demonstrating the visual window feature of the SKG EYE1 eye massager, which allows for unblocked sight during use.

MAINTENANCE

Cleaning

To clean the inner material of the massager, gently wipe it with a damp cloth. If necessary, a mild detergent can be used. Ensure the device is completely dry before storing or using again.

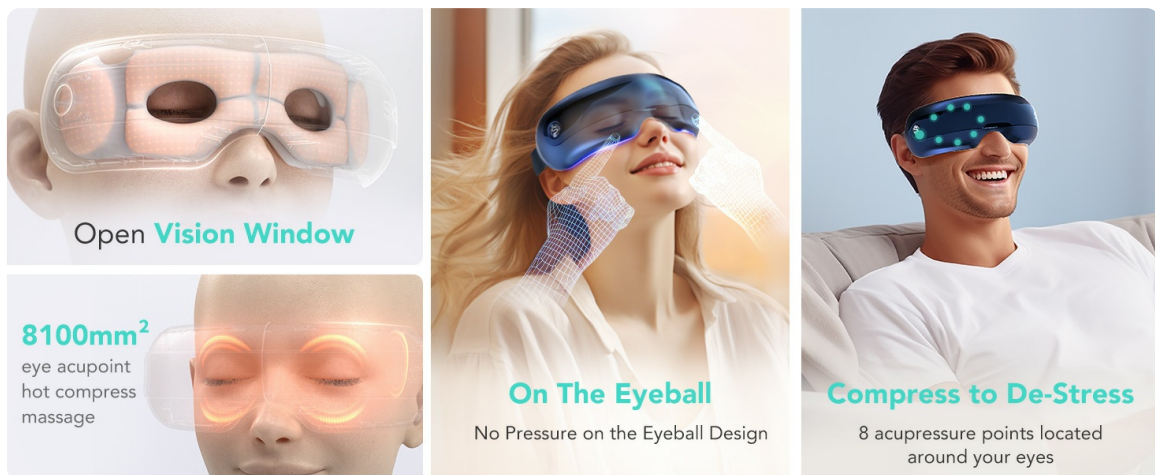


Figure 7: Visual guide on how to clean the inner surface of the SKG EYE1 eye massager using a wet towel.

Storage

Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. The foldable design allows for compact storage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully.
No heat or insufficient heat.	Not in a heat-enabled mode; device malfunction.	Ensure you are in a mode that includes heat therapy. If the issue persists, contact customer support.
Bluetooth not connecting.	Device not in pairing mode; Bluetooth off on audio device; interference.	Ensure the massager is on and your audio device's Bluetooth is enabled. Try restarting both devices.
Music is too loud/quiet or unwanted.	Volume settings; built-in music playing.	Adjust volume on your connected device. Double-click the power button to turn off built-in music.

SPECIFICATIONS

- **Model:** EYE1
- **Brand:** SKG
- **Dimensions:** 8.23 x 4.25 x 2.99 inches (approximate)
- **Weight:** 10.88 ounces (approximate)
- **Battery:** Rechargeable internal battery
- **Heat Temperature:** 104-107°F (40-42°C)
- **Special Features:** Auto Shut-Off, Foldable, Heat, Rechargeable
- **Material:** Plastic
- **Certifications:** CE, RF

WARRANTY AND SUPPORT

SKG products are designed and manufactured to high-quality standards. For warranty information or technical support, please refer to the warranty card included with your product or visit the official SKG website. You may also contact the seller directly for assistance.

Manufacturer: SKG

ASIN: B0CNG99BYV

Item Model Number: EYE1

