COMFYGO Z-4

COMFYGO Z-4 Electric Mobility Scooter User Manual

Model: Z-4

1. INTRODUCTION

This manual provides essential information for the safe operation, assembly, maintenance, and troubleshooting of your COMFYGO Z-4 Electric Mobility Scooter. Please read this manual thoroughly before using the scooter to ensure proper function and safety.

Image: The COMFYGO Z-4 Electric Mobility Scooter, highlighting its 350 lbs load capacity and 3-year frame warranty.

2. SAFETY INFORMATION

The COMFYGO Z-4 scooter is designed with user safety in mind. Adhere to the following guidelines to prevent injury or damage.

- Always operate the scooter on stable, level surfaces.
- Ensure all components are securely fastened before each use.
- The scooter features anti-roll technology for enhanced stability.
- Front and rear brake lights are integrated to improve visibility and safety during operation.
- A speed control device is included to regulate maximum speed, which is set to a maximum of 5 mph.
- Do not exceed the maximum weight capacity of 350 lbs.

PRECISION IN EVERY DIMENSION

Explore the Z-4

Compact, Powerful, Defined by Design

Image: Overview of Z-4 features including speed control and anti-roll technology.

3. PACKAGE CONTENTS

Upon unpacking, verify that all the following components are present:

- Main Frame (approx. 30 lbs)
- Rear Frame (approx. 26 lbs)
- Tiller Bar (approx. 8.5 lbs)
- Lithium-Ion Battery (approx. 8 lbs)
- Super Seat (approx. 17 lbs)
- Front Basket (approx. 2.5 lbs)
- Assembly Manual/Guide
- Charger

4. SETUP AND ASSEMBLY

The Z-4 Electric Mobility Scooter features a detachable frame for ease of transport and storage. Follow these steps for assembly:

1. **Secure the Main Frame:** Place the main frame on a stable surface.
2. **Connect the Rear Frame:** Attach the rear frame to the main frame at its designated connection points.
3. **Lower the Tiller Bar:** Insert the tiller bar into its mounting points and secure it.
4. **Insert Battery:** Place the lithium-ion battery into its compartment.
5. **Install Super Seat:** Mount the super seat onto the main frame.
6. **Attach Basket:** Secure the lightweight front basket to the tiller bar.

STEP-BY-STEP DISASSEMBLY GUIDE

Disassemble the Z-4 Mobility Scooter in a Snap

SECURE
the main
frame
in place

CONNECT
the rear frame
to the main
frame

LOWER
the Tiller Bar
on its designated
mounting points

INSERT
lithium-ion
battery into its
compartment

INSTALL
the Super Seat
onto the main
frame

ATTACH
lightweight
basket on the
mounting points

Main Frame	Rear Frame	Tiller Bar	Battery	Super Seat	Basket
					
30 lbs	26 lbs	8.5 lbs	8 lbs	17 lbs	2.5 lbs

Image: Visual guide for disassembling the Z-4 scooter, illustrating each removable section and its approximate weight.

5. OPERATING INSTRUCTIONS

5.1 Controls

The scooter is operated via controls located on the tiller bar. Familiarize yourself with the throttle, brake levers, and speed adjustment dial before operating.

5.2 Driving

To move forward, gently engage the throttle. Release the throttle to slow down, and use the brake levers for immediate stopping. The scooter has a maximum speed of 5 mph and can handle inclines up to 12 degrees.

5.3 Travel Considerations

The Z-4 scooter is designed for convenient travel. Its foldable and detachable sections allow it to fit into most standard car trunks. The 12AH battery version is authorized for airline and cruise travel. Always confirm specific travel regulations with your airline or cruise line prior to travel.

EFFORTLESS TRAVEL

From Airports to City Streets

AIRLINE APPROVED

CRUISE APPROVED

STANDARD CAR TRUNK FIT

*Acceptable on cruises and airlines with up to 12AH battery version.

Image: The Z-4 scooter demonstrating its suitability for various travel methods, including airline, cruise, and car trunk transport.

6. BATTERY INFORMATION

The Z-4 scooter is powered by a lithium-ion battery. Two battery options are available, offering different ranges:

- **12 AH Battery:** Provides a range of up to 13 miles, suitable for standard usage and shorter trips.
- **30 AH Battery:** Provides a range of up to 30 miles, suitable for longer trips and increased usage between charges.

6.1 Charging the Battery

Connect the charger to the scooter's charging port and then to a standard electrical outlet. The charging indicator light will show the charging status. Disconnect the charger once the battery is fully charged. Avoid overcharging.

POWER YOUR ADVENTURE

Choose Your Range With Different Battery Options

12 AH BATTERY LITHIUM-ION	30 AH BATTERY LITHIUM-ION
  Up to 13 miles*  Suitable for standard usage and shorter trips.	  Up to 30 miles*  Suitable for long trips and increased usage between charges.

Image: Comparison of the 12 AH and 30 AH lithium-ion battery options and their respective ranges.

7. MAINTENANCE

Regular maintenance ensures the longevity and optimal performance of your scooter.

- **Cleaning:** Wipe down the scooter with a damp cloth. Avoid using abrasive cleaners or direct water spray.
- **Battery Care:** Store the battery in a cool, dry place. Charge the battery regularly, even if the scooter is not in frequent use, to maintain its health.
- **Tire Pressure:** Periodically check tire pressure and inflate to the recommended levels.
- **Inspections:** Regularly inspect all nuts, bolts, and connections for tightness.

8. TROUBLESHOOTING

If you encounter issues with your Z-4 scooter, refer to the following common problems and solutions:

- **Scooter does not power on:** Ensure the battery is properly installed and charged. Check all electrical connections.
- **Reduced range:** The battery may not be fully charged, or its capacity may be diminishing over time. Ensure proper charging cycles.
- **Unusual noises:** Inspect for loose components or obstructions in the wheels or motor area.
- **Difficulty steering:** Check tire pressure and ensure the tiller bar is securely fastened.

For persistent issues, contact customer support.

9. SPECIFICATIONS

Feature	Detail
Product Dimensions	42 x 20 x 36 inches (Length x Width x Height)
Item Weight	92 Pounds
Weight Capacity	Up to 350 lbs
Motor Power	250-watt
Max Speed	Up to 5 mph
Max Climbing Degree	Up to 12 degrees
Battery Type	Lithium Ion (1 included)
Battery Range (12 AH)	Up to 13 miles
Battery Range (30 AH)	Up to 30 miles
Material	Metal
Seat Width	22 inches (Adjustable)
Seat Height	17.5 - 19 inches

US BASED COMPANY CALIFORNIA DISPATCH

Image: Detailed dimensions of the Z-4 scooter.

ComfyGO Z-4

Others

- ✓ **Lighter** (8-12 lbs)
- ✓ **Smaller Size**
- ✓ **High Range** (up to 30 miles)

- ✗ **Heavier** (+30 lbs)
- ✗ **Larger Size**
- ✗ **Normal Range** (up to 15 Miles)

MINIMUM WEIGHT
MAXIMUM RANGE

Image: Weight and size comparison of the Z-4 battery with other models.

10. WARRANTY AND SUPPORT

COMFYGO is committed to providing high-standard customer support.

10.1 Warranty Information

- **Metal Frame:** 3-year warranty.
- **Spare Parts:** 1-year warranty.
- **Battery:** 1-year warranty.

10.2 Customer Support

For any inquiries, technical assistance, or warranty claims, please contact our dedicated service team. We aim to design and sell innovative mobility products, restoring mobility and ease to people's lives. Each unit undergoes a three-stage quality control process and a final pre-shipment inspection.


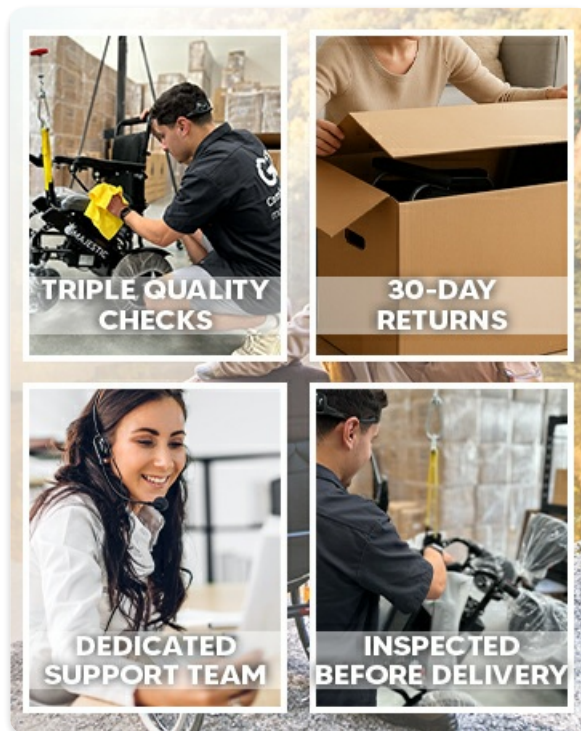
Image: COMFYGO's commitment to quality and customer support, including triple quality checks and a dedicated support team.

10.3 Frequently Asked Questions (FAQs)

- **What is our mission?** We aim to design and sell innovative mobility products, restoring mobility and ease to people's lives.
- **What values do we prioritize?** We value customer needs, trust, passion for products, and compassion for customers, committed to enhancing lives.
- **How do we contribute to society and the environment?** We integrate with organizations to understand societal and environmental needs, support communities and charities, provide employee training, and contribute to the national economy.

Related Documents - Z-4

Model BC-MS310 Instruction Manual	ComfyGo BC-MS310 Foldable Mobility Scooter: Instruction Manual Comprehensive instruction manual for the ComfyGo BC-MS310 foldable mobility scooter, covering product features, operation, maintenance, safety guidelines, and troubleshooting.
---	--

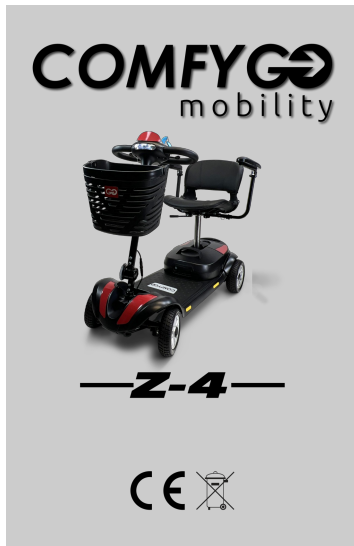
[\[pdf\]](#) User Manual Instructions Warranty

Z 4 EU User manual V1 0 Manuel ComfyGO Mobility Manual EN drive google file d

1A1ISUO9MdQpjl7Y6D77AmJAuAoo8MYk view usp link |||

User Manual / Warranty Table of Contents Preface 2 Product Performance 3 Safety Instructions ... ains the operation, assembly methods, and simple faults solutions. The manual applies to our model: **Z-4** - This manual contains scooter maintenance information and self-checking methods. Please be sure...

lang:en score:33 filesize: 2.12 M page_count: 36 document date: 2022-12-19

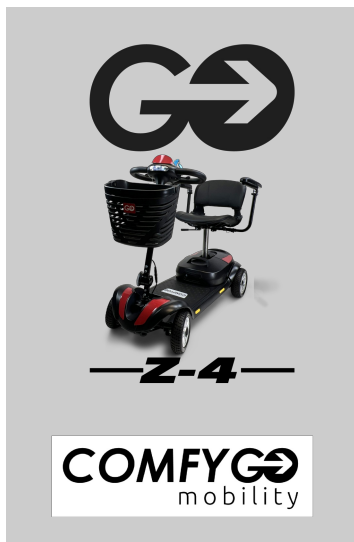
[\[pdf\]](#) User Manual Instructions Warranty

Z 4 EU User manual V1 0 ComfyGo Mobility Scooter Wheels Chairs EN drive google file d 10Q

X65NfalLcs7IWVmlfas3TANda6VLz preview |||

User Manual / Warranty Table of Contents Preface 2 Product Performance 3 Safety Instructions ... ains the operation, assembly methods, and simple faults solutions. The manual applies to our model: **Z-4** - This manual contains scooter maintenance information and self-checking methods. Please be sure...

lang:en score:30 filesize: 1.93 M page_count: 36 document date: 2023-02-28

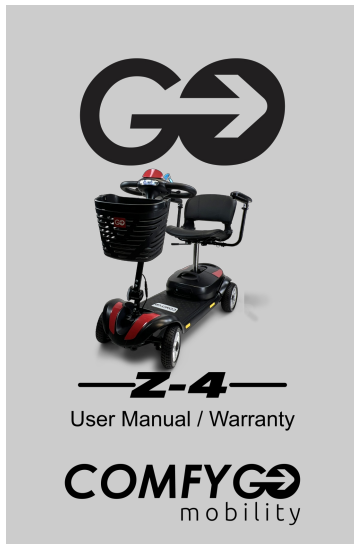[\[pdf\]](#) User Manual Owner's Manual Instructions Warranty

Z 1 User manual V1 Owner s ManualComfyGo 4 Mobility Scooter with Quick Detach Frame —

DailyMobilityShopComfyGo user Manualcdn shopify files 0834 5747 4835 ComfyGo Manual v

1743275385 ManualOwner Ultra Light Electric The Prestige AidComfyGo ManualZ 0cdn 1743275385cdn


0866 0779 0390 1743053227 ||| ||| User Manual / Warranty Table of Contents Preface 2 Product

Performance 3 Safety Instructions 5 Usage and Operation 9 Battery and Charger 14 System Diagnostics

17 Controller 19 Maintenance Areas 21 Tech Support 23 Warranty 24 For all your sp

User Manual / Warranty Table of Contents Preface 2 Product Performance 3 Safety Instructions ... ains the operation, assembly methods, and simple faults solutions. The manual applies to our model: **Z-4** - This manual contains scooter maintenance information and self-checking methods. Please be sure...

lang:en **score:30** filesize: 12.66 M page_count: 36 document date: 2022-04-21

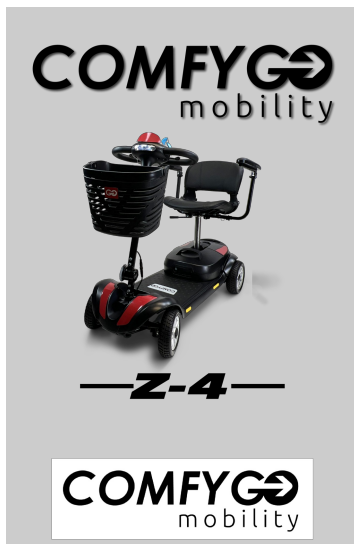[\[pdf\]](#) User Manual

Dutch Download User Manuel Z 4 ComfyGO Mobility Manual NL drive google file d

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lang:nl **score:29** filesize: 1.98 M page_count: 36 document date: 2022-12-19

—Z-4—

User Manual / Warranty

COMFYGO
mobility

[\[pdf\]](#) User Manual Instructions Warranty

Z 1 User manual V1 0 Manual ComfyGo 4 Electric Powered Mobility Scooter Lightweight Detachable F —
Nest v 1686984812 cdn shopify s files 0556 1463 6322 |||

User Manual / Warranty Table of Contents Preface 2 Product Performance 3 Safety
Instructions ... ains the operation, assembly methods, and simple faults solutions. The
manual applies to our model: **Z-4** - This manual contains scooter maintenance
information and self-checking methods. Please be sure...

lang:en **score:28** filesize: 12.66 M page_count: 36 document date: 2022-04-21

—Z-4—

COMFYGO
mobility

[\[pdf\]](#) User Manual

Z 1 User manual V1 0 German Download Manuel 4 ComfyGO Mobility Manual DE drive google file d
1E3jdF4sNLYKKwC5wa5iRwt5yZFYb7 i8 view usp link |||

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lang:de **score:27** filesize: 1.38 M page_count: 36 document date: 2022-12-19

—Z-4—

Manuale d'uso / Garanzia


COMFYGO
mobility

[\[pdf\]](#) User Manual

Z 1 User manual V1 0 Italian Download Manuel 4 ComfyGO Mobility Manual IT drive google file d
1hO68nT9atmQer KxYDQc6tvL 9Ht 5Lv view usp link |||

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lang:it **score:27** filesize: 1.49 M page_count: 36 document date: 2022-12-19

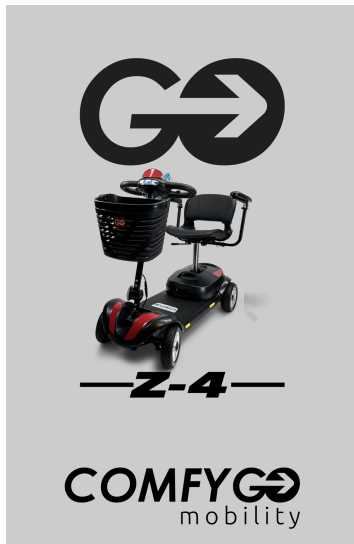
[\[pdf\]](#) User Manual Instructions

Z 1 User manual V1 0 Merve Civelek French Download Manuel 4 ComfyGO Mobility Manual FR drive
google file d 1eS2rlr9LSD3A Q8Vb oXxmQJSJ8XChFXk view usp link |||

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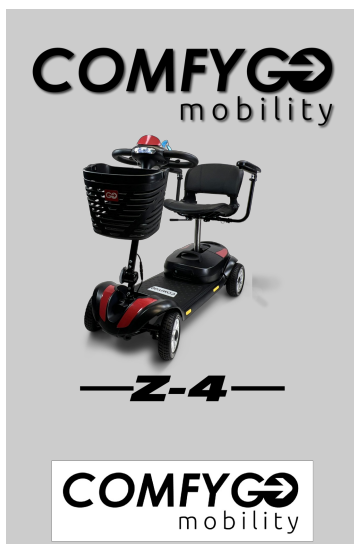lang:fr **score:26** filesize: 1.45 M page_count: 40 document date: 2022-12-19

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