

DiGGiNGSOUND T015

DiGGiNGSOUND T015 Digital Bluetooth Rechargeable Hearing Aids User Manual

Model: T015

INTRODUCTION

Thank you for choosing the DiGGiNGSOUND T015 Digital Bluetooth Rechargeable Hearing Aids. This manual provides detailed instructions for the proper setup, operation, and maintenance of your new hearing aids. Please read this manual thoroughly before using the device to ensure optimal performance and longevity.

PRODUCT OVERVIEW

The DiGGiNGSOUND T015 hearing aids are designed to amplify sound for individuals with mild to moderate-severe hearing loss. They feature Bluetooth connectivity, a rechargeable design, and digital 16-channel sound processing for clear audio. The compact and lightweight design ensures comfortable wear.



Image: The DiGiGiNGSOUND T015 hearing aids shown in their white charging case. The left earbud is marked with 'L' and the right earbud has a red ring, indicating their respective sides.

COMFORTABLE FIT & ULTRA-LIGHT



SUPERIOR
SOUND



DIGITAL
16-CHANNEL



HANDS FREE
CALLING



ONE-TOUCH
CONTROL



Image: An illustration highlighting key features: Superior Sound, Digital 16-Channel processing, Hands-Free Calling, One-Touch Control, and an ultra-light weight of 0.15 oz per earphone.

PACKAGE CONTENTS

Please check the package for the following items:

- DiGiGINGSOUND T015 Hearing Aids (Left and Right)
- Portable Charging Case
- USB Type-C Charging Cable
- Cleaning Brush
- User Manual (this document)

SETUP GUIDE

1. Initial Charging

Before first use, fully charge the hearing aids and the charging case.

1. Place both hearing aids into their respective slots in the charging case. Ensure the L (left) and R (right) indicators align.

2. Connect the provided USB Type-C cable to the charging port on the case and plug the other end into a USB power adapter (not included) or a computer's USB port.
3. The indicator lights on the charging case will show the charging status. A full charge typically takes approximately 2 hours.
4. Once fully charged, disconnect the charging cable. The case can provide 2-3 full charges for both units, offering up to 24 hours of working time.



Image: Details on the portable wireless charging case, showing charging time (≤2 hours), hearing aid usage time (>7 hours), music time (>8 hours), call time (>5 hours), and total working time (24 hours).

2. Wearing the Hearing Aids

Identify the left (L) and right (R) hearing aids. Gently insert each hearing aid into your ear canal until it feels secure and comfortable. The devices are designed to automatically turn on when placed in the ear.

3. Bluetooth Pairing

The hearing aids can be paired with your smartphone or other Bluetooth-enabled devices for calls and music.

1. Ensure the hearing aids are charged and in your ears.
2. Open the Bluetooth settings on your device.
3. Search for available devices and select "DiGGiNGSOUND T015" (or similar name).

4. Once connected, you will hear a confirmation tone.

OPERATING INSTRUCTIONS

Volume Adjustment (Hearing Amplification)

The T015 hearing aids offer 3 levels of hearing amplification. You can adjust the level by a single touch on the earpiece.



Image: A diagram illustrating the three volume levels (Level 1, Level 2, Level 3) and how to cycle through them with a single press on the earpiece. Separate adjustment is possible for each hearing aid.

- **Single Press:** Tap the touch control area on either hearing aid to cycle through Level 1, Level 2, and Level 3. A voice prompt will announce the current level.
- The default volume from the factory is Level 1.
- Each hearing aid can be adjusted independently.

AMPLIFY CLEAR SOUND UP TO 35dB

best choice for mild to moderate-severe hearing loss



Image: Visual representation of the hearing aids' ability to amplify clear sound up to 35dB, suitable for mild to moderate-severe hearing loss.



Image: Examples of the DiGGiNGSOUND T015 hearing aids being used in different environments, such as watching TV, during daily communication, at family gatherings, and while working.

Smart Touch Control (Bluetooth Functions)

The hearing aids feature smart touch controls for managing calls and other functions when connected via Bluetooth.

SMART TOUCH CONTROL

Hearing Amplifier



Answer



Hang up



Reject



Next volume
level



Please don't cover the microphone or it'll whistling.

Image: An illustration demonstrating the smart touch control on the hearing aid, showing icons for Answer, Hang up, Reject calls, and adjusting to the Next volume level. It also warns not to cover the microphone.

- **Answer Call:** Single press the touch control area.
- **Hang Up Call:** Single press the touch control area during a call.
- **Reject Call:** Press and hold the touch control area.
- **Next Volume Level:** Single press to cycle through amplification levels (as described above).

Important: Do not cover the microphone area on the hearing aid, as this may cause whistling or feedback.

MAINTENANCE AND CARE

Proper maintenance ensures the longevity and optimal performance of your hearing aids.

- **Cleaning:** Use the provided cleaning brush to gently remove any earwax or debris from the earpieces and sound outlets. Wipe the hearing aids and charging case with a soft, dry cloth. Do not use water or cleaning solutions directly on the devices.
- **Storage:** When not in use, store the hearing aids in their charging case to protect them from dust and damage.
- **Avoid Moisture:** Keep the hearing aids and charging case dry. Avoid exposure to excessive

moisture, humidity, or extreme temperatures.

- **Battery Care:** Recharge the hearing aids regularly, even if not used daily, to maintain battery health.

TROUBLESHOOTING

If you encounter any issues, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
No sound or weak sound	Low battery; Earwax blockage; Incorrect volume level.	Charge the hearing aids; Clean the earpiece and sound outlet; Adjust to a higher volume level.
Whistling or feedback noise	Improper fit; Microphone covered; Volume too high.	Re-insert the hearing aid for a snug fit; Ensure nothing is covering the microphone; Reduce the volume level.
Cannot connect via Bluetooth	Hearing aids not in pairing mode; Bluetooth off on device; Already connected to another device.	Ensure hearing aids are on; Turn Bluetooth off and on again on your device; Disconnect from other devices.
Charging case not charging	Loose cable connection; Faulty cable/adaptor.	Check USB-C cable connection; Try a different USB-C cable or power adapter.

If the problem persists after trying these solutions, please contact customer support.

SPECIFICATIONS

- **Model:** T015
- **Brand:** DiGGiNGSOUND
- **Color:** White
- **Material:** Silicone
- **Product Dimensions:** 1.38 x 0.7 x 0.04 inches
- **Weight (per earphone):** Approximately 0.15 ounces
- **Amplification:** Up to 35dB
- **Channels:** 16-channel digital sound processing
- **Charging Time:** Approximately 2 hours (for case and aids)
- **Hearing Amplification Time:** Up to 7 hours
- **Music Playing Time:** Up to 8 hours
- **Call Time:** Up to 5 hours
- **Total Working Time (with charging case):** Up to 24 hours
- **Connectivity:** Bluetooth

WARRANTY AND SUPPORT

DiGGiNGSOUND products come with a standard manufacturer's warranty. For specific warranty details, please refer to the warranty card included in your product packaging or visit the official DiGGiNGSOUND website.


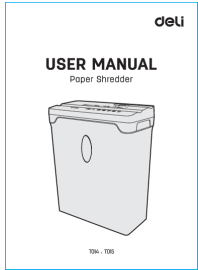



For technical support, troubleshooting assistance, or any questions regarding your DiGGiNGSOUND T015 hearing aids, please contact our customer service team:

- **Email:** support@diggingsound.com
- **Website:** www.diggingsound.com

Please have your model number (T015) and purchase information ready when contacting support.

© 2024 DiGGiNGSOUND. All rights reserved.

Related Documents - T015

	<p>DiggingSound T002 Hearing Aids User Manual</p> <p>User manual for DiggingSound T002 rechargeable hearing aids, covering setup, operation, troubleshooting, and specifications. Learn how to wear, connect, and adjust your hearing aids for optimal performance.</p>
	<p>Deli T014/T015 Paper Shredder User Manual</p> <p>User manual for the Deli T014 and T015 paper shredders, providing essential information on setup, operation, safety precautions, maintenance, product specifications, troubleshooting, and warranty.</p>
	<p>DayBetter WiFi Smart Plug User Manual - T015</p> <p>User manual for the DayBetter T015 WiFi Smart Plug. Learn how to set up, control your appliances remotely via the DayBetter app, integrate with Alexa and Google Home, and find troubleshooting tips.</p>
	<p>Mountain Top EVOM Roll Cover for Pickup Trucks - Features & Specifications</p> <p>Discover the Mountain Top EVOM, a next-generation roll cover for pickup trucks. Featuring advanced security, innovative design, weather resistance, and integrated accessory tracks. Explore compatibility and part numbers for various vehicle models.</p>
	<p>Calefactio Product Catalog 2019: Hydronic Heating and Plumbing Solutions</p> <p>Explore Calefactio's 2019 Product Catalog featuring a comprehensive range of hydronic heating and plumbing solutions, including expansion tanks, air and dirt separators, specialized ASME tanks, glycol make-up packages, condensate neutralizers, pump flanges, ball valves, and flow regulators for residential and commercial applications.</p>