

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Hybodies](#) /

› [Hybodies BDCD-24-AL-1-NEW Heavy Duty Bariatric Commode Instruction Manual](#)

Hybodies BDCD-24-AL-1-NEW

Hybodies Heavy Duty Bariatric Commode Instruction Manual

Model: BDCD-24-AL-1-NEW | Brand: Hybodies

1. INTRODUCTION

Thank you for choosing the Hybodies Heavy Duty Bariatric Commode. This versatile 4-in-1 product is designed to provide comfort, stability, and convenience for individuals requiring assistance with mobility and personal care. It can be used as a bedside commode, a supportive shower chair, an over-toilet chair, or a height extender for existing toilets. This manual provides essential information for safe assembly, operation, maintenance, and troubleshooting.



ACCESSIBLE LIVING, ELEVATED COMFORT



Widened
Design



Multifunctional
Use



Easy
Assembly



Heavy Duty
Aluminium



Senior



Post-op



Pregnancy



Limited Mobility

Image: The Hybodies Heavy Duty Bariatric Commode in its standard configuration.

This commode is especially recommended for the elderly, postoperative individuals, pregnant women, and those with limited mobility, offering a discreet and practical toileting solution that preserves independence and dignity.

2. SAFETY INFORMATION

- **Weight Capacity:** Do not exceed the maximum weight capacity of 350 lbs (158 kg).
- **Stability:** Ensure all legs are evenly adjusted and the commode is stable on a flat, non-slip surface before use.

- **Assembly:** Follow all assembly instructions carefully. Ensure all screws and fasteners are securely tightened.
- **Inspection:** Regularly inspect the commode for any signs of wear, damage, or loose components. Discontinue use if any issues are found and contact customer support.
- **Transferring:** Use caution when transferring to and from the commode. Utilize armrests for support.
- **Cleaning:** Clean the commode regularly with mild soap and water. Avoid abrasive cleaners that may damage the material.

Your browser does not support the video tag.

Video: A load test demonstrating the durability and stability of the Hybodies Heavy Duty Bariatric Commode under significant weight.

3. PACKAGE CONTENTS

Before assembly, ensure all components are present and undamaged. Refer to the included instruction leaflet for a detailed parts list.



Image: All components of the Hybodies Bariatric Commode laid out for assembly.

4. SETUP AND ASSEMBLY

Follow these steps to assemble your Hybodies Bariatric Commode. A detailed step-by-step video is also provided for visual guidance.

1. **Unpack Components:** Carefully remove all parts from the packaging and lay them out on a clean, flat surface.
2. **Assemble Frame:** Connect the main frame components as per the diagrams in the instruction manual. Ensure all locking pins click into place.
3. **Attach Armrests:** Secure the padded armrests to the main frame using the provided screws and tools.
4. **Install Backrest:** Attach the comfortable backrest to the rear of the frame.
5. **Adjust Legs:** Insert the adjustable leg extensions and set them to the desired height. Ensure all four legs are at the same height for stability.
6. **Place Seat and Bucket:** Position the PU seat cushion onto the frame. The commode bucket slides into the designated holder beneath the seat.
7. **Attach Accessories:** Secure the tissue holder and storage bag to the frame as desired.

Your browser does not support the video tag.

Video: A step-by-step tutorial demonstrating the assembly process of the Hybodies Bedside Commode.

5. OPERATING INSTRUCTIONS

5.1 4-in-1 Versatility

The Hybodies commode offers four primary modes of use:

- **Bedside Commode:** Use with the included bucket for convenient toileting near the bed.
- **Shower Chair:** The aluminum frame and waterproof seat make it suitable for use in the shower.
- **Over-Toilet Chair:** Remove the bucket and position the commode over an existing toilet to provide support and raise the seat height.
- **Raised Toilet Seat/Handrail:** Adjust the height to match your toilet for easier sitting and standing, utilizing the sturdy armrests as handrails.

Bedside Commode



Image: The commode configured for bedside use.

Shower Chair



Image: The commode used as a shower chair.



Toilet Raised Seat



Image: The commode positioned as a raised toilet seat.

5.2 Adjustments

The commode features adjustable height and width to accommodate various user needs and environments.

- **Height Adjustment:** Press the push-pins on each leg and slide the inner tube to adjust to one of the 8 available height levels. Ensure all pins are securely engaged before use.
- **Width Adjustment:** The seat width can be adjusted by loosening the side fasteners and extending or retracting the seat frame. Securely tighten fasteners after adjustment.



Image: Visual representation of the commode's adjustable height and width features.

5.3 Thoughtful Features

- **Commode Bucket:** The deep, collapsible bucket with a lid is designed for easy disposal and hygiene.
- **Storage Bag:** A convenient side storage bag is included for personal items or hygiene products.
- **Tissue Holder:** An integrated tissue holder keeps toilet paper within easy reach.
- **Soft Padded Cushion:** The extra-wide PU seat cushion provides enhanced comfort.
- **Comfortable Backrest:** The ergonomic backrest offers additional support and relaxation.
- **Shower Head Clamp:** A clamp is provided to hold a shower head for easy rinsing when used as a shower chair.

Commode Bucket

Dispose of waste easily, promoting cleanliness and hassle-free use.



Image: Detailed view of the commode bucket.

Storage Bag

Store your personal items or hygiene products.



Image: The commode with its attached storage bag.

Tissue Holder

Keep tissues within reach, ensuring convenience and hygiene.



Image: The commode featuring its integrated tissue holder.

6. MAINTENANCE AND CLEANING

- **Regular Cleaning:** Clean all surfaces of the commode regularly with a mild, non-abrasive household cleaner and a soft cloth.

- **Disinfect:** For thorough cleaning, use a disinfectant suitable for bathroom surfaces.
- **Bucket Cleaning:** Empty and clean the commode bucket after each use. Use a brush and disinfectant for optimal hygiene.
- **Drying:** Ensure all parts, especially the frame and seat, are thoroughly dried after cleaning to prevent water spots or corrosion.
- **Storage:** Store the commode in a dry place when not in use.

7. TROUBLESHOOTING


- **Instability:** If the commode feels unstable, check that all leg adjustment pins are fully engaged and that the legs are set to equal heights. Ensure it is on a flat surface.
- **Difficulty Adjusting Height/Width:** Ensure push-pins are fully depressed when adjusting height. For width, loosen fasteners completely before attempting to slide the frame.
- **Bucket Not Fitting:** Ensure the collapsible bucket is fully expanded before placing it in the holder. Check for any obstructions.

8. SPECIFICATIONS

Feature	Detail
Model Number	BD-CD-24-AL-1-NEW
Material	Aluminum
Maximum Weight Recommendation	350 Pounds (158 kg)
Item Weight	11.5 Pounds (5.2 kg)
Product Dimensions (D x W x H)	22"D x 28.1"W x 16.3"H (55.9cm D x 71.4cm W x 41.4cm H)
Seat Depth	22 inches (55.9 cm)
Height Levels	8 adjustable levels

ERGONOMIC AND ADJUSTABLE

Shaped for different individuals
and obese people

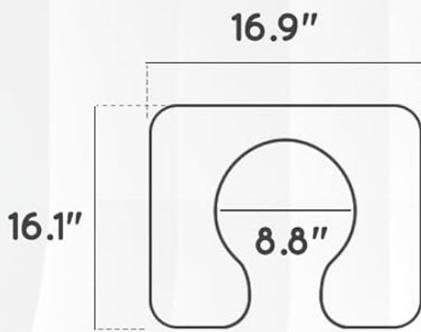

Please confirm the
size careful before purchase.



Image: Detailed dimensions of the Hybodies Bariatric Commode.

9. WARRANTY AND SUPPORT

Hybodies is committed to providing high-quality products. For any questions, concerns, or warranty claims, please contact our customer support team. Refer to your purchase documentation for specific warranty details.

Contact Information: Please refer to the contact details provided with your product packaging or visit the official Hybodies website for support.

© 2026 Hybodies. All rights reserved.