

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

- › [eWheels](#) /
- › [EW-36, Fully Assembled, 45 Mile Range, 350 Lb Capacity, High Performance Long Range Recreational Mobility Scooter for Seniors and Adults, \(Blue\) - User Manual](#)

eWheels EW-36

EW-36 Recreational Mobility Scooter User Manual

Model: EW-36 | Brand: eWheels

INTRODUCTION

Welcome to the user manual for your new eWheels EW-36 Recreational Mobility Scooter. This manual provides essential information for the safe and effective operation, maintenance, and care of your scooter. Please read it thoroughly before using the product.



Image: The eWheels EW-36 Recreational Mobility Scooter in blue, showcasing its design and features from a front-side perspective.

The EW-36 is designed for outdoor adventures, offering high performance with a weight capacity of 350 lbs and rear-wheel drive. It provides an extended range and is equipped with features for safety and convenience.

SETUP AND INITIAL USE

Your EW-36 scooter arrives fully assembled, ensuring a hassle-free setup and immediate use upon delivery. Follow these steps for initial preparation:

1. **Unpack the scooter:** Carefully remove all packaging materials.
2. **Charge the battery:** Ensure the scooter is fully charged before its first use. Connect the charger to the scooter's charging port and a standard electrical outlet. The charging indicator will show when the battery is full.
3. **Adjust the seat and handlebars:** Customize the seat and handlebar positions for your comfort and optimal control. The adjustable saddle seat can accommodate up to two passengers.
4. **Familiarize yourself with controls:** Before riding, locate and understand the function of the twist throttle, adjustable hand brakes, turn signals, high/low beam switch, and rearview mirror.



- 45 Mile Range
- 400 Lb Capacity
- Arrives Fully Assembled

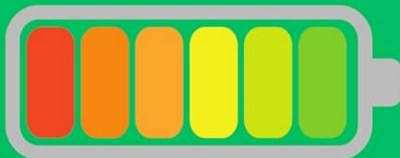


Image: A blue EW-36 scooter with text overlays indicating a 45-mile range, 400 lb capacity (note: product specifications state 350 lbs), and that it arrives fully assembled. A battery charge indicator graphic is also shown.

Comfortable Reclinable Captains Chair



Image: A close-up view of the comfortable, reclinable captain's chair on the EW-36 scooter, emphasizing its ergonomic design.

OPERATING INSTRUCTIONS

The EW-36 is designed for intuitive operation. Always prioritize safety and adhere to local traffic regulations.

1. **Power On:** Press the power button to turn on the scooter.
2. **Acceleration:** Slowly twist the throttle to accelerate. The scooter can reach speeds of up to 15 mph.
3. **Braking:** Use the adjustable hand brakes to slow down or stop the scooter.
4. **Signaling:** Utilize the turn signals for indicating turns and the high/low beam for appropriate lighting.
5. **Storage:** The scooter includes lockable rear and under-seat storage for your convenience.

Easy to Operate



Image: A detailed view of the EW-36 scooter's controls, including the headlight switch, turn signal switch, and a dial for speed adjustment, labeled "Easy to Operate".

Convenient Storage



Image: The open under-seat storage compartment of the EW-36 scooter, demonstrating its "Convenient Storage" feature.

MAINTENANCE

Regular maintenance ensures the longevity and optimal performance of your EW-36 scooter.

- **Battery Care:** Charge the scooter after each use or when the battery level is low to maintain battery health and extend its lifespan. Avoid fully depleting the battery.
- **Tire Pressure:** Periodically check the pressure of the tubeless tires and inflate to the recommended PSI for a comfortable and smooth journey.
- **Cleaning:** Keep the scooter clean by wiping it down with a damp cloth. Avoid using harsh chemicals or direct water spray on electrical components.
- **Brake Inspection:** Regularly inspect the adjustable hand brakes for proper function and wear. Adjust or replace as needed.
- **General Check:** Before each ride, quickly check all lights, signals, and the horn to ensure they are working

correctly.

TROUBLESHOOTING

This section addresses common issues you might encounter with your EW-36 scooter.

Problem	Possible Cause	Solution
Scooter does not power on.	Low battery; power button not pressed; loose battery connection.	Charge the battery; ensure power button is fully pressed; check battery connections.
Reduced range or speed.	Undercharged battery; low tire pressure; heavy load.	Fully charge battery; inflate tires to recommended pressure; reduce load.
Brakes feel loose or ineffective.	Brake cable needs adjustment; brake pads worn.	Adjust brake cables; inspect and replace brake pads if necessary. Consult a technician if unsure.
Unusual noises during operation.	Loose components; motor issue; suspension issue.	Inspect for loose parts and tighten; if noise persists, contact customer support.

SPECIFICATIONS

Key technical specifications for the eWheels EW-36 Recreational Mobility Scooter:

- **Model Name:** EW-36
- **Brand:** eWheels
- **Color:** Blue (as per product listing)
- **Maximum Weight Recommendation:** 350 Pounds
- **Top Speed:** Up to 15 mph
- **Travel Range:** Up to 40 miles on a single charge
- **Assembly:** Fully Assembled
- **Suspension:** Front and Rear
- **Tires:** Tubeless
- **Safety Features:** Turn signals, brake lights, high/low beam, rearview mirror, horn, alarm, rear reflectors.
- **Storage:** Lockable rear and under-seat storage, bag hook.

WARRANTY INFORMATION

Your eWheels EW-36 Electric Mobility Scooter comes backed by a **3-year limited warranty**. This warranty covers potential manufacturing defects or issues that may arise, providing reassurance and support for years to come.







For detailed terms and conditions of the warranty, please refer to the warranty card included with your purchase or visit the official eWheels website.

Please Note: Due to the product's high shipping weight, returns are subject to restocking fees.

LEGAL DISCLAIMER

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a recreational mobility device and should be used in accordance with all applicable laws and regulations.

Related Documents - EW-36

 <p>EW-26 Owner's Manual</p>	<p>EWheels EW-26 Electric Scooter Owner's Manual</p> <p>Comprehensive owner's manual for the EWheels EW-26 electric scooter, covering features, technical specifications, operation, charging, maintenance, safety guidelines, and warranty information.</p>
 <p>EW-19 Owner's Manual</p>	<p>EWheels EW-19 Owner's Manual</p> <p>Comprehensive owner's manual for the EWheels EW-19 recreational mobility scooter, covering intended use, safety guidelines, operation, maintenance, technical specifications, and warranty information.</p>
 <p>EW-20 Owner's Manual</p>	<p>EWheels EW-20 Mobility Scooter Owner's Manual</p> <p>This owner's manual provides comprehensive information for the EWheels EW-20 mobility scooter, including features, technical specifications, operating instructions, charging and maintenance procedures, safety guidelines, and warranty details.</p>
 <p>EW-14 Owner's Manual</p>	<p>EWheels EW-14 Owner's Manual: Operation, Safety, and Maintenance Guide</p> <p>Comprehensive guide for the EWheels EW-14 recreational mobility scooter, covering intended use, safety warnings, operation, technical specifications, battery care, charging, and warranty information.</p>
 <p>EW-11 Owner's Manual</p>	<p>EWheels EW-11 Sport Electric Scooter Owner's Manual</p> <p>Comprehensive owner's manual for the EWheels EW-11 Sport electric scooter, covering intended use, safety warnings, operation, maintenance, specifications, and warranty information.</p>
 <p>EW-72 Owner's Manual</p>	<p>EWheels EW-72 Mobility Scooter Owner's Manual</p> <p>Comprehensive owner's manual for the EWheels EW-72 mobility scooter, covering intended use, safety warnings, operation, maintenance, technical specifications, and warranty information.</p>

