

Orcam MyEye 3 Pro

MyEye 3 Pro - AI Vision Assistant User Manual

Brand: Orcam | Model: MyEye 3 Pro

INTRODUCTION

The Orcam MyEye 3 Pro is a revolutionary wearable device designed to empower individuals with visual impairments, enhancing their independence and confidence in navigating the world. This advanced AI Vision Assistant offers a range of features to help users interact with their environment, read text, and identify objects and faces.

This manual provides essential information on setting up, operating, maintaining, and troubleshooting your MyEye 3 Pro device. Please read it thoroughly to ensure optimal performance and a seamless user experience.

WHAT'S IN THE BOX

Your MyEye 3 Pro package includes the following components:

- OrCam MyEye 3 Pro Device
- Stand and Mount
- Speaker
- USB Data and Power Adapter and Cable
- Eyeglasses Frame

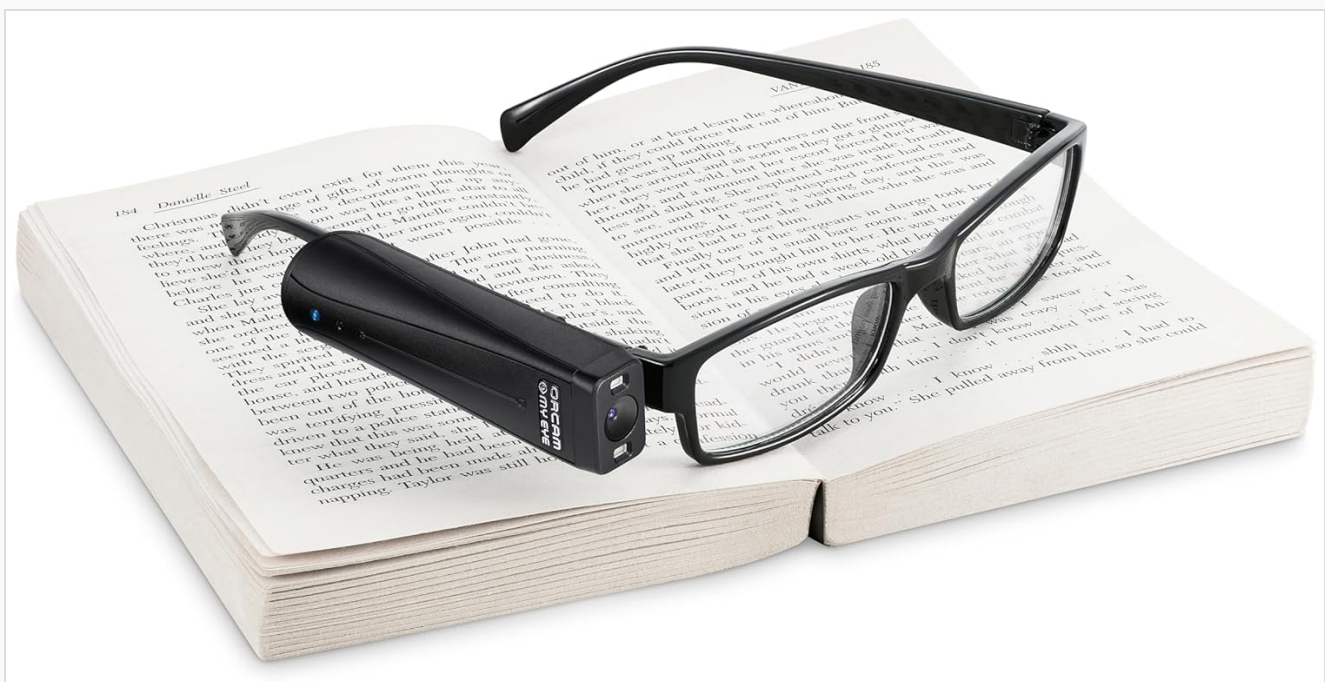


Image: All components included with the Orcam MyEye 3 Pro device.

SETUP

Follow these steps to set up your MyEye 3 Pro for first use:

1. **Charge the Device:** Connect the MyEye 3 Pro to the USB power adapter using the provided cable. Ensure the device is fully charged before initial use. A full charge provides approximately 6 hours of additional use.
2. **Attach to Eyeglasses:** Securely attach the MyEye 3 Pro device to the included eyeglasses frame using the provided mount. Ensure the camera lens is unobstructed and positioned correctly for optimal viewing.
3. **Initial Calibration:** Power on the device. Follow the voice prompts for initial calibration, which may include adjusting the camera angle and performing a brief setup routine to optimize performance for your specific needs.
4. **Pair Speaker (Optional):** If using the external speaker, follow the instructions provided with the speaker to pair it with your MyEye 3 Pro device.



OPERATING YOUR MYEYE 3 PRO

The MyEye 3 Pro offers intuitive controls and advanced features to assist you in various daily tasks.

Reading Text

The device can read text from virtually any surface, including books, newspapers, computer screens, and product labels. Activate reading capabilities using:

- **Voice Commands:** Speak commands such as "Read text" or "Read from here."
- **Hand Gestures:** Use specific hand gestures recognized by the device to initiate reading.
- **Touch Bar:** A physical touch bar on the device allows for quick activation of reading functions.

The scan-to-text capability ensures that the MyEye 3 Pro can read from any point of the text, providing a seamless reading experience.

Interactive AI Assistant (Just Ask)

Utilize the "Just Ask" feature to interact naturally with the AI assistant. You can ask questions related to the text being read or other general topics, receiving instant responses.

Object and Face Recognition

The MyEye 3 Pro can effortlessly recognize faces of loved ones, identify various products, discern money notes, and even recognize colors. This feature provides valuable information about your surroundings, guiding you through different spaces.

Your browser does not support the video tag.

Video: Demonstrates the "What's in front of me" feature of the Orcam MyEye 3 Pro, showcasing its ability to identify objects and provide information about surroundings for navigation.

Advanced Magnification and Translation

Zoom in on text and images for clearer viewing. The device also supports translation features in over 140 languages, expanding accessibility.



Image: A detailed view of the MyEye 3 Pro device.

MAINTENANCE AND CARE

To ensure the longevity and optimal performance of your MyEye 3 Pro, follow these maintenance guidelines:

- **Cleaning:** Gently wipe the device and camera lens with a soft, dry, lint-free cloth. Do not use liquid cleaners or abrasive materials.
- **Storage:** When not in use, store the MyEye 3 Pro in its protective case in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** Avoid fully discharging the battery frequently. For best battery life, charge the device regularly and avoid leaving it uncharged for extended periods.
- **Handle with Care:** The MyEye 3 Pro is a precision electronic device. Avoid dropping it or subjecting it to strong impacts.



Image: MyEye 3 Pro with its protective case, ideal for storage.

TROUBLESHOOTING

If you encounter issues with your MyEye 3 Pro, try the following common troubleshooting steps:

- **Device Not Responding:** Ensure the device is charged. Try restarting the device by holding down the power button.
- **Reading Inaccuracies:** Check that the camera lens is clean and unobstructed. Ensure adequate lighting in your environment. Re-calibrate the device if necessary.
- **Audio Issues:** Verify that the speaker is properly connected or paired. Check the volume settings on the device.

- **Charging Problems:** Ensure the USB cable and power adapter are securely connected and functioning correctly. Try a different power outlet.

If the issue persists after attempting these steps, please refer to the support section for further assistance.

SPECIFICATIONS

Feature	Detail
ASIN	B0D2P6RF84
Manufacturer	Orcam
Date First Available	April 25, 2024
Battery Life	Approximately 6 hours additional use (with enclosed battery)

WARRANTY AND SUPPORT

For detailed warranty information and to register your MyEye 3 Pro device, please visit the official Orcam website or contact their customer support. Keep your proof of purchase for warranty claims.

If you require technical assistance or have questions not covered in this manual, please reach out to Orcam customer support through their official channels. Information regarding optional protection plans may also be available from your retailer.

© 2024 Orcam. All rights reserved.

Documents - Orcam – MyEye 3 Pro



[\[pdf\]](#) Flyer Borchure

Brochures Flyers Mediakit B2B OrCam MyEye 3 Pro Brochure drive google file d

1RLNG9JNT5ApiBqYulpKQc9FLe2fW7WoE view usp link |||

OrCam MyEye gives me the freedom to read my own letters, documents, bank statements and digital tex ... ray OrCam MyEye user A TIME Magazine 2019 Best Invention CES 2022 Innovation Awards Winner OrCam MyEye 3 Pro With the Interactive AI Assistant Step into a revolutionized world of vision assistance ...

lang:en score:41 filesize: 2.22 M page_count: 2 document date: 2024-02-13




[pdf] User Manual

docId C24571F22ED0B5EE385D4A68D8BC8F61 Meny Gantz OrCam MyEye 3 Pro Smart Magnifier Kurzanleitung
reineckerversion de fileadmin user upload 10 Downloads gebrauchsanweisungen elektronische Brillen und mobile
Assistenzsysteme |||

OrCam MyEye 3 Pro Intelligente Lupe Schnellstart-Anleitung Einrichten 1. Schalten Sie die
OrCam MyEye 3 Pro ein. 2. Verbinden Sie die OrCam MyEye 3 Pro ber das mitgelieferte
Micro-USB-Kabel mit seinem magnetischen Anschluss mit Ihrem Computer. 3. Stellen Sie die
OrCam MyEye 3 Pro auf den mitgeliefer...

lang:de score:39 filesize: 83.69 K page_count: 1 document date: 2025-02-20

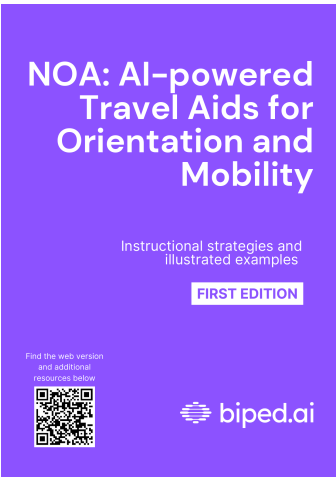


[pdf] Frequently Asked Questions Guide

O M Book DAGcKC1g1Ww BAGWd0ILIVo 0 Claire Pittet biped ai The other most common assistive tool is the use of
a guide dog Guide dogs are carefully trained to navigate obstacles obey directional commands alert for curbs and
steps find crossings doors They also capable making independent decisions ensure their handler s safety including
sometimes disobeying unsafe provide companionship om en o m |||

NOA: AI-powered Travel Aids for Orientation and Mobility Instructional strategies and
illustrated ex ... PT overview. Retrieved November 13, 2024, from
https://openai.com/chatgpt/overview/ 11. OrCam MyEye 3 Pro - Revolutionize Your Vision with
Cutting-Edge AI Technology. n.d. . OrCam Technologie...

lang:en score:23 filesize: 1.93 M page_count: 62 document date: 2025-06-04

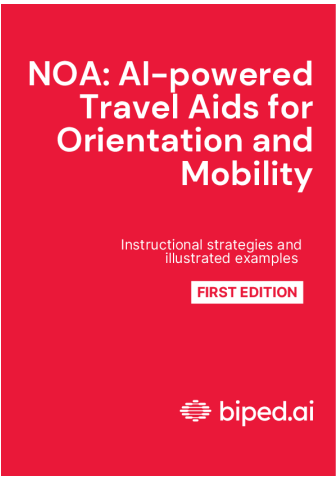

[pdf] Frequently Asked Questions Guide

Copy of LITTLE RED BOOK FROM BIPED DAGXMmxY5KI BAGWd0ILIVo Claire Pittet Mobility Orientation and
Travel Aids for NOA AI powered For people with vision impairments being able to navigate the world on their own is
essential quality life m o b ilit y a n d r ie t io id s lik e g u w h it om guide noa biped ai en |||

NOA: AI-powered Travel Aids for Orientation and Mobility Instructional strategies and
illustrated ex ... PT overview. Retrieved November 13, 2024, from
https://openai.com/chatgpt/overview/ 11. OrCam MyEye 3 Pro - Revolutionize Your Vision with
Cutting-Edge AI Technology. n.d. . OrCam Technologie...

lang:en score:22 filesize: 1.22 M page_count: 46 document date: 2024-11-26

Preliminary Report - Group 1

AudiblEyes
9/3/2024
Anirudh Manika, Austin Ye, Andy Zhu
Page Count: 10

[pdf] Specifications

PrelimReportGroup01 Preliminary Report Group 1 2 bpb us w2 wpmucdn sites wustl edu dist 0 4614 files 2024 10 |||

Preliminary Report - Group 1 AudiblEyes 9/3/2024 Anirudh Manika, Austin Ye, Andy Zhu Page Count: 10 Table of Contents Manika, Ye, Zhu

Background.....1
Need Statement.....

lang:en score:21 filesize: 461.29 K page_count: 13 document date: 0000-00-00



[pdf]

Microsoft PowerPoint Complimentary Services Michael Papp 672542fdea6d26cd181272a1 avrt cdn prod website files 66eb1d4b68d878cfce9a1a75 |||

Rehabilitation Engineering and Assistive Technology Complimentary Services Michael Papp, M.Sc.Eng., ... s For Wheelchairs Turn any wheelchair into a Smart Wheelchair
https://brazemobility.com/ OrCam MyEye 3 Pro Wearable Device Reads text Recognizes faces
Identified products Smart magnifier fea...

lang:en score:19 filesize: 8.91 M page_count: 40 document date: 2024-11-01



[pdf] Guide

8e597baee59ab343aa4d18ec3eb22b63 files ecommercecdns uk 228344 |||

Restoring Independence Since 1996 Low Vision Solutions 95 Freephone: 0800 002 5555
Prices are correc ... y products to automatically recognise e.g. soup Colour detection Beta
orientation features OrCam MyEye 3 Pro - 3,700 - New All the features of the OrCam MyEye 2
Pro with the addition of: AI OCR fo...

lang:en score:14 filesize: 7.97 M page_count: 21 document date: 2024-07-05



[pdf] Guide

d49f5e45894e526069d1dd34aae08de9 files ecommercecdns uk 228344 |||

Restoring Independence Since 1996 Low Vision Solutions 98 Freephone: 0800 002 5555
Prices are correc ... y products to automatically recognise e.g. soup Colour detection Beta
orientation features OrCam MyEye 3 Pro - 3,700 Free Bluetooth earphones RRP 40 All the
features of the OrCam MyEye 2 Pro with t...

lang:en score:14 filesize: 7.97 M page_count: 21 document date: 2025-01-02



VisionAid
Restoring Independence Since 1996

Low Vision Solutions • 100

Freephone:
0800 002 5555

Prices are correct at time of printing but are subject to change without notice. E&OE




[\[pdf\]](#) Guide

76a748a8db53ee4cdee18a11bbef5b58 files ecommerce dns uk 228344 |||

Restoring Independence Since 1996 Low Vision Solutions 100 Freephone: 0800 002 5555

Prices are corre ... y products to automatically recognise e.g. soup Colour detection Beta orientation features OrCam MyEye 3 Pro - 3,700 Free Bluetooth earphones Stand RRP 90

Magnifies any programme s subtitles up ...

lang:en **score:14** filesize: 8.05 M page_count: 21 document date: 2025-02-24