

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Dornvap](#) /

› [The Nebulizer - Nebulizer for Adults Kids, Portable Nebulizer Machine, Handheld Mesh Nebulizer for Home use Travel](#)

## Dornvap B0D3F6ZWCB

# Dornvap Portable Mesh Nebulizer Instruction Manual

Model: B0D3F6ZWCB

## 1. INTRODUCTION AND OVERVIEW

The Dornvap Portable Mesh Nebulizer is designed to deliver a fine mist of medication for effective respiratory symptom alleviation. Its compact, lightweight design and rechargeable battery make it ideal for home, travel, and on-the-go use. This manual provides essential information for the safe and efficient operation and maintenance of your nebulizer.

### Key Features:

- **Advanced Nebulizing Technology:** Delivers fine mist ( $\leq 4\mu\text{m}$  particles) for deep lung penetration.
- **Versatile Use for All Ages:** Features dual-speed functionality with Powerful Mode for adults and Mild Mode for children.
- **Silent Operation:** Operates quietly ( $\leq 30\text{dB}$ ) for discreet use in any environment.
- **Effortless Cleaning:** Equipped with an auto-clean function to prevent clogging.
- **Ultimate Portability:** Compact, lightweight, and rechargeable for convenient use anywhere.

## 2. PRODUCT COMPONENTS

Your Dornvap Portable Mesh Nebulizer package includes the following components:

- Main Nebulizer Unit
- Adult Mask
- Children's Mask
- Mouthpiece
- USB Type-C Charging Cable
- Storage Pouch



Figure 2.1: Main nebulizer unit with included masks and mouthpiece.



Figure 2.2: Various masks and mouthpiece for different users.

### 3. SETUP

1. **Charge the Device:** Ensure the nebulizer is fully charged before first use. Connect the USB Type-C cable to the charging port at the bottom of the device and plug it into a compatible USB power source. The light indicator will show charging status.

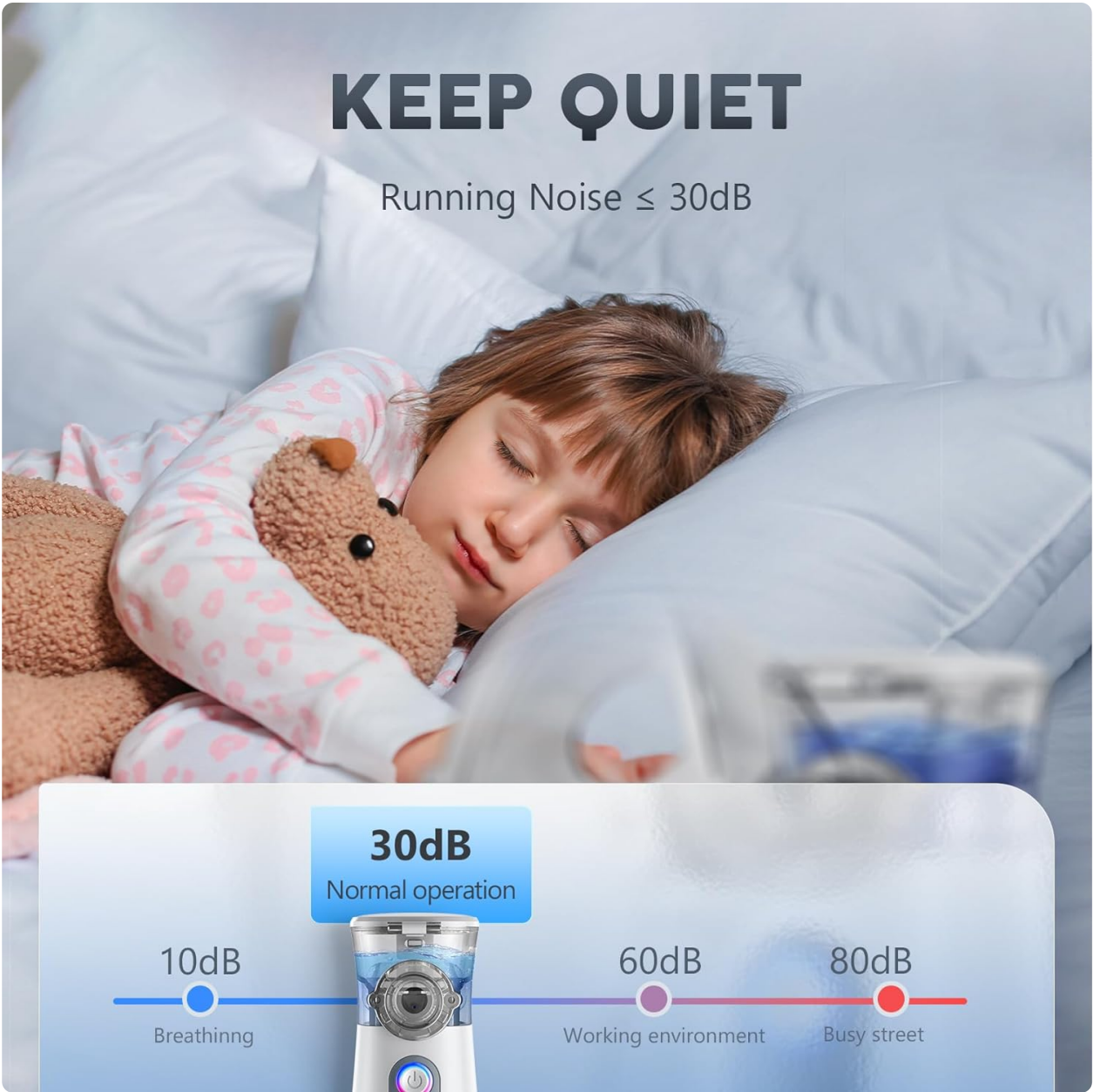



Figure 3.1: Charging port location.

2. **Prepare Medication Chamber:** Open the lid of the medication cup. Ensure the chamber is clean and dry.

**DORNVAP More for You**

 **60+** **COUNTRIES**  
Available in  
60+ Countries

 **120+** **MILLION**  
Products Sold  
Worldwide


 **120+** **MILLION**  
Loyal  
Customers

Figure 3.2: Opening the medication cup lid.

- Add Medication:** Carefully pour the prescribed liquid medication (e.g., normal saline, prescription atomized drugs) into the medication cup. Do not use distilled water, purified water, or oily liquids.



Figure 3.3: Adding medication to the nebulizer.

- Attach Mask/Mouthpiece:** Securely attach the appropriate mask (adult or children's) or mouthpiece to the nebulizer unit.

# Atomised particles

the atomized particles is less than 2um

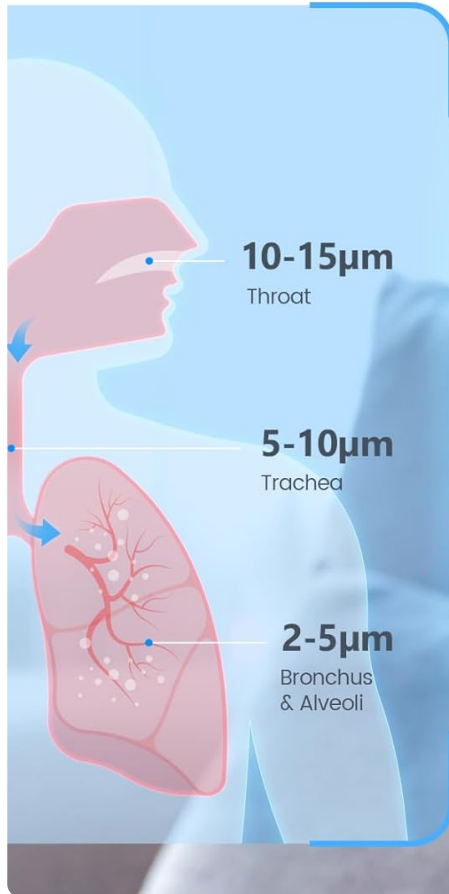


Figure 3.4: Attaching the mask to the nebulizer.

## 4. OPERATING INSTRUCTIONS

1. **Power On:** Press the power button once to turn on the nebulizer. The device will start nebulizing.



Figure 4.1: Powering on the nebulizer.

2. **Select Mode:** The nebulizer offers dual-speed functionality.
- Press once for Powerful Mode (recommended for adults).
  - Press again for Mild Mode (recommended for children).

# TWO WORKING MODES

Long press for 5 seconds , Auto Cleaning Mode On



Figure 4.2: Nebulizer operating modes.

- 3. Begin Treatment:** Place the mask over the nose and mouth, or insert the mouthpiece into the mouth. Breathe normally and deeply to inhale the mist. The silent operation ensures a tranquil treatment experience.

# Relieve Several Diseases

Fever/Bronchitis/Rhinitis/  
Nasal mucus/Asthma



Figure 4.3: Quiet operation for undisturbed use.

4. **Power Off:** Once the medication is fully nebulized or the treatment is complete, press and hold the power button for 5 seconds to turn off the device.

*Note: No official product videos with 'creator\_type: Seller' were found for embedding based on the provided data.*

## 5. MAINTENANCE AND CLEANING

Proper cleaning is crucial for maintaining the nebulizer's performance and hygiene.

1. **Disassemble:** After each use, detach the medication cup from the main unit. Remove the mask or mouthpiece.
2. **Rinse Components:** Rinse the medication cup, mask, and mouthpiece with clean water.
3. **Auto-Clean Function:** For thorough cleaning and to prevent clogging, activate the auto-clean function. Fill the medication chamber with a small amount of vinegar solution (e.g., distilled white vinegar diluted with water), then rinse and dry. Long press the power button for 5 seconds to activate the self-cleaning mode.

4. **Dry Thoroughly:** Allow all components to air dry completely before reassembling or storing.

## 6. TROUBLESHOOTING

If you encounter issues with your nebulizer, please refer to the following common troubleshooting steps:

- **No Mist Production:**
  - Ensure the device is charged.
  - Check if the medication chamber is properly filled and securely attached.
  - Perform the auto-clean function to clear any potential clogs in the mesh.
- **Weak Mist:**
  - Ensure the device is fully charged.
  - Verify that the correct mode (Powerful or Mild) is selected for your needs.
  - Clean the mesh thoroughly using the auto-clean function.
- **Device Not Turning On:**
  - Ensure the battery is charged.
  - Press and hold the power button for a few seconds to ensure proper activation.

If problems persist after attempting these steps, please contact customer support for further assistance.

## 7. SPECIFICATIONS

Feature	Detail
Product Dimensions	129.92 x 1.3 x 3.7 inches
Weight	8.47 ounces
Particle Size	≤4µm
Noise Level	≤30dB
Power Source	Rechargeable Battery (USB Type-C)
First Available Date	May 6, 2024

## 8. WARRANTY AND SUPPORT

For warranty information, technical support, or any inquiries regarding your Dornvap Portable Mesh Nebulizer, please refer to the contact information provided on the product packaging or visit the official Dornvap website. Each item is verified by Transparency with a unique code.

Visit the official Dornvap Store for more information:[Dornvap Store](#)