



Manuals.plus /

› KAUGIC /

› KAUGIC Smart Ear Wax Removal Tool Camera User Manual

KAUGIC Smart Ear Wax Removal Tool Camera

KAUGIC Smart Ear Wax Removal Tool Camera User Manual

Your comprehensive guide to safe and effective ear cleaning with visual assistance.

PRODUCT OVERVIEW

The KAUGIC Smart Ear Wax Removal Tool Camera is an innovative device designed for thorough and safe ear cleaning. It combines a high-definition camera with a water irrigation system, allowing for precise and effective removal of earwax while providing a clear visual of the ear canal.



Image: The KAUGIC Smart Ear Wax Removal Tool Camera, showcasing its main unit, water tank, various nozzles, USB charging cable, and the accompanying smartphone app interface.

Key Features

- **2-in-1 Cleaning & Visualization:** Built-in 1080P HD camera for real-time ear canal viewing via WiFi to your phone.
- **Four Adjustable Water Pressure Modes:** Soft, Normal, Pulse, and Strong (7-24.5 PSI) for customized cleaning.
- **Two Smart Cleaning Modes:** Hands-free or Guided Cleaning via the Soulear App.
- **Long-Lasting Power & Large Water Tank:** 2000mAh battery for up to 5 hours of continuous use and a 200ml water tank.
- **Wide Compatibility:** Works with iOS & Android devices via the Soulear App.

WHAT'S INCLUDED

Your KAUGIC Ear Wax Removal Tool Camera package includes all the necessary components for a complete ear cleaning experience:

- 1 x Ear cleaner with integrated 1080P HD camera
- 10 x Replacement nozzles
- 1 x Ear basin for water collection
- 1 x USB charging cable (Type-C)
- 1 x User manual



Image: A visual representation of all items included in the product package, such as the main device, multiple ear tips, charging cable, and ear basin.

SETUP GUIDE

1. App Installation & Wi-Fi Connection

1. Download the 'Soulear' App from Google Play Store or Apple App Store. You can scan the QR code below for direct access.



Image: Step-by-step guide for downloading the Soulear app via QR code and connecting the device's Wi-Fi.

2. Long press the power switch key on the device to turn it on.
3. On your smartphone, connect to the device's Wi-Fi network (usually named 'Soulear-xxxx').
4. Return to the Soulear App and click 'Connected' to establish the visual link.

2. Device Assembly & Water Filling


5. Unscrew the water tank from the main unit.
6. Fill the 200ml water tank with warm water or a recommended ear cleaning solution. Avoid using hot or cold water to prevent discomfort or dizziness.
7. Securely screw the water tank back onto the main unit.

1080P HD VIDEO IMAGE

Ear cleaner with camera, you
can clearly see the clear effect



Wifi connection

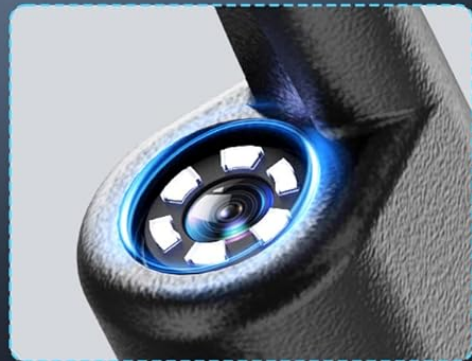

6LED light

Image: The ear cleaner disassembled, showing the main unit, water tank, and charging port, illustrating how to fill the tank.

Your browser does not support the video tag.

Video: An official product video demonstrating the setup and initial use of the Water Powered Ear Cleaner with Camera.

OPERATING INSTRUCTIONS

1. Selecting Pressure Modes

The device offers four adjustable water pressure modes to suit different needs and sensitivities. Press the 'M' button to cycle through the modes:

- **Soft:** Low water pressure, ideal for sensitive ears or first-time users.
- **Normal:** Standard pressure for most users.
- **Pulse:** Alternating water pressure for effective flushing.
- **Strong:** Powerful pressure for stubborn earwax.



Image: A detailed view of the device's control panel, highlighting the four pressure settings and their corresponding icons.

2. Choosing Cleaning Methods

You can choose between two smart cleaning modes:

- **Visual Cleaning:** Use the Soulear App for real-time viewing of your ear canal while cleaning. This method provides enhanced safety and precision.
- **Hands-Free Cleaning:** Clean directly without the app when you are short on time or prefer not to use the visual function.

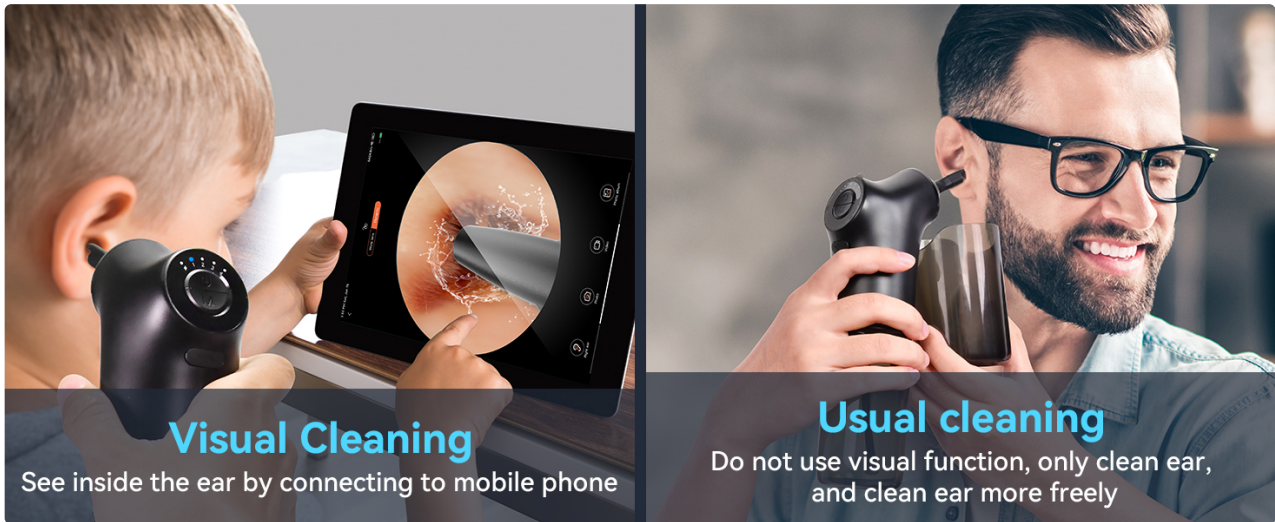


Image: Illustrates the two cleaning methods: Visual Cleaning using the app for real-time view, and Usual Cleaning for quick use without the app.

Your browser does not support the video tag.

Video: An official product video demonstrating the operation and different modes of the Water Powered Ear Cleaner with Camera.

MAINTENANCE

Cleaning the Device

After each use, it is important to clean the device thoroughly to maintain hygiene and prolong its lifespan. Disassemble the water tank and nozzles, rinse them with clean water, and allow them to air dry completely before reassembling or storing.

Replacing Nozzles

The kit includes 10 replacement nozzles to ensure hygienic use for multiple individuals or for regular replacement. Simply pull off the used nozzle and slide on a new one when needed.

TROUBLESHOOTING

If you encounter any issues with your KAUGIC Smart Ear Wax Removal Tool Camera, please refer to the following common solutions:

- **Device Not Turning On:** Ensure the device is fully charged using the provided USB cable.
- **Wi-Fi Connection Issues:** Make sure the device is turned on and your phone's Wi-Fi is connected to the

'Soulear-xxxx' network. Try restarting both the device and your phone.

- **No Water Spray:** Check if the water tank is filled correctly and securely attached. Ensure the nozzle is properly installed.
- **Weak Water Pressure:** Verify that the selected pressure mode is appropriate for your needs. If using a lower setting, try increasing it. Ensure the water tank is not empty.
- **Blurred Camera View:** Ensure the camera lens is clean and free from debris or water droplets.

If problems persist, please contact KAUGIC customer support for further assistance.

PRODUCT SPECIFICATIONS

Specification	Detail
Brand	KAUGIC
Model	Smart Ear Wax Removal Tool Camera
Camera Resolution	1080P HD
Water Pressure Modes	Soft, Normal, Pulse, Strong (7-24.5 PSI)
Battery Capacity	2000mAh
Continuous Use Time	Up to 5 hours
Water Tank Capacity	200ml
Compatibility	iOS & Android
Package Dimensions	10.71 x 6.5 x 3.35 inches
Weight	15.52 ounces
Manufacturer	KAUGIC
Country of Origin	China

WARRANTY & SUPPORT

For any product-related inquiries, technical support, or warranty claims, please contact KAUGIC customer service. Refer to the product packaging or the official KAUGIC website for the most up-to-date contact information.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended for personal hygiene and should be used as directed in this manual.

