

Smart Ring R06

R06 Smart Ring User Manual

Your Guide to Health and Wellness Monitoring

1. INTRODUCTION TO THE R06 SMART RING

The R06 Smart Ring is an advanced wearable device designed to provide comprehensive health and fitness monitoring. It offers continuous tracking of vital metrics such as heart rate, sleep patterns, and activity levels, empowering users with insights into their well-being. Its compact, waterproof design ensures comfort and durability for daily use.

This manual provides detailed instructions on how to set up, operate, maintain, and troubleshoot your R06 Smart Ring to ensure optimal performance and a seamless user experience.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the R06 Smart Ring to prevent injury or damage to the device.

- Medical Disclaimer:** The R06 Smart Ring is not a medical device and is not intended to diagnose, treat, cure, or prevent any disease. Consult a healthcare professional for any health concerns.
- Choking Hazard:** Keep the device and its accessories out of reach of small children.
- Battery Safety:** Do not attempt to replace the battery. Do not expose the device to extreme temperatures, direct sunlight, or fire.
- Water Resistance:** The R06 Smart Ring is IP68 waterproof. While suitable for swimming and showering, avoid prolonged submersion in hot water or exposure to high-pressure water jets.
- Skin Irritation:** If you experience skin irritation, remove the ring and consult a doctor. Ensure the ring is clean and dry before wearing.

3. SETUP GUIDE

3.1 Package Contents

Verify that your package contains the following items:

- R06 Smart Ring
- Charging Cable (USB-A to proprietary connector)
- User Manual (this document)

3.2 Charging the Ring

Before first use, fully charge your R06 Smart Ring.

- Connect the charging cable to a USB power source (e.g., computer USB port, wall adapter 5V/1A).
- Attach the proprietary connector of the charging cable to the charging contacts on the R06 Smart Ring. Ensure a secure connection.
- The ring's indicator light (if present, or app will show status) will indicate charging. A full charge typically takes approximately 1-2 hours.
- Once fully charged, disconnect the ring from the cable.

Note: No specific product images or videos were provided for embedding. Please refer to the physical product for visual guidance on charging.

3.3 App Download and Pairing

The R06 Smart Ring requires a companion mobile application for full functionality and data synchronization.

- Download the App:** Search for the official R06 Smart Ring app on the Apple App Store (for iOS) or Google Play Store (for Android). Alternatively, scan the QR code provided in the packaging (if available).
- Create Account:** Open the app and follow the on-screen instructions to create a new user account or log in if you already have one.
- Enable Bluetooth:** Ensure Bluetooth is enabled on your smartphone.
- Pair the Ring:** In the app, navigate to the device pairing section. Select "Add New Device" or similar. The app will search for nearby R06 Smart Rings. Select your ring from the list to initiate pairing.
- Confirm Pairing:** Follow any prompts on the app to confirm the pairing. Once successful, the app will display a "Connected" status.

3.4 Wearing the Ring Correctly

For accurate data collection, proper placement of the R06 Smart Ring is crucial.

- Wear the ring on your index, middle, or ring finger.
- Ensure the sensors on the inner side of the ring are in firm contact with your skin.
- The ring should fit snugly but comfortably, not too tight to restrict circulation, nor too loose to move freely.

4. OPERATING THE R06 SMART RING

4.1 Automatic Monitoring

The R06 Smart Ring is designed for continuous, automatic monitoring once paired with your device.

- **Heart Rate:** The ring continuously monitors your heart rate throughout the day and night.
- **Sleep Tracking:** It automatically detects sleep onset and wake times, and analyzes sleep stages (light, deep, REM) and sleep quality.
- **Activity Tracking:** The ring tracks your daily steps, distance, and calories burned.

4.2 Data Synchronization

Your R06 Smart Ring automatically synchronizes data with the companion app when within Bluetooth range of your paired smartphone.

- Ensure the app is running in the background or open it periodically to allow data transfer.
- Data is stored on the ring for a limited period (e.g., 7 days). Regular synchronization is recommended to prevent data loss.

4.3 Understanding App Metrics

The companion app provides detailed insights into your health data.

- **Dashboard:** Overview of daily metrics (e.g., current heart rate, sleep score, activity progress).
- **Heart Rate:** View historical heart rate data, including resting heart rate, average, and peak rates.
- **Sleep:** Detailed sleep reports, including sleep stages, duration, efficiency, and disturbances.
- **Activity:** Track daily, weekly, and monthly activity trends. Set goals and monitor progress.
- **Trends:** Analyze long-term trends in your health data to identify patterns and make informed lifestyle adjustments.

5. MAINTENANCE

5.1 Cleaning the Ring

Regular cleaning helps maintain hygiene and sensor accuracy.

- Wipe the ring, especially the inner sensor area, with a soft, damp cloth.
- For stubborn dirt, use a small amount of mild soap and rinse thoroughly with clean water.
- Ensure the ring is completely dry before wearing or charging.
- Do not use abrasive cleaners, solvents, or harsh chemicals.

5.2 Battery Care

To maximize battery lifespan:

- Avoid fully discharging the battery frequently.
- Charge the ring regularly, ideally before the battery level drops below 20%.
- If storing the ring for an extended period, charge it to about 50% and store in a cool, dry place. Recharge every few months.

5.3 Storage

Store your R06 Smart Ring in a clean, dry place away from direct sunlight, extreme temperatures, and corrosive materials.

6. TROUBLESHOOTING

If you encounter issues with your R06 Smart Ring, try the following solutions:

Problem	Possible Solution
Ring not connecting to app	Ensure Bluetooth is enabled on your phone. Make sure the ring is charged. Restart your phone and the app. Try unpairing and re-pairing the device in the app's settings.
Inaccurate readings (HR, Sleep)	Ensure the ring is worn correctly and snugly on your finger. Clean the sensors on the inner side of the ring. Ensure your skin is clean and dry where the sensors make contact.

Problem	Possible Solution
Battery drains quickly	Ensure the ring is fully charged. Avoid extreme temperatures. If the issue persists, contact customer support.
Ring unresponsive	Place the ring on its charger for a few minutes. If there's a reset button (refer to product specifics if available), perform a hard reset. Otherwise, let the battery fully drain and then recharge.

7. SPECIFICATIONS

- **Model:** R06 Smart Ring
- **Functions:** Heart Rate Monitoring, Sleep Tracking, Fitness Tracking (Steps, Distance, Calories)
- **Water Resistance:** IP68 (Suitable for swimming, showering; not for diving or hot water)
- **Battery Life:** Up to 5 days (typical usage)
- **Charging Time:** Approximately 1-2 hours
- **Connectivity:** Bluetooth Low Energy (BLE)
- **Compatibility:** iOS (specific version), Android (specific version) - *Check app store for current requirements.*
- **Material:** (Specific material not provided, typically medical-grade plastic/metal alloy)

8. WARRANTY AND SUPPORT

The R06 Smart Ring comes with a limited warranty against defects in materials and workmanship. The specific warranty period and terms may vary by region and retailer. Please refer to your purchase documentation or the manufacturer's official website for detailed warranty information.

For technical support, troubleshooting assistance, or warranty claims, please contact the customer support team via the contact information provided on the product packaging or the official website. When contacting support, please have your product model (R06) and proof of purchase ready.

Online Resources: Visit the official product website for FAQs, updated manuals, and software downloads.

Related Documents - R06

	Guida Utente Braccialeto Sanitario Intelligente H Band Manuale completo per il braccialeto sanitario intelligente H Band, che copre installazione dell'app, connessione, funzioni, impostazioni, parametri tecnici, considerazioni e domande frequenti.
	Healthy Smart Ring MR6 User Manual and Specifications User manual for the Healthy Smart Ring MR6, detailing product introduction, wearing instructions, charging, app pairing, functions, specifications, frequently asked questions, disclaimers, and FCC compliance information.
	NMJ15 Smart Ring: Comprehensive Health and Fitness Tracking Discover the NMJ15 smart ring, a sophisticated wearable device offering dynamic heart rate monitoring, blood oxygen saturation, pedometer, calorie tracking, multiple exercise modes, sleep analysis, and body temperature measurement. Compatible with Android 8.0+ and iOS 12.0+, it features a long-lasting battery, IP68 waterproof rating, and a stylish design available in silver and rose gold.
	JCRing X6B Smart Health Ring User Manual Comprehensive user manual for the JCRing X6B smart health ring, detailing setup, usage, specifications, troubleshooting, and safety guidelines for monitoring health metrics like heart rate and sleep.
	Kogan Aura Smart Ring User Guide: Setup, Features, and Specifications Comprehensive user guide for the Kogan Aura Smart Ring, detailing setup, operation, app connectivity, health tracking features (activity, sleep, heart rate), specifications, and safety information.



[Smart Health Ring User Manual and Specifications](#)

Comprehensive user guide and technical specifications for the LuckRing Smart Health Ring, detailing its features, functions, usage instructions, and health monitoring capabilities.



[pdf] User Manual Guide

COLMi R02 Smart Ring Manual do usuário anel inteligente R02Manual smartwatch COLMi Smartwatch RingCOLMi User Manualcdn shopify s files 1 0719 4263 5f R03User Mode2024 with Charging Case for Waterproof2024 WaterproofCOLMi ManualSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e ManualSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e00ccfdSee2e04200c12488d96e1984af9 ManualSbfc976df0ec14bd0a7afae04df85722gSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e00ccfdSac0f2ba8ce8d4b89854466a899 ManualSbfc976df0ec14bd0a7afae04df85722gSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e00ccfdS6945481b56544dd1aec0bd198 ManualSbfc976df0ec14bd0a7afae04df85722gSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e00ccfdS6945481b56544dd1aec0bd198 ManualSbfc976df0ec14bd0a7afae04df85722gSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e00ccfdS6945481b56544dd1aec0bd198 ManualSbfc976df0ec14bd0a7afae04df85722gSbfc976df0ec14bd0a7afae04df85722gSee2e04200c12488d96e1984af9e00ccfdS6945481b56544dd1aec0bd198 Overview And Packaging Ring, charging cable, and quick guide. Sensor Blood oxygen red light Charging port Heart rate green light Power On Connect the chargin guide. Sensor Blood oxygen red light Charging port Heart rate green light Power On Connect the charging cable to the ring to booting up the device. Connection ||| rate green light Power On Connect the charging cable to the ring to booting up the device. Connection ||| Quick Guide of Smart Ring Overview And Packaging Ring the ring to booting up the device. Connection ||| Quick Guide of Smart Ring Overview And Packaging Ring, charging cable, and quick guide. Sensor Blood oxygen red light Charging lang:en score:26 filesize: 8.41 M page_count: 1 document date: 2024-02-01


[pdf] User Manual Guide

241030 User manual AMZSA R04AMZSA R04 Smart Ring Men Women With Charging Case Health and Sleep Monitor 5ATM Waterproof Multi spor
Android iOSNew PhoneSe6fbbd8ba7c242cbb41f23f1f3eec1b2Wae01 Se6fbbd8ba7c242cbb41f23f1f3eec1b2Wae01 MISIRUN R05User R09R05 NFC Function To
PhoneSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3eec1b2Wae01 Sa12c4d0c87a14261964
IOSSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3eec1b2WS113dfb62dc9242e2989862b14fc
IOSAMZMAX R02 WaterproofS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbb
R05S82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3e
R05S82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3e
R05S82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3e
R05S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a142619644874171
R05S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a142619644874171
R05S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a142619644874171
R05S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a142619644874171
R08S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a142619644874171
Sba4a176dd28e4be981210494c9582d78Sae01 S5b6fda888a0b490db0daa1d38757d121R R08User R08R08
R08S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3e
Finger R05AMZSA R08S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a142619644874171
Scadb0016b3574eedbc84565029c4a17IQae01 R05New Screen R08R05
R08S749e5376d705406cbde0430885096771QS82c2c93f3caf438280788130b23a7834JSf5c72c231cf84229a2b41be247dee6b53Sa12c4d0c87a1426196448741713c0e5cDSe6fbbd8ba7c242cbb41f23f1f3e
Sd98d498beef04c1eb5b1773e3f689e75pae01 ||| ||| Quick Guide of Smart Ring Overview And Packaging Ring, charging cable, charging case and quick guide. Ser
Quick Guide of Smart Ring Overview And Packaging Ring, charging cable, charging case and quick guide. Sensor Blood Oxygen Red Light Charging Port Heart R
Packaging Ring, charging cable, charging case and quick guide. Sensor Blood Oxygen Red Light Charging Port Heart Rate Green Light Power On Remove the ins
quick guide. Sensor Blood Oxygen Red Light Charging Port Heart Rate Green Light Power On Remove the insulation film first, and put the ring into the ch ||| Quick
Charging Port Heart Rate Green Light Power On Remove the insulation film first, and put the ring into the ch ||| Quick Guide of Smart Ring Overview And Packagin
Remove the insulation film first, and put the ring into the ch ||| Quick Guide of Smart Ring Overview And Packaging Ring, charging cable, charging case and quick g
the ch ||| Quick Guide of Smart Ring Overview And Packaging Ring, charging cable, charging case and quick guide. Sensor Blood Oxygen Red Light Charging Port
Packaging Ring, charging cable, charging case and quick guide. Sensor Blood Oxygen Red Light Charging Port Heart Rate Green Light Power On Remove the ins
quick guide. Sensor Blood Oxygen Red Light Charging Port Heart Rate Green Light Power On Remove the insulation film first, and put the ring into the ch
Quick Guide of Smart Ring Overview And Packaging Ring, charging cable, charging case and quick guide. Sensor Blood Oxygen Red
lang:it score:25 filesize: 2 M page_count: 1 document date: 2024-10-31




[pdf] Specifications Test Report

Measurement Report Kent Test Report2 ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
Shenzhen CTL Testing Technology Co., Ltd. Tel: 86-755-89486194 E-mail: ctl ctl-lab.com TEST REPORT ... r Product Name
lang:en score:21 filesize: 4.21 M page_count: 38 document date: 2024-09-23

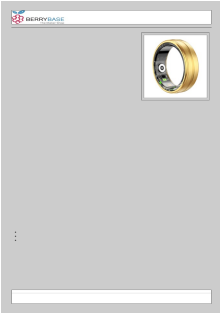


[pdf] User Manual Datasheet

Datasheet for Colmi R06 Smart Ring überwacht Herzfrequenz Schlaf Blutsauerstoff Schrittzähler IP68 gold 10 Der Kalorienzähler Entfernungsmessung Herzfreque
Datasheet for Colmi **R06** Smart Ring, berwacht Herzfrequenz, Schlaf Blutsauerstoff, Schrittzhler, IP68 gold 10 10 Colmi **R06** Smart F
lang:de score:20 filesize: 489.62 K page_count: 4 document date: 2025-05-30

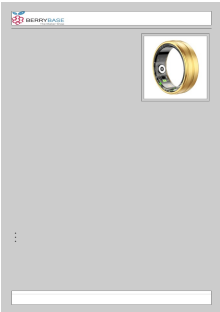


[pdf] User Manual Datasheet

Datasheet for Colmi R06 Smart Ring überwacht Herzfrequenz Schlaf Blutsauerstoff Schrittzähler IP68 gold 10 Die Schlafüberwachung analysiert die Schlafphasen
Datasheet for Colmi **R06** Smart Ring, berwacht Herzfrequenz, Schlaf Blutsauerstoff, Schrittzhler, IP68 gold 10 10 Colmi **R06** Smart F
lang:de score:20 filesize: 3.5 M page_count: 4 document date: 2025-07-12

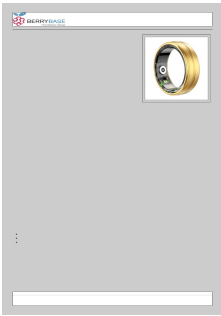


[\[pdf\]](#) User Manual Datasheet

Datasheet for Colmi R06 Smart Ring überwacht Herzfrequenz Schlaf Blutsauerstoff Schrittzähler IP68 gold 11 BerryBase Der Smartring ist die perfekte Lösung für Datasheet for Colmi **R06** Smart Ring, berwacht Herzfrequenz, Schlaf Blutsauerstoff, Schrittzhler, IP68 gold 11 11 Colmi **R06** Smart F
lang:de score:20 filesize: 3.5 M page_count: 4 document date: 2025-07-13

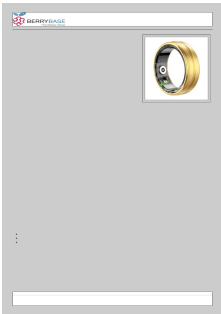

[\[pdf\]](#) User Manual Datasheet

Datasheet for Colmi R06 Smart Ring überwacht Herzfrequenz Schlaf Blutsauerstoff Schrittzähler IP68 gold 10 Mit der QRing App hast du Zugang zu fortschrittliche Datasheet for Colmi **R06** Smart Ring, berwacht Herzfrequenz, Schlaf Blutsauerstoff, Schrittzhler, IP68 gold 10 10 Colmi **R06** Smart F
lang:de score:19 filesize: 3.5 M page_count: 4 document date: 2025-07-20

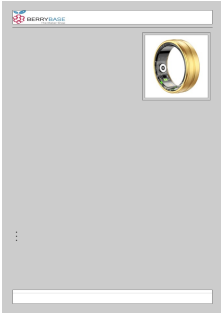

[\[pdf\]](#) User Manual Datasheet

Datasheet for Colmi R06 Smart Ring überwacht Herzfrequenz Schlaf Blutsauerstoff Schrittzähler IP68 gold 10 Der Smartring ist die perfekte Lösung für alle Wert a Datasheet for Colmi **R06** Smart Ring, berwacht Herzfrequenz, Schlaf Blutsauerstoff, Schrittzhler, IP68 gold 10 10 Colmi **R06** Smart F
lang:de score:19 filesize: 3.5 M page_count: 4 document date: 2025-07-13

Shenzhen YaWell intelligent Technology Co., Ltd.
Floor 1-A, Baisha Technology Park, No.3011, Shahe Xi Road, ...
FCC ID: 2AOM3-R06

TABLE OF CONTENTS

1. Introduction	3
1.1 Test Method	3
1.2 Test Item	3
1.3 Test Result	3
1.4 Test Conclusion	3
1.5 Test Date	3
2. Test Method	4
2.1 Test Method	4
2.2 Test Item	4
2.3 Test Result	4
2.4 Test Conclusion	4
2.5 Test Date	4
3. Test Method	5
3.1 Test Method	5
3.2 Test Item	5
3.3 Test Result	5
3.4 Test Conclusion	5
3.5 Test Date	5
4. Test Method	6
4.1 Test Method	6
4.2 Test Item	6
4.3 Test Result	6
4.4 Test Conclusion	6
4.5 Test Date	6

Page 11/20

[\[pdf\]](#) Test Report

TST Test Report1 ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
Shenzhen CTL Testing Technology Co., Ltd. Floor 1-A, Baisha Technology Park, No.3011, Shahe Xi Road, ...
lang:de score:19 filesize: 1.32 M page_count: 19 document date: 2024-09-25




[\[pdf\]](#) Product Photos

Windows External Photos ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
FCC ID: 2AOM3-R06 Page 1 of 4 Report No.: CTL2406283065-WF External Photos of EUT FCC ID: 2AOM3-R06 Page 2 of 4 Repor
lang:de score:18 filesize: 788.1 K page_count: 4 document date: 2024-09-23




[pdf] Label

wang Label Information ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
FCC ID: 2AOM3-**R06** Page 1 of 1 Report No.: CTL2406283065-WF FCC ID Label Information Label location: FCC ID : 2AOM3-**R06** :
lang:de score:18 filesize: 255.74 K page_count: 1 document date: 2024-09-23

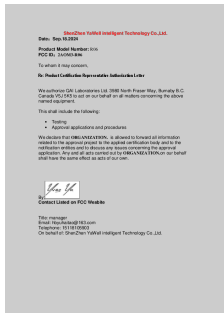

[pdf] Agent Authorization

Agent Authorization Letter Template ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
ShenZhen YaWell intelligent Technology Co.,Ltd. DateSep.18.2024 Product Model Number:**R06** FCC ID2AOM3-**R06** To whom it ma
lang:en score:18 filesize: 75.85 K page_count: 1 document date: 2024-09-23

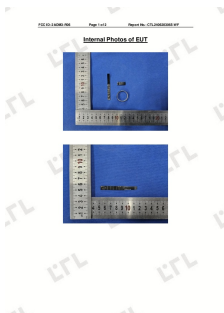


[pdf] Product Photos Teardown

Windows Internal Photos ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
FCC ID: 2AOM3-**R06** Page 1 of 2 Report No.: CTL2406283065-WF Internal Photos of EUT FCC ID: 2AOM3-**R06** Page 2 of 2 Report
lang:en score:18 filesize: 721.59 K page_count: 2 document date: 2024-09-23




[pdf] Product Photos

Windows Test Setup Photos ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
FCC ID: 2AOM3-**R06** Page 1 of 2 Report No.: CTL2406283065-WF Test Setup Photos of the EUT FCC ID: 2AOM3-**R06** Page 2 of 2
lang:en score:18 filesize: 224.85 K page_count: 2 document date: 2024-09-23

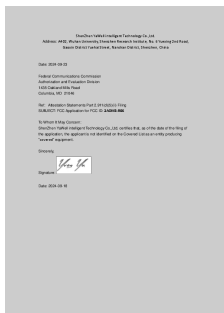

[pdf]

zhouyongbang Attestation Statements2 ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
ShenZhen YaWell intelligent Technology Co.,Ltd. Address: A402, Wuhan University, Shenzhen Research I ... Ref: Attestation Staten
lang:en score:17 filesize: 236.61 K page_count: 1 document date: 2024-09-23

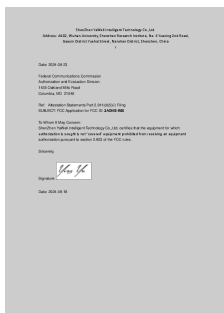

[pdf]

zhouyongbang Attestation Statements1 ShenZhen YaWell intelligent Technology Co Ltd R06 Smart Ring 2AOM3 2AOM3R06 r06
ShenZhen YaWell intelligent Technology Co.,Ltd. Address: A402, Wuhan University, Shenzhen Research I ... 6 Ref: Attestation Stati
lang:en score:17 filesize: 237.65 K page_count: 1 document date: 2024-09-23

[illegible]

[\[pdf\]](#) Decleration of Conformity Test Report

Test Report BLE FCC SHENZHEN ZHISUIXING ELECTRONIC TECHNOLOGY CO LTD R18 Smart Ring 2A3OB 2A3OBR18 r18

Shenzhen Toby Technology Co., Ltd. Report No.: TBR-C-202407-0171-9 Page: 1 of 32 Radio Test Repo ... Under Test EUT EUT No.:
lang:en score:12 filesize: 1.76 M page_count: 32 document date: 2024-08-10

[illegible]

[\[pdf\]](#) Specifications Test Report

admin Test Report SYA0026273 upload sunsky online res certificate |||

Test Report Report No. Applicant Address Manufacturer Address Sample name Sample model Inspection t... in District, Yuehai Stre
lang:en score:11 filesize: 1.36 M page_count: 8 document date: 2024-01-19

[illegible]

[\[pdf\]](#) Decleration of Conformity Test Report

CE Konformitätserklärung CE Fitorbcdn shopify s files 1 0849 8862 5162 Fitorb v 1748687179 |||

EU DECLARATION OF CONFORMITY Product Brand Model Nr. Manufacturer Address : Smart ring : Fitorb : R02, **R06, R06** , R09 :

lang:en score:10 filesize: 88.04 K page_count: 1 document date: 2025-05-31

[illegible]

[\[pdf\]](#) Agent Authorization

TIMCO Engineering Inc admin Agent Authorization Letter Template ShenZhen YaWell intelligent Technology Co Ltd R02 Smart Ring 2AOM3 2AOM3R02 r02

ShenZhen YaWell intelligent Technology Co.,Ltd. A402, Wuhan University, Shenzhen Research Institute, No. 6 Yuexing 2nd Road, (lang:en score:10 filesize: 114.89 K page_count: 1 document date: 2024-01-31