

ScopeAround SA39W/SA39i/M28/MS450 Compatible Ear Specula

ScopeAround Ear Specula Set

PRODUCT INSTRUCTION MANUAL

1. Introduction

Thank you for choosing the ScopeAround Ear Specula Set. This manual provides essential information for the proper and safe use of your new ear specula. These disposable tips are designed to enhance hygiene and safety when using compatible ScopeAround digital otoscopes. Please read this manual thoroughly before use and retain it for future reference.

2. Product Overview

The ScopeAround Ear Specula Set consists of 50 disposable ear specula, specifically designed to fit ScopeAround digital otoscopes with a diameter of 4.3mm and 3.9mm. These include models SA39W, SA39i, M28, and MS450.

Key Features:

- **Enhanced Safety:** Prevents the otoscope camera from entering the ear canal too deeply, protecting the eardrum.
- **Improved Hygiene:** Acts as a barrier to prevent cross-infection between users or during repeated examinations.
- **Optimal Compatibility:** Precisely engineered for a secure fit with specified ScopeAround otoscope models.
- **Convenient Quantity:** A set of 50 pieces ensures a sufficient supply for multiple uses.



Image: The complete set of 50 ScopeAround ear specula.



Image: An ear speculum correctly attached to a ScopeAround digital otoscope.

3. Setup

Before each use, attach a new, clean ear speculum to your compatible ScopeAround digital otoscope.

1. Ensure your ScopeAround otoscope is clean and powered off.
2. Select one ear speculum from the set.
3. Gently push the wider end of the speculum onto the tip of the otoscope camera until it is securely seated. Ensure it is firmly attached and does not wobble.
4. Verify that the speculum does not obstruct the camera lens or the LED lights of the otoscope.



Image: Proper attachment of the ear speculum to the otoscope.

4. Operating

The ear specula are designed to be used in conjunction with your ScopeAround digital otoscope to safely examine the ear canal.

1. With the speculum securely attached, gently insert the tip into the ear canal. The speculum acts as a guide and a protective barrier.
2. Always insert slowly and carefully, observing the live feed from your otoscope. Do not force the speculum into the ear canal.
3. The speculum prevents the otoscope from being inserted too deeply, safeguarding the eardrum from accidental contact.
4. After examination, carefully withdraw the otoscope from the ear.

SAFETY AND PROTECTION

Effective injury prevention protects
the ear canal

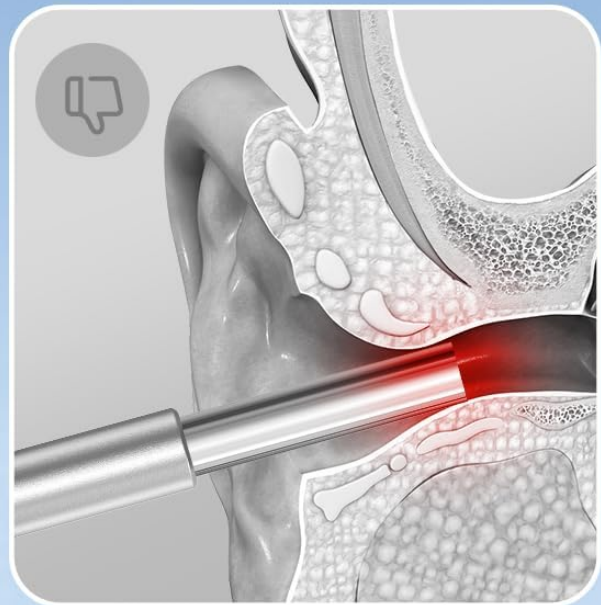
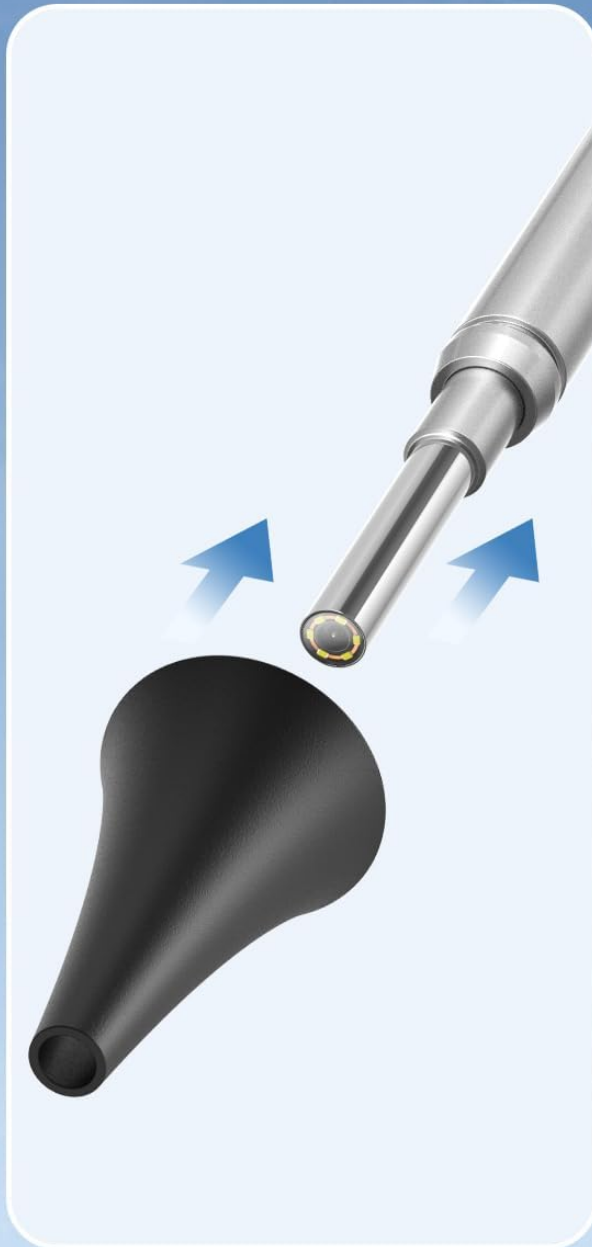


Image: Diagram showing the speculum's role in protecting the ear canal.

Safety and Protection

Effective injury prevention protects the ear canal

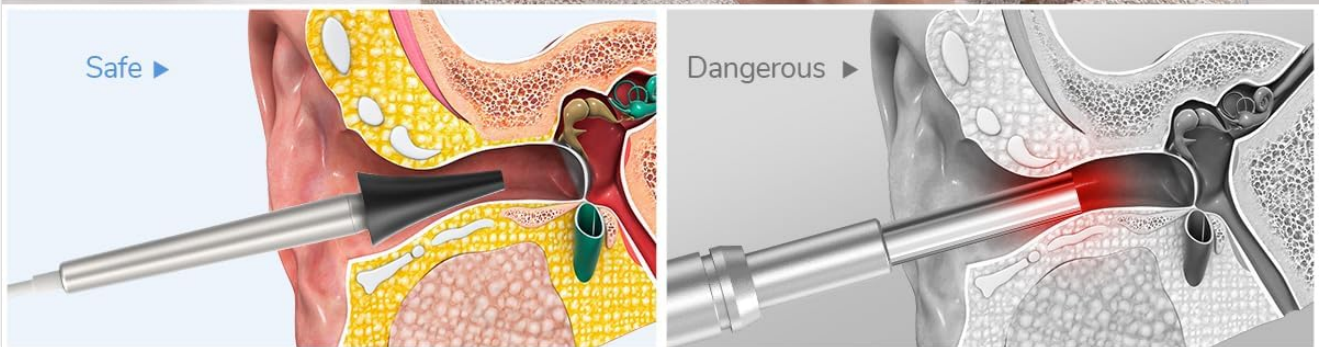


Image: Visual representation of safe usage with the speculum.



Ear Specula Protection

Easy to hurt the eardrum



Image: Comparison of safe (with speculum) vs. dangerous (without speculum) insertion.

5. Maintenance

The ear specula are designed for single-use to ensure maximum hygiene and prevent cross-contamination.

- After each use, remove the used speculum from the otoscope.
- Dispose of the used speculum immediately in a waste bin. Do not reuse specula.
- Clean the tip of your ScopeAround otoscope according to its specific instruction manual after removing the speculum.

6. Troubleshooting

- **Speculum does not fit the otoscope:** Ensure your otoscope is one of the compatible models (SA39W, SA39i, M28, MS450) and has a 4.3mm or 3.9mm diameter tip. If it is, try gently twisting the speculum as you push it on. Do not force it.
- **View is unclear with speculum attached:** Check that the speculum is fully seated and not obstructing the camera lens or LED lights. Ensure the otoscope lens itself is clean.

- **Speculum feels loose:** Ensure it is pushed on firmly. If it continues to feel loose, it may be incompatible or defective. Use a new speculum.

7. Specifications

- **Product Name:** ScopeAround Ear Specula Set
- **Quantity:** 50 pieces
- **Compatibility:** ScopeAround SA39W, SA39i, M28, MS450 digital otoscopes (suitable for otoscopes with a diameter of 4.3mm & 3.9mm)
- **Material:** (Not specified in product data, typically medical-grade plastic)
- **Package Dimensions:** 6.61 x 5.59 x 1.85 inches
- **Package Weight:** 2.08 ounces
- **Manufacturer:** ScopeAround

8. Warranty and Support

For any questions regarding the ScopeAround Ear Specula Set, including warranty information or technical support, please contact ScopeAround customer service directly. As the manufacturer and seller, ScopeAround is best equipped to assist you with product-specific inquiries.

Please refer to the contact information provided with your original otoscope purchase or visit the official ScopeAround website for support details.

© 2025 ScopeAround. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question