



[Manuals.plus](#) /

› [Medacure](#) /

› Medacure Full Body Patient Lift Sling MC-8503XL User Manual

Medacure MC-8503XL

Medacure Full Body Patient Lift Sling MC-8503XL User Manual

Model: MC-8503XL

1. INTRODUCTION

This manual provides essential instructions for the safe and effective use of the Medacure Full Body Patient Lift Sling MC-8503XL. Please read this manual thoroughly before initial use and retain it for future reference. Proper use ensures patient safety and extends the product's lifespan.

2. SAFETY INFORMATION

- **Weight Capacity:** Do not exceed the maximum weight capacity of 850 lbs (385 kg). Ensure the patient's weight is within this limit.
- **Compatibility:** This sling is designed for universal compatibility with most electric and manual patient lifts, including Medacure, Hoyer, ProHeal, Invacare, and Drive brands. Always verify compatibility with your specific lift system before use.
- **Inspection:** Before each use, inspect the sling for any signs of wear, tears, fraying, or damage to the fabric, stitching, or loops. Do not use a damaged sling.
- **Professional Assistance:** Patient transfers should always be performed by trained caregivers or medical professionals. Never attempt to transfer a patient alone if you are not adequately trained or if the patient's condition requires multiple caregivers.
- **Patient Comfort and Safety:** Ensure the patient is properly positioned and comfortable in the sling before lifting. Maintain constant communication with the patient during the transfer process.
- **Secure Attachment:** Always ensure all sling loops are securely attached to the lift's hanger bar or hooks according to the lift manufacturer's instructions.

3. PRODUCT OVERVIEW

The Medacure Full Body Patient Lift Sling MC-8503XL is a durable polyester transfer sling designed for safe and comfortable patient transfers. It features universal compatibility and is available in 2XL and 3XL sizes, both supporting up to 850 lbs.

Key Features:

- **Universal Compatibility:** Works with most major patient lift brands.
- **Full Body Support:** Designed to provide comprehensive support, distributing patient weight evenly.
- **Pressure Point Reduction:** Helps minimize pressure points for enhanced patient comfort during transfers.
- **Color-Coded Straps:** Facilitate various positioning options (sitting, reclining) and ensure correct attachment.
- **High Weight Capacity:** Supports up to 850 lbs.

Universal Compatibility

Works with Medacure and
Most Patient Lift Brands



Image: Universal compatibility with various patient lift systems.


Color Coded Loops



**For Multiple Positions
to Ensure
Comfort & Safety**



Image: Detail of color-coded loops for secure and versatile positioning.

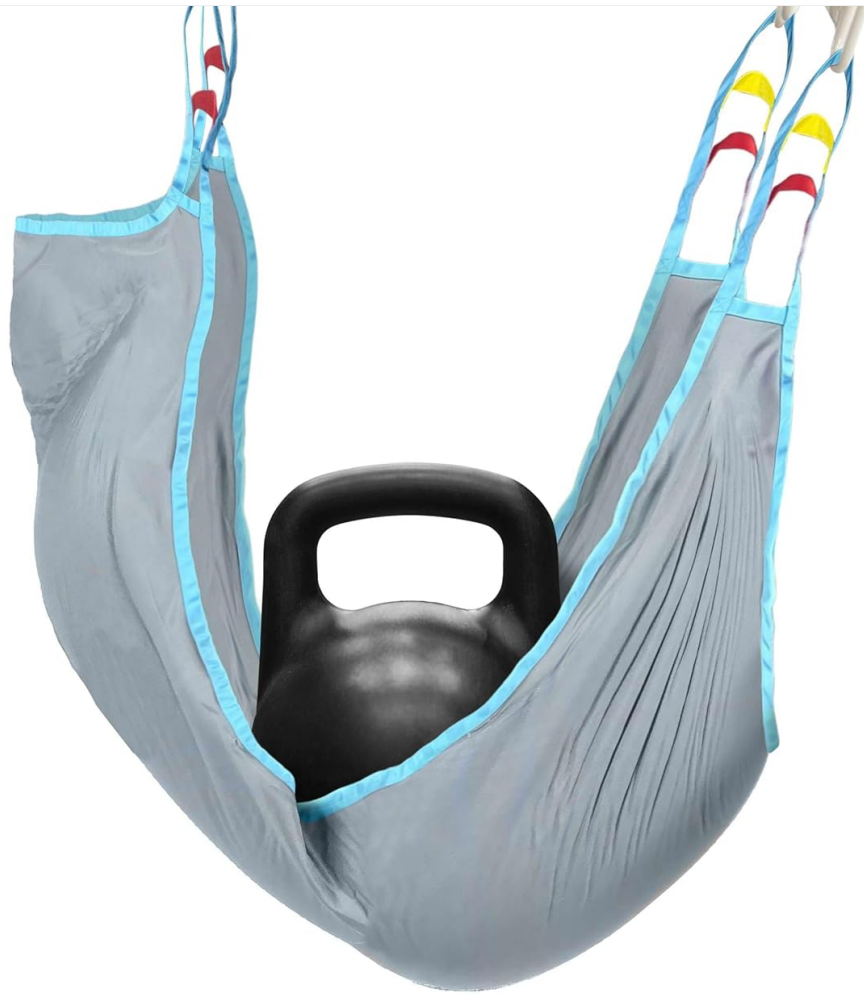

Weight Capacity
850 lbs.

Sizes
2XL and 3XL

Image: The sling demonstrating its 850 lb weight capacity and available sizes.

Full Body Sling

- ✓ Full Body Support
- ✓ Distributes Patient Weight
- ✓ Reduces Pressure Points



Image: The full body sling providing comprehensive support.

4. SETUP

1. **Prepare the Patient Lift:** Ensure your patient lift is in good working condition and positioned correctly for the transfer.
2. **Select Sling Size:** Choose the appropriate sling size (2XL or 3XL) based on the patient's body dimensions to ensure a secure and comfortable fit.
3. **Position the Sling:** Carefully slide the sling under the patient while they are in a supine (lying down) or seated position. Ensure the sling is centered and positioned to provide full body support, with the leg sections correctly placed under the patient's thighs.
4. **Attach to Lift:** Connect the sling's color-coded loops to the corresponding hooks or attachment points on the patient lift's hanger bar. Refer to your patient lift's manual for specific attachment procedures. Ensure all loops are securely fastened.
5. **Verify Connections:** Double-check all connections to ensure they are secure and properly engaged before initiating any lift.

5. OPERATING INSTRUCTIONS

Always follow the operating instructions provided with your specific patient lift system in conjunction with these sling instructions.

1. **Initial Lift:** Slowly raise the patient lift a few inches off the bed or chair. Observe the patient's comfort and the sling's positioning. Adjust if necessary.
2. **Full Lift:** Once satisfied with the positioning, continue to slowly raise the patient to the desired height for transfer.
3. **Transfer:** Carefully maneuver the patient lift to the destination (e.g., wheelchair, commode, bed).
4. **Lowering:** Slowly lower the patient into the new position, ensuring they are stable and comfortable before detaching the sling.
5. **Detaching Sling:** Once the patient is safely seated or lying down, detach the sling loops from the lift's hanger bar.
6. **Positioning for Specific Needs:** The color-coded straps allow for various patient positions. Experiment with different loop combinations to achieve optimal sitting or reclining angles for patient comfort and safety during transfers to a commode/toilet or bed.

6. MAINTENANCE AND CARE

- **Cleaning:** The sling is made of polyester. It can be machine washed with mild detergent in cold water. Do not use bleach.
- **Drying:** Air dry or tumble dry on low heat. Do not iron.
- **Storage:** Store the sling in a clean, dry place away from direct sunlight and extreme temperatures.
- **Regular Inspection:** Periodically inspect the sling for any signs of wear, damage, or degradation. Pay close attention to stitching, seams, and attachment loops. Replace the sling immediately if any damage is found.

7. TROUBLESHOOTING

- **Sling feels insecure:** Ensure all loops are correctly and fully engaged with the lift's hanger bar. Re-check the patient's positioning within the sling.
- **Patient discomfort during lift:** Lower the patient immediately. Re-evaluate the sling's placement and the chosen loop configuration. Adjust as needed to ensure even weight distribution and proper support.
- **Difficulty attaching loops:** Ensure the sling is fully extended and not twisted. Verify that the lift's hanger bar is at an appropriate height for attachment.
- **Sling appears damaged:** Discontinue use immediately. Do not attempt to repair the sling. Contact Medacure customer support or your supplier for replacement.

8. SPECIFICATIONS

Brand	Medacure
Model Name	Full Body Sling 3XL
Item Model Number	MC-8503XL

Material	Polyester & Polyester Blend
Color	Blue
Product Dimensions	11.5"D x 11.5"W x 4.5"H (packaged)
Item Weight	3 Pounds
Weight Capacity	850 lbs
Sizes Available	2XL, 3XL
UPC	850062606456

9. WARRANTY AND SUPPORT

For specific warranty information regarding your Medacure Full Body Patient Lift Sling MC-8503XL, please refer to the documentation provided at the time of purchase or contact Medacure directly. Warranty terms may vary.

For product support, inquiries, or to purchase replacement parts, please visit the official Medacure store or contact their customer service department.

Medacure Store: [Visit Medacure Store on Amazon](#)