

manuals.plus /

› [Medela](#) /

› [Medela Swing Maxi Hands-Free Double Electric Breast Pump User Manual](#)

**Medela 101046709**

# Medela Swing Maxi Hands-Free Double Electric Breast Pump

Model: 101046709

## 1. INTRODUCTION

---

The Medela Swing Maxi Hands-Free Double Electric Breast Pump is designed for mothers seeking an efficient and comfortable pumping experience. Its hands-free design allows for multitasking, making it ideal for busy lifestyles. This pump utilizes advanced technology to optimize milk expression, providing both effectiveness and comfort.

### Key Features

- **More Milk, More Comfort:** Powerful suction combined with anatomically shaped collection cups, featuring a patented 105° angle for increased milk expression and comfort.
- **Ultra-Lightweight:** Each cup weighs just 76g, designed to fit comfortably and discreetly in your bra.
- **3 Parts, Easy to Clean:** Simple assembly and disassembly with only three dishwasher-safe parts per cup.
- **Transparent Material:** Allows for easy monitoring of nipple placement and milk flow.
- **App Connectivity:** Connects with the Medela Family app for real-time tracking, personalized content, and pumping tips.



The ideal  
hands-free solution  
for mums on-the-go.



**Swing Maxi™** Simple. Effective. Hands-Free.

Hands-free

Figure 1: Medela Swing Maxi Hands-Free Double Electric Breast Pump.



## Transparent Material

Check your nipple is placed correctly in the tunnel for effective flow of milk without rubbing while also seeing when and how milk is flowing.



*Figure 2: Comfortable and lightweight design for ease of use.*

## Turning Science into Care



**60+ Years of  
Research**



**Trusted by  
Millions of Mums**



**Swiss Design**

*Figure 3: Discreet design allows for pumping anytime, anywhere.*

## 2. SAFETY INFORMATION

Read all instructions before using the breast pump. Keep this instruction manual for future reference. Do not use the breast pump while bathing or showering. Do not immerse the motor unit in water or other liquids. Always unplug the breast pump immediately after use and before cleaning. Keep out of reach of children.

## 3. COMPONENTS INCLUDED

Your Medela Swing Maxi Hands-Free Double Electric Breast Pump package includes the following items:

- Swing Maxi Hands-Free electric breast pump motor unit
- 2 x Hands-free Collection Cups
- 2 x 21 mm Hands-free breast shields
- 2 x 24 mm Hands-free breast shields

- Tubing
- Lanyard
- Power adaptor and USB-C charging cable

## Medela Family

The Medela Family app is your personal assistant. You can easily track important stats about pregnancy, your baby and breastfeeding, plus get tips and support to help you meet your goals.

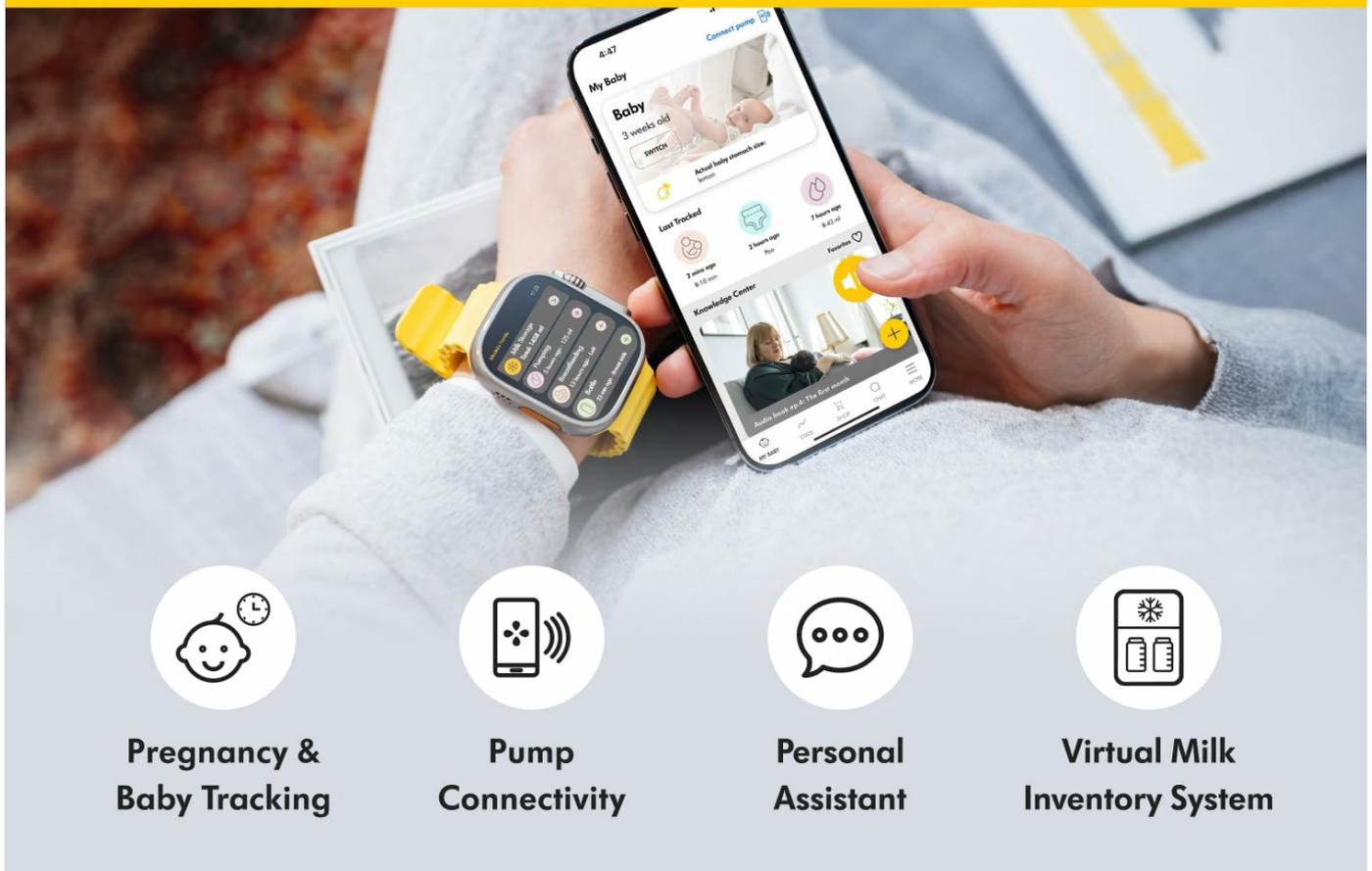


Figure 4: Overview of all included components.

## 4. SETUP AND ASSEMBLY

Before first use, clean and sanitize all parts that come into contact with breast milk. Ensure all parts are completely dry before assembly.

1. Insert the membrane into the breast shield connector.
2. Attach the breast shield to the collection cup.
3. Connect the tubing to the breast shield connector and the motor unit.
4. Ensure all connections are secure to maintain proper suction.



## Intuitive Control

Comes with a user-friendly 4-button system for easy control and a straightforward pumping experience.



*Figure 5: The collection cups consist of only three parts for easy assembly.*

## 5. OPERATION

---

The Medela Swing Maxi features intuitive controls for a personalized pumping experience.

### Controls

- **Power Button:** Turns the pump on/off.
- **Mode Button:** Switches between stimulation and expression modes.
- **Plus (+) Button:** Increases vacuum level.
- **Minus (-) Button:** Decreases vacuum level.

Place the assembled collection cups inside your bra, ensuring your nipple is centered in the breast shield tunnel. Connect the tubing to the motor unit. Press the power button to start pumping. Adjust the vacuum level to your comfort. The transparent design of the cups allows you to visually confirm proper nipple placement and milk flow.



### **3 Easy-to-Clean Parts**

Only 3 parts for easy cleaning and assembly.  
Dishwasher safe.



*Figure 6: Intuitive 4-button control panel for easy operation.*



## Ideal for on the go

The pump remains in your pocket, avoiding added weight on your breast.



Figure 7: Transparent collection cups allow visual monitoring of milk flow.

Your browser does not support the video tag.

Video 1: Official Medela video demonstrating the Swing Maxi Hands-free Breast Pump in use, showcasing its features and ease of operation.

## 6. CLEANING AND MAINTENANCE

---

Proper cleaning and maintenance are essential for hygiene and optimal pump performance. All parts that come into contact with breast milk should be cleaned after each use.

### Cleaning Instructions

1. Disassemble all parts that come into contact with breast milk (breast shields, connectors, membranes, collection cups).
2. Rinse parts with cold water to remove milk residue.
3. Wash parts in warm, soapy water using a mild dish soap.

4. Rinse thoroughly with clean water.
5. Alternatively, all parts are dishwasher-safe (top rack recommended).
6. Air dry on a clean towel or drying rack. Ensure all parts are completely dry before reassembly or storage.
7. The motor unit and tubing should be wiped clean with a damp cloth. Do not immerse them in water.



## Intuitive Control

Comes with a user-friendly 4-button system for easy control and a straightforward pumping experience.



*Figure 8: The three easy-to-clean parts of the collection cups.*

## 7. MEDELA FAMILY APP CONNECTIVITY

Enhance your pumping journey by connecting your Medela Swing Maxi to the Medela Family app. The app provides a personal assistant for tracking, personalized content, smart tools, and pumping tips.

### App Features

- Automatic tracking of pumping sessions in real-time.
- Personalized, science-based content and tips.

- Pregnancy and baby tracking.
- Virtual milk inventory system.



## Discreet Design

Fits comfortably and discreetly in your bra, allowing you to pump anytime, anywhere.



*Figure 9: The Medela Family app provides comprehensive support and tracking.*

## 8. TROUBLESHOOTING

---

If you encounter any issues with your Medela Swing Maxi, please refer to the following common troubleshooting tips:

- **No Suction or Weak Suction:** Check that all parts are correctly assembled and securely connected. Ensure the tubing is not kinked or damaged. Verify that the breast shields are the correct size and form a tight seal against your breast.
- **Pump Not Turning On:** Ensure the pump is fully charged or connected to a power source. Check the power adapter and USB-C cable for any damage.
- **Low Milk Output:** Ensure you are using the correct breast shield size. Try adjusting the vacuum level and cycle speed to find what is most comfortable and effective for you. Ensure you are well-hydrated and relaxed during

pumping sessions.

- **Pain During Pumping:** This may indicate an incorrect breast shield size or improper nipple alignment. Adjust the breast shield and ensure your nipple is centered. Reduce the vacuum level if it is too high. If pain persists, consult a lactation specialist.

For further assistance, please contact Medela customer support.

## 9. TECHNICAL SPECIFICATIONS

---

Feature	Specification
Brand Name	Medela
Model Number	101046709
Operation Mode	Electric
Additional Features	Double Pump, Hands Free Usage, Lightweight, Portable, Wearable
Material Type	Plastic
Color	White
Unit Count	1 Count
UPC	020451467092

## 10. WARRANTY AND SUPPORT

---

Medela products are designed for durability and reliability. For specific warranty details, please refer to the warranty card included in your product packaging or visit the official Medela website. For any questions, technical support, or service inquiries, please contact Medela customer service through their official channels.