

[manuals.plus](#) /

› [BEBIRD](#) /

› [BEBIRD Home 30S Ear Care System - Ear Canal Cleaner with Flexible Scope FHD Camera, UVC LED Light, Visual Ear-Drops Dropper, Wax Removal Tweezers Tool to Clean Ear, 3 Set Cleaning Kits for Adult/Kids](#)

BEBIRD Bebird H30 Smart Cleaning Ear Device

BEBIRD Home 30S Ear Care System

INSTRUCTION MANUAL

1. Introduction

The BEBIRD Home 30S is a state-of-the-art ear care system designed to provide safe, effective, and precise ear cleaning and treatment at home. It integrates a flexible, high-definition camera, a UVC light for hygiene, and a unique visual ear-drops dropper, along with a comprehensive set of tools for various ear care needs. This manual provides essential information for the proper setup, operation, and maintenance of your device.



Figure 1: Overview of the BEBIRD Home 30S Ear Care System and its various components.

2. Setup

2.1 Unboxing and Initial Charge

Carefully remove all components from the packaging. Before first use, ensure the main device is fully charged. Connect the charging cable to the device and a suitable power source. The charging indicator will provide status feedback.

2.2 App Installation and Connection

The BEBIRD Home 30S operates in conjunction with the official Bebird app. Download the latest version of the Bebird app from your device's app store (available for both iOS and Android). Once installed, open the app and follow the on-screen instructions to connect your BEBIRD Home 30S device via Wi-Fi. For stable image transmission, ensure your device's Bluetooth and NFC are turned off during connection, as this has been noted to improve connectivity.

IoT-Enabled Ear Care System in Home



IoT-Enabled Home 30S will gradually unveil & refine more features via the latest version of Bebird app.

✓ Sharing Mode ✓ Ear canal assessment



4 using modes in bebird App for anyone in home



Figure 2: The Bebird app interface, showcasing various modes for different users and ear care needs.

3. Operating Instructions

3.1 Flexible Scope and FHD Camera

The BEBIRD Home 30S features a flexible scope with a 3.1mm Full-HD camera aperture, providing a clear, real-time view of the ear canal. This allows for precise navigation and identification of earwax or other foreign objects. Gently insert the scope into the ear canal while observing the live feed on the Bebird app. Short press the designated button on the device or within the app to take photos, and long press to record videos.

World's First Flexible & FHD In-ear Camera



Figure 3: The flexible scope and FHD camera provide a detailed view of the ear canal, enabling precise and safe operation.



Video 1: Demonstration of the flexible scope and FHD camera in action, showing its ability to navigate and provide clear visuals within the ear canal.

3.2 UVC Light for Ear Hygiene


The device incorporates a UVC light for enhanced hygiene. After the ear canal illuminated LED light is turned off, the UVC light automatically activates for 9 seconds, suppressing 99.9% of unseen microorganisms on the earscope tip. This feature ensures a clean tip for each use. Avoid direct eye contact with the UVC light.

Innovative UVC Light Auto-On Ear Hygiene

9-second disinfection for clean, hygienic use every time

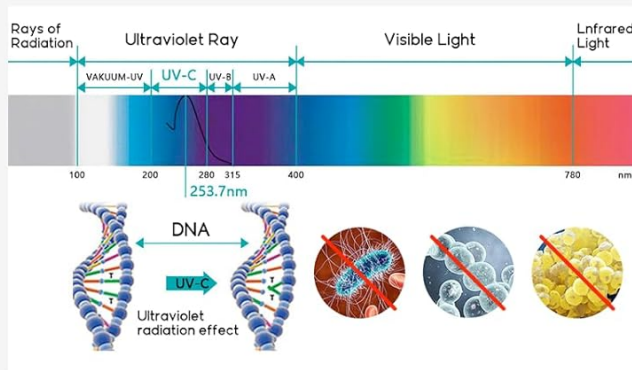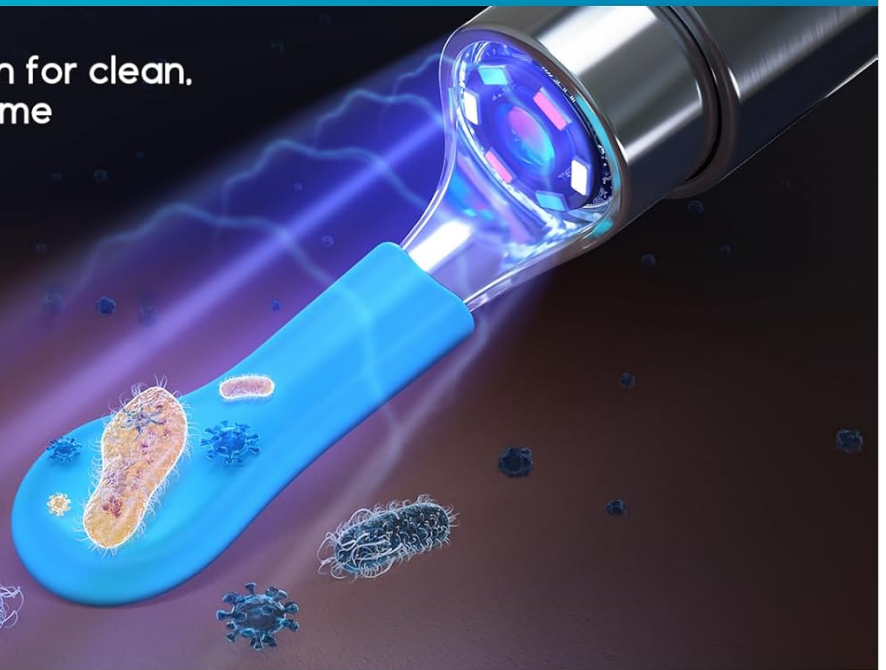

Before Sterilization

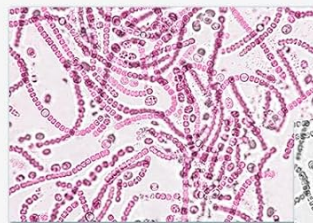

After sterilization



Upon turning off the device, the UVC light automatically activates, eradicating 99% of bacteria in 9 seconds. Avoid direct eye contact!

Figure 4: The UVC light feature automatically disinfects the earscope tip after use, promoting hygiene.

3.3 Visual Ear Medicine Dropper

The BEBIRD Home 30S includes a visual dropper feature, allowing for precise application of ear drops or oil. This is particularly useful for treating wax buildup, ear canal infections, or other ear conditions where targeted liquid delivery is necessary. Use the live camera feed to guide the dropper to the exact location requiring treatment.

Visual & Precise Dropper 🦻 Ear Treatment 🦻

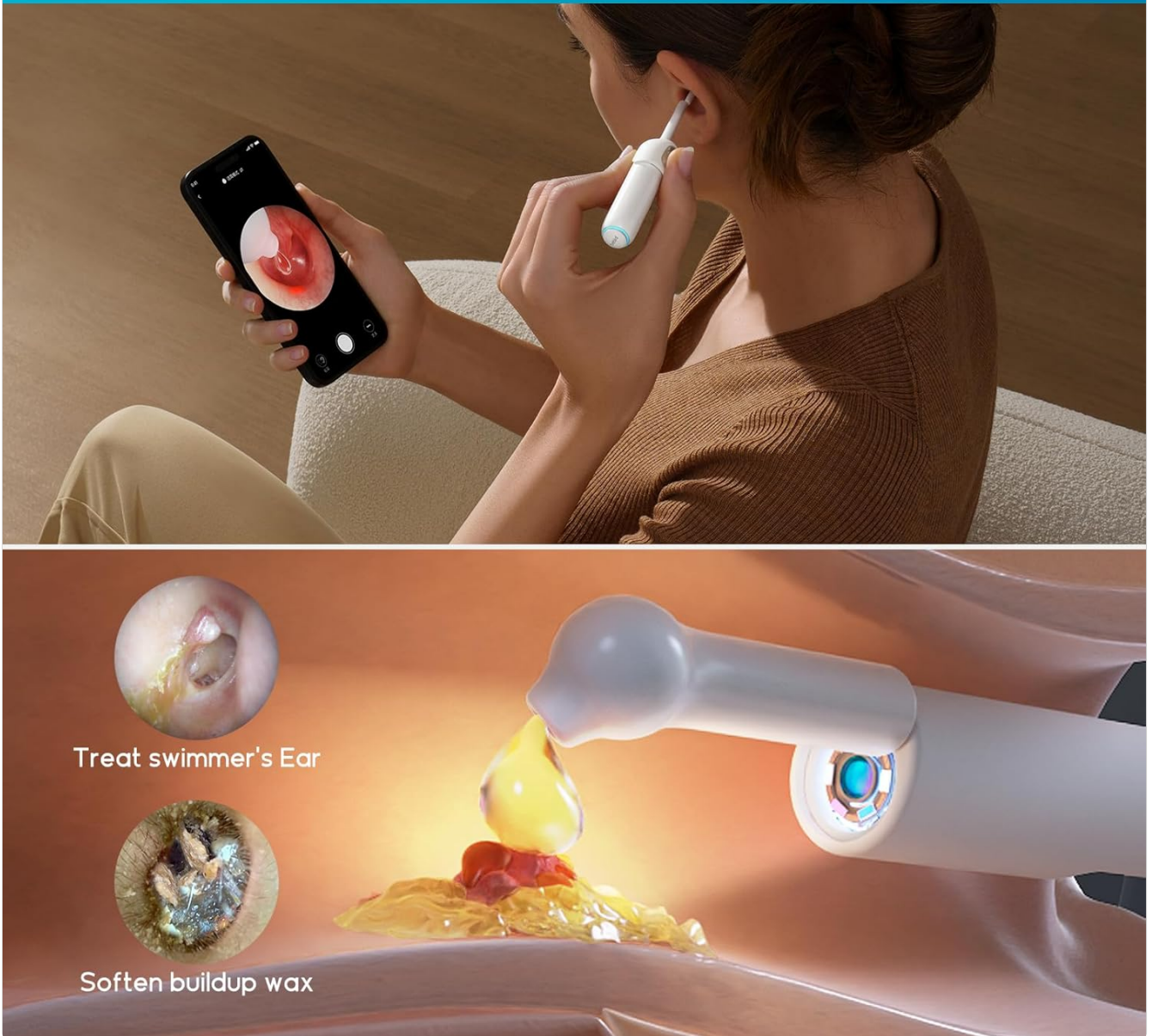


Figure 5: The visual ear medicine dropper allows for accurate and targeted application of liquids within the ear canal.

3.4 Comprehensive Cleaning Tool Set

The system comes with over 30 attachments, including various tips and tools designed for different types of earwax and cleaning scenarios. These are housed in magnetically-attached bases for organized storage. The set includes:

- **Grabbing Tweezers:** For firm wax or foreign objects.
- **Sponge Picks:** Ideal for runny or oily wax.
- **Gel-Stick Ends:** For dry wax debris.
- **Spoon/Scoop:** Suitable for dry soft wax.
- **Rubber Digger:** Effective for both dry and wet wax.

Ultra-soft silicone ear tips are included for the safest and gentlest care, especially for newborn infants.



Figure 6: A selection of the comprehensive ear canal cleaning tools included with the BEBIRD Home 30S.

3.5 Safety Guidelines

Always use the device with caution. Ensure a clear view of the ear canal via the app before attempting any removal. For children, adult supervision is required. The flexible design and soft silicone tips are intended to minimize risk, but gentle handling is paramount. If you experience any discomfort or pain, discontinue use immediately.

Silicone Material & Safe Design Safe for kids

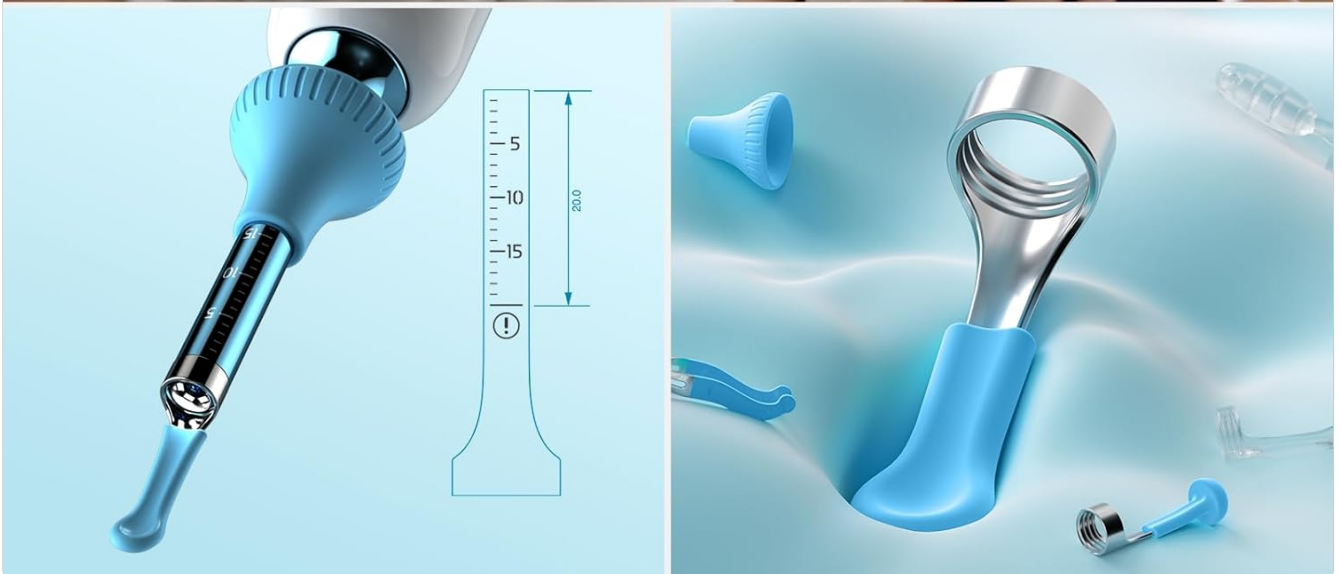


Figure 7: The BEBIRD Home 30S is designed with safety in mind, including soft silicone tips suitable for children under adult supervision.

4. Maintenance

4.1 Cleaning the Device and Attachments

After each use, clean the earscope tip and any used attachments thoroughly. Use a soft, damp cloth to wipe the main device. The UVC light provides automatic disinfection for the tip, but physical cleaning is also recommended. Ensure all parts are dry before storing.

4.2 Storage

Store the BEBIRD Home 30S and its accessories in the provided storage base in a cool, dry place, away from direct sunlight and extreme temperatures. Keep out of reach of children.

4.3 Battery Care

The device is equipped with a long-lasting battery, providing up to 80 minutes of working time. To prolong battery life, avoid fully discharging the device frequently. Charge it regularly, even if not in constant use.

5. Troubleshooting

Problem	Possible Cause	Solution
Device not connecting to app	Bluetooth/NFC interference, incorrect Wi-Fi settings, app issues.	Turn off Bluetooth and NFC on your phone. Ensure you are connected to the device's specific Wi-Fi network. Restart the app and the device.
Image transmission unstable/laggy	Weak Wi-Fi signal, intelligence sensor chip processing.	Ensure a strong Wi-Fi connection between the device and your phone. Move closer to the device. The intelligence sensor chip is designed for stable transmission, but environmental factors can affect it.
Device not turning on	Low battery, power button issue.	Charge the device fully. If the issue persists, contact customer support.
Attachments not fitting securely	Incorrect attachment, debris in connection point.	Ensure you are using the correct attachment for the specific task. Clean the connection point on both the device and the attachment.

6. Specifications

Model Number: Bebird H30 Smart Cleaning Ear Device

Camera Resolution: Full HD (10 Megapixel)

Scope Diameter: 3.1mm

Battery Life: Up to 80 minutes working time

Disinfection: Built-in UVC light (99.9% effective in 9 seconds)

Connectivity: Wi-Fi enabled, IoT-enabled

Compatibility: iOS and Android devices

Dimensions: 6.54 x 5.79 x 2.72 inches

Weight: 1.61 Pounds

Manufacturer: Heifeng Zhizao (Shenzhen) Technology Co., Ltd




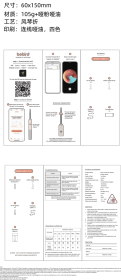
7. Warranty and Support



The BEBIRD Home 30S comes with a **1-year warranty** from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use. Please retain your proof of purchase for warranty claims.

For any problems encountered during use, technical assistance, or warranty inquiries, please contact BEBIRD customer support. Refer to the product packaging or the official Bebird website for the most current contact information.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

Related Documents - Bebird H30 Smart Cleaning Ear Device

	<p>bebird Home 30S Smart Ear Cleaning Device User Manual</p> <p>Comprehensive user manual for the bebird Home 30S Smart Ear Cleaning Device, covering setup, app connection, accessory usage, troubleshooting, and safety guidelines.</p>
	<p>bebird R3 Smart Visual Ear-Clean Rod User Manual</p> <p>Comprehensive user manual for the bebird R3 Smart Visual Ear-Clean Rod, detailing setup, usage, troubleshooting, and specifications. Learn how to connect, replace accessories, and safely use this innovative ear cleaning device.</p>
	<p>Bebird C3 Smart Visual Ear Cleaning Rod User Guide</p> <p>This document provides instructions and specifications for the Bebird C3 Smart Visual Ear Cleaning Rod, including setup, usage, troubleshooting, and product details.</p>
	<p>Bebird K3 Visual Ear Cleaner User Manual and Specifications</p> <p>Detailed user guide for the Bebird K3 and K3 Pro visual ear cleaners, covering app setup, Wi-Fi connection, accessories, troubleshooting, technical specifications, safety warnings, and FCC compliance.</p>

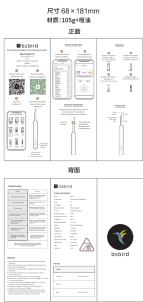	<p>Bebird C3 Smart Visual Ear Cleaning Rod - User Manual and Specifications</p> <p>Comprehensive guide for the Bebird C3 smart visual ear cleaning rod, covering setup, usage, troubleshooting, product specifications, and warranty information. Learn how to connect, use accessories, and maintain your device.</p>
	<p>Bebird Note3 Smart Visual Ear Cleaning Rod User Manual and Specifications</p> <p>Comprehensive guide to the Bebird Note3 smart visual ear cleaning rod, covering setup, usage, troubleshooting, and product specifications. Learn how to use the device for safe and effective ear cleaning.</p>