



Manuals.plus /

- › QUMOX /
- › QUMOX Electric Ear Wax Removal Kit E30 User Manual

QUMOX E30

QUMOX Electric Ear Wax Removal Kit E30 User Manual

Model: E30

1. INTRODUCTION

This manual provides instructions for the safe and effective use of your QUMOX Electric Ear Wax Removal Kit, Model E30. This device is designed to gently clean ear canals by irrigating them with water, helping to remove earwax, dirt, and debris. Please read this manual thoroughly before first use and retain it for future reference.



Image 1.1: QUMOX Electric Ear Wax Removal Kit in operation.

2. SAFETY INFORMATION

Important: Always consult a medical professional before using this device if you have any pre-existing ear conditions, ear pain, discharge, recent ear surgery, or a perforated eardrum. Do not use this device on children without medical advice.

- Use only lukewarm water. Water that is too hot or too cold can cause dizziness or discomfort.
- Do not insert the ear tip too deeply into the ear canal.
- Discontinue use immediately if you experience pain, dizziness, or hearing changes.
- Keep the device out of reach of children.
- Ensure the device is fully charged before use to avoid interruptions.
- Use only the provided disposable ear tips to maintain hygiene and prevent infection.

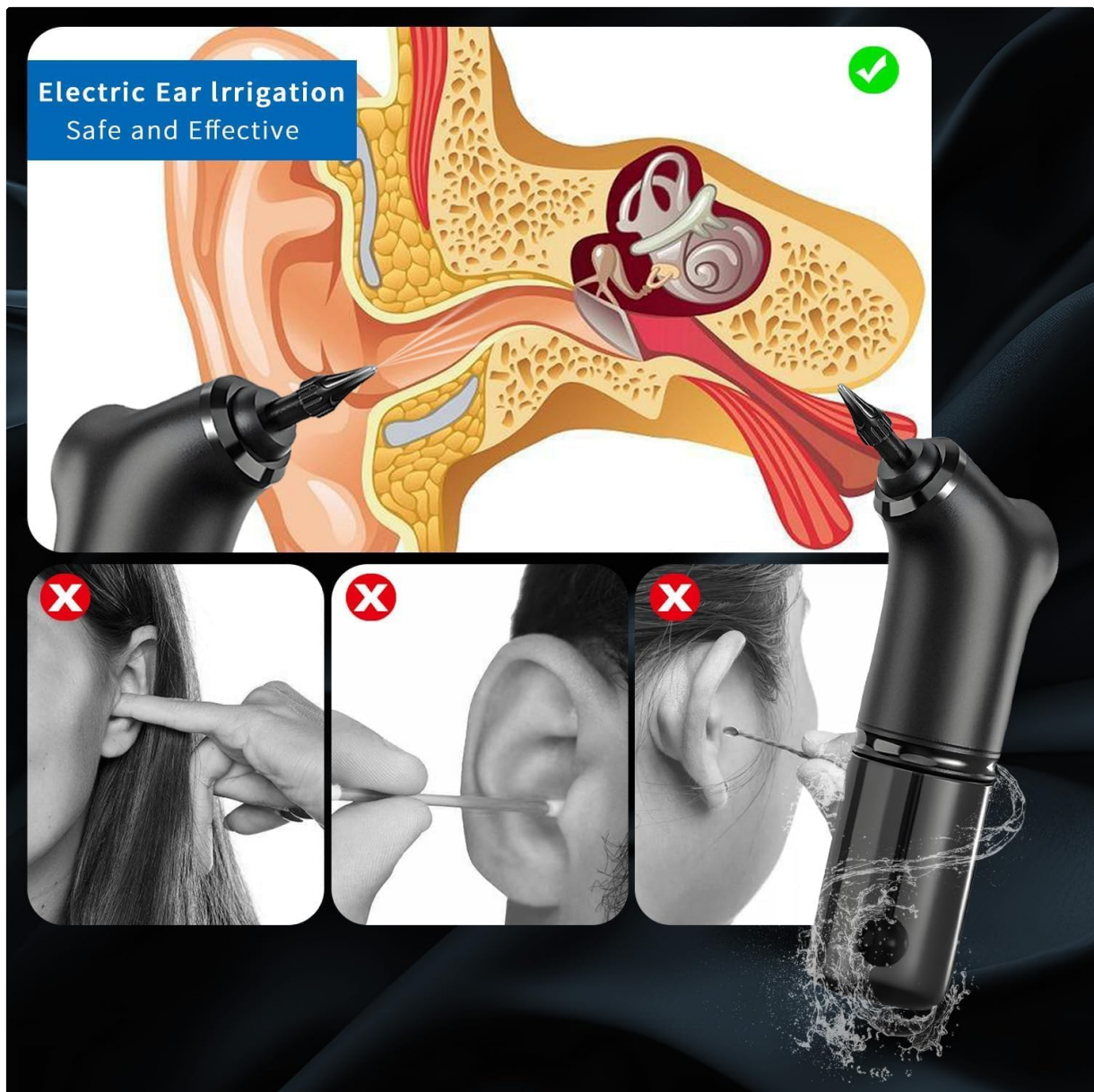


Image 2.1: Illustration of safe ear cleaning practices with the device versus unsafe methods.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- Electric Ear Wax Cleaning Kit (Main Unit)
- Ear Wax Remover with Gravity Straw
- Water Bottle
- Ear Catch Basin
- 10 Disposable Soft Ear Tips
- Type-C Charging Cable
- Instruction Manual (this document)



Image 3.1: All components included in the QUMOX Ear Wax Removal Kit.

4. PRODUCT FEATURES

- **Gentle and Effective Cleaning:** Designed to remove earwax, dirt, and debris without harming sensitive ear canal skin.
- **Hygienic Design:** Includes 10 disposable soft ear tips for sanitary use and reduced risk of infection.
- **Four Adjustable Pressure Settings:** Offers Normal, Soft, Pulse, and Strong modes (7 to 24.5 PSI) for customized cleaning.
- **Long-lasting Battery:** Equipped with a 2000 mAh battery, providing up to 7 hours of use on a full charge.
- **Type-C Charging:** For quick and efficient recharging.

5. SETUP

1. **Charge the Device:** Before first use, fully charge the device using the provided Type-C charging cable. Connect the cable to the device's charging port and a suitable USB power adapter (not included). The indicator light will show charging status.
2. **Prepare Water:** Fill the water bottle with lukewarm water. Ensure the water temperature is comfortable and not too hot or cold.
3. **Attach Ear Tip:** Select a new, disposable soft ear tip and securely attach it to the nozzle of the main unit.
4. **Assemble:** Attach the water bottle to the main unit, ensuring it is securely fastened to prevent leaks.

6. OPERATING INSTRUCTIONS

1. **Positioning:** Lean your head over a sink or use the provided ear catch basin to collect the water and dislodged earwax.
2. **Power On:** Press and hold the power switch until the indicator light blinks.

3. **Select Mode:** Press the "M" button to cycle through the four pressure settings: Normal, Soft, Pulse, and Strong. Choose the setting that feels most comfortable and effective for your needs.
4. **Begin Irrigation:** Gently insert the ear tip into the ear canal. Do not force it. Press the power switch again to start the water flow.
5. **Clean Ear:** Allow the water to irrigate the ear canal. Move the device slightly to ensure thorough cleaning. The water will flow out into the sink or catch basin.
6. **Repeat if Necessary:** If needed, repeat the process for the same ear or switch to the other ear.
7. **Power Off:** Press and hold the power switch to turn off the device after use.



Image 6.1: Device controls and pressure settings.

7. MAINTENANCE AND CLEANING

1. **Discard Ear Tip:** After each use, remove and discard the used disposable ear tip.
2. **Empty Water Bottle:** Empty any remaining water from the bottle.
3. **Rinse Components:** Rinse the water bottle and the main unit's nozzle with clean water.
4. **Wipe Down:** Wipe the exterior of the device with a clean, damp cloth. Do not submerge the main unit in water.

5. **Dry Thoroughly:** Ensure all components are completely dry before storage.

8. STORAGE

Store the QUMOX Electric Ear Wax Removal Kit in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it out of reach of children.

9. TROUBLESHOOTING

| Problem | Possible Cause | Solution |
|-----------------------|---|--|
| Device not turning on | Low battery; Device not properly switched on. | Charge the device fully; Press and hold the power switch until the indicator light blinks. |
| Weak water pressure | Low battery; Water bottle not securely attached; Clogged ear tip. | Charge the device; Ensure the water bottle is tightly screwed on; Replace the ear tip. |
| Water leakage | Water bottle not securely attached; Damaged ear tip. | Ensure the water bottle is tightly screwed on; Replace the ear tip. |
| No water flow | Water bottle empty; Ear tip clogged. | Refill the water bottle; Replace the ear tip. |

10. SPECIFICATIONS

- **Model Number:** E30
- **Battery:** 1 Lithium Ion battery (included), 2000 mAh
- **Charging Port:** Type-C
- **Pressure Settings:** 4 levels (7 to 24.5 PSI)
- **Dimensions:** Approximately 211.6 mm (8.33 inches) height x 56.4 mm (2.22 inches) width
- **Weight:** Approximately 1.04 Pounds





Image 10.1: Device dimensions.

11. WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided at the point of purchase or visit the official QUMOX website. Keep your purchase receipt as proof of purchase.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.