

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [HOMLAND](#) /
- › [HOMLAND Foldable Rollator Walker with Seat, Model R2-RD, User Manual](#)

## HOMLAND R2-RD

# HOMLAND Foldable Rollator Walker with Seat, Model R2-RD

User Manual

## INTRODUCTION

---

This manual provides comprehensive instructions for the safe assembly, operation, maintenance, and troubleshooting of your HOMLAND Foldable Rollator Walker with Seat, Model R2-RD. Designed for seniors, this lightweight rolling walker features 8-inch wheels, adjustable handles and seat height, and a 350 lb weight capacity, ensuring enhanced mobility and comfort. Please read this manual thoroughly before using the rollator to ensure safe and correct operation.



Image: The HOMLAND Foldable Rollator Walker with Seat, featuring a red frame, black seat, and 8-inch wheels.

## **SAFETY INFORMATION**

---

Always prioritize safety when using your rollator. Adhere to the following guidelines:

- Always ensure the frame is fully open and locked before use.
- Confirm both brakes are engaged when sitting on the seat or standing still.
- Regularly check that all wheels are secure and functioning correctly.
- Adjust the handles to your proper height to maintain stability and prevent falls.
- Do not exceed the maximum weight capacity of 350 pounds.
- Avoid using the rollator on uneven or slippery surfaces without extreme caution.

- Ensure the pinch-proof folding points are clear of hands or clothing during folding and unfolding.

## Reliable Braking System

Features three modes (free, deceleration, stop) for ultimate control.



Image: Diagram highlighting the rollator's safety features, including ergonomic hand grips, reliable braking system, multi-reflector design, and pinch-proof design.

## ASSEMBLY INSTRUCTIONS

The HOMLAND Rollator Walker is designed for tool-free assembly. Follow these steps to set up your walker:

1. **Unfold the Frame:** Carefully pull the sides of the rollator apart until the frame is fully extended and locks into place. Ensure the locking mechanism is engaged to prevent accidental folding.
2. **Attach Handles:** Insert the adjustable handle tubes into the main frame. Align the holes and secure them with the provided knobs.
3. **Adjust Handle Height:** Press the adjustment buttons on the handle tubes and slide the handles up or down to your desired height. Release the buttons to lock the handles in place. The ideal height is when your wrists are level with the handgrips when standing upright with arms relaxed.
4. **Verify Stability:** Before first use, ensure all connections are secure and the rollator stands firmly on all four wheels.

**HOMLAND**

## R2 Ultimate Rollator

Strength, Comfort, and Safety Combined!



*Experience the difference with our rollator - where mobility meets comfort and safety!*

Image: A person performing tool-free assembly on the HOMLAND rollator, securing a handle with a knob.

## Enjoy a Wide Breathable Backrest



## Extra-wide Memory Foam Seat

that prioritizes comfort and safety, making every outing a pleasure.

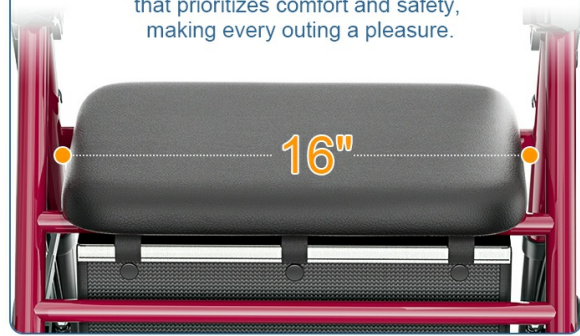


Image: Diagram illustrating the dual-adjustable height for both handles (33"-41") and seat (20"-23"), designed for users from 4'7" to 6'6".

## OPERATING INSTRUCTIONS

Using your HOMLAND rollator effectively involves understanding its features and proper techniques.

### Walking with the Rollator

1. **Positioning:** Stand upright with the rollator slightly ahead of you. Keep your body within the frame of the walker for optimal balance.
2. **Movement:** Push the rollator forward a short distance, then step forward with one foot, followed by the other. Maintain a steady pace.
3. **Braking:** The rollator features a reliable braking system with three modes: free, deceleration, and stop.
  - To **slow down or brake temporarily**, gently squeeze the brake levers on the handles.
  - To **lock the brakes** for sitting or standing still, push the brake levers down until they click into the locked position. Pull them up to release.



### Ergonomic Hand Grips

Soft grips reduce hand fatigue, allowing easy use by individuals with limited strength.



Image: Close-up view of the ergonomic hand grips and brake levers, showing how to engage the brake (squeeze) and lock (push down) functions.

### Using the Seat

1. **Engage Brakes:** Always lock both brakes before attempting to sit on the rollator's seat.

2. **Sit Down:** Turn around and slowly lower yourself onto the extra-wide memory foam seat. Use the handles for support if needed.
3. **Rest:** The breathable backrest provides comfort during rest.
4. **Stand Up:** Use the handles to push yourself up, then release the brakes before resuming walking.

## Under Seat Storage with Side Pocket



Image: Detailed view of the rollator's breathable backrest and extra-wide memory foam seat, emphasizing comfort and support.

### Folding for Storage and Transport

The rollator features a quick-fold design for easy storage and transport.

1. **Clear Storage:** Remove any items from the under-seat storage bag.
2. **Lift Seat:** Lift the seat upwards.
3. **Collapse Frame:** Pull the folding strap or mechanism (usually located under the seat) to collapse the frame inwards. Ensure your hands are clear of the pinch-proof folding points.
4. **Secure:** The rollator will fold into a compact size, suitable for placing in a car trunk or storing in a small space.

## Product Details

### HOMLAND Walkers for Seniors

Max Weight Capacity	350 lb	Handles Height	33"-41"
Product Dimension	28"(L) x 25"(W) x 33-41"(H)	Folding Dimension	34"(L)*25"(W)*10"(H)
Product Weight	20.4 lb	Wheels	8" PU Wheels
Seat (width*depth*height)	16"(W)*12(D)*20"-23"(H)	Seat Height	20-23"
Material of frame	Carbon steel	Material of Seat	PU+Memory foam
Accessories	Storage bag	Material of handles	PVC
Applicable height	4'7"-6'6"	Applicable scenarios	Outdoor & Indoor



US Based Hotline Support



Specialized Technical Support



30-Day Hassle-Free Return

Image: A sequence showing the quick-fold design of the rollator and a person placing the folded rollator into a car trunk.

### MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your rollator.

- **Cleaning:** Wipe down the frame and seat with a damp cloth and mild soap. Avoid abrasive cleaners.
- **Wheel Inspection:** Periodically check the 8-inch puncture-proof wheels for wear and tear. Ensure they rotate freely and are securely attached.
- **Brake Check:** Test the brakes regularly to ensure they engage and lock effectively. If brakes feel loose or ineffective, consult the troubleshooting section or contact support.
- **Fastener Check:** Ensure all screws, bolts, and adjustment knobs are tight and secure.
- **Storage:** Store the rollator in a dry place, away from extreme temperatures.

## TROUBLESHOOTING

Here are solutions to common issues you might encounter:

Problem	Possible Cause	Solution
Rollator does not fold or unfold smoothly.	Obstruction in folding mechanism; mechanism not fully engaged.	Check for any debris or objects obstructing the folding points. Ensure the folding strap is pulled firmly.
Brakes are not engaging or locking properly.	Brake cable too loose; worn brake pads.	Inspect brake cables for damage. If loose, they may need adjustment by a qualified technician. If pads are worn, contact support for replacements.
Wheels are noisy or do not turn freely.	Debris in wheel mechanism; loose wheel axle.	Clean wheels and axles of any dirt or hair. Ensure wheel axles are tightened. If problem persists, contact support.
Handles feel wobbly after adjustment.	Adjustment knobs not fully tightened; handles not locked into height positions.	Ensure the height adjustment buttons have fully popped out into the desired hole. Tighten the securing knobs firmly.

## SPECIFICATIONS

Detailed specifications for the HOMLAND Foldable Rollator Walker, Model R2-RD:



Image: A table summarizing key product details and specifications for the HOMLAND Rollator Walker.

Feature	Detail
---------	--------

Feature	Detail
Model Number	R2-RD
Max Weight Capacity	350 lb
Product Dimensions (L x W x H)	27.5"D x 21.7"W x 33.41"H (Unfolded)
Product Weight	20.4 Pounds
Handle Height Adjustment	33" - 41"
Seat Height Adjustment	20" - 23"
Seat Dimensions (W x D)	16"W x 12"D
Wheel Size	8-inch PU Wheels
Material (Frame)	Carbon Steel (Reinforced)
Material (Seat)	PU + Memory Foam
Applicable User Height	4'7" - 6'6"
Color	Red

## Foldable and Convenient

Its quick-fold design allows  
for easy transport and storage



Image: Visual representation of the rollator's 350 lbs weight capacity, highlighting its robust construction.

## Tool-Free Assembly



Image: Comparison showing the HOMLAND rollator's 8-inch large wheels versus smaller 6-inch wheels found on other models, emphasizing enhanced mobility.

### WARRANTY AND SUPPORT

HOMLAND is committed to providing quality products and excellent customer service.

- **Lifetime Warranty:** The metal frame of your HOMLAND rollator is covered by a lifetime warranty.
- **Parts Warranty:** Non-metal parts are covered by a 1-year warranty, with an option for a 1-year extension.
- **Customer Support:** For any questions, concerns, or warranty claims, please contact our US-BASED customer support.
- **Additional Resources:** Access to user manuals and installation videos is available through HOMLAND customer support channels.