

## HOMLAND C4

# HOMLAND Raised Toilet Seat with Handles

## MODEL: C4

Instruction Manual

### 1. Introduction

The HOMLAND Raised Toilet Seat with Handles, Model C4, is designed to provide assistance for individuals who require a higher toilet seating position and additional support when sitting down or standing up. This product features adjustable height and width, a padded seat for comfort, and sturdy construction to enhance safety and independence in the bathroom.

### 2. Safety Information

- **Weight Capacity:** This raised toilet seat is designed to safely support users up to 400 lbs (181 kg). Do not exceed this weight limit.
- **Stability:** Ensure all legs are evenly adjusted and the unit is stable before each use. The non-slip rubber feet must be in full contact with the floor.
- **Assembly:** Follow all assembly instructions carefully. Incorrect assembly can lead to instability and potential injury.
- **Maintenance:** Regularly inspect the product for any signs of wear, damage, or loose components. Discontinue use if any issues are found and contact support.
- **Placement:** Position the unit centrally over the toilet. Ensure sufficient space around the unit for safe entry and exit.
- **Medical Advice:** Consult with a healthcare professional or physical therapist to determine if this product is suitable for your specific needs.

# Stable & Safety



Image: The HOMLAND Raised Toilet Seat features dual-crossbar support for enhanced stability, safely supporting up to 400 lbs. The aluminum alloy construction uses 2mm thick tubing for durability.

## 3. Package Contents

Please verify that all components are present before beginning assembly:

- Raised Toilet Seat Frame (with pre-attached padded seat and armrests)
- Organizer Pouch (1)
- Toilet Paper Holder (1)
- Assembly Hardware (snap fasteners, screws)

# 4-Step Raised Toilet Seat Fit Check

Confirm Before Buying for Perfect Fit




## Toilet Fit

Recommended Toilet Size  
Height  $\leq 21.3"$  | Width  $\leq 24"$



## Bathroom Fit

Recommended Bathroom  
Width  $\geq 28"$

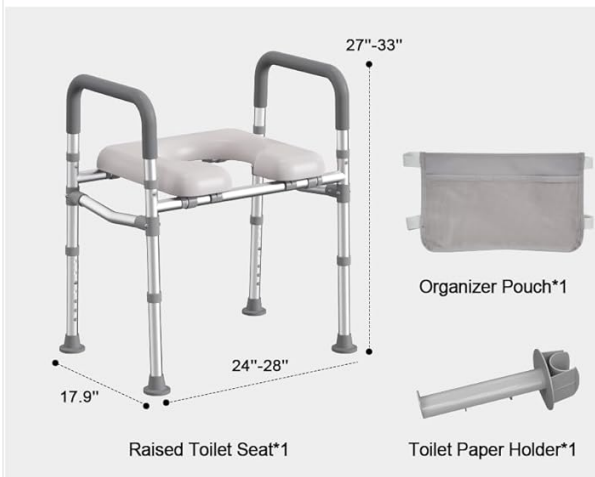

## What's in the Box?

## IT DOES NOT INCLUDE THE FOLLOWING FEATURES:

- Backrest: No backrest included.
- Splash Protection: No anti-splash guard.

## Toilet Seat Riser Only

Image: Included components: the raised toilet seat, an organizer pouch, and a toilet paper holder. Ensure all parts are present before assembly.

## 4. Assembly Instructions

Assembly typically takes less than 10 minutes and requires no additional tools. Follow these steps:

1. Unpack all components and lay them out. Identify the main frame pieces, often labeled 'A' and 'B'. Ensure 'A' parts connect only to 'A' parts, and 'B' parts only to 'B' parts.
2. Connect the transverse tubes to the front and rear connecting tubes. Secure them using the C-type clips provided, ensuring they snap firmly into the holes.
3. Insert the front and rear ends of the handrails into their respective connecting tubes. Ensure the spring pins engage securely into the adjustment holes.
4. Attach the padded seat to the assembled frame. Ensure it is securely fastened.
5. Adjust the height and width as needed (refer to Section 5 for detailed instructions).
6. Attach the organizer pouch and toilet paper holder to the frame as desired.

Your browser does not support the video tag.

Video: This video demonstrates the assembly process for the HOMLAND Raised Toilet Seat. It shows how to connect the frame components and attach the seat and accessories.

## 5. Adjustments (Operating)

---

### 5.1. Height Adjustment

The HOMLAND Raised Toilet Seat offers 7 levels of height adjustment. The seat height can be adjusted from 18.1 inches to 24 inches, and the handle height from 27 inches to 33 inches. To adjust the height:

1. Locate the spring pins on each leg.
2. Press the spring pin inward and slide the inner leg tube up or down to the desired height setting.
3. Ensure the spring pin fully engages into one of the adjustment holes.
4. Repeat for all four legs, ensuring they are set to the same height for stability.

### 5.2. Width Adjustment

The product width is adjustable from 24 inches to 28 inches to accommodate various toilet sizes and user preferences. To adjust the width:

1. Locate the width adjustment mechanisms on the frame, typically near the base of the armrests.
2. Loosen any securing mechanisms (e.g., knobs or clips) that hold the width in place.
3. Slide the side frames inward or outward to achieve the desired width.
4. Secure the adjustment mechanisms firmly to lock the width in position.

# 6 Adjustable Width



Fit Most  
Toilets



Square Toilet



V Shape Toilet



Elongated Toilet

Image: The toilet seat's width can be adjusted between 24 and 28 inches, fitting various toilet shapes including square, V-shape, and elongated designs.

# 7 Adjustable Height

Suitable for most people of different height



Recommended Toilet Height  $\leq 21.3"$

Image: The seat height is adjustable from 18.1 to 24 inches, and handle height from 27 to 33 inches, with 7 distinct height settings.

## 6. Features

- **Stable & Durable Design:** Constructed from aluminum alloy with a welding reinforced design and dual-crossbar support, ensuring stability and a weight capacity of 400 lbs.
- **Comfortable Padded Seat:** Features an enlarged 16.5"x16.5"x1.6" soft padded seat for comfort during prolonged use. The seat is waterproof for easy cleaning.
- **Adjustable Height & Width:** Seat height adjusts from 18.1" to 24", handle height from 27" to 33", and product width from 24" to 28", accommodating various users and toilets.
- **Non-Slip Armrests & Rubber Feet:** Handles are equipped with soft gray foam grips for secure support. Each leg has non-slip rubber feet for stability and safety.
- **Convenience Features:** Includes a toilet paper holder and a storage bag for added utility.
- **FSA/HSA Eligible:** This product qualifies for FSA/HSA eligibility.

# Comfort Padded Seat

Thicker padded make you more comfortable



Image: The HOMLAND Raised Toilet Seat features a comfortable padded seat and is FSA/HSA eligible.

# More Convenience Features



Image: Additional features include anti-slip padded armrests, an organizer pouch, and a convenient toilet paper holder.

## 7. Care and Maintenance

- **Cleaning:** The padded seat is waterproof and can be easily wiped clean with a damp cloth and mild soap. Avoid abrasive cleaners that may damage the surface.
- **Frame Cleaning:** Clean the aluminum frame with a damp cloth. Dry thoroughly after cleaning.
- **Inspection:** Periodically check all screws, clips, and adjustment pins to ensure they are secure. Verify that the rubber feet are not worn or damaged.
- **Storage:** Store the unit in a dry place away from direct sunlight when not in use for extended periods.

## 8. Troubleshooting

- **Instability:** If the unit feels unstable, ensure all four legs are adjusted to the same height and that the spring pins are fully engaged. Check that the width adjustment mechanisms are securely tightened. Verify that the non-slip rubber feet are clean and making full contact with the floor.
- **Difficulty Adjusting Height/Width:** Ensure no debris is obstructing the adjustment holes or sliding

mechanisms. Apply gentle pressure while pressing the spring pins or loosening adjustment knobs.

- **Missing Parts:** If any parts are missing upon arrival, refer to the 'Package Contents' section and contact customer support immediately.

## 9. Specifications

---

|                          |                                             |
|--------------------------|---------------------------------------------|
| Model Number             | C4                                          |
| Material                 | Aluminum Alloy Frame, Padded Seat           |
| Weight Capacity          | 400 lbs (181 kg)                            |
| Adjustable Seat Height   | 18.1" - 24" (46 cm - 61 cm)                 |
| Adjustable Handle Height | 27" - 33" (68.5 cm - 84 cm)                 |
| Adjustable Product Width | 24" - 28" (61 cm - 71 cm)                   |
| Padded Seat Dimensions   | 16.5" x 16.5" x 1.6" (42 cm x 42 cm x 4 cm) |
| Product Weight           | Approximately 6.6 lbs (3 kg)                |
| Color                    | Gray                                        |

## 10. Warranty and Support

---

HOMLAND is committed to providing quality products and customer satisfaction. If you encounter any issues with your HOMLAND Raised Toilet Seat, please do not hesitate to contact our US-based support team. We aim to respond to all inquiries within 24 hours.

For support, please refer to the contact information provided with your purchase or visit the official HOMLAND website.