

manuals.plus /

› [CANMOV](#) /

› [CANMOV Power Lift Recliner Chair User Manual - Model CD0276](#)

## CANMOV CD0276

# CANMOV Power Lift Recliner Chair User Manual

Model: CD0276

## INTRODUCTION

---

Thank you for choosing the CANMOV Power Lift Recliner Chair. This manual provides essential information for the safe assembly, operation, and maintenance of your new recliner. Please read these instructions thoroughly before use and retain them for future reference.

# Oversized Design

The seating size is up to 25 inches wide,  
spacious and comfortable



Figure 1: The recliner features an oversized design for enhanced comfort.



## IMPORTANT SAFETY INFORMATION

- **Weight Limit:** Do not exceed the maximum weight capacity of 350 lbs (158 kg).
- **Adult Use:** This product is designed for adult use. Children should not operate the chair without adult supervision.
- **Clearance:** Ensure adequate clearance around the chair during operation to prevent obstruction or injury.
- **Power Cord:** Keep the power cord away from heat sources, sharp edges, and wet areas. Do not operate the chair with a damaged cord or plug.
- **Medical Conditions:** If you have any health concerns, consult your doctor before using the massage and heat functions.
- **USB Ports:** The USB ports are designed for low-power devices only, such as smartphones and tablets. Do not attempt to charge high-power devices.
- **Assembly:** Ensure all components are correctly assembled and secured before using the chair.

## PACKAGE CONTENTS

Your CANMOV Power Lift Recliner Chair is shipped in three separate boxes. Components may arrive at different times. Please verify all parts are present before beginning assembly.

# PACKAGE

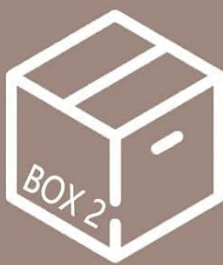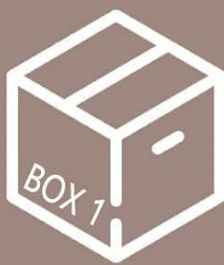

The item comes in 3 BOXES and may not arrive at the same time.

Figure 2: The recliner is delivered in three separate boxes.

Expected components include:

- Chair Base
- Chair Backrest
- Side Arms (Left and Right)
- Power Adapter and Cord
- Remote Control (for massage/heat)
- Assembly Hardware (bolts, washers, etc.)
- Instruction Manual (this document)



## SETUP AND ASSEMBLY

---

Assembly typically requires two people. Follow these steps carefully:

1. **Unpack Components:** Carefully remove all parts from the boxes and place them on a soft, clean surface to prevent damage.
2. **Attach Backrest:** Slide the backrest onto the chair base until it clicks securely into place. Ensure both sides are fully engaged.
3. **Connect Wiring:** Locate the electrical cables extending from the backrest and the chair base. Connect them according to the color-coded or labeled connectors. Ensure all connections are firm.
4. **Attach Side Arms:** Slide the side arms onto the designated slots on the chair base. Secure them with the provided hardware, if applicable.
5. **Connect Power:** Plug the power adapter into the chair's power input, then plug the power cord into a standard electrical outlet.
6. **Initial Test:** Before sitting, perform a brief test of the lift and recline functions using the side control panel to ensure proper operation.



Figure 3: View of the chair's base mechanism, illustrating the sturdy construction.

## OPERATING INSTRUCTIONS

---

### Power Lift and Recline Functions

The chair features a dual-motor system for independent control of the backrest and footrest, as well as a power lift function.

# Lift&Recline Control



Figure 4: Control panel and various chair positions (Reclined, TV, Up).



Figure 5: Detailed view of the integrated control buttons and USB port.

- **Recline:** Use the "Recline" button on the side control panel to adjust the backrest and footrest simultaneously. The chair reclines between 110 and 150 degrees.
- **Lift:** Press the "Lift" button to gently raise the chair, assisting you in standing up.
- **Adjustable Headrest:** Use the "Power Headrest" button to adjust the angle of the headrest independently for optimal neck support.
- **Tilt:** The "Tilt" button allows for slight adjustments to the overall chair angle.

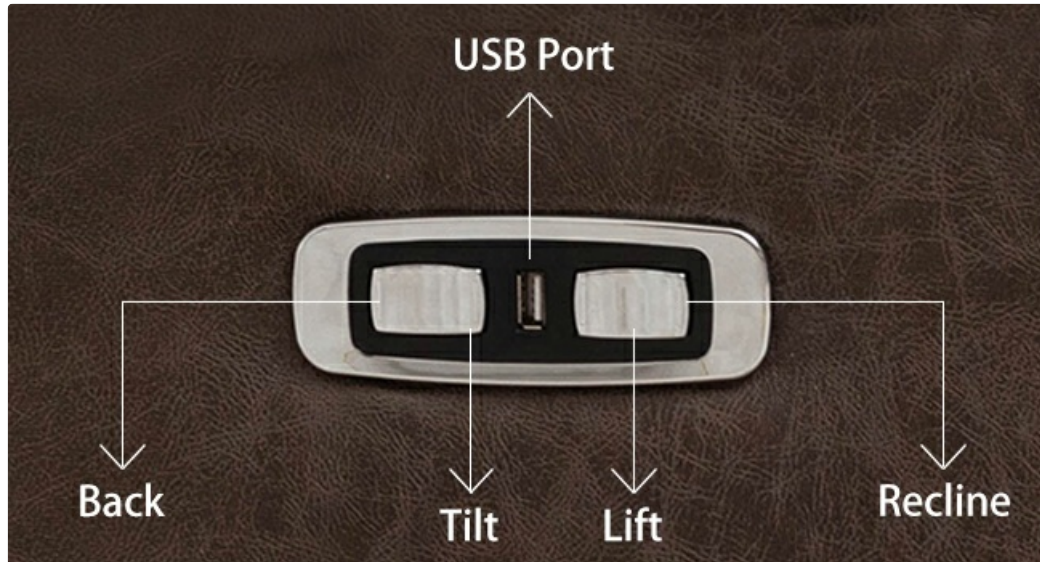


Figure 6: Illustrates various functional positions: Reading, Watching TV, Napping, and Stand-up assist.

### Massage and Heat Functions

The chair is equipped with a remote control for massage and lumbar heat functions.



# MASSAGE & HEATING



Figure 7: Overview of massage and heating features with remote control.

- **Massage Points:** The chair has 4 powerful massage parts (back, lumbar, thighs, legs) with 8 vibrating nodes.
- **Massage Modes:** Select from 5 different massage modes using the remote control. Each massage point can be operated individually.
- **Lumbar Heat:** Activate the lumbar heat function via the remote control to promote blood circulation.
- **Timer Function:** Set the massage duration for 15, 30, or 60 minutes using the remote control.

## USB Charging Ports and Cup Holders

- **USB Ports:** Located on the side control panel, the USB outlet allows charging of low-power devices like smartphones and tablets.
- **Concealed Cup Holders:** Two concealed cup holders are integrated into the armrests. Push to open and close.
- **Side Pockets:** Convenient side pockets are available for storing small items like remote controls or magazines.





Figure 8: Concealed cup holder in the armrest.



Figure 9: Integrated USB port, control buttons, and side storage pocket.



## CARE AND MAINTENANCE

- **Cleaning Faux Leather:** Wipe the faux leather surface with a soft, damp cloth. For stubborn stains, use a mild soap solution and wipe dry immediately. Avoid harsh chemicals or abrasive cleaners.
- **Spills:** Blot spills immediately with a clean, dry cloth. Do not rub.
- **Dusting:** Regularly dust the chair with a soft, dry cloth.
- **Mechanism:** Periodically check all bolts and connections to ensure they are tight. Do not lubricate the electrical components.
- **Power Cord:** Inspect the power cord for any signs of damage. If damaged, discontinue use and contact customer support.



Figure 10: Detail of the faux leather material and stitching.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Chair does not operate (no power).	Power cord unplugged or loose connection. No power at outlet.	Check if the power cord is securely plugged into the chair and the wall outlet. Test the outlet with another device.
Lift or recline function is not smooth.	Obstruction in mechanism. Loose connections.	Ensure there are no objects obstructing the chair's mechanism. Check all assembly connections.
Massage or heat not working.	Remote control battery low or not connected. Power issue.	Ensure the remote control is properly connected and functioning. Verify the main power to the chair.
USB port not charging devices.	Device requires too much power. Loose connection.	Ensure the device is low-power (e.g., smartphone, tablet). Check the USB cable and connection.

## PRODUCT SPECIFICATIONS

**Brand:** CANMOV

**Model:** CD0276

**Color:** Smoke Gray

**Material:** Faux Leather, Foam Fill

**Product Dimensions (Overall):** 32"D x 41"W x 41"H (81.3 cm D x 104.1 cm W x 104.1 cm H)

**Seat Depth:** 24 inches (61 cm)

**Backrest Length:** 25 inches (63.5 cm)

**Weight Capacity:** 350 lbs (158 kg)

**Item Weight:** 142 pounds (64.4 kg)

**Reclining Positions:** 3 (adjustable from 110° to 150°)

**Special Features:** Power Lift, Dual Motor, Adjustable Headrest, 8-Point Vibration Massage, Lumbar Heat, USB Charging Port, 2 Concealed Cup Holders, Side Pockets.

**Power Source:** Electric

**UPC:** 768795929530

## Product Dimensions



Figure 11: Detailed product dimensions and weight capacity.

## WARRANTY AND CUSTOMER SUPPORT



**Warranty Information:** The manufacturer typically offers a 1-year warranty on this product. Please retain your proof of purchase for warranty claims.

**Customer Support:** For any questions, concerns, or assistance with assembly, operation, or warranty claims, please contact CANMOV customer support. The manufacturer aims to respond to all inquiries within 24 hours.

Visit the official CANMOV store for more information:[CANMOV Store](#)

