

VITCOCO DQ-FRs3

VITCOCO Visual Ear Wax Removal Tool

Model: DQ-FRs3

IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using this device. Failure to follow safety guidelines may result in injury or damage to the device.

- Do not insert the device too deeply into the ear canal.
- Always use the ear protection cover when operating the tweezers.
- Supervise children and pets during use.
- Do not use if you have an ear infection, ear drum perforation, or any other ear condition. Consult a medical professional if unsure.
- Clean the ear spoons and lenses thoroughly after each use with water or an alcohol cotton pad. Do not immerse the entire device in water.
- Keep out of reach of children.

PACKAGE CONTENTS

- VITCOCO Visual Ear Wax Removal Tool (DQ-FRs3)
- Silicone Ear Spoons (various shapes)
- Ear Protection Cover
- USB Charging Cable
- User Manual

PRODUCT OVERVIEW

The VITCOCO Visual Ear Wax Removal Tool is an advanced ear cleaning device featuring an integrated HD camera and dual-mode functionality (tweezer and ear pick). It provides a clear view of the ear canal, enabling precise and safe ear wax removal.



Image: The VITCOCO Visual Ear Wax Removal Tool, showing the main device, various ear spoons, a USB charging cable, and a smartphone displaying the ear canal view.

Key Features:

- **Upgraded Tweezer Mode:** Designed for sensitive ear canals with a 360° ring button for easy operation and precise control (0.05mm error).
- **3-in-1 Visual Otoscope:** Combines ear tweezers, ear pick, and a visual camera for comprehensive ear canal inspection and cleaning.
- **HD In-Ear Camera:** Features a 360° wide-angle lens, 3.2mm diameter, 6-axis gyroscope for stable image transmission, and 1080p resolution.
- **Wide Compatibility:** Compatible with Android and iOS devices via Wi-Fi connection.
- **Long Battery Life:** 60-minute charging time, 70 minutes of use, and up to 30 days standby time.

SETUP GUIDE

Step 1: Charging the Device

Before first use, fully charge the device using the provided USB charging cable. Connect the cable to the device and a compatible USB power source. The indicator light will show charging status.

Step 2: App Installation

Download the companion application to your smartphone or tablet. Search for "Soulear" in the Apple App Store (for iOS devices) or Google Play Store (for Android devices).

Alternatively, scan the QR code below to directly access the app download page:



Image: A visual guide showing the four steps: 1. Search for 'Soulear' app, 2. Long press button to turn on device, 3. Connect to device's Wi-Fi, 4. Return to app to start using.

For direct download, visit: <https://qr19.cn/A0kiJ4>

Step 3: Connecting to the Device

1. Long press the power button on the device to turn it on. The indicator light will flash.
2. On your smartphone/tablet, go to Wi-Fi settings.
3. Locate and connect to the Wi-Fi network named "Soulear-xxxx" (where xxxx is a unique identifier).
4. Return to the "Soulear" app. The device should now be connected, and you will see a live view from the camera.

OPERATING INSTRUCTIONS

Basic Operation

- **Power On/Off:** Long press the power button on the device. The indicator light will illuminate when on.
- **Mode Switching:** The device supports both Tweezer and Ear Pick modes. Ensure the correct attachment is securely fitted before use.

Ear wax cleaning mode



Tweezers Mode



Ear Pick Mode



Image: Illustration of the ear wax cleaning modes, demonstrating how the tweezers and ear pick are used within the ear canal.

Using the Tweezer Mode

Attach the tweezer head to the device. Ensure the ear protection cover is in place. Use the 360° ring switch to control the opening and closing of the tweezers. Carefully guide the tweezers into the ear canal while observing the live feed on your connected device. Gently grasp and remove larger pieces of earwax or foreign objects. Avoid excessive force.

High Precision of Tweezers Otosopes



Image: A detailed view of the high-precision tweezers and a user demonstrating its use with a smartphone for visual guidance.

Using the Ear Pick Mode

Attach a suitable silicone ear spoon to the device. Carefully insert the ear pick into the ear canal, using the live camera feed for guidance. Gently scoop out earwax. Do not scrape or push wax deeper into the ear canal.


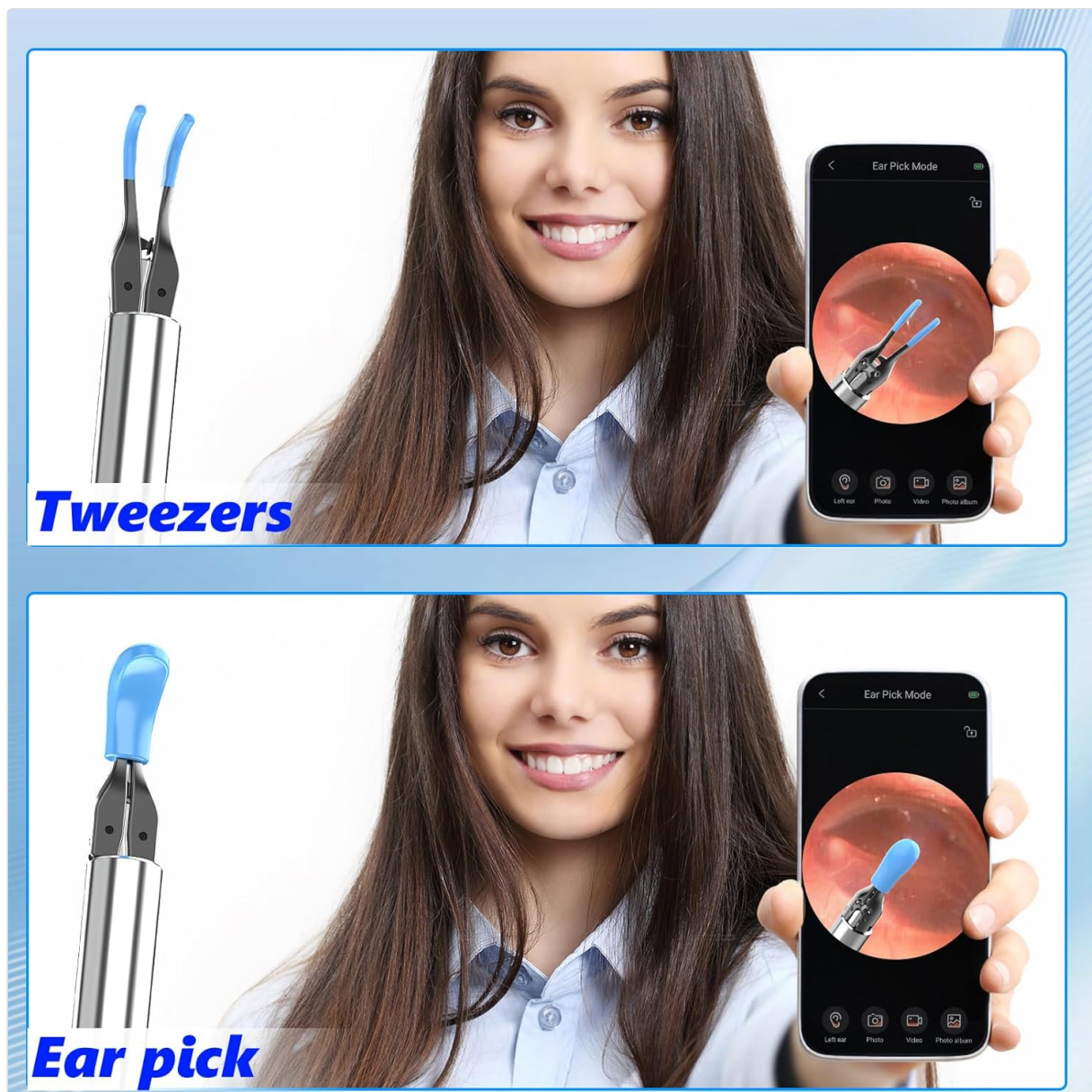


Image: A visual comparison between the tweezer attachment and the ear pick attachment, highlighting their distinct designs.

Capturing Images and Videos

While using the app, you can capture still images or record videos of the ear canal. Look for the camera and video icons within the app interface. These functions allow you to document the condition of the ear or track progress over time.

High Definition Resolution Provides Clear Images

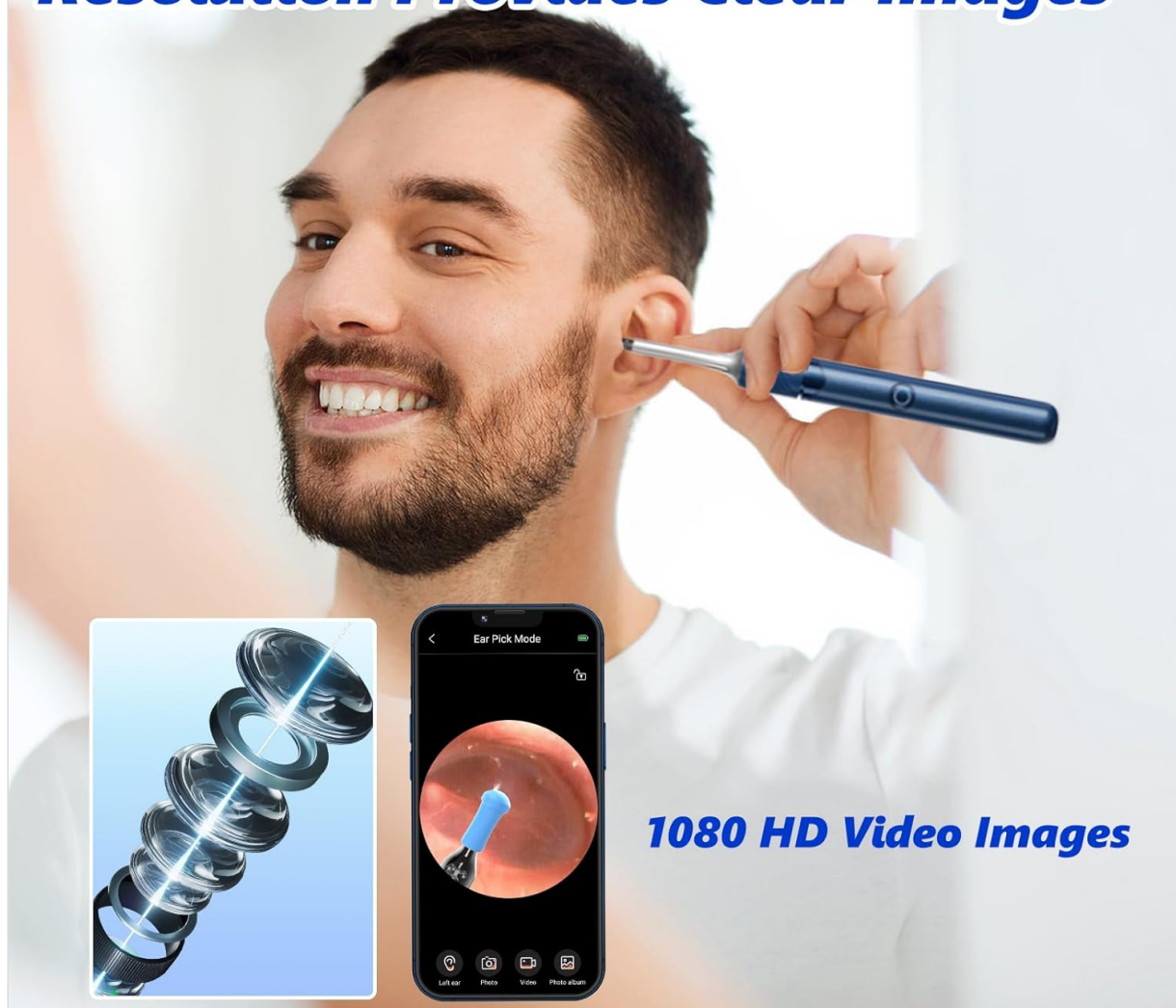


1080 HD Video Images

Image: A person using the device, emphasizing the high-definition camera's ability to provide clear images of the ear canal on a connected smartphone.

Wide Applications

Beyond ear wax removal, the visual otoscope can be used for inspecting other areas:

- **Ears:** For ear wax removal and general ear canal inspection.
- **Nose:** To inspect nasal passages.
- **Mouth:** For oral cavity examination.
- **Pets:** Suitable for inspecting the ears of pets.

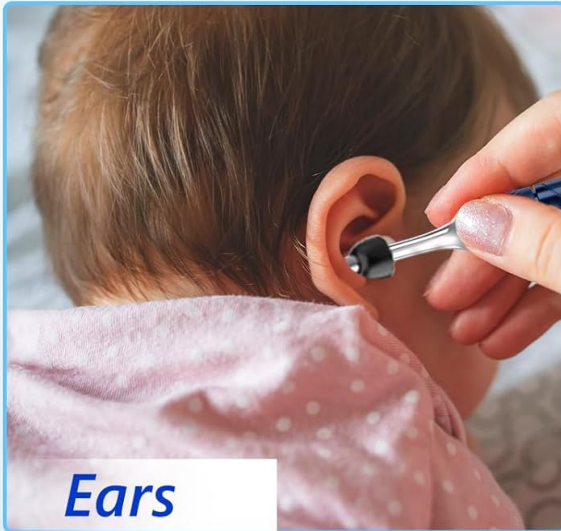

Ears



Nose

Wide applications



Mouth



Pets

Image: Demonstrates the versatility of the tool for inspecting ears, nose, mouth, and even pets' ears.

MAINTENANCE AND CLEANING

Proper maintenance ensures the longevity and hygiene of your device.

- After each use, remove the ear spoon or tweezer attachment.
- Rinse the ear spoons and the lens tip with water or wipe them thoroughly with an alcohol cotton pad.
- **DO NOT** immerse the entire device or the main lens body in water.
- Store the device in a clean, dry place away from direct sunlight.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not power on.	Low battery.	Charge the device fully using the provided USB cable.

Problem	Possible Cause	Solution
Cannot connect to Wi-Fi.	Device not turned on or Wi-Fi not enabled. Incorrect Wi-Fi selected.	Ensure the device is powered on. Go to your phone's Wi-Fi settings and select the "Soulear-xxxx" network. Restart both the device and your phone if necessary.
No image in the app.	Device not connected to app. App permissions not granted.	Verify Wi-Fi connection. Check app permissions (camera, storage) in your phone's settings. Restart the app.
Image is blurry or unclear.	Lens is dirty.	Clean the camera lens gently with an alcohol cotton pad.

SPECIFICATIONS

Feature	Detail
Model	DQ-FRs3
Brand	VITCOCO
Camera Resolution	1080p HD
Lens Diameter	3.2 mm
Gyroscope	6-axis
Connectivity	Wi-Fi (2.4GHz)
Compatibility	iOS, Android
Charging Time	Approx. 60 minutes
Usage Time	Approx. 70 minutes
Standby Time	Up to 30 days
Dimensions	18.6 x 7.9 x 1.9 cm
Weight	80 g

WARRANTY AND SUPPORT

This VITCOCO product is covered by a standard manufacturer's warranty against defects in materials and workmanship. Please refer to the product packaging or the official VITCOCO website for specific warranty terms and duration.

For technical support, troubleshooting assistance, or warranty claims, please contact VITCOCO customer service through their official channels. Contact information can typically be found on the product packaging or the brand's official website.

Related Documents - DQ-FRs3

	<p>VITCOCO S3 Intelligent Visual Ear Picker User Manual and Guide</p> <p>Comprehensive user manual for the VITCOCO S3 Intelligent Visual Ear Picker (Otoscope). Learn how to set up the Endscope Tool app, connect via Wi-Fi, use the device features, and troubleshoot common issues. Includes product specifications, precautions, and contact information.</p>
	<p>CASIO PQ-10, 11, 12, DQ-541, 581, 582, 712 User Manual</p> <p>User manual for CASIO PQ-10, 11, 12, DQ-541, 581, 582, 712 alarm watches. Learn how to set time, use the alarm and snooze functions, replace batteries, and view specifications.</p>
	<p>БУДИЛЬНИК CASIO PQ-10, 11, 12, DQ-541, 581, 582, 712 РУКОВОДСТВО ПОЛЬЗОВАТЕЛЯ</p> <p>Руководство пользователя для электронных будильников CASIO моделей PQ-10, PQ-11, PQ-12, DQ-541, DQ-581, DQ-582 и DQ-712. Описание настройки времени, использования и установки будильника, функции Snooze, подсветки дисплея, замены батареи и спецификаций.</p>
	<p>Casio PQ/DQ Series Alarm Clocks: Quick Start Guide</p> <p>Concise user guide for Casio PQ-10, PQ-11, PQ-12, DQ-541, DQ-581, DQ-582, DQ-712 alarm clocks. Learn how to set time, alarms, and use snooze functions.</p>
	<p>Будильник Casio PQ-10, 11, 12, DQ-541, 581, 582, 712: Руководство пользователя</p> <p>Полное руководство пользователя для будильников Casio моделей PQ-10, PQ-11, PQ-12, DQ-541, DQ-581, DQ-582, DQ-712. Инструкции по установке времени, использованию будильника, подсветке дисплея, замене батареи и техническим характеристикам.</p>