Manuals+

Q & A | Deep Search | Upload

VITCOCO VITCOCO Rechargeable Ear Wax Removal Kit

VITCOCO Rechargeable Ear Wax Removal Kit User Manual

Model: B0DFH7BV6P | Brand: VITCOCO

1. Introduction

The VITCOCO Rechargeable Ear Wax Removal Kit is designed to provide a safe, effective, and convenient solution for ear hygiene. This device utilizes a gentle water irrigation system with multiple cleaning modes to help remove earwax buildup, promoting better ear health and hearing clarity. Its rechargeable battery and waterproof design ensure ease of use and maintenance.

This manual provides comprehensive instructions for the proper setup, operation, maintenance, and troubleshooting of your VITCOCO Ear Wax Removal Kit. Please read this manual thoroughly before first use to ensure safe and optimal performance.



Image: The VITCOCO Ear Wax Removal Kit, including the main device, ear basin, various ear tips, and cleaning tools.

2. SAFETY INFORMATION

Please observe the following safety precautions to prevent injury or damage to the device:

- Do not use if you have ear infections, perforated eardrums, ear tubes, or any other ear injury or condition. Consult a medical professional if unsure.
- Always use lukewarm water (approximately 37°C ± 2°C / 98.6°F ± 3.6°F). Water that is too hot or too cold can cause dizziness or discomfort.
- Do not insert the nozzle too deeply into the ear canal. The flared design of the nozzles is intended to prevent overinsertion.
- If you experience pain, dizziness, or any unusual sensation during use, stop immediately and consult a medical professional.
- Keep out of reach of children. Adult supervision is required if used by or on children.
- Ensure the charging port is dry before charging.
- This device is for personal use and should not be shared without proper sterilization of the ear tips.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is not a medical device and should not replace professional medical advice or treatment.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- Ear cleaner device (1)
- Ear basin (1)
- Three-hole earplugs (5)
- Five-hole earplugs (5)
- Type-C charging cable (1)
- Instruction manual (1)
- Ear spoons (7)
- Ear spoon brush (1)



Image: The comprehensive ear cleaning spoon set and the main device, illustrating the various accessories included in the kit.

4. Product Overview and Features

The VITCOCO Ear Wax Removal Kit is equipped with advanced features for effective and comfortable ear cleaning:

- 4 Ear Cleaning Modes: Choose from Normal, Gentle, Pulse, and Powerful modes to suit your needs.
- Safe and Effective Nozzles: Includes 5 five-hole nozzles for a gentler, wider rinse and 5 three-hole nozzles for a more concentrated stream. Flared design prevents over-insertion.
- Long Battery Life: 2000 mAh battery provides up to 12 hours of continuous use in mild mode.
- IPX7 Waterproof: The entire housing is waterproof, allowing for safe use in the shower and easy cleaning.
- 200ml Water Tank: Sufficient capacity for a thorough cleaning session without frequent refills.
- Automatic Shut-off: Device switches off automatically after 10 minutes of inactivity to conserve battery.



Image: Detailed view of the device's control panel, highlighting the power button, mode button, and indicators for the four distinct cleaning modes: Normal, Soft, Pulse, and Strong.

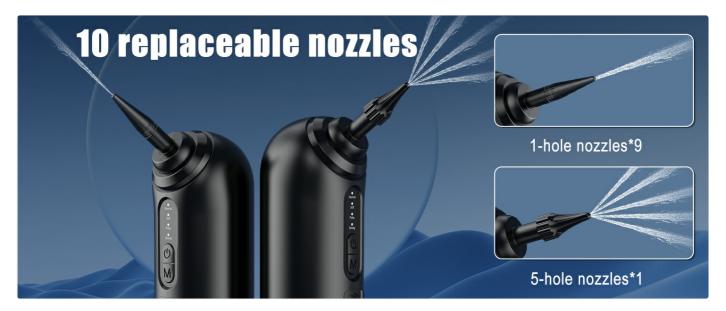


Image: Illustration of the two types of replaceable nozzles included: the 1-hole nozzle for concentrated stream and the 5-hole nozzle for a broader, gentler spray.

5. SETUP INSTRUCTIONS

Follow these steps to prepare your device for use:

- 1. **Charge the Device:** Before first use, fully charge the device using the provided Type-C charging cable. The red LED will flash when the battery is low.
- 2. **Check Waterproof Insert:** Ensure the waterproof insert for the charging port is securely in place before filling with water.
- 3. **Prepare Water:** Disassemble the water tank from the main unit. Fill the water tank with lukewarm water (approximately 37°C ± 2°C / 98.6°F ± 3.6°F). Do not use hot or cold water.
- 4. **Install Nozzle:** Select your desired ear tip (three-hole or five-hole) and insert it onto the device's nozzle. Rotate the nozzle to ensure it is securely installed.
- 5. Assemble the Unit: Reattach the filled water tank to the main device, ensuring it is securely fastened.



Image: A visual guide demonstrating the six key steps for setting up and using the ear wax removal kit, from checking the waterproof insert to gentle cleaning.

6. OPERATING INSTRUCTIONS

Once the device is set up, follow these steps for ear cleaning:

- 1. Initial Use: Press and hold the power button to turn on the device.
- 2. Select Mode: Short press the mode button (M) to cycle through the four cleaning modes:
 - Normal: Gentle yet effective for everyday cleaning.
 - Soft: Deep cleaning for a refreshed feel, suitable for sensitive ears.
 - Pulse: Alternating wash for a unique and thorough clean.
 - Strong: Powerful wash to tackle stubborn earwax buildup.
- 3. **Positioning:** Hold the ear basin under your ear to catch the water and dislodged earwax. Position the nozzle at the entrance of your ear canal, ensuring it does not go too deep.
- 4. **Gentle Cleaning:** Press the power button to start the water flow. Move the nozzle gently around the ear canal entrance. Do not force the nozzle or water into the ear.
- 5. Pause/Stop: Short press the power button to pause the water flow. Long press to turn off the device.
- 6. **Repeat if Necessary:** If needed, refill the water tank and repeat the process for the same ear or switch to the other ear.

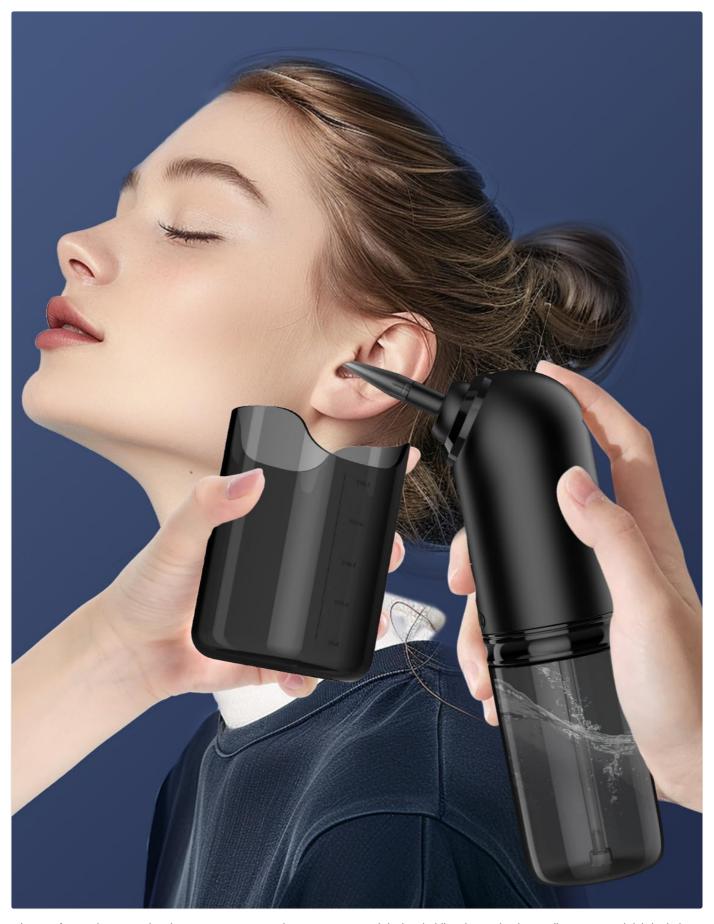


Image: A user demonstrating the correct way to use the ear wax removal device, holding the ear basin to collect water and debris during the cleaning process.

7. MAINTENANCE AND CLEANING

Proper maintenance ensures the longevity and hygiene of your device:

· After Each Use:

- Empty any remaining water from the tank.
- Rinse the water tank and ear tips thoroughly with clean water.
- · Clean the ear basin.
- Use the provided ear spoon brush to clean the ear spoons after use.
- **Device Cleaning:** The main unit is IPX7 waterproof and can be rinsed under running water. Ensure the charging port cover is sealed. Do not submerge the device for extended periods.
- **Drying:** Allow all components to air dry completely before reassembling or storing.
- Storage: Store the device and accessories in a clean, dry place away from direct sunlight and extreme temperatures.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully using the Type-C cable.
No water flow or weak water flow.	Water tank empty or not properly assembled. Nozzle clogged.	Refill water tank and ensure it's securely attached. Check nozzle for blockages and clean if necessary.
Red LED flashing.	Low battery indicator.	Recharge the device.
Device automatically shuts off.	Inactivity for 10 minutes.	This is a normal power-saving feature. Turn the device back on to continue use.

9. SPECIFICATIONS

• Model: VITCOCO Ear Wax Removal Kit (ASIN: B0DFH7BV6P)

• Battery Capacity: 2000 mAh Lithium Ion

• Continuous Use (Mild Mode): Up to 12 hours

• Water Tank Capacity: 200 ml

• Waterproof Rating: IPX7

• Charging Interface: Type-C

• Dimensions: Approximately 6.65 x 6.46 x 3.35 inches (Package Dimensions)

• Weight: Approximately 1.01 Pounds (Package Weight)

• Manufacturer: VITCOCO



Image: Close-up view of the device's transparent water tank, highlighting its 200ml capacity, designed for extended use without frequent refills.



Image: An illustration depicting the internal 2000mAh battery, emphasizing its long-lasting power for up to 5 hours of continuous use on a single charge.

10. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or contact VITCOCO directly through their official channels. Protection plans may be available for purchase separately.

Manufacturer: VITCOCO

Date First Available: July 18, 2025

For further assistance, please visit the VITCOCO Store on Amazon.