



[Manuals.plus](#) /

> [mivo](#) /

> Mivo GM-901T Rechargeable Hearing Aids User Manual

## mivo GM-901T

# Mivo GM-901T Rechargeable Hearing Aids User Manual

Model: GM-901T

## INTRODUCTION

---

This manual provides detailed instructions for the safe and effective use of your Mivo GM-901T Rechargeable Hearing Aids. These devices are designed to enhance sound clarity and reduce background noise, offering an improved listening experience. Please read this manual thoroughly before using your hearing aids and retain it for future reference.

## IMPORTANT SAFETY INFORMATION

---

Always read and follow all safety instructions. Failure to do so may result in injury or damage to the device.

- Consult a medical professional if you experience sudden hearing loss, dizziness, or ear pain.
- Keep the device out of reach of children and pets to prevent accidental ingestion of small parts.
- Do not expose the hearing aids to extreme temperatures, humidity, or direct sunlight.
- Avoid dropping the device or subjecting it to strong impacts.
- Do not attempt to repair or modify the hearing aids yourself. Contact customer support for assistance.
- Remove hearing aids before showering, swimming, or using hairspray/perfume.

## WHAT'S IN THE BOX

---

Carefully unpack your Mivo GM-901T hearing aids and ensure all items are present:

- Mivo GM-901T Hearing Aids (Pair)
- Rechargeable Charging Case
- USB Charging Cable
- Multiple Sizes of Silicone Ear Tips
- Cleaning Brush

- Small Screwdriver (for adjustments, if applicable)
- User Manual (this document)



Image: The Mivo GM-901T hearing aids, showing their compact size and the accompanying charging case.

## PRODUCT FEATURES

---

- **Exceptional Clarity:** Designed to enhance speech and everyday sounds for improved communication.
- **User-Friendly Design:** Simple controls and comfortable fit for ease of use.
- **Advanced Noise Reduction:** Automatic noise reduction feature to minimize background distractions.
- **Rechargeable Convenience:** Long-lasting performance with a quick 2-hour charge time.
- **Discreet and Stylish:** Compact design that blends seamlessly for confident wear.
- **16-Channel Sound Processing:** Advanced digital chip technology for precise sound management.
- **Nano Titanium Diaphragm:** Filters noise and restores sound quality for clearer audio.



Image: A visual representation of the 16-channel sound processing chip, highlighting its advanced digital technology.

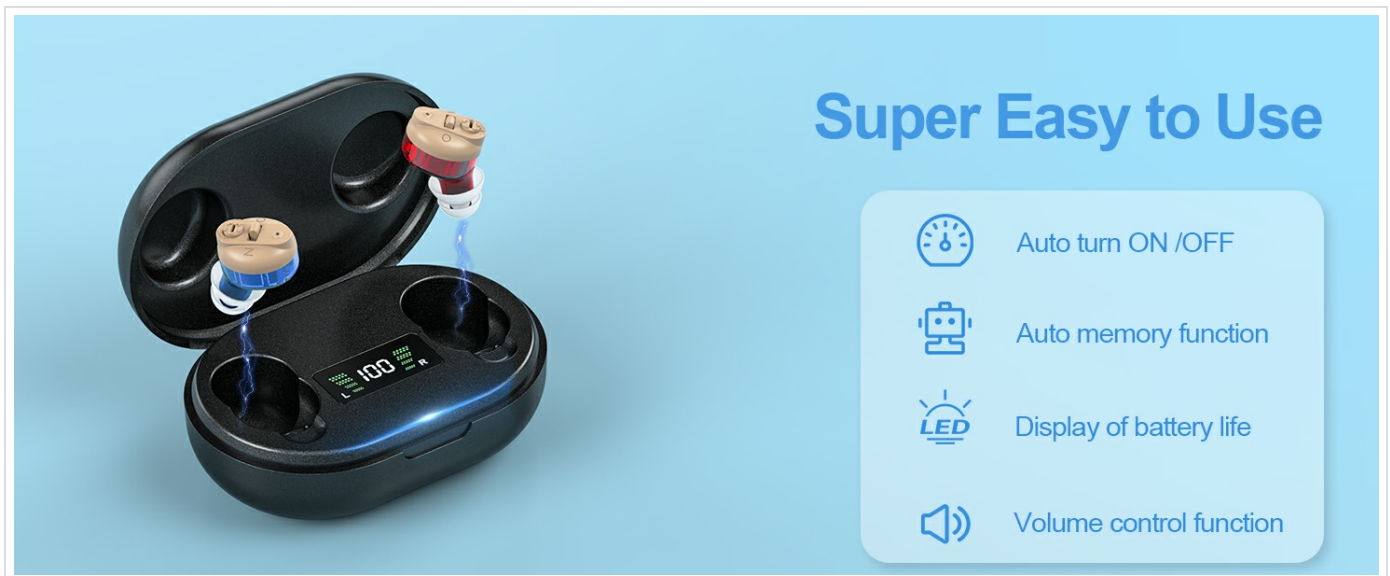


Image: An exploded view of the hearing aid's internal speaker components, illustrating the nano titanium diaphragm for sound quality restoration.

## SETUP

### 1. Charging the Hearing Aids

Before first use, fully charge your hearing aids and the charging case.

1. Place the hearing aids into their respective slots in the charging case. Ensure they are correctly seated.
2. Connect the USB charging cable to the charging case and a suitable power adapter (not included).
3. The digital display on the charging case will indicate the charging status. A full charge typically takes about 2 hours.
4. Once fully charged, the hearing aids provide approximately 30 hours of single running time, and the case holds enough power for 4 full recharges, offering 5-6 days of total use.

# Recharge Anytime, Anywhere

## Compact and Convenient

**2 hours**  
Charging time

**30 hours**  
Single running time

**5-6 days**  
Total running time

**4 recharges**  
Charging compartment



Image: The Mivo hearing aids placed in their charging case, with the digital display showing battery charge levels for both the case and individual hearing aids.

## 2. Selecting and Fitting Ear Tips

Choosing the correct ear tip size is crucial for comfort and optimal sound quality.

1. Experiment with the different sizes of silicone ear tips provided.
2. Select a tip that fits snugly in your ear canal without causing discomfort. A good fit helps prevent feedback (whistling) and ensures sound is directed effectively.
3. Gently attach the chosen ear tip to the hearing aid.

## The Perfect Gift for Better Hearing

Start your hearing journey today—for yourself, a parent, or a friend. Every moment counts!



Image: A senior man wearing a Mivo hearing aid, with an illustration of various ear tip sizes to ensure a comfortable and effective fit for different users.

### 3. Inserting the Hearing Aids

1. Hold the hearing aid with the ear tip pointing towards your ear canal.
2. Gently insert the hearing aid into your ear canal until it feels secure and comfortable.
3. Ensure the hearing aid is flush with your outer ear.

# Super Comfortable Sound Experience

8 sizes, suitable for all people, more comfortable



Image: A close-up view demonstrating the correct insertion of a Mivo hearing aid into the ear, showing its discreet placement.

## OPERATING INSTRUCTIONS

---



Image: A detailed diagram illustrating the key components and controls of the Mivo hearing aid, including the power switch, volume adjustment, microphone, and ear tip.

## 1. Power On/Off

- **Power On:** The hearing aids automatically turn on when removed from the charging case.
- **Power Off:** Place the hearing aids back into the charging case to automatically turn them off and begin charging.
- Alternatively, use the designated Power Switch on the device for manual control.

## 2. Volume Adjustment

Adjust the volume to a comfortable level using the "Turning Volume" control located on the hearing aid.

- Rotate the volume wheel gently to increase or decrease the sound amplification.
- Start with a lower volume and gradually increase it until sounds are clear and comfortable, avoiding excessively loud settings.

## 3. Noise Reduction

The Mivo GM-901T features automatic noise reduction to enhance speech clarity in various environments.

- The device is designed to automatically detect and minimize background noise, allowing you to focus on conversations.
- This feature is particularly effective in environments with consistent background sounds.


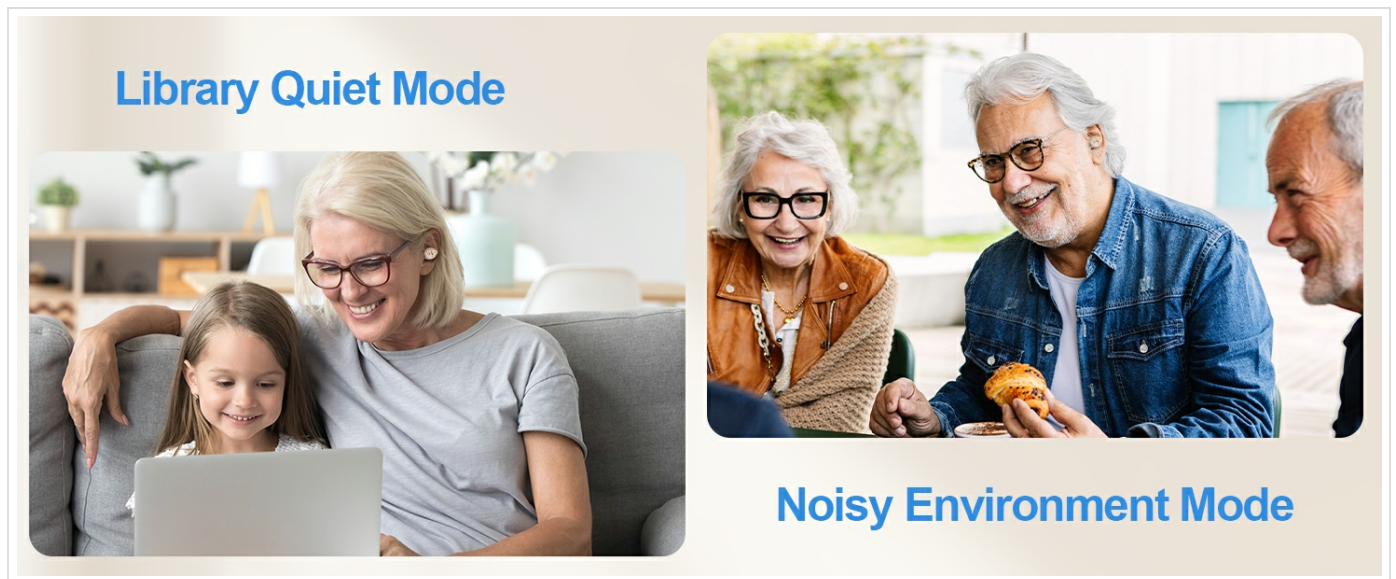


Image: An illustration depicting a social gathering, demonstrating the effect of automatic noise cancellation in reducing background distractions for clearer hearing.

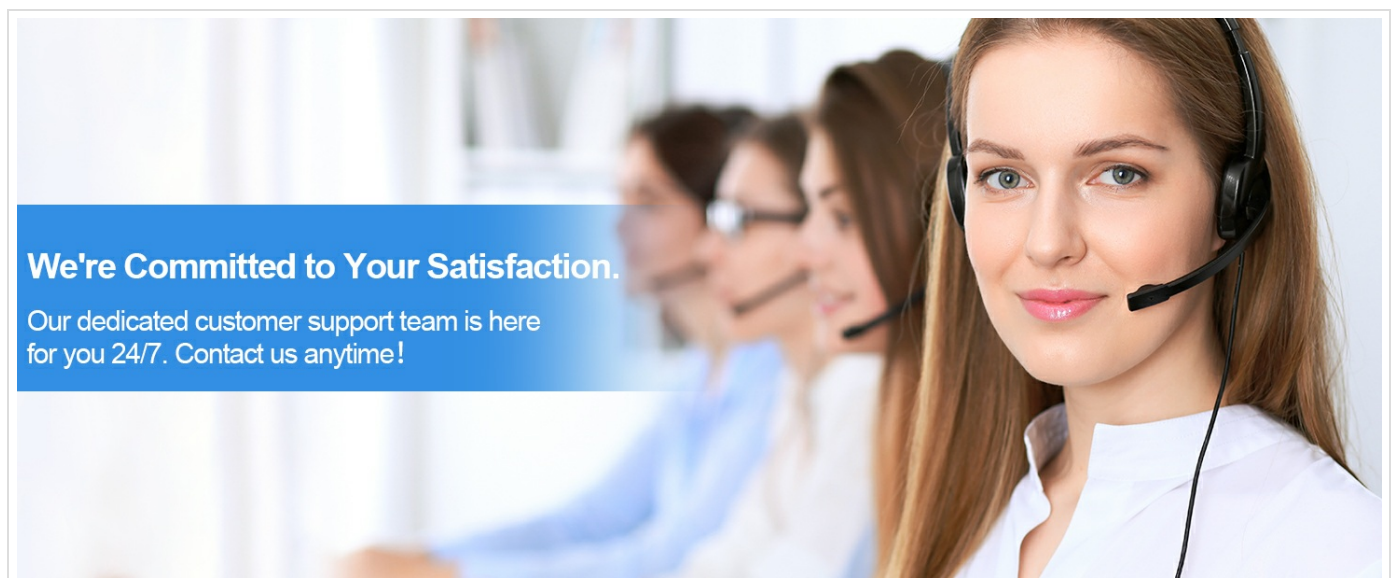

## We're Committed to Your Satisfaction.

Our dedicated customer support team is here for you 24/7. Contact us anytime!

Image: Two comparative scenes illustrating the hearing aids' performance in different acoustic environments: a quiet library setting and a noisy outdoor gathering.

## MAINTENANCE

---

### 1. Cleaning Your Hearing Aids

Regular cleaning is essential for optimal performance and hygiene.

- Use the provided cleaning brush to gently remove any earwax or debris from the ear tips and sound ports daily.
- Wipe the body of the hearing aids with a soft, dry cloth. Do not use water or cleaning solutions directly on the device.
- Ensure the charging contacts on both the hearing aids and the charging case are clean and free of debris.

## 2. Storage

- When not in use, store your hearing aids in the charging case to protect them and keep them charged.
- Store the charging case in a cool, dry place, away from direct sunlight and extreme temperatures.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
No sound or weak sound	Low battery Ear tip blocked by earwax Volume too low Improper fit	Charge the hearing aids. Clean the ear tip and sound port. Increase the volume. Adjust the hearing aid for a secure fit, try a different ear tip size.
Whistling or feedback noise	Improper fit or loose ear tip Volume too high Earwax buildup in ear canal	Re-insert the hearing aid, ensure a snug fit. Try a different ear tip size. Lower the volume. Consult a doctor for earwax removal if necessary.
Distorted sound or static	Low battery Moisture in the device Device malfunction	Charge the hearing aids. Ensure the device is dry. Contact customer support.
Hearing aid does not charge	Charging cable or adapter issue Dirty charging contacts Charging case low on power	Try a different USB cable or power adapter. Clean the charging contacts on both the hearing aids and the case. Charge the charging case.

## SPECIFICATIONS

- **Model Name:** GM-901T
- **Brand:** Mivo
- **Style:** In-Ear
- **Material:** Plastic
- **Color:** Black
- **Item Weight:** 5.6 ounces (total package)
- **Package Dimensions:** 4.92 x 3.78 x 1.77 inches
- **Batteries:** 1 Lithium Ion battery required (included)
- **Charging Time:** Approximately 2 hours
- **Single Running Time:** Approximately 30 hours
- **Total Running Time (with case recharges):** 5-6 days (up to 150 hours)
- **Manufacturer:** HONG YU KANG TECHNOLOGY (GUANGDONG) CO., LTD.

## WARRANTY AND SUPPORT

Your Mivo GM-901T Hearing Aids come with a **2-year warranty** from the date of purchase, covering manufacturing

defects.

For technical assistance, troubleshooting, or warranty claims, please contact Mivo customer support. Our dedicated team is available to assist you.

Contact information can typically be found on the product packaging or the official Mivo website.

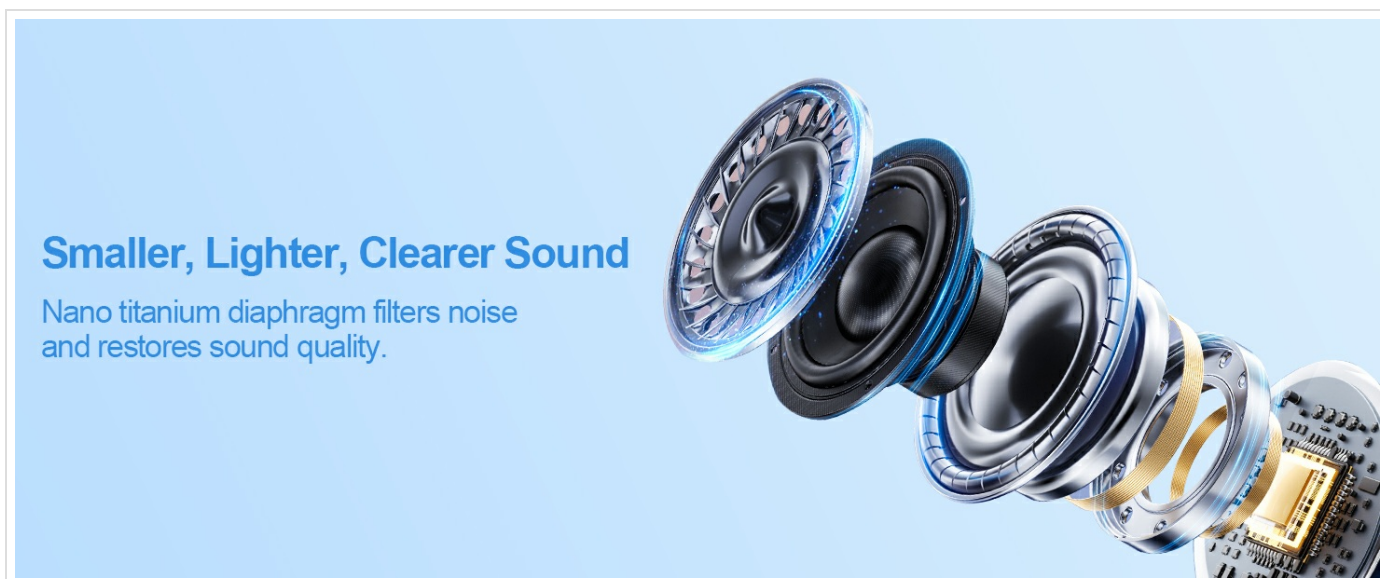


Image: A senior man wearing a Mivo hearing aid, illustrating the lightweight and comfortable design for extended wear.